

Instructions for use

Perseus A500



WARNING

To properly use this medical device, read and comply with these instructions for use.

Anesthesia workstation Software 2.0n

This page has been left blank intentionally.

Contents

Information about this document	5	Getting started	79
Applicability of these instructions for use	5	Turning on Perseus	80
Typographical conventions	5	Checking the device configuration	
Use of terms	5	Checking the operational readiness	
Illustrations	5	Emergency start-up	
Further documents Trademarks	6 6	Operation	
Safety information definitions	6	Safety information	
User group requirements	7	Starting the therapy	
For your safety and that of your patients	8	Adjusting the therapy	
		Special forms of therapy	
General safety information	9	Maneuvers (optional)	
Product-specific safety information	12	Using fields with special functions	
Application	14	Organizing the screen display	
Intended use	15	Displaying additional data	
Indications/Contraindications	15	Setting the sound volume	
Further information on application	16	Adjusting the alarms	
r drifter information on application	10	Changing the patient data	
Overview	17	Exporting data	
Hardware	18	Other settings	
Functional scope	33	Ending the therapy	
Gas flow plan	36	Change of patient	110
Operating concept	39	Tests	120
Screen	40	Information on the availability of individual	
	40	device functions	
Color concept		Testing the device	
Selecting and setting Remote control for the ceiling-mounted version	43	Performing the tests	
(combination with Dräger ceiling supply units)	44	Completing the checklist	126
Safety sensor	45	Ending operation	135
Assembly and preparation	46	Turning off the device	136
Perseus as a ceiling-mounted version	47	Alarms	138
Mounting of accessories	49	Displaying alarms	139
Before first operation	52	Response to alarms	140
Intrahospital transport	56	Adopting alarm settings when changing the	
Establishing the gas supply	57	ventilation mode	143
Ensuring the gas scavenging	62	Alarm delay and alarm escalation	
Preparation for a day of operations / after		Activation of alarms after breath detection	
reprocessing	63	Intelligent alarm behavior	
Selecting and connecting patient-specific		5	• •
accessories	67		
Connecting or exchanging consumables	76		

Configuration	148	Ven
Device settings	149	Brea
Setting the date and time	149	Exte
Defining the start settings	150	Ane
Transferring device configurations	171	Mea
Activating software options	172	Disp
Overview of configurable screen contents	173	Ope Dev
Troubleshooting	180	Rele
Leakage	181	Diag
Power supply failure	182	Dev
Emergency O2 delivery (electronically		EM
controlled gas mixer)	183	Emi
Failure of the gas supply	184	Con
Failure of fresh-gas delivery (electronically		Оре
controlled gas mixture)	186	
Ventilator failure	186	Prir
Gas measurement failure	187	Des
Flow measurement failure	187	Min
Screen fault/user interface failure	188	Influ
Complete failure	188	on c
Problems with the anesthetic gas receiving		Infir
system (AGS)	189	Sch
Problems with cylinder pressure reducers	190	Con
Problems with the Vapor View option	190	Anr
Support request	190	
Alarm – Cause – Remedy	191	Abb Syn
Service	214	Pro
Safety information	215	Ove
Transport for maintenance purposes	216	_
Definition of service terminology	216	Pas
Inspection	217	Con
Maintenance	219	A50
Repair	220	Inde
Disposal	221	
Safety information		
Disposing of the device		
Disposing of accessories	222	
Technical data	223	
General information	224	
Ambient conditions	224	
Fresh-gas delivery		
Fresh-gas delivery (electronically controlled		
gas mixer)	225	
Fresh-gas delivery (mechanically controlled		
gas mixer)	226	

Ventilator	231 231 237 238 244 247 249 250 251
Principles of operation	257
Description of the ventilation modes	258268269
Infinity ID accessory support	269
Schematic diagram of alarm tones	271
Combination lock on the drawer (optional)	272
Annex	274
-,	275 277 280 282
Password	285
Configuration password for Perseus A500 Software 2.0n	285
Index	287

Information about this document

Applicability of these instructions for use

This document applies to products with software version 2.03 or higher.

Typographical conventions

- 1 Consecutive numbers indicate steps of action, with the numbering restarting with "1" for each new sequence of actions.
- Bullet points indicate individual actions or different options for action.
- Dashes indicate the listing of data, options, or objects.
- (A) Letters in parentheses refer to elements in the related illustration

- A Letters in illustrations denote elements referred to in the text.
- > The greater-than symbol indicates the navigation path in a dialog window.

Bold, italicized text indicates labels on the device and texts that are displayed on the screen.

Use of terms

Dräger uses the term "accessories" not only for accessories in the sense of IEC 60601-1, but also for consumables, removable parts, and attached parts.

The product "Perseus A500" is also referred to as "Perseus".

Illustrations

Illustrations of products and screen content in this document may differ from the actual products depending on configuration and design.

Further documents

The reprocessing of this product is described in the separate reprocessing instructions supplied with the product.

Trademarks

Trademarks owned by Dräger

Trademark
Perseus [®]
Infinity [®]
D-Vapor [®]
Drägersorb [®]
MEDIBUS [®]
ServiceConnect [®]
WaterLock [®]

The following web page provides a list of the countries in which the trademarks are registered: www.draeger.com/trademarks

Trademarks owned by third-party manufacturers

Trademark	Trademark owner
Selectatec [®]	Datex-Ohmeda

Trademark used under license	Trademark owner
BIPAP	Respironics

Safety information definitions

WARNING

A WARNING statement provides important information about a potentially hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

A CAUTION statement provides important information about a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient or in damage to the medical device or other property.

NOTE

A NOTE provides additional information intended to avoid inconvenience during operation.

User group requirements

The term "user group" describes the personnel responsible who have been assigned by the operating organization to perform a particular task on a product.

Where product specific knowledge or tools are required, the service activities must be carried out by specialized service personnel. The specialized service personnel was trained by Dräger for these service activities on this product.

Duties of the operating organization

The operating organization must ensure the following:

- Every user group has the required qualifications (e.g., has undergone specialist training or acquired specialist knowledge through experience).
- Every user group has been trained to perform the task.
- Every user group has read and understood the relevant chapters in this document.

User groups

Clinical users

This user group operates the product in accordance with the intended use.

Users have medical specialist knowledge in the field of anesthesia. Users have knowledge of device monitoring and perioperative care.

Reprocessing personnel

This user group carries out the necessary activities to reprocess the product.

Reprocessing personnel has specialist knowledge in the reprocessing of medical devices.

Service personnel

This user group installs the product and performs the service activities.

Service personnel has specialist knowledge in electrical and mechanical engineering and experience in the servicing of medical devices.

For your safety and that of your patients

General safety information	٤
Strictly follow these instructions for use	ç
Service	ć
Accessories	ć
Patient safety	10
Patient monitoring	10
Electromagnetic compatibility (EMC)	11
Training	11
Keeping the instructions for use	11
Product-specific safety information	12
Essential performance	13
Handling Infinity ID components	13

General safety information

The following WARNING and CAUTION statements apply to general operation of the medical device.

WARNING and CAUTION statements specific to its subsystems or particular features appear in the respective sections of these instructions for use or in the instructions for use of any other product being used with this device.

Strictly follow these instructions for use

WARNING

Risk of incorrect operation and incorrect use

Any use of the medical device requires full understanding and strict observation of all sections of these instructions for use. The medical device must only be used for the purpose specified under "Intended use" (see page 15) and in conjunction with appropriate patient monitoring (see page 10). Strictly observe all WARNING and CAUTION statements throughout these instructions for use and all statements on medical device labels.

Failure to observe these safety information statements constitutes a use of the medical device that is inconsistent with its intended use.

Service

WARNING

Risk if service is not performed regularly

If service is not performed regularly, malfunctions may occur, which can result in personal injury and property damage.

Perform the service in accordance with the chapter "Service".

Accessories

WARNING

Risk due to incompatible accessories

The use of incompatible accessories may adversely affect the functional integrity of the product. Personal injury and property damage may occur as a consequence.

Use only compatible accessories. The accessories that are compatible with this product are listed in the list of accessories supplied with the product.

WARNING

Risk due to faulty accessories

The use of faulty accessories may compromise the functional integrity of the product. Personal injury and property damage may result.

Use only accessories that are in good working order.

WARNING

Risk due to misoperation, incorrect use, or incorrect reprocessing

Strictly observe the instructions for use of all accessory parts, e.g.:

- Water traps
- Flow sensors
- CLIC adapter
- CLIC absorber
- Soda lime
- Breathing hoses
- Masks
- Filter
- Bronchial suction
- Vaporizer
- Manual resuscitator
- AGSS terminal unit

Patient safety

The design of the medical device, the accompanying documentation, and the labeling on the medical device are based on the assumption that the purchase and the use of the medical device are restricted to persons familiar with the most important inherent characteristics of the medical device.

Instructions and WARNING and CAUTION statements are therefore largely limited to the specifics of the Dräger medical device.

The instructions for use do not contain any information on the following points:

- Risks that are obvious to users
- Consequences of obvious improper use of the medical device
- Potentially negative effects on patients with different underlying diseases

Medical device modification or misuse can be dangerous.

Patient monitoring

The user of the medical device is responsible for choosing suitable monitoring that provides appropriate information about medical device performance and the patient's condition.

Patient safety may be achieved by a wide variety of means ranging from electronic surveillance of medical device performance and patient condition to simple, direct observation of clinical signs.

The responsibility for selecting the best level of patient monitoring lies solely with the user of the medical device.

Electromagnetic compatibility (EMC)

Medical electrical equipment is subject to special precautionary measures concerning electromagnetic compatibility. During installation and before initial operation, follow the information in section: "EMC declaration" (page 251).

This device can be affected by other electrical devices.

WARNING

Risk due to electrostatic discharge

Malfunctions that endanger the patient may occur if no protective measures against electrostatic discharge are employed in the following situations:

- When touching the pins of connectors that carry the ESD warning symbol.
- When establishing connections with these connectors.

To prevent malfunctions, observe the following measures and train the relevant personnel:

- Observe the ESD protective measures. Such measures may include wearing antistatic clothing and shoes, touching a potential equalization pin before and while making the connection, or using electrically insulating and antistatic gloves.
- Observe the requirements for the electromagnetic environment. Observe the following section: "Electromagnetic environment" (page 251).

WARNING

Risk due to electromagnetic disturbance

Wireless communication devices (e.g., cellular phones) and medical electrical equipment (e.g., defibrillators, electrosurgical devices) emit electromagnetic radiation. When such devices are operated too close to this device or its cables, the functional integrity of this device may be compromised by electromagnetic disturbances. As a result, the patient could be put at risk.

- Maintain a distance of at least 0.3 m (1.0 ft) between this device and wireless communication devices, to ensure that the essential performance of this device is fulfilled.
- Maintain an adequate distance between this device and other medical electrical equipment.

Training

User training is offered by the responsible Dräger organization, see www.draeger.com.

Keeping the instructions for use

WARNING

Risk of operating errors

The instructions for use must be kept in an accessible location for users.

Product-specific safety information

WARNING

Risk due to device failure

Device failure can compromise the correct therapy functionality of the device.

To ensure immediate remedial action in case of device failure, the device may only be operated under permanent supervision of users. Always have a manual resuscitator ready.

WARNING

Dräger recommends that the user remains in the vicinity of the anesthesia machine, i.e. within a distance of up to 4 meters (12 feet). This facilitates fast recognition and response in the event of an alarm.

WARNING

Risk due to modifications

Modifications to the product may lead to malfunctions and unforeseen risks. This may result in injury to the patient or the user or in property damage.

Do not modify this product.

WARNING

Risk of fire

The flow sensor can ignite medications or other substances based on highly flammable substances.

- Do not nebulize medications or other substances that are easily flammable or spray them into the device.
- Do not use substances containing alcohol.
- Do not allow flammable or explosive substances to enter the breathing system or the breathing circuit.
- Do not use cyclopropane or ether.

WARNING

Risk due to failure of flow measurement

Deposits that were not removed during reprocessing can damage the measuring wires in the flow sensor or cause a fire.

- Before inserting the flow sensor check for visible damage, soiling, and particles.
 Repeat this check regularly.
- Replace flow sensors when damaged, soiled, or not particle-free.

CAUTION

Risk of crushing

Movable device parts or attached components may cause crushing due to clamping. Pay special attention to edges, movable parts, and corners when working with the following components:

- Column cover
- Breathing system cover
- Drawers
- Extensible writing tray
- Swivel arms for mounted devices
- Accessories such as gas cylinders, vaporizers, CLIC absorber, and CLIC adapter

WARNING

Risk due to electromagnetic fields

Although the medical device does not exceed the applicable limiting values for electromagnetic fields, such radiation can interfere with the functioning of pacemakers.

All wearers of pacemakers should maintain a distance of at least 25 cm (10 in) between pacemaker and medical device.

Essential performance

General

- Supply of the anesthesia workstation with O2: If the O2 supply (central gas supply or gas cylinder) fails, an alarm is issued.
- Supply of the patient with adequately oxygenated breathing gas:
 If the breathing gas contains insufficient levels of O2, an alarm is issued.
- Patients are not supplied with excessively high anesthetic gas concentrations:
 If excessively high anesthetic gas concentrations are delivered, an alarm is issued.
- Monitoring of the airway pressure:
 Alarms are issued depending on the set alarm limits.

Gas measurement

- Breathing gas monitoring:
 Measurement of the gas composition with ISO accuracy.
- Monitoring of the breathing gas concentrations:
 If a set alarm limit is exceeded or the gas measurement fails, alarms will be issued.

Handling Infinity ID components

Ownership or purchase of this medical device with RFID technology only includes the right to use the medical device and RFID technology in conjunction with products approved by Dräger and in strict compliance with these instructions for use. No intellectual property rights or any rights to the use of the medical device or RFID technology are hereby granted, either explicitly or implicitly, which are contrary to the above-mentioned conditions.

Application

Intended use	15
Indications/Contraindications	15
Indications Contraindications	
Further information on application	16
Environments of use	16

Intended use

The Perseus anesthesia workstation is intended for use in anesthetizing adults, children, and neonates and can be used for automatic and manual ventilation, pressure-supported spontaneous breathing, and spontaneous breathing.

Perseus is equipped with airway monitoring, gas measurement and device monitoring, O2 insufflation, and an anesthetic gas receiving system.

Anesthesia is achieved through a mixture of pure oxygen and Air (medical compressed air) or pure oxygen and nitrous oxide, with the addition of volatile anesthetic agents.

Ventilation is accomplished on the patient through a laryngeal mask airway, a full-face mask, or an endotracheal tube.

The integrated breathing system can be used with partial rebreathing (low-flow or minimum-flow).

A non-rebreathing system such as the Kuhn or Medec Water System may be used at the external fresh-gas outlet (optional).

Indications/Contraindications

Indications

Perseus is specified for inhalational anesthesia and/or patient ventilation in accordance with the intended use during surgical or diagnostic interventions

Contraindications

- Perseus applies medications such as oxygen, nitrous oxide, or volatile anesthetic agents, among others. For contraindications to the applied medications, strictly follow the instructions for use of the medication.
- Do not use soda lime based on potassium hydroxide. Otherwise, there is a risk of CO formation.
- Only use pelletized soda lime. Otherwise, there is a risk of faulty measurement or faulty delivery as well as progressive damage to the breathing system due to dust.
- For patients suspected of malignant hyperthermia: Do not use any volatile anesthetic agent or Perseus with residual concentrations of these gases above 5 ppm.

 Do not perform low-flow anesthesia on patients with ketoacidosis or patients under the influence of alcohol. Otherwise, there is a risk of acetone accumulation in the patient.

The user is responsible for setting the gas delivery and ventilation according to the individual patient status. Patient status must be continually monitored for any potential changes.

Further information on application

Environments of use

Perseus is designed for use in rooms in which therapeutic or diagnostic interventions can be performed under permanent supervision of users.

WARNING

Risk of explosion and fire

This medical device is neither approved nor certified for use in areas where oxygen concentrations greater than 25 Vol%, combustible or explosive gas mixtures are likely to occur.

WARNING

Risk of device malfunctions and/or patient injury and user injury

Magnetic fields can negatively influence the correct functioning of the medical device and therefore endanger the patient or user.

Do not use the medical device in rooms where devices for magnetic field applications are used (e.g., magnetic resonance imaging).

Also, do not use Perseus in the following environments:


- Outside buildings
- On intensive care units
- During patient transport
- In vehicles, airplanes, or helicopters

Overview

Hardware	18
Front (trolley version)Front (ceiling-mounted version)	18 19 19
Plug-in connection with Vapor View option for up to 3 vaporizers	20 21 22 23 24 26 27 28
Auxiliary power sockets (trolley version) Gas mixing unit (electronically controlled) Gas mixing unit (mechanically controlled)	30 31 32
Functional scope	33
Device versions, options, and accessories Gas delivery Ventilation modes Additional operation modes Monitoring Display on the screen Protocoling Gas supply Gas scavenging Data exchange, interfaces	33 33 33 34 34 34 34 35
Gas flow plan	36
Breathing systemGas supply (electronically controlled gas mixer)	36 37
Gas supply (mechanically controlled gas mixer)	38

Hardware

Front (trolley version)

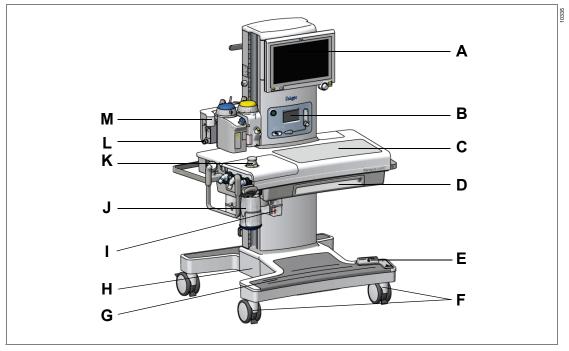


- A Screen with integrated workplace illumination
- **B** Gas mixing unit (electronically controlled gas mixer shown here)
- C Work surface
- D Drawer or pull-out writing tray with lockable compartment (optional)
- E Central brake
- F Castors with central brake
- **G** Footrest
- **H** Trolley
- Anesthetic gas receiving system with flow indicator
- J CO₂ absorber
- K APL valve
- **L** External fresh-gas outlet (optional)

M Plug-in connector with Vapor View option for up to 3 vaporizers (the version for 2 vaporizers is shown here)

Front (ceiling-mounted version)

Differences from the trolley version:

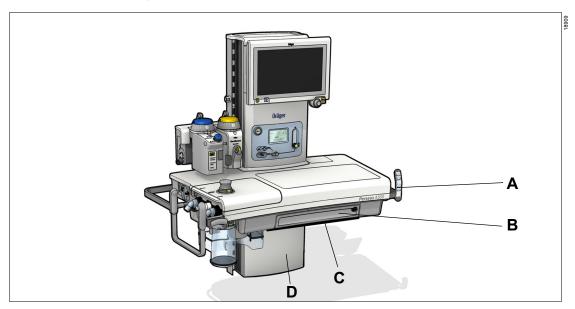


- A Remote control for controlling and positioning the ceiling supply unit (functions dependent on the type of the ceiling supply unit)
- **B** Lockable drawer (optional)
- C Safety sensor

D Device column without trolley

Although the remote control and the safety sensor are mounted on Perseus, from an electrotechnical point of view they are part of the connected ceiling supply unit. The safety sensor functions with any ceiling supply unit.


Screen

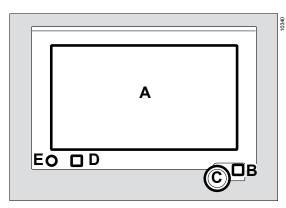

- A Touchscreen to call up functions or dialog windows
- **C** Rotary knob with background illumination to select and confirm settings
- D Key of for switching the workplace illumination on and off and dimming the illumination intensity in 3 steps (dark, medium, and bright)
- **E** Key ⊕ for turning the device on or off


Plug-in connection with Vapor View option for up to 3 vaporizers

The version for 2 vaporizers is shown here:

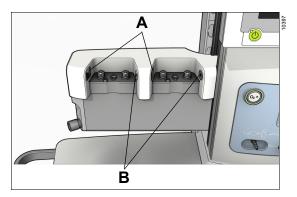

- A Sensor unit
- **B** Illumination unit

Breathing system

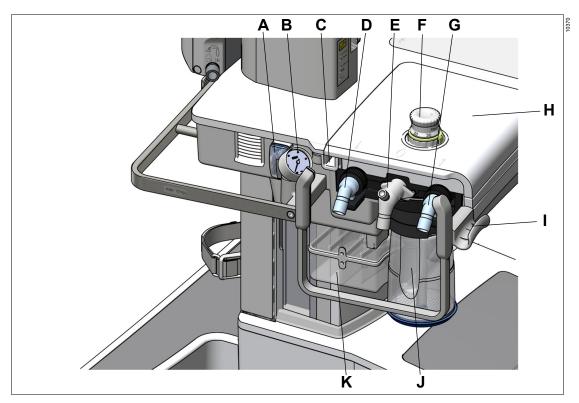

- A Water trap with connection for sample line
- **B** Pressure gauge (optional) for indicating the pressure in the internal breathing system
 The pressure in the internal breathing system is also indicated on the status display, see (G) on page 31 and (H) on page 32.

WARNING


Pressure gauge indication possibly inaccurate

Always compare the pressures indicated on the pressure gauge with those on the status display. Only use the pressure gauge as the primary source of information if the pressure indicator on the status display has failed.

- C Collecting channel
- **D** Inspiratory port

- **E** Bag elbow with circuit plug, e.g., to seal the Y-piece during an automatic test
- F APL valve
- **G** Expiratory port
- **H** Breathing system cover
- I Holder, e.g., for breathing bag hose
- J CO₂ absorber
- **K** Anesthetic gas receiving system


Side view from left

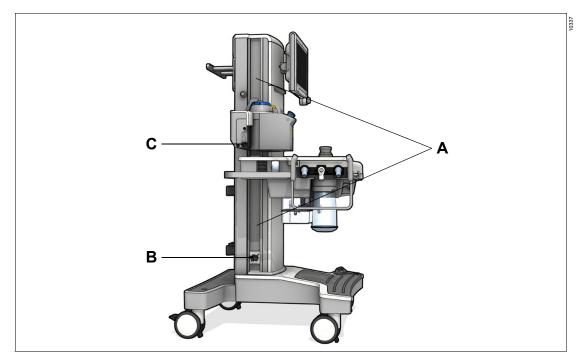

- **A** GCX rail for mounting additional workstation components
- **B** Strain relief for AGS hose, adjustable height
- **C** External fresh-gas outlet (optional)

Device column

The following illustration shows the left side of the device:

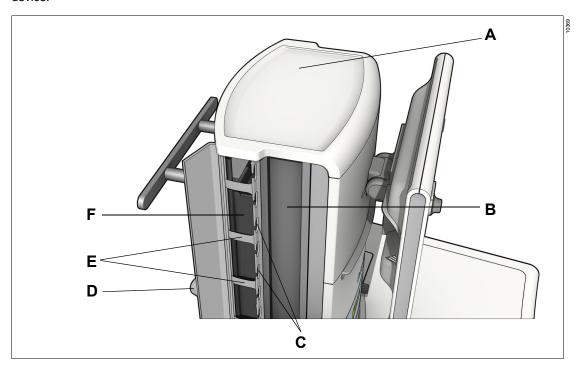


- A Column cover
- **B** GCX rail
- **C** Recesses for cables which lead e.g., into the device arms
- **D** Screw for closing the cable channel lid
- E Tabs to hold the cable
- F Cable channel

On the ceiling-mounted version, there is a compartment located under the column cover for stowing accessories, e.g., a cable.


Rear (trolley version)

Version with screw connections for gas cylinders

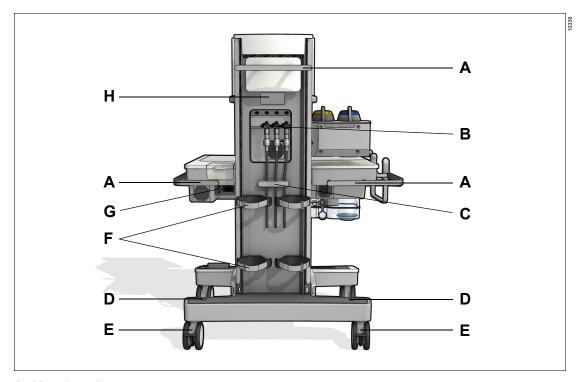

- A Mounting rail
- **B** Gas inlets
- C Strain relief for compressed gas hose (available for the trolley version only)
- **D** Connection for optional pole (38 mm)
- E Castor with castor brake
- Gas cylinder holders (optional) with hook-and-loop strap (available for the trolley version only)
- **G** Interfaces
- H Rating plate


Version with hanger yokes for gas cylinders with pin-index connections (optional)

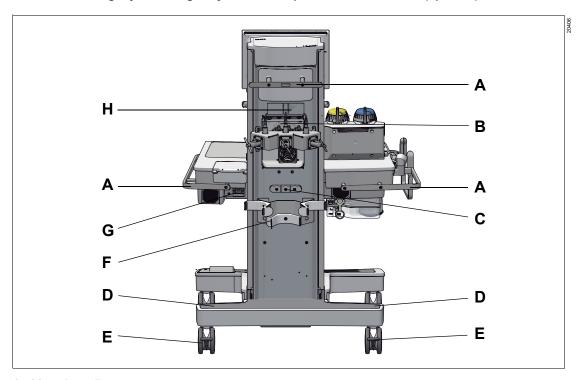

- **A** Mounting rail
- **B** Gas inlets
- C Strain relief for compressed gas hose (only available for the trolley version shown here or in combination with the flexibility trolley)
- **D** Connection for optional pole (38 mm)
- E Castor with castor brake
- F Gas cylinder holders (optional) with hook-andloop strap (available for the trolley version only)
- **G** Interfaces
- **H** Rating plate


Rear (ceiling-mounted version)

Differences from the trolley version:

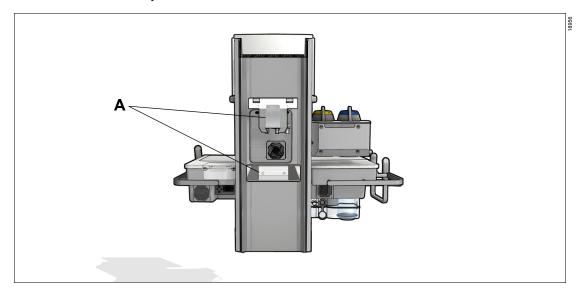


A Mounting interface for docking to the ceiling supply unit

Interfaces

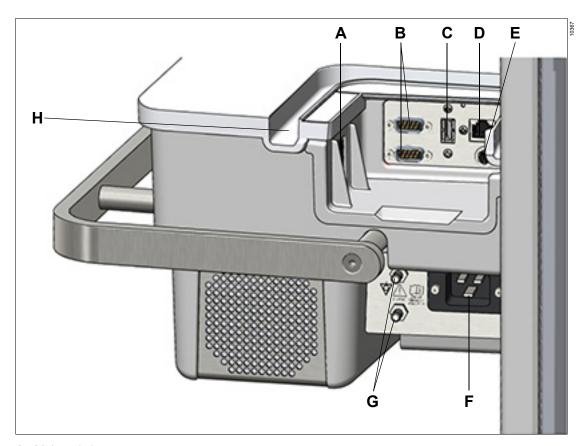


- A Main switch
- **B** Serial interfaces (COM 1 and COM 2)
- **C** USB interface
- **D** Network interface
- E Interface for workstation light (optional)
- **F** IEC connector (socket for power cable)
- G Potential equalization pin
- H Collecting channel

Gas inlets

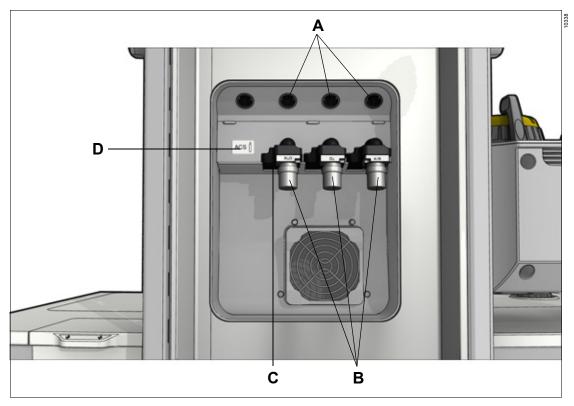

- A Connections for gas pressure measurement for gas cylinders (optional)
- **B** Connections for central gas supply (N2O optional)
- **C** Connections for external gas cylinders (optional)
- **D** Advanced Cylinder Support label (if present)


Hanger yokes (optional) for gas cylinders with pin-index connections

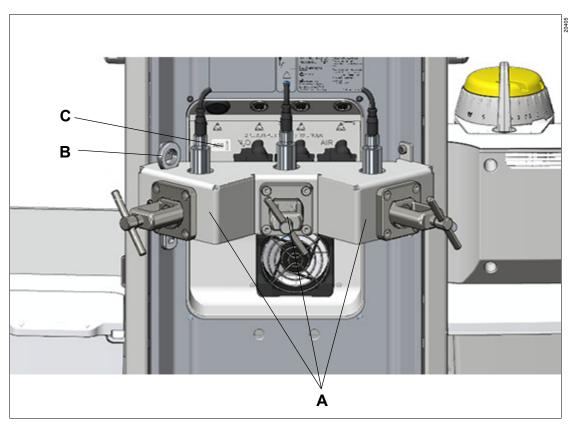

- A Hanger yokes (optional) for gas cylinders with pin-index connections
- **B** Wrench for opening and closing the gas cylinder valves
- **C** Advanced Cylinder Support label (if present)


Auxiliary power sockets (trolley version)

View with column cover removed:

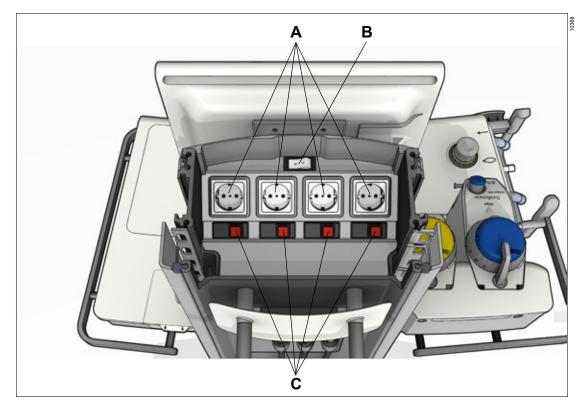

- A Auxiliary power sockets, depending on equipment variant
- **B** Isolation transformer switch (optional)
- **C** Circuit breakers or fuses, depending on equipment variant

Gas mixing unit (electronically controlled)

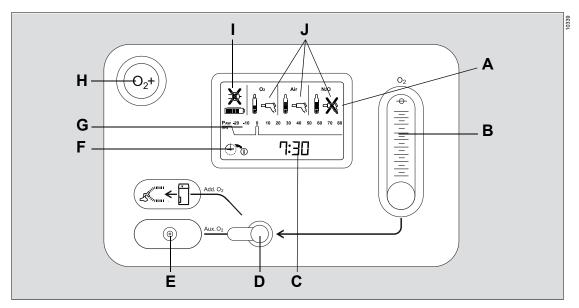

- A Status display
- **B** O2 flowmeter (for O2 insufflation *Aux. O2* and emergency O2 delivery *Add. O2*)
- C Current time or time for Auto On
- D O2 switch (for switching between O2 insufflation Aux. O2 and emergency O2 delivery Add. O2)
- E Outlet for O2 insufflation, e.g., for nasal cannula
- F Symbol for programmed Auto On
- **G** Display of pressure in the internal breathing system, see page 21
- H O2+ key (O2 flush)
- Symbols for mains power supply and power supply from internal battery
- J Symbols for gas supply (O2, Air, N2O) through central supply and gas cylinders

WARNING


Risk of mix-up due to deviating arrangement

In certain countries, the arrangement of the gases on the status display may differ from the arrangement of the flow tubes depicted on the screen.

Always pay attention to the respective labeling.

Explanation of the symbols which may be displayed, see page 277.

Gas mixing unit (mechanically controlled)

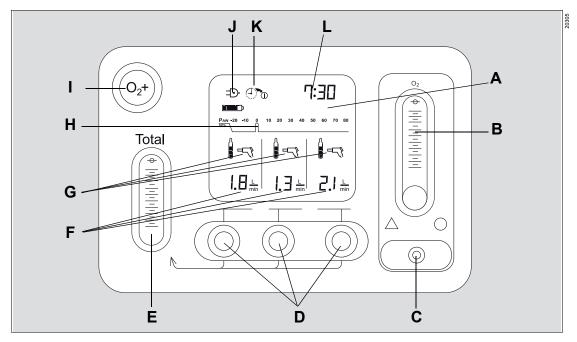

- A Status display
- B O2 flowmeter (for O2 insufflation Aux. O2)
- C Outlet for O2 insufflation, e.g., for nasal cannula
- **D** Flow control valves for fresh gas (O₂, Air, N₂O)
- **E** Total flow tube for fresh gas
- F Display of the set fresh-gas flows
- **G** Symbols for gas supply (O2, Air, N2O) from central supply and gas cylinders
- **H** Display of pressure in the internal breathing system, see page 21
- I O2+ key (O2 flush)
- J Symbols for mains power supply and power supply from internal battery
- K Symbol for programmed Auto On
- L Current time or time for **Auto On**

Explanation of the symbols which may be displayed, see page 277.

Functional scope

Device versions, options, and accessories

Some functions are optional and differ from the individual device configuration. Not all device versions or options are available worldwide.

Perseus is intended for use with the options and accessories listed in the associated list of accessories.

Gas delivery

Perseus can deliver mixtures of medical gases to which an anesthetic agent is added by means of a vaporizer:

Usable gas mixtures (electronically controlled gas mixer)

- O2/Air
- O2/N2O (optional)

Usable gas mixtures (mechanically controlled gas mixer)

O2/Air/N2O

Usable anesthetic agents

- Halothane
- Enflurane
- Isoflurane
- Sevoflurane
- Desflurane

Ventilation modes

- Manual / Spontaneous
 On selection: CPAP (optional)
- CPAP / PSV (optional)
- PC CMV
- PC BIPAP
- PC BIPAP / PS (optional)
- PC APRV (optional)
- VC CMV
- VC SIMV
- VC SIMV / PS (optional)
- VC CMV / AutoFlow
- VC SIMV / AutoFlow
- VC SIMV / PS / AutoFlow (optional)

For a detailed description of the ventilation modes and the additional settings, see page 258.

Additional operation modes

- External fresh-gas outlet (optional)
- Pause (optional)
- Monitoring

Monitoring

Perseus can monitor the following:

- Airway pressure
- Minute volume
- Inspiratory tidal volume
- Inspiratory anesthetic gas concentration
- Inspiratory O2 concentration
- Inspiratory and expiratory CO₂ concentrations
- Apnea (pressure, flow, and CO₂)
- Occurrence of anesthetic gas mixtures
- Degree of filling of the breathing bag

Display on the screen

Perseus can display the following:

- Waveforms
- Graphical trends
- Numeric trends
- Loops
- Alarm history
- Logbook
- Numeric parameters
- Prediction of the FiO₂ concentration (optional)
- Prediction of the anesthetic gas concentration (optional)
- Pre-configured lists for measured values and set values

Protocoling

Among other things, Perseus can save the following data in a logbook:

- Measured values
- Set values and their changes
- Patient data
- Ventilation modes
- Events (e.g., alarms, confirmed alarms, switchon time and switch-off time)
- Test results
- Gas consumption and anesthetic agent consumption

Gas supply

Gas supply variants

Gas	Central gas supply	Gas cylinder connection with
O2	Yes	permanently mounted Dräger pressure reducer (optional) or third-party pressure reducer (via country- specific NIST or DISS connection)
Air		permanently mounted
N2O	Yes (optional)	Dräger pressure reducer (optional)

Gas scavenging

Gas scavenging is conducted through the integrated anesthetic gas receiving system (AGS). The particle filter it contains filters the ambient air. This ensures the display accuracy of the flow indicator.

Data exchange, interfaces

Serial interfaces

Two serial interfaces, COM 1 and COM 2, are provided for data communication using the Dräger MEDIBUS data protocol.

USB interface

After a Dräger USB flash drive is connected, a USB interface enables, e.g., the following actions:

- Saving screen content as screenshot
- Saving and loading device configurations
- Saving system test results or protocols (partly optional) as text file

Network interface

If a corresponding service contract has been obtained, the Dräger RemoteService function can be executed via a network connection and the hospital network.

Perseus can be connected to the Dräger ServiceConnect Gateway or a DrägerService computer.

If the connected network offers an NTP service, the time on the device can be synchronized with the time on the NTP server.

For more information, see "Connections to IT networks" on page 253.


Infinity ID accessory support

- Exchange interval monitoring
- Anti-interchange security for breathing hoses

For more information, see page 269.

Gas flow plan


Breathing system

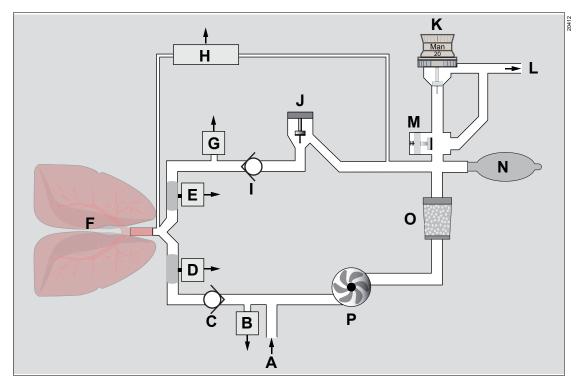

- A Gas supply from gas mixer
- **B** Inspiratory pressure measurement
- C Inspiratory valve
- D Inspiratory flow sensor
- **E** Expiratory flow sensor
- F Patient
- **G** Expiratory pressure measurement
- H Patient gas measurement module
- I Expiratory valve
- J PEEP/Pmax valve
- K APL valve
- L Anesthetic gas receiving system
- M Changeover between automatic ventilation and Manual / Spontaneous


- N Breathing bag
- O CO₂ absorber
- P Blower module TurboVent 2

Gas supply (electronically controlled gas mixer)

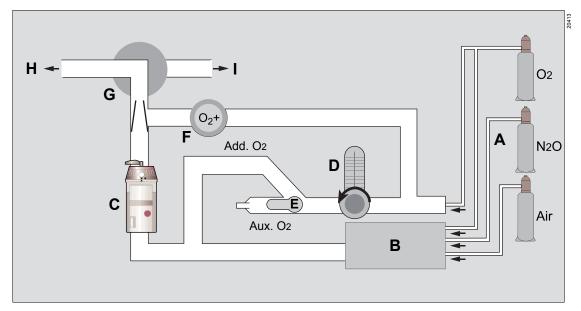

- A Gas supply (central supply or gas cylinders)
- **B** Gas mixer
- **C** Vaporizer
- D O2 flowmeter
- E O2 switch
- F O2 flush
- **G** Switch-over valve
- **H** Breathing system
- I External fresh-gas outlet (optional)

Gas supply (mechanically controlled gas mixer)

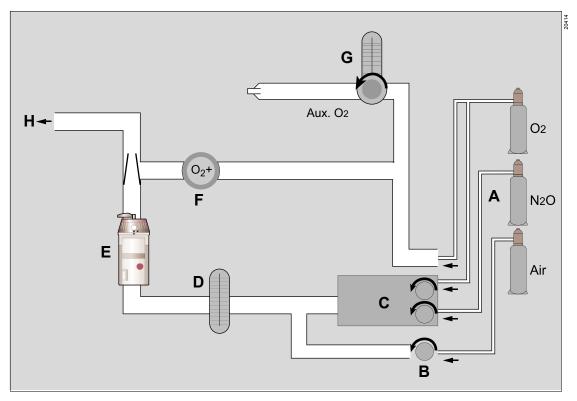

- **A** Gas supply (central supply or gas cylinders)
- **B** Flow control valves
- C Minimum O2 delivery
- **D** Total flow tube
- E Vaporizer
- F O2 flush
- **G** O₂ flowmeter
- H Breathing system


Operating concept

Screen	40
Main screen Therapy bar Dialog windows Quick setup window	40 40 41 41
Color concept	42
Colors of the control elements	42 42
gases	43
Day and night colors	43
Selecting and setting	43
Setting of parameters	43 43 44
Remote control for the ceiling-mounted version (combination with Dräger ceiling supply units)	44
Releasing the locking brakes	44 45 45
Safety sensor	45

Screen

Main screen

The main screen displays the most important information regarding anesthesia and ventilation.

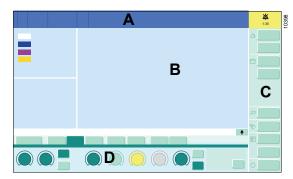

A Header bar

The header bar contains the following fields:

- Patient category
- Patient data
- System information (date, time, device name)
- Alarms, messages, and instructions for the user, see page 191
- Information regarding temporarily deactivated alarms

B Monitoring area

The following information is displayed in the monitoring area:

- Gas measurement
- Waveforms
- Parameter fields
- Loops (Pressure-Volume and Flow-Volume)
- Mini-trends
- Virtual flow tubes
- Prediction for anesthetic agents or FiO2 (optional)

Information regarding configurable fields, see page 173.

C Main menu bar

The main menu bar contains permanently assigned buttons to open dialog windows and activate functions.

These buttons are assigned to various groups. Touching a button opens the corresponding dialog window with the same name or activates the corresponding function.

For more information, see page 282.

D Therapy bar

The ventilation settings can be adjusted in the therapy bar.

Electronically controlled gas mixer:


- Tabs for selecting ventilation mode
- Therapy controls for ventilation parameters
- Therapy controls for fresh-gas delivery

Mechanically controlled gas mixer:

- Tabs for selecting ventilation mode
- Therapy controls for ventilation parameters

Therapy bar

The illustration shows the expanded therapy bar for the electronically controlled gas mixer:

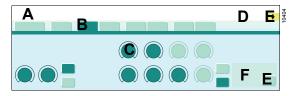

- A Name of active ventilation mode
- B Tabs
- C Therapy controls
- D Message field for information
- **E** Buttons to expand and shrink the therapy bar.
- F Field with additional information:
 - Additional and calculated set values
 - Spontaneous breathing activity of the patient

Start values

Arrows ▼ on the scales of the therapy controls mark the values resulting from the patient data and start settings. The start values can be configured, see page 150.

Linked therapy controls

Certain parameters can be linked to other parameters. If one parameter is changed, the linked parameter is also selected and changed. Among other things, this applies to the adjustment of ventilation pressures, ventilation times or during electronically controlled fresh-gas delivery.

Example: The device can be configured so that a change to the PEEP setting automatically causes a change to Pinsp; as a result, the difference between PEEP and Pinsp and therefore the tidal volume remain constant

Linking therapy controls, see page 162.


Additional information

Some settable parameters may be limited or be mutually restricted so that certain combinations of therapy settings are not possible, e.g., Ti 6.9 s at RR 100 /min.

If a condition is reached in which a parameter cannot be changed any more, Perseus displays a corresponding message in the message field (D).

Dialog windows


Dialog windows consist of one or several pages which are displayed by touching the corresponding horizontal or vertical tab.

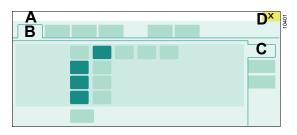

- A Dialog window title
- **B** Horizontal tab to open a page
- C Vertical tab to open subordinate structures
- **D** Button for closing the dialog window

Quick setup window

The **Quick setup** window is a context-sensitive dialog window. Depending on the selected parameter field or waveform, it contains various setting possibilities, e.g., for limits, scale, or content.

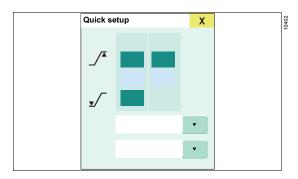

This window can be opened by touching the corresponding parameter field in the monitoring area. In the event of an alarm, the window can be opened automatically, see page 155.

Color concept


Colors of the control elements

The availability of functions and settings is indicated by certain colors of the therapy controls, of the therapy bar, and in dialog windows.

Therapy controls and buttons

Color	Example	Meaning
Dark green		Available element: function activated
Yellow		Selected ele- ment: not yet con- firmed with rotary knob
Light green		Available element: function not activated
Dark gray		Operating element: currently not available, function activated
Gray		Unavailable element

Rotary knob

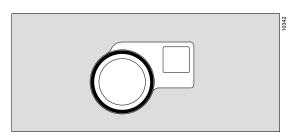

The rotary knob lights with different colors.

Color	Meaning
Blue	Therapy in progress
Orange	A selected function or setting must be confirmed.
Flashing orange	Functions or settings, which are still not confirmed, will be reset within the next 5 seconds.

Waveforms and parameters

Waveforms for mechanical breaths are displayed in the colors specified in the start settings, see page 151.

In the flow waveform, spontaneous breathing and pressure support are displayed in a light brown color.

Measured values whose specified accuracy cannot be maintained are displayed in dark gray.

Color coding for anesthetic agents and medical gases

Standardized color coding complying with ISO 5359 / ISO 32 / ISO 5360 is used to identify anesthetic agents and medical gases.

The colors for O₂, Air, and N₂O are adapted to locally applicable standards.

Day and night colors

There are 3 color schemes that can be selected:

- Day light
- Day dark
- Night

Setting the color schemes, see page 108.


Selecting and setting

Setting of parameters

Each of these settings requires confirmation by pressing the rotary knob.

1 Select

Touch the operation element (A). The color turns yellow. For therapy controls, the unit of the parameter to be set is displayed.

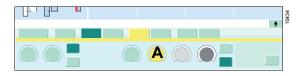

2 Set

Turn the rotary knob to set the value. For some therapy controls, faster turning raises the increment width.

3 Confirm

Press the rotary knob to confirm the value. The color of the operating element changes to green.

In the following chapters of this document, this sequence of action is represented by simplified explanations:

"Set the value."

Or

"Touch the button."

Canceling the setting procedure or the change procedure

If a change to a parameter should be canceled (color is still yellow), the following options exist to retain the previous setting:

- Touch the changed parameter again. This resets the selection of and the change to the parameter.
- Select another parameter. This selection resets the change to the other parameter.
- Do not press the rotary knob. After 15 seconds, the change is reset and signal tones sound during the last 5 seconds (timeout).

Activation of buttons

Some buttons are immediately active without additional confirmation. The color immediately turns to dark green.

Examples:

- Selecting a view
- Deactivating the CO₂ alarms.

Operating the flow control valves

The flow control valves of the mechanically controlled gas mixer and the O₂ flowmeter are operated as follows:

Opening the flow control valve

• Turn the flow control valve counterclockwise.

Closing the flow control valve

 Turn the flow control valve clockwise to the end stop.

In the subsequent chapters of this document, the following is represented by simplified explanations:

- "Open the flow control valve."
- "Close the flow control valve."

Remote control for the ceiling-mounted version (combination with Dräger ceiling supply units)

The remote control controls the following functions of the ceiling supply unit:

- Releasing the locking brakes for repositioning
- Height adjustment

The keys on the control panel on the media column of the respective ceiling supply unit can also be used for this purpose.

WARNING


Risk of operating errors and incorrect use Strictly follow the instructions for use for the ceiling supply unit.

WARNING

Risk of injury or material damage

When the medical device is docked to the ceiling supply unit, unintentional movement of the device may trap persons or objects and, in the case of persons, cause crushing injuries.

- Take particular care when moving the medical device.
- Prevent accidental activation of the remote control, e.g., by objects.

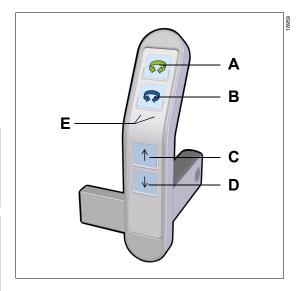

Releasing the locking brakes

Function of the keys on double-arm ceiling supply units:

- A Key (green) to release the brake on the ceiling bearing
- **B** Key (blue) to release the brake on the intermediate bearing

Function of the keys on single-arm ceiling supply units:

- **A** Key (green) has no function.
- B Key (blue) to release the brake on the ceiling bearing
- Keep key (A) or (B) pressed while simultaneously moving the arm system and Perseus.

Height adjustment

Prerequisite: Remote control is connected to the ceiling supply unit.

- C Up arrow key for lifting the ceiling supply unit
- D Down arrow key for lowering the ceiling supply unit
- Keep key (C) or (D) pressed and then release it when the ceiling supply unit has reached the desired position.

Status display

E LEDs for status indication

Indication	Meaning
Green LED lights.	The remote control is connected to the ceiling supply unit and is ready for use.

Depending on the software of the ceiling supply unit, the following status conditions may be indicated:

Indication	Meaning
Green LED flashing.	The control system of the ceiling supply unit has a malfunction.
Red LED lights.	The mechanical connection between Perseus and the ceiling supply unit is not working properly.

Safety sensor

The ceiling-mounted version of Perseus is fitted with a safety sensor, see page 19.

The safety sensor is triggered by resistance from below, e.g., as a result of obstacles such as objects or persons. When the safety sensor is triggered, lowering of the ceiling supply unit is stopped.

To continue lowering the ceiling supply unit after the safety sensor has triggered, perform the following steps:

- Briefly lift the anesthesia machine using the ceiling supply unit.
- **2** Make sure that no obstacles are present.

Assembly and preparation

Perseus as a ceiling-mounted version	47
Prerequisites and notes	47
Combination with Dräger ceiling supply units	47
Combination with other ceiling supply units Flexibility trolley for ceiling-mounted version	48
	48
Mounting of accessories	49
Safety information	49
Mounting on the sides of the device column Mounting on the mounting rails	49 50
Special features of the ceiling-mounted version	50
	50
Before first operation	52
Establishing the mains power supply	52
Charging the battery	53
Connecting additional devices to the auxiliary power sockets	53
Establishing potential equalization	54
Connecting devices or computers to the data	•
interfaces	54
Intrahospital transport	56
Increasing the tipping stability during transport	56
Parking the medical device	57
Visual inspection after transport	57
Establishing the gas supply	57
Central gas supply	58
Connecting the gas cylinders Fitting the vaporizers	58 61
	-
Ensuring the gas scavenging	62
Preparation for a day of operations / after	
reprocessing	63
Assembling the breathing systemInserting the breathing system	63 66
Mounting the flexible arm (optional) or the bag	00
elbow	67
Selecting and connecting patient-specific	
accessories	67
Safety information	67
Fitting the breathing circuit and the filters Breathing bag	68 71
Di Cati mig Dag	, ,

Observing resistance and compliance	
Connecting or exchanging consumables	76
Single-use CO2 absorber	76
Reusable CO2 absorber	76
Water trap	77
Connecting the sample line	

Perseus as a ceiling-mounted version

Prerequisites and notes

Docking and undocking

Perseus may only be docked to a ceiling supply unit by service personnel. When docking and undocking, observe the instructions for use of the ceiling supply unit and the service cart.

WARNING

Risk of injury or damage to the device

Safety functions on the Dräger ceiling supply unit can be temporarily disabled with the keyswitch.

After the device has been docked or undocked by service personnel, make sure that the key is not left in the keyswitch but instead is kept in a secure place.

Establishing the cable connection between Perseus and the ceiling supply unit

WARNING

Risk of injury

If the cable connection between Perseus and the ceiling supply unit is not established, the safety sensor will be inoperative.

Before using Perseus, make sure that the cable connection has been established.

- Establish the cable connection between Perseus and the ceiling supply unit.
- Observe the instructions for use of the ceiling supply unit.

Combination with Dräger ceiling supply units


Keyswitch

The ceiling supply unit must be set up for operation with Perseus. Located on the media column there is a keyswitch with the following setting positions:


- A Operation of Perseus or Atlan
- **B** Docking and undocking of Perseus or Atlan
- C Docking, undocking, and operation of other Dräger devices

The keyswitch can be marked with 2 different product labels.

Product label 1:

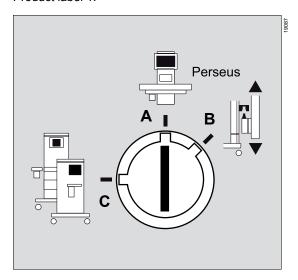

Product label 2:

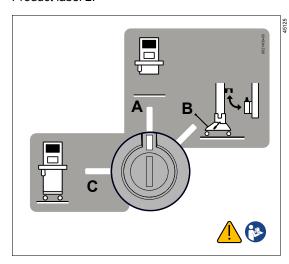

 Make sure that the keyswitch on the ceiling supply unit is in position (A). Contact service personnel if necessary.

WARNING

Risk of injury

If the ceiling supply unit has not been set up for operation with Perseus or the keyswitch is in an incorrect setting position, the safety sensor will be inoperative.

Only operate Perseus with ceiling supply units that have been set up for Perseus and have the keyswitch in the correct setting position. If necessary, contact service personnel for clarification and specialized service personnel if the keyswitch is faulty.

Combination with other ceiling supply units

Perseus can be used with ceiling supply units from other manufacturers. In this case, a special kit (see list of accessories) is required. The ceiling supply unit must be prepared before Perseus can be mounted.

Dräger will provide the manufacturer of the ceiling supply unit with the document "Perseus A500 anesthesia workstation - ceiling-mounted version - interface description for ceiling supply units". The instructions in this document must be followed during construction.

Flexibility trolley for ceiling-mounted version

The ceiling-mounted version of Perseus can be mounted on a flexibility trolley and can then be used as a trolley version.

When Perseus is mounted on the flexibility trolley, the remote control is inoperative and the safety sensor is disabled.

Information on attaching accessories to Perseus on the flexibility trolley can be found in the chapter "Mounting of accessories".

Mounting of accessories

Information on mounting accessories is described in the assembly instructions.

Safety information

CAUTION

Risk of tipping over

When support arms are used, the end stops must be in perfect condition. Otherwise there is a risk that the support arms will swing out of control.

Check the functional integrity of the support arm end stops after the following activities:

- After fitting accessories
- After transporting the device

Mounting on the sides of the device column

Depending on the arm length, the following maximum weights are permissible:

Maximum weight			
Arm length	Version with plug-in connections for 2 vaporizers	Version with plug-in connections for 3 vaporizers	
215 mm	20 kg	13.5 kg	
(8.5 in)	(44 lbs)	(29.8 lbs)	
300 mm	15 kg	8.5 kg	
(11.8 in)	(33 lbs)	(18.7 lbs)	
400 mm	15 kg	8.5 kg	
(15.7 in)	(33 lbs)	(18.7 lbs)	
570 mm	15 kg	8.5 kg	
(22.4 in)	(33 lbs)	(18.7 lbs)	

Note the following information in addition:

- If more than one arm is installed on a device side, the weights or arm lengths must be reduced. Do not exceed the maximum tilting moment.
- Maximum weight on the right side of the device: 35 kg (77 lbs)
- Maximum weight on the left side of the device:
 - Versions with plug-in connections for 2 vaporizers: 25 kg (55 lbs)
 - Versions with plug-in connections for 3 vaporizers: 15 kg (33 lbs)
- Maximum installation height: 1400 mm (55.1 in)

The maximum tilting moment is calculated based on the following equation:

Variants with plug-in connections for 2 vaporizers:

One attached P= (A * B) <8500 mm*kg component: (<738 in*lbs)

Multiple attached P1 + P2 + P... <8500 mm*kg

components: (<738 in*lbs)

Versions with plug-in connections for 3 vaporizers:


One attached P= (A * B) <5000 mm*kg component: (<434 in*lbs)

Multiple attached P1 + P2 + P... <5000 mm*kg

components: (<434 in*lbs)

Explanation:

- A: Distance (including arm length) to outer edge of the attached component
- B: Arm weight (B1) + weight of attached component (B2)
- P: Product of distance and weight

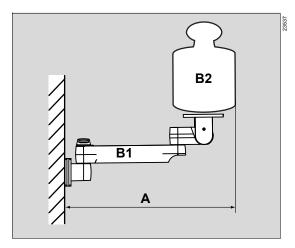

Example of a permissible configuration:

Maximum length	Arm weight	Weight of the attached component
400 mm	4 kg	5 kg
(15.7 in)	(8.8 lbs)	(11 lbs)
300 mm	3 kg	12 kg
(11.8 in)	(6.6 lbs)	(26.5 lbs)

P1= 400 mm * (4 kg + 5 kg) = 3600 mm*kg (P1=15.7 in * (8.8 lbs + 11 lbs) = 310.9 in*lbs)

P2= 300 mm * (3 kg + 12 kg) = 4500 mm*kg (P2= 11.8 in * (6.6 lbs + 26.5 lbs) = 390.6 in*lbs)

3600 mm*kg + 4500 mm*kg = 8100 mm*kg (310.9 in*lbs + 390.6 in*lbs = 701.5 in*lbs)

8100 mm*kg < 8500 mm*kg (701.5 in*lbs < 738 in*lbs)

If a combination of arms on one device side exceeds the value of 8500 mm*kg (738 in*lbs), tipping stability as per IEC 60601-1 is no longer given. Check the tilting stability.

Mounting on the mounting rails

Depending on the position of the mounting rail, the following weights are permissible:

Position	Maximum weight
On the side of the working surface	10 kg (22 lbs)
At the top on the rear of the device	2.5 kg (5.5 lbs)

WARNING

Risk of tipping over

If the weight of the mounted accessories exceeds the permissible maximum weight, the medical device may tip over.

Observe the maximum weight per arm.

WARNING

Risk of tipping over

If the weight of the accessories is unevenly distributed about the device, the medical device may tip over.

Distribute the weight evenly.

Special features of the ceiling-mounted version

WARNING

Risk of tipping over

If the ceiling-mounted version has been placed on the flexibility trolley or the service cart, the information for fitting the accessories still applies as for the ceiling-mounted version.

- Pay attention to the maximum weight of the accessories.
- Pay attention to the distribution of the weight.

CAUTION

Risk of injury or material damage

If an accessory is fitted to Perseus which projects below the lowest point of the anesthesia workstation, the accessory may strike objects or persons when the ceiling supply unit is lowered.

Only fit accessories that maintain a distance of at least 12 cm (4.7 in) from the floor when the anesthesia workstation is lowered to its lowest point.

The total weight of all accessory parts that are fastened to standard rails or are placed on the shelves must not exceed 20 kg (44 lbs). Additionally, all product labels regarding the maximum load that are attached to standard rails, shelves, and ceiling supply units must be observed.

CAUTION

Risk of damage to the device

If the maximum load capacity of the ceiling supply unit is exceeded, the load cannot be supported.

Refer to the instructions for use of the ceiling supply unit and observe the maximum load capacity.

WARNING

Risk of injury or material damage

If Perseus is not properly mounted on the flexibility trolley or is being dismounted from it, the trolley may tip over or Perseus may fall off.

The mechanical coupling and uncoupling of Perseus and the flexibility trolley must only be performed by specialized service personnel.

WARNING

Risk of injury or material damage

If Perseus is mounted on a flexibility trolley and additional devices have been attached, it is possible that the maximum load capacity will be exceeded when Perseus is coupled to the ceiling supply unit again.

- Observe the maximum load capacity of the ceiling supply unit.
- Use only those combinations of patient monitoring and accessories that are permissible for the ceiling-mounted version of Perseus.
- Observe the maximum weights and tilting moments.

Before first operation

NOTE

Prior to initial commissioning, the breathing system and blower module TurboVent 2 must be reprocessed. Carry out the reprocessing as described in the reprocessing instructions supplied with the product.

Establishing the mains power supply

The mains voltage must correspond to the voltage range indicated by the rating plate on the rear of the device.

WARNING

Risk due to incorrect mains voltage or missing protective ground


If the device is connected to a power socket with incorrect mains voltage or a power socket without a protective ground, an electric shock may occur.

Connect the device only to power sockets with correct mains voltage and a protective ground.

NOTE

The mains plug must be readily accessible so that the power supply to Perseus can be interrupted quickly in the event of a device failure.

- Connect the power cable to the device and secure it.
- 2 Lay the power cable in the cable conduit on the right side of the device.

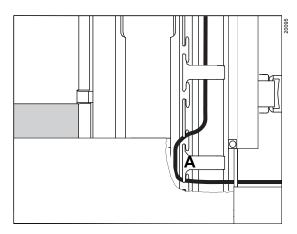

3 Pass the power cable around the lower lug (A) for strain relief.

WARNING

Risk of electrical overload

If the device is connected to additional power socket strips, this may lead to increased leakage current. The leakage current may exceed the permissible values.

- Do not connect the device to additional power socket strips.
- Do not connect additional power socket strips to the auxiliary power sockets under the column covering.
- **4** Plug the power cable into the mains power socket on the wall.
- 5 Check the displays for mains voltage and battery on the status display.
- Set the main switch to position *I*, see page 27.

Charging the battery

The internal batteries are automatically charged by the mains power supply.

WARNING

Risk of device malfunction

If the batteries are not sufficiently charged, it may not be possible to maintain operation for long enough if the mains power supply fails.

Before first operation or after storage, charge the batteries for at least 8 hours.

WARNING

Risk due to reduced power supply from the internal battery

Batteries are wear parts. The capacity of the battery diminishes with the period of use.

Check the functional capability of the battery by performing regular inspections.

WARNING

Risk of patient injury

If Perseus is operated or connected to the mains power supply at ambient temperatures above 35 °C (95 °F), the battery cannot be charged properly. The power supply out of the battery may be limited.

Do not expose the device to temperatures above 35 °C (95 °F) on a permanent basis.

Connecting additional devices to the auxiliary power sockets

Perseus is equipped with auxiliary power sockets and circuit breakers or fuses under the column covering, see page 30.

- 1 Connect the power cable of the additional device to an auxiliary power socket.
- 2 For Perseus with isolation transformer: Switch on the power for the isolation transformer and the auxiliary power sockets using the isolation transformer switch.

Make sure that the maximum power consumption of additional devices does not exceed permissible values

WARNING

Risk of fire

Sources of ignition such as electrosurgical equipment, laser surgical equipment, sparks, damaged cables, or damaged plugs can cause fires in combination with oxygen or nitrous oxide. As a result, the patient or user may be at risk.

- Maintain a distance of at least 200 mm (7.9 in) between electrical connections and components which conduct oxygen and nitrous oxide.
- Cables and connections must be sufficiently insulated and must not be damaged. Check cables for damage daily.
- Do not allow the device to come into contact with sources of ignition.

WARNING

Risk of electric shock

The connection of devices to auxiliary power sockets can lead to an increased leakage current. If the protective ground of one of these devices fails, the leakage current may rise above the permissible values. Only connect devices with the permission of the respective manufacturer. Have the leakage current checked by service personnel.

If the permissible value is exceeded, use a mains power socket on a wall instead of the auxiliary power socket of the device.

WARNING

Risk of electric shock and of device malfunction

Any connected devices or device combinations not complying with the requirements in these instructions for use may compromise correct functioning of the medical device.

- Do not connect high-frequency surgery equipment to the auxiliary power sockets of the anesthesia machine.
- Before using the medical device, refer to and strictly comply with the instructions for use of all connected devices and device combinations.

WARNING

Risk of electrical shock due to penetration of fluids

When the device is connected to the mains power supply, the auxiliary power sockets underneath the column cover receive current.

- Reattach the column cover correctly after additional devices have been connected to the auxiliary power sockets.
- Only use the device with the column cover closed.

WARNING

Risk of fire

Components such as power supply units which are connected to the auxiliary power sockets on Perseus become warm during operation. If these components are placed under the column cover, they may overheat.

Do not place components that become warm under the column cover.

Establishing potential equalization

Differences in electrical potential between devices can be reduced by potential equalization.

Potential equalization does not replace the protective ground connection.

During operation, the potential equalization connectors must be readily accessible and must be removable without tools.

Connecting the potential equalization cable

- 1 Connect the potential equalization cable to the potential equalization pin on the device.
- 2 Connect the potential equalization cable to a potential equalization connector of the hospital (e.g., wall, ceiling supply unit, operating table).

Connecting devices or computers to the data interfaces

This device is equipped with data interfaces such as RS232, LAN, and USB. These interfaces can be used to set up a network in accordance with IEC 60601-1.

WARNING

Risk of electric shock

If USB devices with their own power supply are used, the patient leakage current may increase.

- Only use USB devices that do not have their own power supply.
- Do not connect USB devices by means of charger cables.

WARNING

Risk of electric shock

Only connect devices to the serial interfaces, or devices and networks to the network interface, which have a maximum nominal voltage of 24 Vdc and meet one of the following standards:

- IEC 60950-1 / IEC 62368-1: ungrounded SELV circuits
- IEC 60601-1 (as of 2nd edition): touchable secondary circuits

WARNING

Risk of electric shock

Connecting devices to the serial data interfaces can lead to an increased leakage current. If the protective ground connector of one of these devices fails, the patient leakage current may rise above the permissible values.

- Have the leakage current checked by service personnel.
- If the permissible value is exceeded, disconnect the devices from the serial interfaces.
- Do not touch the connectors of the interfaces and the patient simultaneously.

WARNING

Risk due to impermissible data connections

The impermissible use of data interfaces can result in new hazards.

- Only make connections to data interfaces with permission from the responsible organization (IT representative and the hospital equipment officer).
- Follow the information in the chapter "Connections to IT networks".

Establishing a data connection

• Connect the device to a network or a computer.

See page 27 for a picture of the interfaces.

Only use the cables from the list of accessories.

For further information on configuring the respective interface, see page 168.

Intrahospital transport

Transport comprises any movement of a medical device that does not serve solely for positioning.

WARNING

Risk of damage to the pressure reducers

If the device collides with an obstacle during transport, the pressure reducers may be damaged.

- Align the pressure reducers so that they are protected from collisions during transport.
- Before transporting, close the valves on the gas cylinders.

WARNING

Risk of tipping over during transport

The medical device may tip over if handled incorrectly.

- The medical device may only be moved by people who have the physical ability to do so.
- Dräger recommends that the medical device always be transported by two people. This improves maneuverability.
 When transporting over inclines, around corners, or over thresholds (e.g., through doors or in elevators), make sure that the medical device does not bump against anything. Do not pull the medical device over hoses, cables, or other obstacles lying on the floor.
- Do not activate the central brake while the medical device is being moved.
- To push the medical device, hold on to the mounting rails on the work surface.
- Do not lean against the medical device.

WARNING

Risk of injury

The device may tip over when it is transported over inclines.

Always move the device using two people.

CAUTION

Risk of tipping over

When support arms are used, the end stops must be in perfect condition. Otherwise there is a risk that the support arms will swing out of control.

Check the functional integrity of the support arm end stops after the following activities:

- After fitting accessories
- After transporting the device

The Perseus ceiling-mounted version can be transported on the service cart for servicing purposes.

Observe the following section: "Transport for maintenance purposes" (page 216)

Increasing the tipping stability during transport

- 1 Carefully fold the holding arm with any mounted equipment against the device, e.g., patient monitoring, data management systems, and syringe pumps. When folded in, these components should not project beyond the mounting rails if possible.
- 2 Remove loose objects from the attached arms and the shelves.
- 3 Remove the heavy objects of more than 8 kg from the mounting rails, e.g., the vaporizers.
- 4 Clear the writing tray and slide it completely into the device.
- 5 Position the optional flexible arm for the breathing bag close to the device.
- 6 Push the optional drawers in.
- 7 Lock the lockable optional drawers.

Parking the medical device

When parking, always engage the brakes (central brake for front castors, individual wheel brakes at rear), especially on inclined surfaces.

Visual inspection after transport

- 1 Visually check the medical device for damage, particularly the hoses and cables.
- 2 Any damage must be repaired by service personnel before using the device.

Establishing the gas supply

WARNING

Risk due to gas supply failure

All gas supplies (central gas supply, gas cylinders) must be correctly connected since otherwise the backup system (gas cylinders) will not be available if gas supply fails.

- Make sure that all compressed gas hoses are correctly connected to the rear side of the device.
- After connecting the gas supply, check for correct function.
- Even when the anesthesia machine is connected to the central gas supply, the gas cylinders should remain at the device with valves closed as backup.

WARNING

Risk of device malfunction

Gas supply (central gas supply or gas cylinders): To avoid damage to devices connected to the gas supply, only use medical gases. In particular, observe the national and international standards regarding the use of medical gases.

WARNING

Risk due to failure of central gas supply

If the central gas supply fails, additional devices such as such devices which are connected to the gas outlets on the device, will no longer be supplied with gas. The device will be supplied from the gas cylinders.

When using such additional devices, always monitor the gas supply independently of the main device.

WARNING

Danger to the patient and user

The device may be damaged if the strain relief of the compressed gas hoses is not used correctly.

Use the strain relief of the compressed gas hoses correctly.

WARNING

Risk of fire

Ignition sources in combination with oxygen can lead to fires.

Do not position oxygen sources in the vicinity of ignition sources, e.g., electrical connectors.

Central gas supply

- Screw the compressed gas hoses for the central gas supply to the gas inlets on the rear of the device by hand, see page 28.
- 2 The following applies to the trolley version and to the ceiling-mounted version on the flexibility trolley: Insert the compressed gas hoses into the strain
 - relief and screw the strain relief tightly into place, see page 24.
- **3** Connect the compressed gas hoses to the terminal units.
- 4 Check if all gas supplies are correctly connected. Check if the gases to be supplied are available by observing the status display (see page 31).

Connecting the gas cylinders

WARNING

Risk due to gas supply failure

Using a pressure reducer without the required pressure sensor in place of a Dräger pressure reducer or connecting a central supply hose to the connection for the external gas cylinders will prevent the availability of the gas cylinders from being monitored during the system test and while the system is in use. The back-up functionality may be compromised.

If the gas cylinder monitoring function is not available, the user must employ suitable pressure monitoring conforming to ISO 80601-2-13 which allows the pressure indicator to be read from the user's operating location at the front. Do not connect central supply hoses to the connectors for external gas cylinders.

WARNING

Risk of injury or damage to the device

Defective pressure reducers can cause the connected device to malfunction.

- Service the pressure reducer regularly.
 Follow the information in the chapter
 "Maintenance".
- If signs of wear are visible, replace the pressure reducer.


Connections of gas cylinders and pressure reducers must be undamaged as well as free of dust, particles, and grease. Otherwise, there is risk of fire.

When handling pressure reducers, observe the relevant national laws and regulations.

Special features of the ceiling-mounted version

There are no gas cylinder holders on the rear of the device.

Connecting the gas cylinders on variants with screw connections

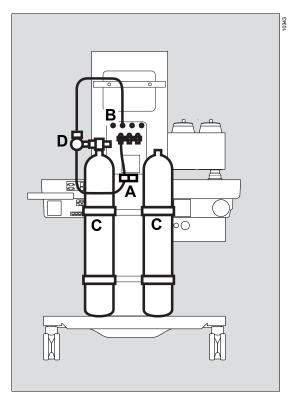

- Insert the compressed gas hoses into the strain relief (A) and screw the strain relief tightly into place.
- 2 Connect the pressure measurement lines to the connectors (B) above the gas inlets.
- 3 Set the full gas cylinders (C) into the gas cylinder holders and secure them with the hook-and-loop straps.

WARNING

Risk of injury or damage to the device

Pressure reducers have an internal relief valve. If a defect occurs, gas can escape into the ambient air.

Do not block or cover the relief valve.

WARNING


Risk of damage to the device

When connecting the pressure reducers, ensure that they do not protrude beyond the device.


4 Tightly screw the pressure reducers (D) to the gas cylinder valves. The connections must fit each other directly; do not use transition pieces.

Connecting the gas cylinders on variants with pin-index connections

Gas cylinder with pin-index connection:

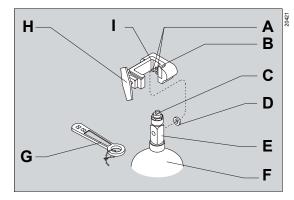

Rear of device:

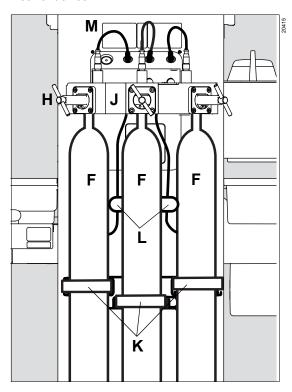

Before first use:

- Insert the compressed gas hoses into the strain relief (L) and screw the strain relief tightly into place.
- 2 Connect the pressure measurement lines to the connections (M).
- 3 Two fixing positions are possible for the gas cylinder holder (K). Adjust the position of the gas cylinder holder to the size of the gas cylinder in use. Contact service personnel to do this.
- 4 Remove the protection cap from the head of the cylinder.

When changing gas cylinders:

- 5 Remove the old sealing washer (D).
- 6 Insert a new sealing washer (D) at the cylinder holder (I).
- 7 Make sure that both pin-index pins (A) are present below the gas inlet (B).
- 8 Align the gas cylinder (F) so that the pin-index holes on the head of the cylinder (E) are pointing towards the pin-index pins (A) on the cylinder holder (I).
- 9 Insert the cylinder head (E) of the gas cylinder (F) from below into the cylinder holder (I) of the hanger yoke (J).
- **10** Allow the pin-index pins (A) to engage in the pin-index holes.
- 11 Turn the handle (H) on the cylinder holder (I) clockwise. The tip of the threaded retaining pin will then be turned into the visible recess on the cylinder head.
 - Make sure that the gas cylinder (F) is suspended vertically.
- **12** Tighten the handle (H) on the cylinder holder (I).
- **13** Secure the gas cylinders (F) with hook-and-loop straps (K).

If required, the gas cylinder valve (C) can be opened with the supplied wrench (G).

Handling O₂ gas cylinders

WARNING

Risk of explosion

When pressurized, O2 is self-igniting in combination with oil or grease.

Do not oil or grease the gas cylinder valve or the pressure reducer of the O2 cylinder. Do not touch with oily or greasy fingers.

Gas cylinder valves must only be opened and closed slowly. Do not use any tools with the screw connection variants.

Have service personnel repair any leaky or stiff gas cylinder valves.

Fitting the vaporizers

Depending on the configuration, Perseus can be operated with vaporizers which have a Dräger Auto Exclusion plug-in adapter or a Selectatec plug-in adapter. Dräger recommends using only vaporizers that have been tested and are listed in the list of accessories.

WARNING

Risk due to incorrect anesthetic agent delivery

If the vaporizer is filled with the wrong anesthetic agent or if it is not filled sufficiently, incorrect anesthetic gas concentrations or concentrations that are too low can occur as a result.

- Follow the instructions for use for the vaporizer.
- Compare the color coding and labeling on the vaporizer used with the anesthetic agent bottle and the anesthetic agent indicated on the screen.

The vaporizers used must conform to the ISO 8835-4 or ISO 80601-2-13 standard. If the internal patient gas measurement module fails, an independent gas measurement system complying with ISO 80601-2-55 must be used.


WARNING

Risk due to improperly mounted vaporizers

Incorrectly mounted vaporizers can cause leakage. This can cause the fresh-gas delivery to be too low or contaminate the ambient air. Patient and user can be endangered.

- Make sure that the connected vaporizers are hanging vertically and are not jammed by objects leaning on the device.
- When using D-Vapor vaporizers, make sure that the power cable is not pinched.
- After mounting the vaporizers, perform a leakage test.

The following illustration shows the Dräger Vapor 3000.

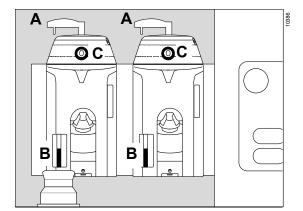

- Set the vaporizer on the plug-in connector evenly and firmly.
- 2 Turn the locking lever (A) clockwise. The lever is in the locked position when it is pointing left.
- 3 Check the vaporizer filling level through the sight glass (B). If necessary, fill the vaporizer.
- **4** Turn the control dial to position **0**; the key (C) latches into place.
- 5 Check the locking: Set the control dial on a vaporizer to a position other than 0 and make sure the other vaporizer remains locked in its 0 position. Repeat the test with the other vaporizer.
- **6** Set both control dials to position **0**.

When connecting D-Vapor 3000 vaporizers

- Connect the power cable to an auxiliary power socket beneath the column covering of Perseus.
- **2** If needed, establish potential equalization.
- 3 If needed, stow the cable in the cable conduit.

Vapor View option

When combined with the Dräger Vapor 3000 or D-Vapor 3000, the Vapor View option enables the following functions:

- Illumination of control dial and sight glass on vaporizer
- Reading the anesthetic agent type
- Reading the control dial position
- Monitoring the filling level
- Prediction of the anesthetic gas concentration

Ensuring the gas scavenging


CAUTION

Risk of ambient air contamination


If the anesthetic gas receiving system is not connected to the disposal system, contamination of the ambient air with anesthetic gas may result. Ensure the following:

- The anesthetic gas receiving system is correctly connected to the disposal system.
- The anesthetic gas scavenging system is functioning properly. Check the flow indicator.

Perseus is equipped with an integrated anesthetic gas receiving system (AGS).

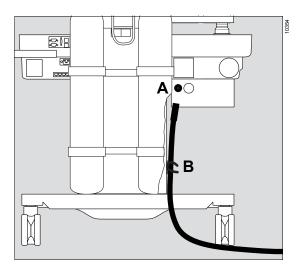

- Connect the scavenging hose to the nozzle (A) on the receiving system.
- 2 Secure the scavenging hose with the strain relief (B).

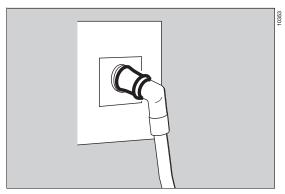

3 Connect the probe of the scavenging hose to the terminal unit of the scavenging system.

As an option, the integrated anesthetic gas receiving system can be operated in combination with a control valve. Observe the assembly instructions for the control valve.

Preparation for a day of operations / after reprocessing

Assembling the breathing system

WARNING

Risk of insufficient anesthetic gas concentrations

If the component connections of the breathing system are not leak-tight enough, ambient air may be added to the anesthetic gas mixture.

Make sure that all components of the breathing system are connected tightly.


- 1 Check all parts for damage or wear:
 - TurboVent 2 blower module
 - API valve
 - Upper housing of the breathing system
 - Lower housing of the breathing system
 - Valve cages and valve plates
 - Flow sensors
 - Ports
 - Incident flow meshes in the inspiratory limb of the lower part of the breathing system housing and in the expiratory port
 - Seals and sealing rings

- **2** The following parts must be free from deposits:
 - Incident flow meshes
 - Non-return valves (red and blue)


If necessary, remove deposits on the valve plates of the non-return valves with a soft cloth.

If the incident flow meshes are damaged, contact DrägerService.

Installing the TurboVent 2 blower module

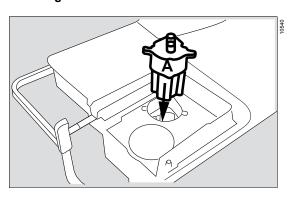


1 Insert the TurboVent 2 blower module (A) into the breathing system mount.

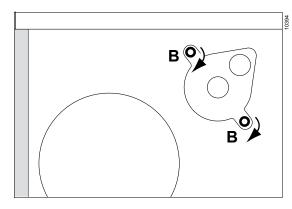

2 Screw in the quick release screws (B) tightly with a clockwise rotation of 90°, e.g., with a coin.

Inserting the non-return valves


There are recesses in the valve cages of the non-return valves. These recesses are arranged differently on the red valve cage (inspiratory valve) than on the blue valve cage (expiratory valve). When inserting the valve cages, these recesses align with the corresponding lugs on the lower housing of the breathing system.

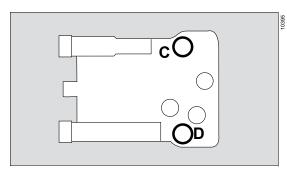

Perform the following steps for non-return valves (C) and (D):

- 1 Align the valves correctly.
- 2 Insert the valves.

Fitting the APL valve

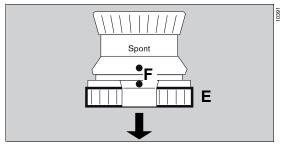

- 1 Set the APL valve vertically on the upper housing of the breathing system.
- 2 Align the APL valve correctly and screw it tight with the knurled nut (E).

When the breathing system is installed, points (F) must point towards the user.


WARNING

Risk of an incorrectly set pressure limitation

If the lower dot marked on the APL valve is not aligned correctly, APL valve settings may be made incorrectly or read off incorrectly.

When fitting, take care that the lower dot is facing the user during operation.

Mounting the upper part of the breathing system housing

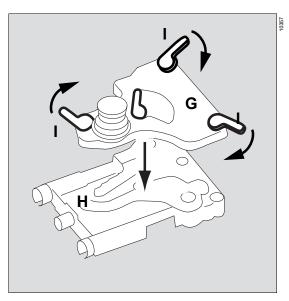

- 1 Place the upper part (G) of the breathing system on the lower part (H).
- 2 Pay attention to the correct alignment of the three levers (I), see illustration above.
- 3 Turn the 3 levers (I) by about 120° clockwise.


Inserting the flow sensors and the ports

WARNING

Risk of fire

High temperatures arise in the flow sensors, for example, in operation or while calibrating during a system test. Due to the high temperatures, residual vapors of flammable disinfectants (e.g., alcohols) and residues that were not removed during reprocessing may ignite.

- Ensure particle-free cleaning and disinfection.
- After disinfection, allow the flow sensor to air for at least 30 minutes.
- Before inserting the flow sensor, check for visible damage and soiling, such as residual mucus, medication aerosols, and particles.
- Replace flow sensors when damaged, soiled, or not particle-free.

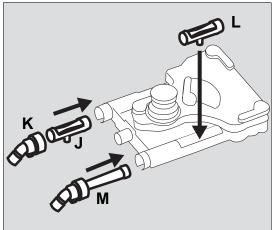


- 1 Insert flow sensors (J) and (L).
- 2 Insert the inspiratory port (K) and screw it tight with the knurled nut.
- 3 Insert in the expiratory port (M) and screw it tight with the knurled nut.


Mounting the CLIC adapter (optional)

The CLIC adapter allows the following disposable CO₂ absorbers to be used:


- Infinity ID CLIC Absorber 800+
- CLIC Absorber 800+
- Infinity ID CLIC Absorber Free
- CLIC Absorber Free

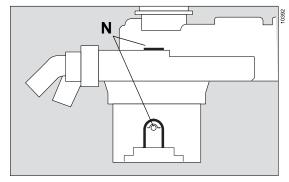

- 1 Remove the cleaning plate from the adapter.
- 2 Close the adapter so that it engages.
- 3 Pay attention to the correct orientation of the adapter: The two bypass symbols -O- (N) must align in a vertical line.
- 4 Screw on the CLIC adapter.

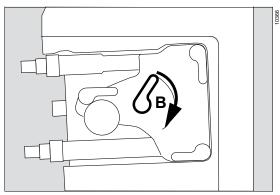


2 Lock the inserted breathing system: Turn the middle lever (B) clockwise about 120°.

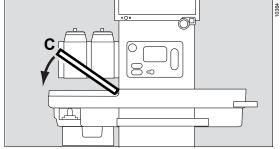

3 Set the breathing system cover (C) into place and click it into position.

Inserting the breathing system

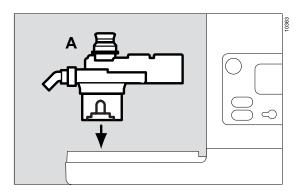


 Insert the assembled breathing system (A) vertically into the breathing system mount.

Mounting the flexible arm (optional) or the bag elbow

The breathing bag can be mounted either on the flexible arm or, using the bag elbow and a breathing hose, mounted directly on the breathing system.

Mounting the flexible arm

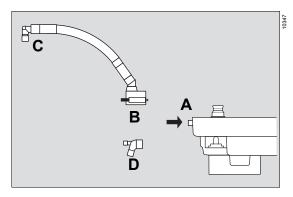

- 1 Attach the attachment piece of the arm (B) to the connection piece (A) on the breathing system and screw down tightly with the two knurled screws. Check that the arm is fixed securely!
- 2 Attach the elbow (C) to the end of the flexible
- 3 Align the flexible arm so that collisions with other mounted assemblies are prevented.

Mounting the bag elbow

 Attach the bag elbow (D) directly to the connection piece on the breathing system.

Selecting and connecting patient-specific accessories

Safety information

WARNING

Risk due to leakage in coaxial breathing circuit

Leakage in the inner hose of a coaxial breathing circuit may result in CO2 rebreathing or inadequate gas exchange. The device can only detect such leakage if a separate test with a coaxial test adapter is performed.

- Use a coaxial test adapter to check the inner hose for leakage, see page 125.
 Next, perform a leakage test on the entire breathing circuit.
- Monitor the measured gas concentrations during ventilation.

NOTE

This device is made without natural rubber latex.

To minimize the risk of exposure to latex, use breathing bags and breathing hoses that are not made with natural rubber latex.

Fitting the breathing circuit and the filters

WARNING

Risk due to particles and dust

In order to protect the patient from particles and dust, a filter must be used between the inspiratory limb of the breathing system and the patient.

Use a patient-side filter or a filter at the inspiratory port.

WARNING

Risk of infection

The breathing system may be contaminated with infectious agents/pathogens. This may have the following causes:

- No bacteria filters have been used at the Y-piece or at the expiratory port.
- The breathing system is being used for the first time.

Perform the following steps:

- Reprocess the breathing system before first use.
- Reprocess the breathing system if necessary.
- Carry out the reprocessing as described in the reprocessing instructions supplied with the product.
- To prevent future contamination, use bacteria filters close to the patient.

Perseus can be used with Infinity ID breathing hoses or conventional breathing hoses. If no leakage test has yet been performed after switching on the system, hose compliance and hose resistance will automatically be adopted when Infinity ID breathing hoses are connected.

WARNING

Risk due to occluded components in the breathing circuit

If filters, hoses, or endotracheal tubes are occluded and sample gas is being taken from between the patient and the occluded component, the sample gas flow may lead immediately to negative pressure in the lungs.

Ensure the following when ventilating pediatric patients and neonates:

- When fine pored filters are used, do not connect the sample line between the tube and the filter.
- When no filters but only an HME is in use, set the alarm limits for MV low and Paw high to values that are suitable for this patient so that an incipient occlusion can be detected.

NOTE

When applying tidal volumes in the range of the maximum or minimum values indicated for each patient category, use the smaller breathing bag and the smaller breathing circuit.


NOTE

In order to make use of the functionality of the Infinity ID breathing circuit, dispense with the inspiratory and expiratory microbial filters and fit the Y-piece with a filter instead. In cases which preclude use of a microbial filter at the Y-piece, the ID functionality of the Infinity ID breathing hoses cannot be used.

 Select suitable accessories for the respective patient category.

	Adults		Pediatric patients	Neonates
Tidal volume	>700 mL	301 to 700 mL	50 to 300 mL	<50 mL
Breathing bag	3 L	2 L	1 L	0.5 L
Breathing circuit	Adults		Pediatric	Neonates (or pediatric)
Filter	Filter or HMEF			Use filters with low resistance and compliance.

2 Assemble the breathing circuit and connect it to the Y-piece and to the inspiratory and expiratory ports on the breathing system, see page 70, "Table with recommended hose configurations".

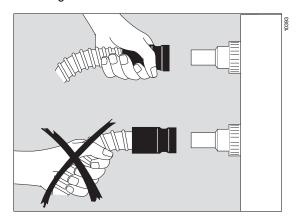


When attaching or removing the breathing hoses, always hold them at the connection sleeve and not at the spiral ribbing.

Table with recommended hose configurations

Filter or HMEF with connector for sample line between Y-piece and patient:


Pediatric patients

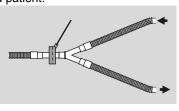

Or


Adults

Filter at inspiratory port and expiratory port, connector for sample line at Y-piece:

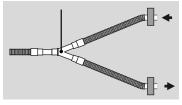

Neonates


Filter at inspiratory port and expiratory port. connector for sample line as close to patient as possible:

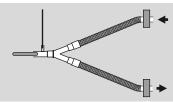


Lateral connectors for the sample line support CO₂ measurement and help to flush the dead space between Y-piece and hose adapter. If using an expiratory filter is not possible (e.g., due to an intrinsic PEEP due to air trapping), hygienically reprocess the device after use with this patient.

Carry out the reprocessing as described in the reprocessing instructions supplied with the product.

When fine pored filters (A) are used, do not connect the sample line between the tube and the filter:

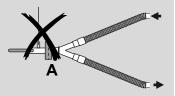


When no filters, but only an HME (B) is in use, set the alarm limits for MV low and Paw high to suitable values:

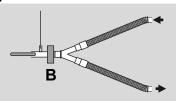

Breathing bag

Attaching the breathing bag

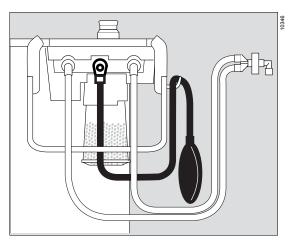

 Connect the breathing bag to the breathing bag hose using the connection nozzle. Attach the breathing bag hose to the bag elbow.

WARNING


Risk due to pinched breathing bag

If the breathing bag is pinched, excessive airway pressures or a lack of fresh gas may result.

When attaching and aligning, pay attention to the following:

- -The breathing bag is visible.
- -The breathing bag is not pinched.
- -The breathing bag can inflate freely.
- 2 Hang the breathing bag over the hose holder so that it hangs vertically downwards. To allow the bag to unfold freely, place it to the right over the hose holder so that it is not obstructed by cables or breathing hoses.

Attaching the breathing bag to the flexible arm (optional)

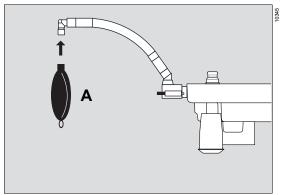

1 Attach the breathing bag (A) to the elbow.

Observing resistance and compliance

CAUTION

Risk due to accessory components in breathing circuit

When using additional components or hose configurations which deviate from the standard breathing circuit, the inspiratory and expiratory resistance values may be increased beyond standard requirements.

- When such configurations are used, the user must pay special attention to the measured values.
- After adding or replacing components such as breathing hoses, vaporizers, or soda lime, perform a leak test.

CAUTION

Risk due to misleading data

Exchanging breathing hoses, filters, vaporizers, or soda lime can change the determined leakage values or compliance values of the anesthesia machine and thus affect therapy.

- Perform a leakage test after replacing breathing hoses, particularly flex hoses, vaporizers or soda lime.
- Perform a leakage test after changing the length of flex hoses.

CAUTION

Risk due to changed hose lengths

Changed hose lengths can change resistance and compliance. Especially for neonates, this may cause increased or decreased ventilation volumes.

For neonates in particular, do not use flex hoses.

Accessories or accessory components such as filters affect the dead space, compliance, or resistance of the breathing circuit.

In addition, changes in resistance and compliance arise over time as a result of moisture in the breathing gas or deposits of secretions.

Calculating the resistance of the breathing system and connected accessories

The sum of the resistances in the inspiratory limb must not be less than –6.0 hPa (cmH2O). The sum of the resistances in the expiratory limb must not exceed 6.0 hPa (cmH2O).

Only include resistance data in calculations that were taken under the same flow conditions:

Patient category	Flow
Adults	30 L/min
Pediatric patients	15 L/min
Neonates	2.5 L/min

The following formula is used to calculate the resistance (R):

RInspiration =

RBreathing system_insp - RInsp hose - RBreathing bag hose

- RInsp filter (port) - RInsp filter (Y-piece)

RExpiration =

RBreathing system_exp + RExp hose + RExp filter (port) + RExp filter (Y-piece)

If necessary, take into consideration additional parts such as water traps or additional hoses. The specifications for the resistance of the breathing system can be found on page 249. The specifications for all other accessories can be found in the respective instructions for use.

In these instructions for use, the specifications for the resistance in the inspiratory limb are regarded as negative values. The resistance values given in the instructions for use for the accessories must therefore be subtracted from the inspiratory resistance of the breathing system.

Example of calculation: Adult hose with filter on Y-piece, no filters on the ports.

	Resistance [hPa (cmH2O)]			
	Inspiratory at 30 L/min		Expiratory at 30 L/min	
Breathing system with reus- able CO ₂ absorber and MX50115 soda lime dust filter	RBreathing_system_insp	-0.4	RBreathing_system_exp	2.6
Breathing hose	- Rinsp_hose	0.5	+ RExp_hose	0.5
Breathing bag hose	- RBreathing_bag_hose	0.3		
Filter at inspiratory port	- RInsp_filter(port)	0		
Filter at expiratory port			+ RExp_filter(port)	0
Filter at Y-piece	- RInsp_filter(Y-piece)	2	+ RExp_filter(Y-piece)	2
Result	Rinspiration	-3.2	RExpiration	5.1

Since Rinspiration is greater than –6 hPa (cmH2O) and Rexpiration is less than 6 hPa (cmH2O), this configuration may be used.

Depending on the breathing circuit in use, the connected accessories, and the resistance of the patient's airways, air trapping (incomplete expiration) may occur with some ventilation settings.

The effects are, for example, a reduced minute volume in pressure-controlled ventilation or higher mean airway and peak pressures in volume-controlled ventilation.

Air trapping can be prevented on the anesthesia machine by the following measures:

- Adjusting the respiratory rate and inspiratory time
- Changing the configuration of breathing hoses and accessories that carry gases during the expiratory time

It is the responsibility of the user of the medical device to select the most suitable remedial measure.

Connecting a non-rebreathing system (optional)

This connection is only possible with the *External fresh-gas outlet* option.

WARNING

Insufficient gas supply to the patient

Non-rebreathing systems are only suitable and intended for manual ventilation or spontaneous breathing and may only be connected to the external fresh-gas outlet.

When using a non-rebreathing system, ensure an adequate gas monitoring.

WARNING

Risk of excessively high airway pressure

Without a pressure-relief valve or breathing bag, airway pressure may become too high.

Only connect breathing systems with breathing bags or pressure-relief valves which comply with ISO 8835-2.

Strictly follow the instructions for use for the non-rebreathing system and the transfer hose.


To prevent contaminating the ambient air, connect the gas outlet of the breathing system with the inlet on the AGS. Use transfer hose with overpressure valve complying with ISO 8835-3.

WARNING

Risk of faulty gas delivery

O2 and CO2 and any anesthetic gases must also be monitored for non-rebreathing systems.


The sample line must be connected to the elbow and the water trap on the anesthesia machine.

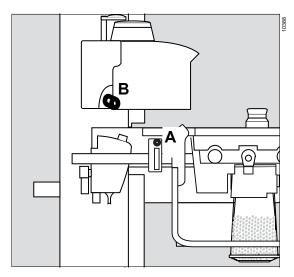

- 1 Tightly screw the sample line to the elbow of the non-rebreathing system and water trap (A). For elbows without port for sample line:
 - Fit the T-piece with T-piece filter directly to the elbow and tightly screw the sample line to the T-piece filter. The part numbers for the T-piece and the T-piece filter are listed in the list of accessories.

Or

- If necessary, use the gas sampling port of the filter on the Y-piece.
 Ensure the correct course of the sample line. Do not use adapters.
- 2 Connect the fresh-gas hose of the nonrebreathing system to the external fresh-gas outlet (B).

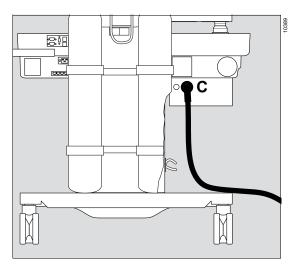

- 3 Remove the sealing plug from the inlet nozzle (C) of the AGS.
- 4 Use the transfer hose to connect the nonrebreathing system with the inlet nozzle on the AGS (C).

CAUTION

Risk due to open AGS inlet nozzle

After using a non-rebreathing system, reinsert the sealing plug into the inlet nozzle to prevent contamination of the ambient air with anesthetic gases.

Connecting or exchanging consumables

WARNING

Risk of high inspiratory CO₂ values

Use of soda lime over prolonged periods can increase the inspiratory CO₂ values.

Check the color of the soda lime regularly and exchange as needed, especially if the inspiratory CO₂ value increases unexpectedly.

Single-use CO₂ absorber

 Connect or replace the CLIC absorber or the Infinity ID CLIC absorber.

Strictly observe the instructions for use of the single-use CO₂ absorber.

WARNING

Risk of insufficient ventilation

If the CO₂ absorber is not correctly locked into place, system leakage may occur.

After mounting and replacing, make sure the CO₂ absorber is firmly locked into place.

Reusable CO₂ absorber

As an alternative to single-use CO2 absorbers, the reusable CO2 absorber may also be used.


CAUTION

Risk of chemical burns


Soda lime is caustic and is a strong irritant for eyes, skin, and airway.

Handle this absorption material carefully and do not spill it.

Dismounting and emptying

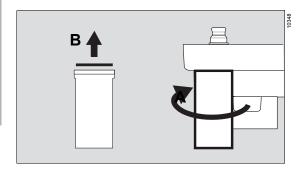


- 1 Turn the CO₂ absorber (A) clockwise and remove it from below
- 2 Remove and dispose of the disposable dust filter (B).
- 3 Empty used soda lime and dispose of according to the instructions for use.

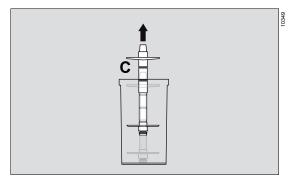

4 If it is necessary to clean the absorber insert (C), remove the absorber insert from the absorber container. Leave the inner and outer sealing rings on the absorber insert.

Filling and mounting

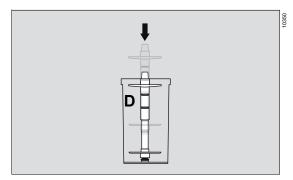

- After any cleaning, push the absorber insert back into the absorber container (D) completely.
- 2 Fill CO2 absorber with fresh soda lime to the upper mark.


Recommendation: Only use Drägersorb 800 Plus or Drägersorb Free.

WARNING

Risk of hypoventilation and incorrect gas measurement

Reuse of the disposable dust filter can increase the filter resistance and impair the ventilation function of the device.

If the reusable CO₂ absorber is used, always use a disposable dust filter. Replace the disposable dust filter with every change of soda lime.

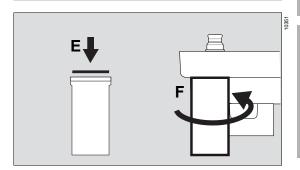

- 3 Insert a new disposable dust filter (E). Only use dust filters from the list of accessories. Only use undamaged filters.
- 4 Insert CO2 absorber (F) on the breathing system from below and then rotate it counterclockwise as far as it will go.

Strictly observe the instructions for use for Drägersorb 800 Plus or Drägersorb Free soda lime

Water trap

 Empty or replace the water trap according to its instructions for use.

WARNING

Risk due to full water trap

When the water trap is full, the gas measurement can fail and insufficient ventilation may occur.

Check the water level in the water trap regularly and empty it as necessary.

CAUTION

Risk due to gas measurement failure

Contaminants, damage, or overfilling of the water trap can influence gas measurement.

- Observe the instructions for use of the water trap.
- Always operate the device with a mounted water trap.

WARNING

Risk of infection

The water trap may contain infectious fluid.

- Proceed carefully when emptying and take protective measures if necessary.
- Follow the infection prevention policies and reprocessing regulations of the healthcare facility.

WARNING

Risk of incorrect measured gas values

Blocked water traps or blocked sample lines prevent correct gas measurement. As a result, incorrect gas values may be displayed.

- Use only Dräger sample lines.
- Replace the O-rings of the water trap mount every 2 years.

CAUTION

Risk of misleading data

Silicon can get into the measuring cuvette and disrupt measurement.

Do not spray the O-rings of the water trap holder with silicon spray.

WARNING

Risk of fire

Silicon or aerosol residues in the water trap can cause fires.

Do not spray the O-rings of the water trap holder with silicon spray.

Connecting the sample line

WARNING

Danger due to incorrectly connected sample line


If the sample line is connected to the incorrect connections, e.g., infusion pump connections, liquids instead of sample gas can be sucked in. As a result, the gas measurement will not display correct values.

When attaching the sample line, make sure not to attach it incorrectly.

WARNING

Risk due to leakage from improperly connected or damaged sample line

Connect the sample line correctly, otherwise faulty gas measurements may result.

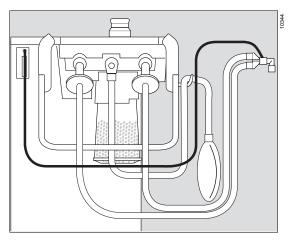

 Screw on the sample line at the Y-piece, HME filter, or hose adapter and at the water trap.
 Ensure the correct course of the sample line.
 Do not use adapters.

WARNING

Risk of incorrect measured gas values

Blocked water traps or blocked sample lines prevent correct gas measurement. As a result, incorrect gas values may be displayed.

- Use only Dräger sample lines.
- Replace the O-rings of the water trap mount every 2 years.

Getting started

Turning on Perseus	80
Checking the device configuration	81
Checking the operational readiness	81
Emergency start-up	82
Electronically controlled gas mixer Mechanically controlled gas mixer	

Turning on Perseus

Prerequisite: The device has been reprocessed and assembled to be ready for operation in accordance with the reprocessing instructions (see page 46).

CAUTION

Risk of device malfunction

Condensation may form when the device is moved from a cold storage location to a warm environment.

To prevent condensate formation and resulting failures of electrical components, do not turn on the device after abrupt temperature changes for 1 to 2 hours

WARNING

Risk of explosion and fire

Do not set the device into operation if oxygen leakage is suspected in the medical device or its vicinity.

Stop all oxygen supplies and contact service personnel.

WARNING

Risk due to unbraked device

An unbraked device may accidentally move during operation.


Prevent this by actuating the central brake or the castor brakes and check their function.

WARNING

Risk of device malfunction

Some safety systems are only checked during start-up.

Restart Perseus at least once per month to maintain proper functionality.

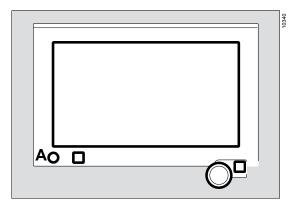

- 1 Connect the device to the mains power supply.
- 2 Set the main switch to position I.
- 3 Turn on Perseus: Press the () (A) button.

If the battery charge is sufficient, Perseus will also start without the power plug being plugged in as long as the main switch is not set to position $\boldsymbol{\theta}$.

When the starting procedure is completed, the **Standby** page is displayed.

Checking the device configuration

Perseus can be customized to suit the requirements of the user. Settings for the following features are possible:


Start settings for the ventilation

- Alarm limits
- General device behavior

For further information about the configuration, see page 148.

Checking the operational readiness

On the **Standby** page, the readiness for operation resulting from the system test is indicated by color (A).

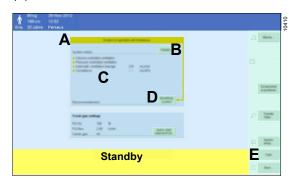

Color	Meaning
Green	System is fully operational.
Yellow	System is operational with limitations. There are functional restrictions. Take further measures to ensure patient safety (e.g., external monitoring).
Red	System is not operational. Call service personnel.
Gray	The system has not been tested.

If the system is not fully operational, the most important irregularities (C) are displayed and a specific test (D) is recommended as a remedy where possible.

Additionally, the current system leakage is displayed in area (C).


To view details regarding the status of the system, touch the *Details...* button (B) or the *Tests...* button (E), see page 121.

Dräger recommends performing the system test every 24 hours. Otherwise, it will not be possible to ensure that the device is functional.

Emergency start-up

Electronically controlled gas mixer


1 Adjust the APL valve.

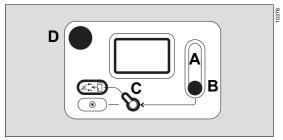

- 2 Set O2 switch (C) to the Add. O2 position.
- 3 Open the flow control valve (B) and set the desired O2 flow. If necessary, press the O2+ (D) key to quickly fill the breathing bag.
- 4 Monitor the set flow on the O₂ flowmeter (A).
- **5** Set the anesthetic gas concentration at the vaporizer.
- 6 Manually ventilate the patient.
- 7 Switch on the device.
- 8 As soon as the **Standby** page is displayed, start the therapy, see page 87.
- 9 Set the O₂ switch (C) to **Aux. O₂** to stop the increased fresh-gas flow.
- 10 Close the flow control valve (B).

Mechanically controlled gas mixer

1 Adjust the APL valve.

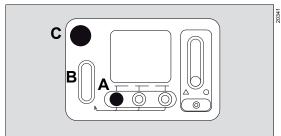

- 2 Open the flow control valve (A) and set the desired O2 flow. If required, press the O2+ key (C) to fill the breathing bag quickly.
- 3 Monitor the set flow on the total flow tube (B).
- **4** Set the anesthetic gas concentration at the vaporizer.
- **5** Manually ventilate the patient.
- 6 Switch on the device.
- 7 As soon as the **Standby** page is displayed, start the therapy, see page 87.

Operation

Safety information	85	Customizing the current view	
Starting the therapy	87	Using loops Displaying mini-trends	
Loading the patient data	87	Displaying alarm limits and units of	101
Checking the patient data	88	measurement	107
Setting and starting the therapy	88	Adjusting the sweep speed and the scale	107
Useful tips	89	Changing the color scheme and the screen	
Adjusting the therapy	90	brightness	108
Setting the APL valve	90	Displaying additional data	108
Adjusting the fresh-gas flow	91	Viewing current measured values	108
Using the O2 flush	91	Logbook	
Using the vaporizer	91	Trends	
Changing the ventilation mode	92	Displaying installed options	
Synchronizing the breaths	92	Displaying an overview of the accessories and	
Ventilating pediatric patients and neonates	92	consumptions	109
Special forms of therapy	93	Setting the sound volume	110
Using non-rebreathing systems	93	Setting alarm tone volume and breathing	
Activating or deactivating the CBM mode	93	sound volume	110
Only gas concentration measurement/use of		Adjusting the alarms	110
Monitoring	94		
Pausing the therapy/using Pause (optional)	94	Setting the alarm limits	
Maneuvers (optional)	95	Activating or deactivating CO2 alarms Automatic xMAC monitoring	
Inspiration hold, expiration hold	95	Deactivation of automatic xMAC monitoring	
One-step recruitment	95 96	-	
Multi-step recruitment	97	Changing the patient data	113
Alarm behavior during the maneuvers	99	From a setting and a to	440
·		Exporting data	
Using fields with special functions	100	General information	
Breathing gas measurement and xMAC		Exporting the screen contents	
display (MAC multiples)	100	Exporting trends and data	113
Prediction of the anesthetic gas concentration (optional)	101	Other settings	114
Prediction of the inspiratory O2 concentration		Switching the breathing system warmer on or	
(optional)	102	off	
Econometer (optional)	103	Creating additional logbook entries	
Stopwatch		Resetting user-specific settings	115
Countdown		Ending the therapy	116
Volumeter		Switching to the Standby mode	116
Low-flow wizard (optional)	104	Using O2 insufflation	
Organizing the screen display			
Available views	105		
Changing the current view	105		

Operation

Change of patient	118
Cleaning and disinfecting the workstation	118
Checking or replacing consumables	118
Checking the device	118
Flushing the breathing system	118

Safety information

WARNING

Risk of strangulation

Negligent placement of hoses, cables, sample lines, and similar device components can endanger the patient.

Use particular caution when establishing connections to the patient.

CAUTION

Risk due to falsified gas measurement values and failure of the gas measurement

Aerosols can damage the patient gas measurement module and the membrane of the water trap. Propellants can falsify the gas measurement values.

Do not use medication nebulizers.

WARNING

Risk of faulty gas delivery

O2 and CO2 and anesthetic gases must be monitored every time breathing systems are used.

Ensure adequate gas monitoring.

WARNING

Risk of misinterpretation

Misdiagnosis or misinterpretation of the measured values or other parameters can endanger the patient.

Do not make therapeutic decisions based solely on individual measured values and monitoring parameters. Therapeutic decisions must be made solely by the user.

CAUTION

Risk of injury

Looking directly into the LEDs of the workplace illumination and of the illumination unit on the vaporizer plug-in connection can cause damage to the retina.

Do not look directly into the LEDs. Make sure that the patient is not dazzled by the LEDs.

CAUTION

Risk due to unsuitable illumination

If illumination without neutral colors is used during the medical examination of the patient, this may result in, e.g., misinterpretation of skin coloring.

- Do not use the device's workplace light during examinations.
- For examinations, use an examination light conforming to IEC 60601-2-41.

WARNING

Risk of misinterpretation

Measured gas values and waveforms such as the CO2 waveform are determined on the basis of the composition of the sample gas. The composition of the sample gas is affected by many factors and their interactions, especially in patients with low body weight. This may result in biased measurements or waveforms.

The following factors affect the sample gas measurement:

- Dead space
- Patient resistance, patient compliance
- Type of surgical procedure
- Gas sampling site
- Breathing circuit, filter, sample line, tube
- Ventilation settings and the resulting ventilation
- Leakage
- Spontaneous breathing
- Cardiogenic oscillations

Observe the following:

- Do not make therapeutic decisions based solely on individual measured values or parameters.
- If possible, minimize the effects of the factors described above, e.g., take the sample gas from a gas sampling site close to the patient, minimize leakage, adjust the ventilation settings.

CAUTION

Danger due to inaccurate gas measurement values

Due to the technical characteristics of gas measurements, the measured gas values might be inaccurate at high respiratory rates and at certain I:E ratios.

Observe the technical data.

WARNING

Risk of infection

Disease-causing germs can be transmitted between patients via the circuit plug. If a used Y-piece or filter is fitted to the circuit plug, and then later a reprocessed component is fitted (e.g., during a leakage test), the new component may become contaminated.

Do not fit hoses that have already been used and that have attached filters or Y-pieces to the circuit plug. Instead, hang them over the handle on the left-hand side of the device.

WARNING

Risk due to occluded components in the breathing circuit

If filters, hoses, or endotracheal tubes are occluded and sample gas is being taken from between the patient and the occluded component, the sample gas flow may lead immediately to negative pressure in the lungs.

Ensure the following when ventilating pediatric patients and neonates:

- When fine pored filters are used, do not connect the sample line between the tube and the filter.
- When no filters but only an HME is in use, set the alarm limits for MV low and Paw high to values that are suitable for this patient so that an incipient occlusion can be detected.

CAUTION

Danger from an unintentional therapy start

If Perseus is in the Standby mode and the flow control valves are open, ventilation starts in the Man / Spon mode.

Take care that this does not happen unintentionally.

WARNING


Risk of fire

Components such as power supply units that heat up are unable to cool down in enclosed storage locations and may cause a fire.

Do not keep components that heat up in the drawers or under the column cover.

Starting the therapy

Prerequisite: The device is in Standby mode.

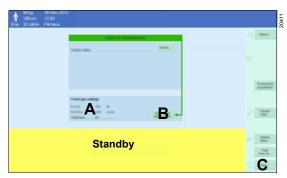

There are 2 possibilities for starting the therapy:

Quick start with Manual / Spontaneous mode

 To perform a quick start with the displayed fresh-gas settings (A), touch the Quick start Man/Spon button (B).

Or

Normal start with customized settings

- 1 Touch the **Start...** button (C).
- 2 Adjust the patient data and ventilation settings.

Loading the patient data



There are two possibilities for loading patient data:

- A Defining a new case
- **B** Continuing a case

Depending on the selection, the patient data in area (C), e.g., weight or age, are filled in automatically.

Defining a new case

Select to start the therapy for a new patient:

 Depending on the patient category needed, touch the New adult, New pediatric, or New neonate button (A).

The ventilation parameters and alarm limits are set to the configured start settings, see page 150. The set value for *Ti* is set automatically based on *RR* so that the *I:E* ratio set in the start settings are the result, see 158.

Continuing a case

Select to start the therapy with the settings of the last case:

• Touch the **Continue case** button (B).

The ventilation parameters and alarm limits are adopted from the preceding case.

After turning on the device, no previous case is available. Perseus then starts with the configured start settings.

Checking the patient data

WARNING

Risk due to incorrect settings

Different standard alarm limits or therapy settings might be configured for medical devices within the same area. The user must observe the following:

- Make sure that the set values are suitable for the patient.
- Make sure that the alarm system is neither rendered useless by setting extreme values for the alarm limits nor deactivated by switching off the alarms.
- Check the start settings for alarms and alarm settings each time the ventilation mode is changed.
- Only switch off alarms if the safety of the patient is not affected.

CAUTION

Risk due to incorrect setting for patient age

Incorrectly setting the patient age can lead to inappropriate xMAC values and therefore to an inappropriate anesthetic gas delivery.

Always set patient age correctly.



Adjust patient data (A) if necessary.

When these data are adjusted, appropriate therapy settings are suggested e.g., for tidal volume, respiratory rate, alarm limits. For more information, see page 269.

Setting and starting the therapy



1 Select a ventilation mode (A) as needed.

The following ventilation modes are available:

- Man / Spon
- PSV (optional)
- PC
- APRV (optional)
- VC-AF
- VC

The following operation modes are also available:

- Ext. FGO (optional)
- Pause

For additional information about the ventilation modes, see page 258.

2 Set the fresh-gas delivery.

Perseus is equipped with a minimum O2 delivery that ensures that a minimum quantity of oxygen is delivered, see page 268.

Electronically controlled gas mixer:

Select the carrier gas (B).
 Set the O2 concentration (D) and fresh-gas flow (C).

To prevent the soda lime from drying out at an increased rate, the electronically controlled gas mixer limits the setting range for the fresh-gas flow *FG flow*. This restriction depends on the set tidal volume *VT* and the set respiratory rate *RR*.

Mechanically controlled gas mixer:

WARNING

Risk of incorrect gas delivery

Use the flow values on the status display as the primary source of information. Always compare these values with the values on the total flow tube.

Only use the total flow tube as the primary source of information in the following cases:

- The status display has failed.
- The values on the status display do not match the values on the total flow tube.
 - Use the flow control valves to adjust the fresh-gas flows.
 Use the total flow tube to check the total flow set, see page 32.

CAUTION

Risk of patient injury

Unsuitable soda lime can result in disintegration products from the anesthesia gases.

Use suitable soda lime such as Drägersorb Free.

WARNING

Risk of patient injury

Under certain conditions, acetone can accumulate in the body in patients under anesthesia. As a result, the patient could be put at risk.

Do not perform low-flow anesthesia on patients with ketoacidosis or patients under the influence of alcohol.

CAUTION

Risk of patient injury

The use of minimum-flow or low-flow settings can lead to accumulation of metabolic by-products in the breathing system.

When using minimum-flow setings or low-flow settings follow the recommendations of professional societies (e.g., regular flushing of the breathing system).

- **3** Adjust the ventilation settings (E).
- **4** To start therapy, press the rotary knob. A signal tone is emitted when therapy starts.

Useful tips

Quick start with Man / Spon mode

		Mechanically con- trolled gas mixer	
•	In the Standby mode, touch Quick	•	Open the flow control valves for fresh
	start Man/Spon.		gas.

Opening the start dialog in Standby mode

Touch the screen in the monitoring area.

Or

Squeeze the breathing bag.

Starting when time is limited

When time is limited, it is possible to bypass the adjustment of the patient data and the ventilation settings. Start the therapy as follows:

- 1 Touch the screen.
- 2 Check the displayed start values.
- **3** To start therapy, press the rotary knob.
- **4** Adjust the patient data and ventilation settings as soon as possible.

Adjusting the therapy

Setting the APL valve

The pressure limitation set with the APL valve only takes effect during manual ventilation or spontaneous breathing.

WARNING


Risk of excessively high airway pressures

If the ventilator fails, the device switches into the *Man / Spon* ventilation mode.

The APL valve should also be set to a pressure limitation value suitable for the patient when using automatic ventilation modes. If the ventilator fails, ventilate the patient manually.


 The selection between manual ventilation (Man) and spontaneous breathing (Spont) is made at the APL pressure limitation valve, see page 18.

Manual ventilation

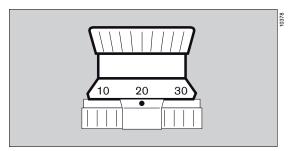

 Set the APL valve to the desired maximum airway pressure.

The patient can be manually ventilated with the breathing bag. The pressure is limited by the set value.

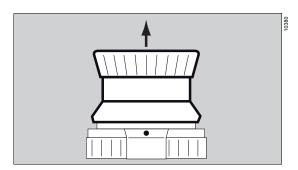

In the *Manual / Spontaneous* mode, lifting the valve relieves pressure from the breathing system.

Spontaneous breathing



• Turn the APL valve counterclockwise as far as it will go.

The points are vertically aligned and the valve is raised

The pressure limitation is canceled and the valve is open for free spontaneous breathing.

Adjusting the fresh-gas flow

When a very high fresh-gas flow is reduced to a very low value, the gas concentration in the entire breathing system is optimized and changed as smoothly as possible.

Electronically controlled gas mixer

The setting range for the fresh-gas flow is temporarily restricted to values greater than 3 L/min. After the new fresh-gas flow is confirmed, it is not altered abruptly but, rather, continuously over a few seconds until the new target value is reached. After the reduced flow is reached, the full setting range is available once more, so that values below 3 L/min can also be set.

Mechanically controlled gas mixer

After a sharp reduction of the fresh-gas flow, the internal gas recirculation is temporarily ramped up without significantly altering the ventilation. This may result in visible effects in the displayed waveforms.

This is indicated in the waveforms by the messages *Preparing low flow phase* and *Measurement temporarily inaccurate*.

Using the O₂ flush

The O2 flush is used for flushing and quickly filling the breathing system and breathing bag with oxygen. The vaporizers are bypassed for this.

 Press the O2+ key. O2 flows for as long as the key is held down.

The gas concentration can change abruptly when the O₂ flush is used.

Using the O₂ flush has the following effects:

- In manual ventilation, pressing the O2 flush results in a rapid rise in pressure to the APL level.
- In automatic ventilation, permanently pressing the O2 flush may result in a slight rise of the PEEP level. However, this rise has no effect on the peak pressure.

Using the vaporizer

WARNING

Risk due to incorrect anesthetic agent delivery

If the vaporizer is filled with the wrong anesthetic agent or if it is not filled sufficiently, incorrect anesthetic gas concentrations or concentrations that are too low can occur as a result.

- Follow the instructions for use for the vaporizer.
- Compare the color coding and labeling on the vaporizer used with the anesthetic agent bottle and the anesthetic agent indicated on the screen.
- Operate the vaporizer according to its instructions for use.

When the Vapor View option is installed, the control dial and sight glass on the Dräger Vapor 3000 and D-Vapor 3000 are illuminated:

Control dial position	Illumination
T to < 0	Off
0	Medium
>0	Bright


WARNING

Risk due to faulty Vapor View option

If the Vapor View option is faulty, the following information could be displayed incorrectly on the screen: Anesthetic gas type, filling level, control dial position, prediction of anesthetic gas concentration

- Do not rely solely on the display on the screen, but always pay attention to the settings and values at the vaporizer as well.
- Do not make therapy decisions based solely on a displayed prediction.

Changing the ventilation mode

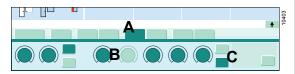

 In the therapy bar, touch the tab (A) of the new ventilation mode.


When the ventilation mode is changed, the start settings are adopted from the parameters of the previous ventilation mode and the patient data, or they are sensibly derived.

In addition, the alarm settings are adjusted to appropriate values, see page 143.

- 2 Adjust the therapy as needed with the therapy controls (B) or buttons (C).
- 3 Activate the ventilation mode by pressing the rotary knob. A signal tone is emitted when the mode is changed.

Synchronizing the breaths

Turning on the synchronization activates the set pressure support, for example, see page 258.

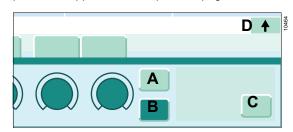

- 1 Turn the synchronization on or off with the buttons **Sync. on** (A) or **Sync. off** (B).
- 2 As needed, the expanded therapy bar can be displayed using the *More* buttons (C) or (D); here, additional parameters (*Trigger*, Δ*Psupp*, etc.) can be adapted to the patient.

WARNING

Risk of insufficient ventilation

In ventilation modes in which breaths are to be triggered only by the patient (e.g., *PSV*), adverse settings or sensor failure can lead to insufficient ventilation.

Set the respiratory rate to a suitable value so that a minimum ventilation of the patient is maintained.

WARNING

Risk of insufficient ventilation

The display of the spontaneous minute volume *MVspon* shows the volume that results from the patient's own breathing and machine support. If machine support is triggered by small tidal volumes of the patient, a large part of *MVspon* is achieved by machine support and not by the patient's own breathing. In this case, *MVspon* shows a high value although the actual spontaneous minute volume is very low.

Do not base therapy decisions solely on the value displayed for *MVspon*.

Ventilating pediatric patients and neonates

For tidal volumes below 300 ml:

 Use suitable ventilation accessories, see chapter "Selecting and connecting patientspecific accessories" starting on page 67.

For tidal volumes below 20 mL or when using unblocked tubes:

 Use pressure-controlled ventilation, see table "Ventilation modes and effective parameters" starting on page 260.

Special forms of therapy

Using non-rebreathing systems

Only available with the option *External fresh-gas* outlet.

Prerequisite: The non-rebreathing system is connected, see page 74.

CAUTION

Risk of gas contamination

The extracted sample gas is also returned to the internal breathing system during operation with an external fresh-gas outlet if the Perseus sample gas measurement is used.

Use breathing circuit to close the breathing system or set the APL valve to spontaneous breathing. Flush the breathing system each time patients or anesthetic gas are changed!

Conducting the fresh gas to the external outlet

- 1 Start *Ext. FGO* operating mode.
- 2 Adjust the fresh-gas delivery. Set vaporizer as required.

Strictly observe the instructions for use of the non-rebreathing system.

Activating or deactivating the CBM mode

The CBM mode allows patient monitoring without unnecessary alarms during extracorporeal oxygenation of the patient by a heart-lung machine.

Properties of CBM mode:

- All gas concentrations are measured independent of the respiratory phase.
- The CO₂ apnea and pressure apnea alarms are inactive.


The CBM mode can be used in all active ventilation modes.

When ventilation modes are changed, the CBM mode remains active. Changing to the **Standby** mode deactivates the CBM mode.

Deactivating the CBM mode activates the apnea alarms.

Activating

- 1 Open the *Alarms* dialog window.
- 2 Touch the Settings tab (A).

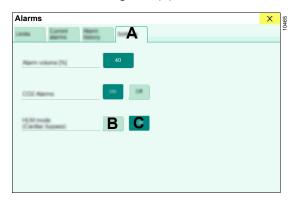

3 For Cardiac bypass mode (CBM), touch the button On (B).

Deactivating

- For Cardiac bypass mode (CBM), touch the button Off (C).
- In the main menu bar, touch the Exit CBM button.

Only gas concentration measurement/use of *Monitoring*

No ventilation in **Monitoring** operation mode. Gas delivery is also stopped with electronically controlled gas mixers. The gas concentration measurement remains active and waiting for respiratory phases.

Use this operating mode for regional anesthesia, for example.

Prerequisite: The device is in **Standby** mode or **Man / Spon** mode.

Activating

Start *Monitoring* operating mode.

Deactivating

• After **Standby** or change to ventilation mode.

For more information, see page 162.

Pausing the therapy/using Pause (optional)

In **Pause** mode, ventilation is stopped. Gas delivery is also stopped with electronically controlled gas mixers. The gas concentration measurement remains active and waiting for respiratory phases.

This mode is equipped with an adjustable therapy control *Timer*, depending on the patient category. When the set time has elapsed, this emits an alarm as reminder that ventilation should be resumed. There is no automatic resumption of ventilation. Setting the therapy control *Timer* to "Off" deactivates the alarm. The total elapsed time is also displayed.

Use this operating mode for regional anesthesia or short pauses in therapy, such as disconnection or intubation, to reduce contamination of the ambient air with anesthetic gases due to an open Y-piece.

Activating

- 1 Start **Pause** operating mode.
- 2 Adjust the therapy control *Timer* if necessary.

Returning to the previous mode

- 1 Touch the **Resume ventilation** button.
- 2 Confirm the ventilation mode.

For more information, see page 155 and page 162.

Maneuvers (optional)

Perseus has various maneuvers for lung recruitment. During a maneuver, various data concerning the lung mechanics are displayed so that the user can assess the progress of the maneuver.


The reminder function (option) reminds the user to perform a maneuver. A reminder is issued after the first change to a ventilation mode with medium or high breathing support and also at settable intervals after the ending of a maneuver. The **Consider recruitment** message is displayed in the waveforms for flow and pressure.

The use of lung recruitment maneuvers is the sole responsibility of the user.

Dräger recommends that the patient's hemodynamics always be monitored while the maneuvers are being performed.

Available maneuvers

- Insp./Exp. hold
- One-step recruitment
- Multi-step recruitment
- 1 Open the **Procedures** dialog window.
- 2 Touch the tab for the required maneuver.

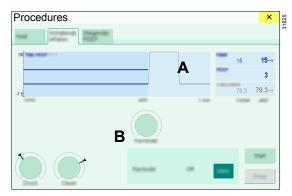


Relevant waveforms, trends, or measured values are displayed in area (A) of the dialog window. The parameters for the particular maneuver are displayed and set in area (B).

Inspiration hold, expiration hold

Perseus provides functions whereby a breath can be initiated or extended, or expiration can be extended.

This can be useful in situations where the lungs of the patient are not supposed to move, e.g., during use of imaging techniques.

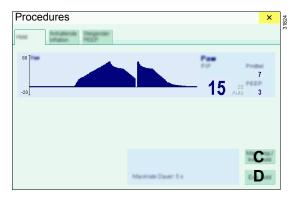

Manual inspiration/inspiration hold

This maneuver is available in the volume-controlled modes, the pressure-controlled modes, and in the *PSV* mode, and offers the following options:

- In the expiratory phase between 2 mandatory breaths, a breath can be manually triggered and held. The ventilation pattern of the manually triggered breath corresponds to the ventilation pattern for the active ventilation mode.
- A mandatory breath can be extended.

Manually triggering a breath

Touch the *Man. insp./Insp. hold* button (C) briefly.

Manually extending a breath

 Touch the *Man. insp./Insp. hold* button (C) and keep it pressed for the desired time.

Perseus will trigger an extended breath or will extend an already triggered automatic breath.

The breath is ended automatically:

- After a maximum of 40 seconds in the Adult patient category
- After a maximum of 30 seconds in the **Ped** patient category
- After a maximum of 5 seconds in the **Neo** patient category

Holding or extending the expiration

 Touch the Exp. hold button (D) and keep it pressed for the desired time.


Perseus will extend the expiration and delay the next breath.

The expiration is ended automatically:

- After a maximum of 40 seconds in the Adult patient category
- After a maximum of 30 seconds in the *Ped* patient category
- After a maximum of 5 seconds in the **Neo** patient category

One-step recruitment

This maneuver applies a set pressure for a specific duration and thus enables, e.g., an extended breath with adjustable pressure.

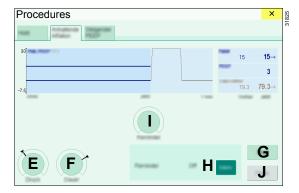

- 1 Using the therapy controls, set the **Pressure** (E) and the **Duration** (F).
- 2 If necessary, touch the *More* button (H) and set a reminder with the therapy control (I).
- 3 Touch the **Start** button (G) and confirm.


The rise in pressure from the **PEEP** level to the set pressure level and the decrease in pressure at the end of the maneuver take place at 20 hPa/s.

The maneuver ends automatically after the time has elapsed. To cancel the maneuver, touch the *Cancel* button (J) and confirm. If required, the *PEEP* can be adjusted before confirming.

The *Paw high* alarm limit is checked at the start of the maneuver. If the alarm limit is too low, it is changed to a value 5 hPa above the set pressure. After the maneuver, the alarm limit is reset to its original value.

Multi-step recruitment

This maneuver applies a sequence of pressure-controlled mandatory breaths at variable pressure levels. The inspiratory pressure and expiratory pressure are increased in steps and then reduced again to the starting levels. For parameters such as *RR* and *Ti*, the values from the previous ventilation mode are used. If there is no value available for the *Slope* parameter, a value of 0.2 seconds is used.

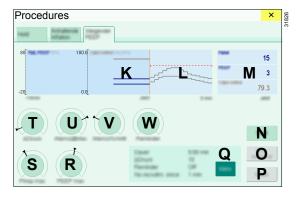

- K Trends of the measured values of pressure, compliance, or tidal volume
- L Predicted trends of inspiratory pressure and expiratory pressure during the maneuver ("what if..." function)

The exact progress of the maneuver is influenced by the following values:

- Parameters entered in this dialog
- Values set for RR and PEEP in the therapy har
- M Values of pressure, compliance, or tidal volume
- Using the *Pinsp max* (S) and *PEEP max* (R) therapy controls, set the pressures that are to be reached during the maneuver.

2 If required, touch the *More* button (Q) to display the following settings:

T ΔPressure

Pressure difference between *Pinsp* and *PEEP* with which the ventilation is performed during the maneuver. After the set maximum pressures for *PEEP max* or *Pinsp max* is reached, the pressure difference is adapted accordingly so that the set values for both *PEEP max* and *Pinsp max* can be reached.

U Breaths@Max

Number of breaths at the **Pinsp max** pressure level

V Breaths/Step

Number of breaths after which the next pressure

step begins

W Reminder

Time after which a reminder for a new maneuver is issued

3 Touch the **Start** button (N) and confirm.

During the maneuver, the inspiratory pressure is automatically increased in steps to the set maximum *Pinsp max* value and the expiratory pressure is increased to the maximum *PEEP max* value. Both pressures are held at the highest pressure level for a certain number of breaths *Breaths@Max* and then reduced again in steps. The level of the pressure rise and the pressure decrease is dependent on the selected patient category. The appropriate settings are made in the system setup.

The *Paw high* alarm limit is checked at the start of the maneuver. If the alarm limit is too low, it is changed to a value 5 hPa above the set pressure. After the maneuver, the alarm limit is reset to its original value.

If the pressure is not to be increased any further during the maneuver, touch the **Descend now** button (O) and confirm. The pressure will then be reduced again in steps. The duration of the maneuver is reduced as a result.

To cancel the maneuver, touch the *Cancel* button (P) and confirm. If required, the *PEEP* can be adjusted before confirming. This can be useful, e.g., when the desired effect has already been achieved during the maneuver. The ventilation will now continue immediately with the previous settings and the adjusted *PEEP*.

Due to the technically related delay between the set and the measured airway pressures, it may happen that the current and the last measured pressures differ by up to one step level.

The cursor can be used to display the previous measured values. To do this, touch screen area (K) and move the cursor with the rotary knob. To facilitate reading the data, screen area (K) will temporarily not be updated. To show the current measured values again, touch area (L).

Alarm behavior during the maneuvers

During the maneuvers, some alarms are adapted as follows:

	Maneuver				
Alarm	Insp./Ex	rp. hold	One step respect	Multi atau na amit	
	Man. insp./Insp. hold	Exp. hold	One-step recruit- ment	Multi-step recruit- ment	
Pressure alarms:					
Airway pressure high			Alarm limit is adjusted to Pressure + 5	Alarm limit is adjusted to Pinsp max + 5	
Airway press. continuously high	Stopped ¹⁾		When set to Auto , the alarm limit is automatically adjusted to Pressure +3	When set to Auto , the alarm limit is automatically adjusted to half way between PEEP max and Pinsp max. Example: PEEP max = 20 hPa Pinsp max = 40 hPa Adjusted alarm limit = 30 hPa	
Airway pressure not achieved	Stopped ¹⁾	Stopped ¹⁾	Stopped ¹⁾		
Volume alarms:			•		
Inspiratory tidal volume high	Suppressed in VC and VC - AF		Suppressed in VC and VC - AF	Suppressed in VC and VC - AF	
Minute volume low	linute volume low Displayed no sooner than 90 seconds after the start of the maneuver.				
Apnea alarms:					
Apnea (no pres- sure)					
Apnea (no flow)	Stopped ¹⁾				
Apnea (no CO ₂)					
Apnea					

¹⁾ Existing alarms are preserved. No new alarms will be issued.

Some technical alarms stop the maneuver automatically, e.g., alarms caused by a sensor failure.

Using fields with special functions

Breathing gas measurement and xMAC display (MAC multiples)



The MAC value is a simple navigation aid for anesthetic agent delivery.

Perseus displays the measured inspiratory and expiratory values for *O2*, *N2O*, and anesthetic gases, and the *xMAC* in the monitoring area. Measured nitrous oxide concentration or anesthetic gas concentration is only displayed when it is not zero.

The **xMAC** is the MAC multiple calculated from the current expiratory measured values and the age-dependent MAC values. If no respiratory phases are detected, expiratory values and **xMAC** cannot be displayed.

The integrated MAC algorithm is based on the MAC values shown in the following table. These values are guiding values only. The binding values are specified on the package information leaflet of the anesthetic agent.

The MAC values depend on the age of the patient. The values specified in the table apply to a patient age of 40 years.

	1 MAC corresponds to: (in 100 % O2)
Halothane	0.77 Vol%
Enflurane	1.7 Vol%
Isoflurane	1.15 Vol%
Desflurane	6.0 Vol%
Sevoflurane	2.1 Vol%
N2O	105 Vol%

The age-corrected MAC values are calculated according to the equation of W.W. Mapleson (British Journal of Anaesthesia 1996, pp. 179-185).

The equation applies to patients older than 1 year. $MACage corrected = MAC^{1} \times 10^{(-0.00269 \times (age -40))}$

^{1) 40} years

For gas mixtures, the respective multiples for **N2O** and anesthetic agents are added according to the following equation.

Example

exp. Iso. = 0.65 Vol%; exp. **N2O** = 69 %; Age = 32 years

MACage corrected for Iso.: $MAC^{1)} = 1.21 \text{ Vol}\%$ MACage corrected for **N2O**: $MAC^{1)} = 110 \text{ Vol}\%$

xMAC = 0.54 + 0.63 = 1.2

The influence of other drugs (opioids or intravenous hypnotics) is not considered in the **xMAC** calculation.

1) 32 years

Prediction of the anesthetic gas concentration (optional)

Perseus offers the option of displaying predictive concentrations of anesthetic agent in the breathing gas. This requires the Vapor View option and Dräger Vapor 3000 or D-Vapor 3000.

The anesthetic gas prediction is appropriate for patients with the following data only:

Height: 150 to 200 cm

(59.1 to 78.7 in)

Weight: 40 to 140 kg

(88 to 308 lbs)

Age: 18 to 90 years


The anesthetic gas prediction is not suitable for the following patients:

- Alcohol dependents
- Greatly overweight patients
- Patients with ASA ≥IV
- Patients with severe circulatory disorders or a cardiopulmonary bypass

These restrictions in the patient data result from the scientific models all calculations are based on.

Using the anesthetic gas prediction

Prerequisite: The **Agent prediction** parameter field must be configured in the monitoring area and respiratory phases must be detected.

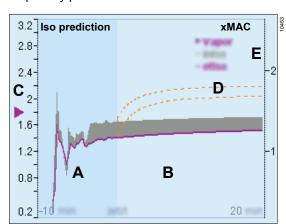

The following is displayed:

- A Section for the trend of inspiratory and expiratory anesthetic gas concentration
- **B** Section for the prediction
- C Current vaporizer setting
- **D** Prediction curve ("what if..." function)
- E xMAC scale

WARNING

Risk of incorrect therapy settings

The anesthetic gas prediction is based on mathematical models and does not provide individual patient values.

Do not make therapy decisions based solely on the displayed anesthetic gas prediction.

The anesthetic gas concentration of the currently delivered anesthetic agent (or of the measured primary anesthetic agent if delivery is switched off) is displayed as a color diagram. The corresponding expiratory concentration is shown by the color coding for the anesthetic agent.


The previously measured concentrations are displayed in the left-hand section (A) and the predicted concentrations in the right-hand section (B). The predicted concentrations are dependent on the vaporizer setting, the set fresh-gas flow, and on various measurements.

Electronically controlled gas mixer only: During the setting procedure, 2 dotted lines (D) are displayed for the fresh-gas flow (the so-called "what if..." function). When the value is altered with the therapy control, the lines move and show a preview of the predicted concentration.

Prediction of the inspiratory O2 concentration (optional)

Perseus offers the option of displaying the predicted inspiratory O2 concentration in the breathing gas. This requires the O2 prediction option. The *FiO2 prediction* parameter field must be configured in the monitoring area.

The FiO2 prediction is suitable only for patients with a body weight of more than 30 kg (66 lbs).

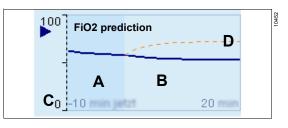

The following is displayed:

- A Section for the trend
- **B** Section for the prediction
- C Current O2 fresh-gas concentration
- **D** Prediction curve ("what if..." function)

WARNING

Risk due to incorrect O₂ prediction

Therapy decisions should not be made solely on the basis of a displayed O₂ prediction.


The inspiratory O₂ concentration is displayed as a line.

The previously measured concentration is displayed in the left-hand section (A) and the predicted concentration in the right-hand section (B). The predicted concentration is dependent on the set O₂ flow and the measured patient values.

Electronically controlled gas mixer only: During the setting procedure, a dotted line (D) is displayed for the O2 flow (the so-called "what if..." function). When the value is altered with the therapy control, the line moves and shows a preview of the predicted concentration.

Econometer (optional)

During operation, Perseus monitors the breathing bag for sufficient filling.

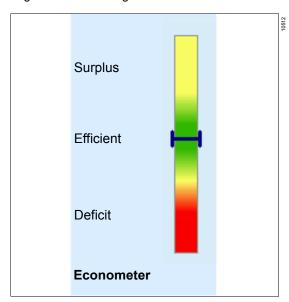

The bar graph indicates the qualitative fresh-gas supply.

Range	Color	Meaning	
Surplus	Yellow	Indication of a possibility to save fresh gas and, therefore, volatile anes- thetics	
Efficient	Green	 No action necessary 	
		 Breathing bag suffi- ciently filled 	
		 Sufficient reserve capacity available 	
Deficit	Red	 Insufficient fresh-gas supply 	
		 Fill the breathing bag immediately 	

An insufficiently filled breathing bag can trigger the *Fresh gas low or leakage* or *Emergency air inlet activated* alarms, for example.

Stopwatch



Using the stopwatch

- 1 Touch the **Start** button (A) to start.
- 2 Touch the **Stop** button (A) to stop.
- 3 To reset the stopwatch to zero, touch the Reset button (A).

Countdown



Setting the countdown

- 1 Touch the **Set** button (B) or the parameter field.
- 2 Set the countdown time.


Using the countdown

The countdown always starts with the last time set.

- 1 Touch the **Start** button (A) to start.
- 2 Touch the **Stop** button (A) to stop.
- To reset the countdown to zero, touch the **Reset** button (A).

Volumeter

The volumeter can be used to monitor and assess the ventilation during spontaneous breathing, manual ventilation, or mechanical ventilation.

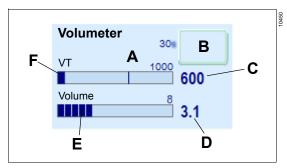

The bar graph indicates the inspiratory and expiratory tidal volume.

At the end of the inspiration, the delivered tidal volume is displayed as a bar (A).

At the end of the expiration, the difference between the inspiratory and the expiratory tidal volume (F) is displayed.

The expiratory tidal volume is displayed next to the bar graph (C).


Using the volumeter (minute volume measurement)

- 1 Touch the Start button (B) to start. The bar graph displays the individual measured breaths in segments (E). The summed volume (D) is displayed next to the bar graph.
- 2 Touch the **Stop** button (B) to stop.
- 3 To reset the volumeter and time display to zero, touch the **Reset** button (B).

The volumeter stops automatically after 60 seconds. The measured values are displayed for 4 minutes and then deleted.

Low-flow wizard (optional)

The low-flow wizard displays bar graphs for the required fresh-gas flow and the total flow. Both bar graphs are to the same scale. The required fresh-gas flow is calculated on the basis of balancing the gas volumes.

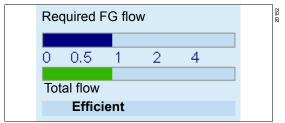

An evaluation of the total flow is displayed below the bar graph:

Evaluation	Color	Meaning
Too high	Yellow	The fresh-gas flow is possibly too high. If the fresh-gas flow can be reduced, both fresh gas and anesthetic agent will be saved.
Efficient	Green	No action is necessary.
Too low	Red	The fresh-gas flow is too low. • Check the fresh-gas flow.
		 Check the position of the breathing bag.
Refill bag	Red	The fresh-gas flow is too low.
		 Refill the breathing bag immediately, e.g., with O2 flush.


An insufficient fresh-gas flow may trigger the **Fresh gas low or leakage** or **Emergency air inlet activated** alarms, for example.

Organizing the screen display

Available views

Perseus offers two view settings for the monitoring area:

Standard view

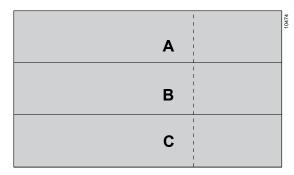

The three waveforms (A), (B), and (C) are displayed along with their associated parameter fields.

Expert view




In addition to the standard view, the four additional parameter fields (D), (E), (F), and (G) are displayed.

Changing the current view

Two buttons are available for changing the current view:

- Views... button
- İ̀d View button

Changing with the Views... button

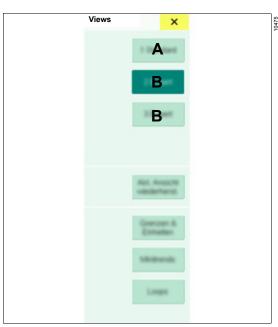

- Touch the Views... button. The Views dialog window will be opened.
- **2** Touch the button for the desired view:
 - A Opens the standard view
 - B Opens the expert view

The views can be renamed, see page 150.

Changing with the View button

• Touch the lill View button.

The screen displays the second view iii.

• Touch the lil View button.

The screen displays the third view iii.


• Touch the lil View button.

The screen displays the first view iii.

Customizing the current view

The monitoring area can be customized during operation:

1 Touch a waveform or parameter field. The **Quick setup** window opens.

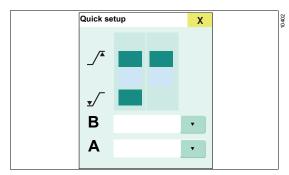

2 For *Content* (A), select the desired content from the list.

For a list of the possible screen content, see page 173.

3 For Scale (B), select the desired scale.

WARNING

Risk of insufficient monitoring

National and medical regulations require certain parameters to be displayed.

Always consider the relevant regulations when configuring the screen layout.

Restoring the current view


The changes in the current view can be reset to the saved standard.

- 1 Open the *Views* dialog window.
- 2 Touch the **Restore current view** button.

Using loops

The following loops are available:

- Pressure-Volume loop
- Flow-Volume loop

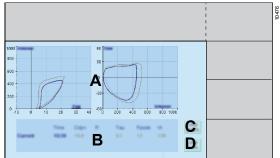


- 1 Open the *Views* dialog window.
- 2 Touch the **Loops** button.

The following is displayed:

- A Current loop and the last 5 loops
- **B** Parameters:
 - Cdyn
 - R
 - TC

The area (E) can be configured so that the Pressure-Volume loop is displayed.

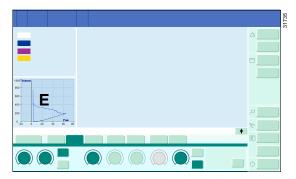

Displaying or deleting reference loops

 To save a reference loop, touch the Save ref. button (C).


Or

 To delete the reference loop, touch the Delete ref. button (D).

These buttons are alternatively displayed in the **Quick setup** window.

Displaying mini-trends

Mini-trends (B) can be displayed for the waveforms (A).

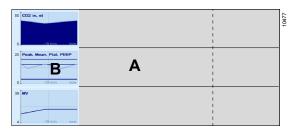


- 1 Open the *Views* dialog window.
- 2 Touch the *Mini-trends* button.

Larger and more detailed graphic and numeric trend displays, see page 109.

Displaying alarm limits and units of measurement


The alarm limits and the units of measurement can also be displayed in the waveform and parameter fields.

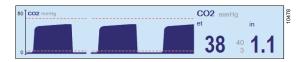

- 1 Open the *Views* dialog window.
- 2 Touch the Limits & units button.

Adjusting the sweep speed and the scale

- 1 Open the **System setup** dialog window.
- 2 Touch the *Screen layout* > *Waveforms* tab (A).

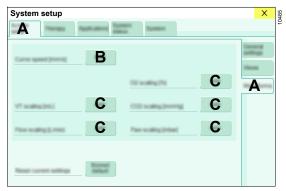

Setting the sweep speed


• Touch button (B). Set the sweep speed.

Adjusting the waveform scale

 To change the waveform scale, touch one of the buttons (C) and select the value.

Changing the color scheme and the screen brightness

- 1 Open the **System setup** dialog window.
- 2 Touch the **Screen layout > General** tab (A).

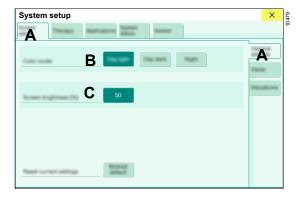

- **3** To change the color scheme, touch one of the buttons at (B).
- 4 Set the screen brightness at (C).


Displaying additional data

Viewing current measured values

During operation, measured values for ventilation, gases, and the device are available in tabular overviews.

- 1 Open the *Trends/Data* dialog window.
- 2 Touch the **Values** tab (A).


The vertical tabs (B) display different combinations of parameters.

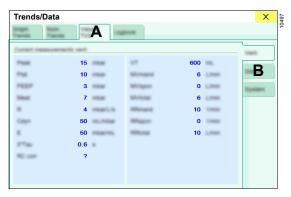

Logbook

The logbook can save up to a maximum of 20000 entries. Logbook data are displayed in table form.

- 1 Open the *Trends/Data* dialog window.
- 2 Touch the Logbook tab (A).

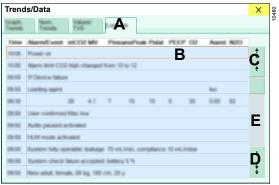

Use the rotary knob or the arrow buttons to scroll the cursor (B) up (C) or down (D) in the logbook. To scroll quickly, touch the gray area (E).

The entries in the logbook cannot be deleted and are retained even after the device has been switched off and on again or following a power supply failure. When the storage limit is reached, the oldest entries are overwritten.


Creating entries and associated settings, see page 115.

Trends

Trends are displayed in the form of a graphic or a table.

- 1 Open the *Trends/Data* dialog window.
- 2 Touch the Graphical trends (A) or Tabular trends tab (B).

The following illustration shows the graphical trend:

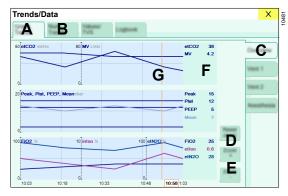

The vertical tabs (C) display different combinations of parameters.

Zooming

In both trend displays, the displayed time period can be enlarged or diminished.

Changing the time period:

- Touch one of the following buttons (D) according to choice:
 - Zoom +
 - Zoom -

Displaying the standard time period and the current point in time:

Touch the Reset zoom button (E).

Moving the cursor

The exact measured values for a specific point in time can be displayed numerically in area (F). To do so, move the cursor.

- Use the rotary knob to move the cursor (G). Or
- Touch the corresponding area on the screen.

Displaying installed options

Listing of the additionally installed software options.

- Open the System setup dialog window.
- 2 Touch the *Licenses/Options* tab.

Displaying an overview of the accessories and consumptions

- 1 Open the **System setup** dialog window.
- 2 Touch the System status tab.

Vertical tab	Overview	
Accessories	Accessories (when Dräger Infinity ID accessories are used) and information as to when the accessory must be replaced.	
Supply	Status display of the connected gas and power supplies.	
Consumption	Gas consumption During operation:	
	for the current caseIn Standby:	
	 for the last case 	
	 since the last reset 	

In *Standby* > *System setup* > *System status* > *Consumption*, the gas consumption levels can be reset to zero, see page 170.

Setting the sound volume

Setting alarm tone volume and breathing sound volume


WARNING

Risk of operation error

The acoustic alarm signals might not be heard if functions such as "Breathing sound" are used or when operating in a noisy environment.

Always set the alarm tone to be sufficiently loud.

- 1 Open the **System setup** dialog window.
- 2 Touch the **System > Sound volume** tab (A).

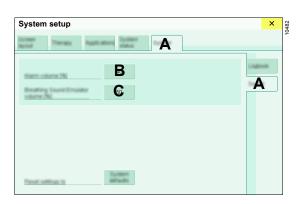

Alarm tone volume

 Set the desired value for Alarm tone volume (B).

Breathing sound volume (optional)

 Set the desired value for Breathing sound volume (C).

Adjusting the alarms


Setting the alarm limits

For a current case, the alarm limits can be set in two ways:

- Setting via the Quick setup window
- Setting via the Alarms dialog window

Setting via the Quick setup window

1 Touch the respective waveform or parameter field. The *Quick setup* window opens.

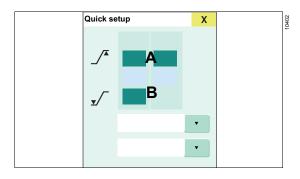


- 2 Set the upper alarm limit (A).
- 3 Set the lower alarm limit (B).

Setting via the Alarms dialog window

In the *Alarms* dialog window, the alarm limits can be set either manually or automatically.

- 1 Open the *Alarms* dialog window.
- 2 Touch the *Limits* tab (C).

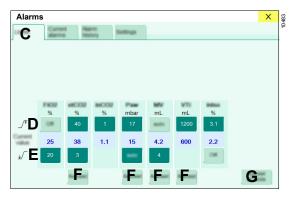

Manual setting

- 1 Set the upper alarm limits (D).
- **2** Set the lower alarm limits (E).

Automatic setting

Alarm limits can be adapted to current measured or setting values.

- 1 Touch the *Autoset limits* button in the main menu bar.
 - The Autoset all button (G) is selected.
- 2 Confirm to adapt the alarm limits for all parameters.

Ör

To adjust the alarm limits for an individual parameter, touch one of the buttons *Autoset* (F) and confirm.

As an alternative, open the dialog window via *Alarms* > *Limits* (B).

Only use the automatic setting when measured or setting values are stable to prevent artifacts of the adjustment algorithm.

The lower alarm limit for the **xMAC** level is also adjusted during automatic setting, see page 112.


Configuration and algorithm, see page 154.

Activating or deactivating CO₂ alarms

The CO₂ monitoring (affects the alarms for *inCO*₂, *etCO*₂, and *CO*₂ apnea) can be activated or deactivated.

Deactivation is indicated in the header bar and in the parameter field by the 🖄 symbol.

- 1 Open the *Alarms* dialog window.
- 2 Touch the Settings tab (A).

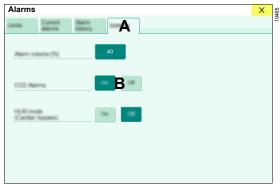

3 For CO2 alarms, touch the button (B): On: Alarms are activated. Off: Alarms are deactivated

Or


- Use the CO2 alarms off button in the main menu bar to activate or deactivate the alarms. This button is only visible in the following ventilation modes:
 - Manual / Spontaneous
 - External fresh-gas outlet
 - Pause

The alarm system is immediately activated when the CO₂ monitoring is activated.

Automatic xMAC monitoring

Perseus provides monitoring of the **xMAC** level. The monitoring is automatically activated as soon as the following conditions are met:

- Anesthetic gas is administered.
- The inspiratory xMAC value is greater than the expiratory xMAC value.
- The expiratory xMAC value reaches approximately 0.3.

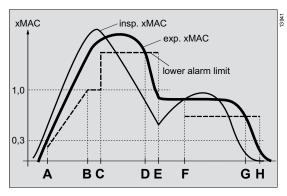

If the **xMAC** value rises, the lower alarm limit for the **xMAC** level will be automatically adjusted to the anesthetic gas concentration after activation (A). The lower alarm limit (B) can thus reach a maximum value of 1.0. The lower alarm limit can be recalculated by touching the **Autoset** button (C). In special anesthesia situations, the **xMAC low** alarm can be adjusted in this way and consequently also exceed the value of 1.0 if necessary.

If the expiratory **xMAC** value falls below the alarm limit (D or G), Perseus generates the **xMAC** low alarm with low priority. If the alarm is not acknowledged with the **ALARM RESET** button, the priority is raised to medium priority after 60 seconds.

Deactivation of automatic *xMAC* monitoring

If the **xMAC low** alarm (D or G) is acknowledged with the **ALARM RESET** button (E or H), the monitoring is deactivated. This prevents renewed alarms as a result of the anesthetic gas concentration continuing to fall at the end of anesthesia (H). If the anesthesia is continued (E), the monitoring will be automatically reactivated as soon as the inspiratory **xMAC** value rises above the expiratory **xMAC** value (F).

In CBM mode the lower alarm limit is adjusted both upwards and downwards so that no alarm is issued during this time. Similarly, the value is not limited to 1.0 during this time.

Changing the patient data

Patient data can be changed during operation.

1 Open the Patient dialog window.



2 Make the desired changes.

Changes influence, e.g., therapy suggestions, as is recognizable by the position of the arrow ▼ at the therapy controls.

The current therapy settings remain unaffected.

When patient category is changed, age, weight, and height are adapted as needed so they remain within the described limits, see page 157.

Exporting data

General information

Prerequisite: USB flash drive is connected to the USB interface.

During a saving process, the button is displayed as activated (dark green).

The data are stored in the "Draeger\ExportData" directory.

Exporting the screen contents

The screen contents can be exported to a USB flash drive as a screenshot.

 Touch the Export screenshot button in the main menu bar.


The screenshot will be saved as a ".png" file.

Exporting trends and data

In the **Standby** mode, the following data can be exported to a USB flash drive:

- System test results
- Logbook Selection from the following time periods is possible:
 - Last case
 - Today
 - All
- Alarm history
- Trends

- 1 Open the *Trends/Data* dialog window.
- 2 Touch the Export tab (A).

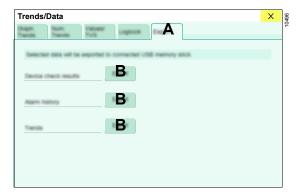

3 Touch the desired button (B).


The data will be saved as a ".txt" file.

Other settings

Switching the breathing system warmer on or off

The breathing system warmer should only be switched off in special situations (e.g., for intentional reduction of the body temperature of the patient).

- 1 Open the **System setup** dialog window.
- 2 Touch the *Therapy* tab (A).

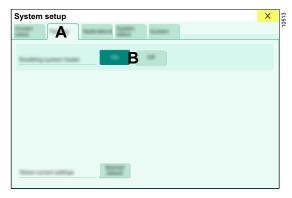

3 To switch the warmer on or off, touch the corresponding button (B).

When switching to **Standby** mode, the warmer is reset to the value configured in **System setup**.

CAUTION

Risk due to faulty or switched-off breathing system warmer

Increased condensation or accumulation of water in the breathing system and the hoses may occur.

Increase fresh-gas flow as required. Remove condensate from hoses, water traps, and breathing system regularly. Have service personnel repair the faulty breathing system warmer.

CAUTION

Risk of inaccurate measured values


The accuracy of the flow measurement may be compromised if the breathing system warmer is switched off.

When the breathing system warmer is switched off, do not make therapy decisions solely on the basis of the displayed values for flow and volume.

Creating additional logbook entries

The following events can create logbook entries with measured values for the parameters *etCO*2, *MV*, *Pmean*, *PIP*, *Pplat*, *PEEP*, *FiO*2, exp. concentration of the primary anesthetic gas, and *etN*2*O*:

- Adjustable time intervals
- Alarms with high or medium priorities
- 1 Open the **System setup** dialog window.
- 2 Touch the **System** > **Logbook** tab (A).

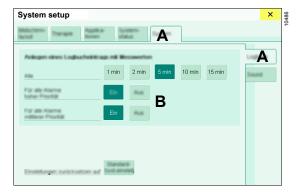

3 Touch the corresponding button (B).

Resetting user-specific settings

Changes made in the **System setup** dialog window during operation can be reset to the start settings.

- 1 Open the **System setup** dialog window.
- 2 Open the corresponding dialog window.
- 3 Touch the System defaults button and confirm.

Ending the therapy

Switching to the Standby mode

 In the main menu bar, touch the Standby... button.

Or

Press the () key below the screen.

2 In the **Standby** dialog window, confirm the automatically selected **Standby** button.

With mechanically controlled gas mixer:

3 Close the flow control valves.

Using O₂ insufflation

WARNING

Risk of fire

In combination with oxygen or nitrous oxide, ignition sources such as electrosurgical devices and laser surgical devices can cause fires.

- Prevent leakage, e.g., at endotracheal tubes, laryngeal mask airways, face masks, Y-piece, breathing system including hoses, filters, and breathing bag, at the external fresh-gas outlet, and at the outlet for O2 insufflation.
- Use only intact and leak-free hoses at the outlet for O₂ insufflation.
- Before beginning laser surgery or electrosurgery, flush with sufficient air (<25 % O2), and flush beneath the surgical drapes as well.
- Close the flow control valve on the O2 flowmeter to the end stop.
- When O2 outlets are in use (e.g., for insufflation), do not use any ignition sources in the immediate vicinity.
- Do not position oxygen sources in the vicinity of ignition sources, e.g., electrical connectors.

WARNING

Risk of a device fire

If there is a fire on the patient side, the device may ignite.

If there is a fire, disconnect the device from all oxygen-carrying supplies.


WARNING

Risk due to overpressure


When the patient is connected to the outlet for O2 insufflation without a relief valve, increased pressure may be applied to the patient.

Only connect the patient in a way that allows excess gas to escape (e.g., through a relief valve).

Electronically controlled gas mixer:

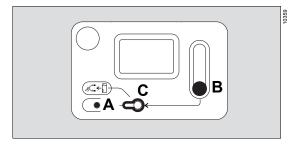

Mechanically controlled gas mixer:

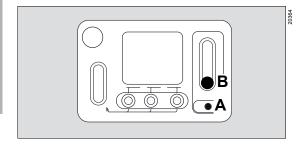

Prerequisite:

- The appropriate accessory is connected to the outlet for O₂ insufflation (A).
- With electronically controlled gas mixer: The O2 switch is horizontal in the *Aux. O2* position (C).

Beginning O₂ insufflation

 Open the flow control valve (B) of the O2 flowmeter.

Ending O2 insufflation

 Close the flow control valve (B) of the O2 flowmeter.

Change of patient

Cleaning and disinfecting the workstation

 Clean and disinfect the workstation according to hospital hygiene regulations. Follow the reprocessing instructions supplied with the product.

Checking or replacing consumables

Prerequisite: Device is in Standby mode

Sample gas measurement

WARNING

Risk of infection

Used sample lines and water traps may be infectious due to the breathing gases that passed through them.

- Replace the sample line regularly in the following situations:
- If the sample line is connected to the filter on the Y-piece, replace it daily.
- If there is no filter fitted to the Y-piece and the sample line is connected directly to the Y-piece, replace the sample line after every patient.
- Remove the sample line from the water trap.
- Initially leave the water trap fitted to prevent infectious fluid from spurting out.
 Remove the water trap only after surface disinfection.
- Empty or replace the water trap according to its instructions for use.
- Check the water trap of the sample gas measurement; empty or replace as needed.
- 2 When no filter is used, replace the sample line and dispose of used sample lines.

Vaporizer filling level

 Check the vaporizer filling level at the sight glass. If necessary, fill the vaporizer.

CO₂ absorber

 Check the coloration of the soda lime and replace it if necessary, see page 76.

Breathing hoses and filters

- Replace the hoses and filters according to hospital hygiene regulations.
- 2 Select and connect a suitable breathing circuit and filters, see page 67.

Checking the device

Prerequisite: Device is in Standby mode

- 1 Perform the leakage test, see page 122.
- 2 Flush the breathing system if necessary.

Flushing the breathing system

With the *Flush* function, excess moisture and residual anesthetic agent present in the breathing system and breathing circuit can be reduced by flushing with ambient air.

Prerequisite: Device is in Standby mode

- 1 Touch the *Flush* button. Alternatively, flushing can also be started via the *Tests...* button in the main menu bar.
- 2 Prepare the device according to the displayed checklist.
- 3 Touch the Start button.
- 4 After the end of flushing, the device switches to **Standby** mode.

If necessary, end the flushing prematurely with *Cancel*.

WARNING

Risk due to residual concentrations of anesthetic agent

Even after flushing the breathing system and breathing circuit, residual concentrations of anesthetic agent might not be sufficiently reduced.

For patients suspected of having malignant hyperthermia, do not rely solely on the flushing function.

For further information and recommendations for therapy settings for patients with suspected malignant hyperthermia, contact the responsible national Dräger organization.

CAUTION

Risk of impaired ventilation due to faulty drying function

The breathing system and breathing hoses may still retain traces of moisture even after the flushing and drying functions are performed.


Check the breathing system and breathing hoses for condensate.

Tests

Information on the availability of individual device functions	121
Testing the device	121
Available test types	122
Performing the tests	123
Special characteristics of Dräger coaxial hoses	
System test and leakage test Leakage assistant	123
Completing the checklist	126
Prerequisites	126
Vaporizers	127
Gas flow	
Breathing circuit	
Valves	
Components	131
Accessories	132

Information on the availability of individual device functions

 In the Standby mode, touch the Details... button or the Tests... button.

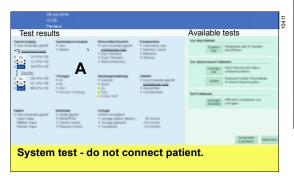

In the *Test results* list (A), the following information is displayed in color:

- Influence of the individual device functions on the availability of the entire device
- Latest detected leakage values

Color	Meaning
Green	Successfully tested, fully available
Yellow	A non-critical fault has been detected. The device can be operated with restricted function.
Red	A serious fault has been detected. Operation is not possible or is forbidden.
Gray	Function not tested.

WARNING

Risk due to faulty sensors

If the sensors for flow, O2, N2O, CO2, or anesthetic agent are not operational, the patient will not be adequately monitored.

Before using the medical device, ensure a suitable substitute monitoring.

Testing the device

CAUTION

Risk due to device malfunction and/or patient injury

If the system test is canceled, it is possible that some malfunctions might not be detected. Thus, more attention is required during operation.

Perform the system test every day. If the system test is canceled during execution, perform it again as soon as possible.

WARNING

Risk of patient injury

During the system test, the system is pressurized.

To prevent patient injury, do not perform the system test on the medical device if a patient is connected.

Available test types

Test	Type and duration	Perform the test	Description	
System test Auto-		– Daily	Initialization:	
	matic, about 8 min	 After service activities 	 Checking of components which frequently cause operational restriction: e.g., high leakage, incorrect setting of APL valve. 	
			Test:	
			Calibration of all valves and sensors	
			Test of all device functions	
			 Includes the leakage test 	
			 The test can be programmed, see page 136 	
Leakage test	Auto- matic,	 After filling the CO₂ absorber 	 Calculation of leakage, system compliance, and system resistance 	
	approx. 2 min	 After changing the hose configuration (e.g., changed hoses, changed lengths of extendable hoses etc.) 	 Calibration of valves and flow sensors, if required. In this case, the test is extended by approx. 3 minutes. 	
		 After replacing the breathing system 		
		 After replacing the flow sensors 		
Leakage assistant	Manual, as required	After occurrence of prob- lems with leakage during automatic tests	sure and the leakage to support the manual check. Changes are immedi-	
		 If leakage in breathing system and breathing circuit is suspected 	ately visible.Leakage test of a connected vaporizer	
		 If leakage at vaporizer is suspected 	 Leakage test of the inner hose of a connected Dräger coaxial breathing 	
		If coaxial breathing circuits are used	circuit. Use the coaxial test adapter for this. For more information, see page 124.	

Performing the tests

Special characteristics of Dräger coaxial hoses

When using coaxial breathing hoses, always use the leakage assistant to check the inner hose of the coaxial breathing circuit for leakage prior to the system test or leakage test, see page 125.

System test and leakage test

Prerequisite:

- Electronically controlled mixer: Central O2 supply or central Air supply is connected.
- Mechanically controlled mixer: Central O2 supply is connected.
- The anesthetic gas receiving system is correctly connected.


The two tests consist of a checklist followed by an automatic test.

The checklist can be presented in tabular form or as a walk-through mode.

Button (D) can be used to switch from the tabular checklist to the walk-through mode.

Whether tests always start in the walk-through mode can be specified in the system setup, see page 166.


1 Touch the button for the required test.

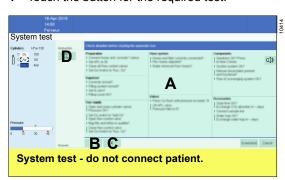

2 Complete the checklist (A), see page 126.

3 If all components are ready for operation, touch the button (B). The automatic test will start. If a component is not ready for operation, touch the button (C).

The walk-through mode will start.

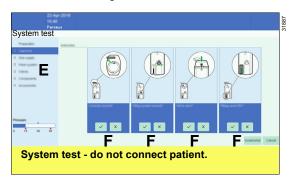


The components (E) are polled sequentially.

Buttons (F) are used to document whether the test passed.

Button	Meaning	
\checkmark	Test passed	
×	Test failed	

The automatic test starts after all the tests in the walk-through mode are complete.

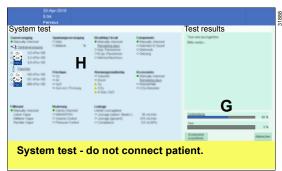


The test progress is displayed in area (G) while the automatic test is running. All test results are displayed in area (H).

After the test, the final test result is displayed on the *Standby* page, see page 81.

Test interruption due to irregularities

If an irregularity is detected during the automatic test, the following occurs:

- The test is interrupted.
- An acoustic signal sounds.
- Information on the cause and remedy are displayed.

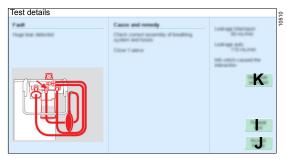

Remedying the cause:

- 1 Remedy the cause of the interruption. If there is leakage, the leakage assistant (K) can be used to support troubleshooting.
- 2 Touch the *Repeat* button (I) and repeat the test of the component.

Accepting the irregularity:

 Touch the Accept button (J) and continue the test.


Accepted irregularities prevent the total result of the system test from indicating "fully operational" and are protocoled in the logbook.

Leakage assistant

Prerequisite when vaporizers are connected:

- Vaporizer is vertically level and securely mounted on the plug-in adapter.
- Filling inlet is closed.

In this test, a continuous pressure is generated and the current leakage value is displayed.

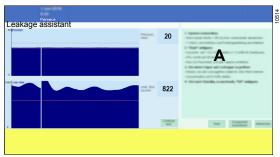

If the leakage value changes as a result of changes made to the device (e.g., loosening or readjusting hose connections), this can help locate the cause of the leakage.

The displayed leakage value may differ from the value that was determined in the leakage test. The reasons for this are the different measuring methods and different pneumatic ranges.

- 1 Touch the button for the test (*Leakage* assistant).
- 2 Follow the instructions (A).

Checking a vaporizer for leakage

Prerequisites:


- Gas mixer:
 - With electronically controlled gas mixer:
 The flow control valve of the O2 flowmeter is closed.
 - With mechanically controlled gas mixer:
 All flow control valves are closed.
- The vaporizers are closed.
- The breathing circuit is correctly connected.
- The APL valve is set to 30.

Perform the test:

- Touch the Start button.
 - Wait until the leakage value is stable.
 Memorize the value.
 - Set the anesthetic gas concentration on the vaporizer to >0.2 Vol% (>2 Vol% for desflurane).
 - Wait until the leakage value is stable.
 Memorize the value.
 - If the two values differ from one another by more than 50 mL/min, check the vaporizer and the vaporizer interface for leakage.
 - Set the anesthetic gas concentration to 0 Vol%.
- 2 Touch the OK button to switch back to the display of available tests.

Checking a coaxial breathing hose for leakage

The leakage assistant can be used to identify leakage in the inner hose of a coaxial breathing circuit. When Dräger coaxial breathing circuits are used, the inner hose must be connected directly to the expiratory port using a test adapter.

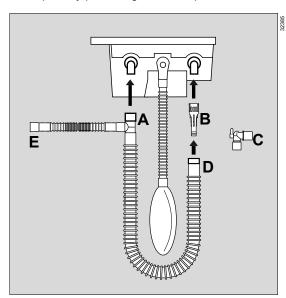

Prerequisites:

- Gas mixer:
 - With electronically controlled gas mixer:
 The flow control valve of the O2 flowmeter is closed.
 - With mechanically controlled gas mixer:
 All flow control valves are closed.
- The vaporizers are closed.
- The APL valve is set to 30.

Perform the test:

- Connect the connector (A) to the inspiratory port.
- 2 Connect the coaxial test adapter (B) to the expiratory port.
- 3 Remove the elbow (C) from the hose. Connect the connector (D) to the coaxial test adapter (B).
- 4 Touch the Start button.
 - Wait until the leakage value is stable.
 If the leakage value exceeds
 500 mL/min, use another breathing circuit.
- 5 Touch the OK button to switch back to the display of available tests.
- 6 Fully connect the breathing circuit.
 - Remove the coaxial test adapter.
 - Plug in the elbow again.
 - Connect the connector (E) (expiration) to the expiratory port.
 - Perform the system test or the leakage test.


Completing the checklist

This section describes how the checklist is processed, using as an example a device with an electronic gas mixer and factory defaults.


The scope of the checklist can vary due to differing system settings.

The instructions on the screen take precedence.


Prerequisites

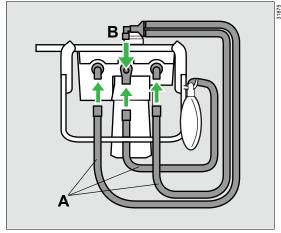


- 1 Connect the hoses (A).
- 2 Seal the Y-piece (B).

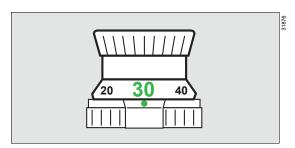

3 Set the APL valve to 30.

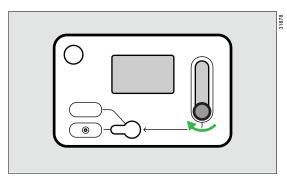


4 Close all the flow control valves.

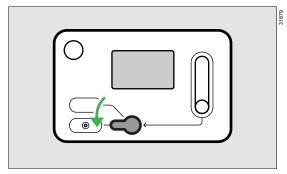


5 Set the O2 switch to Aux. O2.

Vaporizers


For each vaporizer, check:

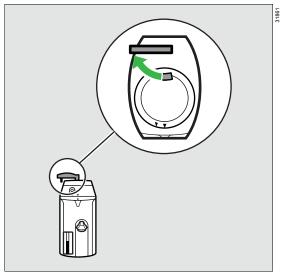


1 The locking lever points left, indicating the vaporizer is locked.

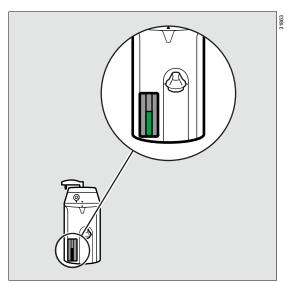

2 Check the filling level in the sight glass. Refill anesthetic agent if required. When using a Dräger Vapor 3000 or D-Vapor 3000 with the Vapor View option, a yellow triangle on the screen indicates that the filling level has fallen below the refill mark.

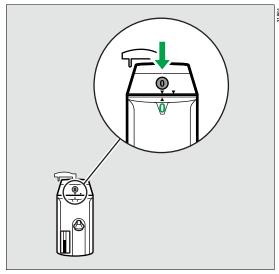

3 The control dial is set to position *0* and the key is locked in placed.

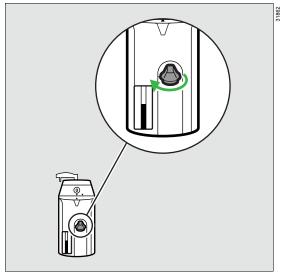

4 The filling inlet is closed.

Gas flow

WARNING

Risk due to impermissible gas supply pressure

Impermissible gas supply pressures may cause incorrect gas composition.

Check the supply pressures of the central gas supply and of the gas cylinders before operation.

WARNING

Risk due to impermissible gas supply

Using oxygen supplies with less than 100 % O2 may cause incorrect gas composition.

When using O₂, only use 100 % O₂.

WARNING


Danger when using O2 concentrators

The following effects may occur:

- Discrepancies between the set value and the actual value for fresh-gas flow and O2 concentration in the fresh gas
- Inaccurate measured values for volume, anesthetic agent consumption, econometer (optional), and low-flow wizard
- Inaccuracy in the FiO₂ prediction
- Inaccurate volume delivery in volumecontrolled ventilation modes
- Accumulation of argon in low-flow operation and minimal-flow operation

Do not use any O2 concentrators.

Gas cylinders

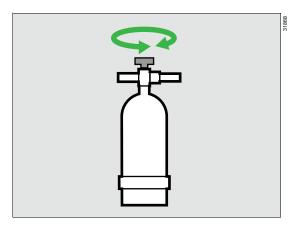


- Open the gas cylinder valves slowly. Check that the displayed pressures are sufficient.
 - When using pressure reducers without electronic pressure measurement, read the pressure from the pressure gauge.
- 2 Close the gas cylinder valves.


On devices that are equipped with Advanced Cylinder Support, the gas cylinder valves can remain open during operation. These devices are identified by an appropriate label near the gas inlets, see page 28.

Checking the emergency O₂ delivery (with electronically controlled gas mixer)


Prerequisite: Y-piece is sealed.

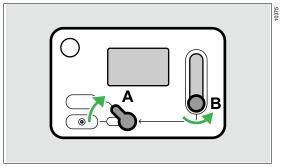

- 1 Set O2 switch (A) to the Add. O2 position.
- 2 Open the flow control valve (B) of the O2 flowmeter and set the desired O2 flow.

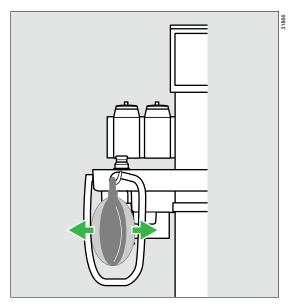


The breathing bag fills and the inflow of gas is audible.

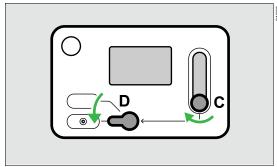

- 3 Close the flow control valve (C).
- 4 Reset O2 switch (D) to the Aux. O2 position.

Checking the flow control valves (mechanically controlled gas mixer)

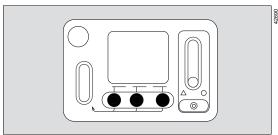


- 1 Open the flow control valve for O2 and set a flow of approximately 2 L/min.
 - The breathing bag inflates.
 - The inflow of gas is audible.
 - The total flow tube shows the total freshgas flow.
- 2 Repeat the settings for the Air and N2O flow control valves.
- **3** Close all the flow control valves.

Breathing circuit

Prerequisite:


- Breathing system is complete and locked.
- Breathing system cover is in place.

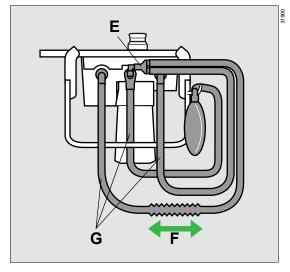

- **1** Hoses (G) and filters, e.g., at the Y-piece (E), are connected properly.
- 2 Extendable hoses (F) are extended to the foreseen application length.

NOTE


Do not change the length of flex hoses after the test is done.

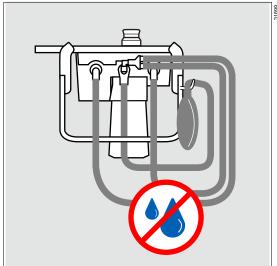


3 Remove the water from the hoses.

Valves


Checking the O₂ flush

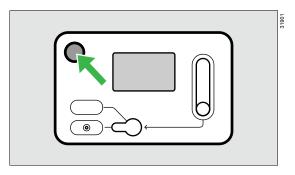

1 Press the *O2+* key and keep it pressed until the pressure exceeds 15 hPa.

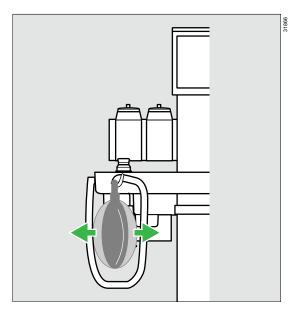

The breathing bag fills and the inflow of gas is audible.

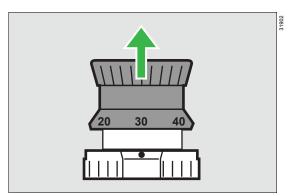


2 Lift the APL valve.

The pressure is released.

Components

Loudspeakers

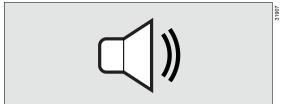


 Touch the button and wait for 2 acoustic signals.

NOTE


If the acoustic signals are not emitted, contact service personnel.

Bronchial suction

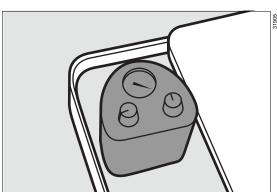


 Check the functional integrity of the bronchial suction.

Manual resuscitator

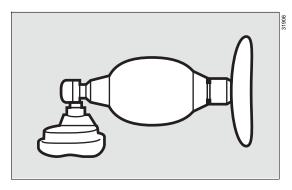

- 1 Make sure there is a manual resuscitator on the device.
- Check the functional integrity of the manual resuscitator.

Anesthetic gas receiving system

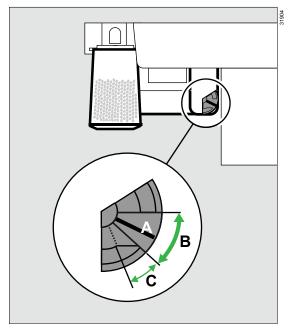


 Have the flow for the anesthetic gas scavenging system set so that the red flow indicator (A) floats in range (B) ("normal range").

When the flow indicator is floating in range (C) ("restricted range"), certain fresh-gas flow rates should not be exceeded, see "Anesthetic gas receiving system" in chapter "Technical data". This can prevent contamination of the ambient air.

Accessories

Soda lime

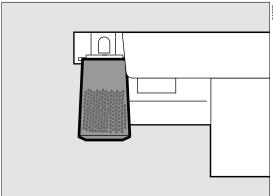

Make sure that the soda lime is not discolored. Change the soda lime if it is discolored or when its maximum period of use has been reached, see page 76.

When Infinity ID functionality is configured:

Absorbers of the type Infinity ID CLIC
 Absorber will automatically be detected,
 and the exchange date will automatically be set.

When Infinity ID functionality is not configured:

- Absorbers will not be detected, e.g., reusable CO₂ absorbers
- Update the exchange date manually: Touch the *Reset* button after the soda lime has been replaced.

WARNING

Risk due to soda lime drying out

The soda lime loses moisture. If the moisture level falls below the minimum moisture level, the following adverse reactions occur regardless of the type of soda lime and the inhalational anesthetic agent used:


Decreased CO2 absorption, increased generation of heat in the CO2 absorber resulting in increased breathing gas temperature, formation of CO, absorption and/or degradation of the inhalational anesthetic agent.

- Do not use unnecessarily high fresh-gas flows.
- Only use the O2 flush if necessary.
- With electronically controlled gas mixer:
 Only use the emergency O2 delivery if necessary.
- With mechanically controlled gas mixer:
 Do not leave the flow control valves open unnecessarily long.

NOTE


Only perform leakage tests with CLIC absorber locked into place because this affects the system compliance values.

Sample line

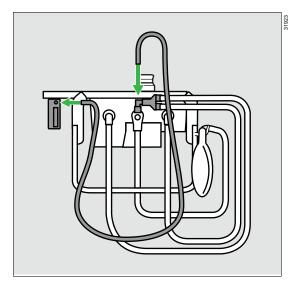

Connect the sample line.

Water trap

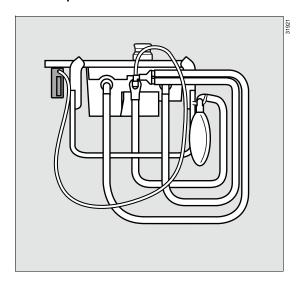

- 1 Check the water level in the water trap.
- 2 Pay attention to the period of use of the water trap and replace if necessary.

When Infinity ID functionality is configured:

 Water traps of type Infinity ID WaterLock 2 will automatically be detected and the exchange date will automatically be set.

When Infinity ID functionality is not configured:

- Water traps will not be detected.
- Update the exchange date manually: Touch the *Reset* button after a new water trap has been installed

Ending operation


Turning off the device	136
End of operation	136
Programming the system test (Auto On option)	
Flushing and drying the breathing system	
(Flush & Dry function)	137
Storing the device	137
Disconnecting the mains power supply	137

Turning off the device

Dräger recommends shutting down Perseus during longer periods of non-use such as overnight or on weekends. This can lower power consumption and prolong the life span of the medical device without negatively influencing device availability.

When shutting down the device, several functions are available to the user.

Function	
Auto On (optional)	Prepares the device for automatic start-up, including system test.
Flush & Dry	Flushes and dries the breathing system.

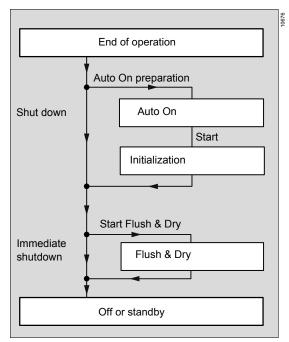

For more information, see page 165.

End of operation

- Make sure that all flow control valves are closed
- 2 Press the () button.
- 3 Follow the instructions on the screen.

Programming the system test (*Auto On* option)

Using the **Auto On** function, the device can be programmed so that it performs an automatic system test and is tested at a specific time.

During the initialization, the *Auto On* function checks the components that frequently cause irregularities: e.g., high leakage, incorrect setting of APL valve. After that, the device remains in the *Standby* mode until the automatic system test, or it is shut down.

- **1** Press the () button.
- 2 Touch the **Auto On preparation** button.
- 3 Prepare the device in accordance with the instructions on the screen. For details on the test steps, see page 126.
- 4 Set the desired day and time when the device should next be ready for operation.
- 5 Select whether the device should remain in **Standby** mode or shut down.
- 6 Touch the Start button.
- 7 Take note of any messages that appear.

After the Auto On initialization is complete, the time remaining until the system test is displayed in the **Standby** mode. When the device is being shut down, the corresponding symbol is shown in the status display, see page 31.

For a final test of the operational readiness of the device, it is only necessary to perform some brief manual checks after the automatic system test.

Flushing and drying the breathing system (*Flush & Dry* function)

Dräger recommends always flushing the breathing system before switching off to prevent damage or failure due to residual anesthetic agent and residual moisture.

During this process, first ambient air is used for flushing because this does not dry out the soda lime. Then a small amount of Air or O2 from the central gas supply or the gas cylinders is used to dissipate residual moisture.

- 1 Prepare the device in accordance with the instructions on the screen.
- 2 Touch the Start Flush & Dry button.
- 3 If necessary, cancel the flushing with *Cancel* to return to the *Standby* mode.

After the flushing procedure, Perseus shuts down or remains in the **Standby** mode (**Auto On** option).

Storing the device

To store Perseus

- 1 Set the main switch to position **0**.
- Disconnect from mains power supply if necessary.

Disconnecting the mains power supply

Pull out the mains plug.

On the ceiling-mounted version, the corresponding cable is labeled with **Perseus**.

Alarms

Displaying alarms	139
Optical alarm signals	
Response to alarms	140
Displaying information on alarms Suppressing the alarm tone Downgrading and acknowledging alarm	140 141
messages	141
Opening the alarm history	
Adopting alarm settings when changing the ventilation mode	143
Activating the alarms related to volume	
Alarm delay and alarm escalation	144
Activation of alarms after breath detection.	147
Intelligent alarm behavior	147
Combined alarms	
Suppressed alarms	


Displaying alarms

Alarms are signaled optically and acoustically.

Optical alarm signals

On the screen

In the event of an alarm, the system displays the relevant alarm message in the alarm message field (A). The parameter field (B) for the alarm-generating parameter flashes.

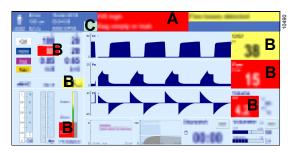

Up to 8 alarms can be displayed simultaneously in the alarm message field in the header bar (A). If more alarms occur, the *All alarms* button (C) is displayed in the header bar. Touching this button opens the dialog window *Alarms* > *Current alarms* with information about all active alarms, see page 140.

On Dräger Vapor 3000

When the Vapor View option is installed, alarms concerning the vaporizer are signaled by blinking illumination on the Dräger Vapor 3000 if their cause can be remedied by adjusting or refilling the vaporizer, for example:

- Inspiratory xMAC high
- Filling level of vaporizer low
- Vaporizer open
- xMAC low

Acoustic alarm signals

It always is the alarm with the highest priority that is acoustically signaled. The signal is emitted until either the cause of the alarm is remedied or the Audio paused key is pressed.

Depending on the overall alarm situation, it is possible that the 10-tone sequence for the high alarm priority is only issued as a 5-tone sequence due to the coinciding occurrence of alarms.

Regardless of the set alarm volume, the **No O2 delivery** alarm is issued at maximum volume.

Alarm priorities

Perseus assigns the appropriate priority to each alarm message.

The background color of the alarm message field indicates the alarm priority of the active alarms. The parameter field of the parameter triggering the alarm flashes in the color matching the alarm priority.


Background color	Alarm priority	Marking	Meaning
Red	High	!!!	Immediate action is necessary in order to avert imminent danger.
Yellow	Medium	!!	Prompt action is necessary in order to avert a danger.
Cyan	Low	!	Attention is necessary, but a delayed response is sufficient.

Response to alarms

Displaying information on alarms

Touch the alarm in the header bar.
 Or
 Open the *Alarms* dialog window and touch the

Open the *Alarms* dialog window and touch the *Current alarms* tab (A).

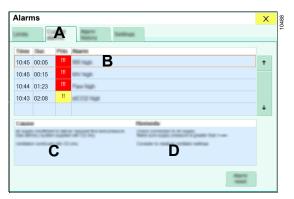


2 In the list (B), touch the corresponding alarm or select it with the rotary knob. 3 Refer to the information under *Cause* (C) and *Remedy* (D) to remedy the error.

A list of all possible alarms can be found in chapter "Troubleshooting", see page 191

Suppressing the alarm tone

The alarm tone can be suppressed for a maximum of 2 minutes.

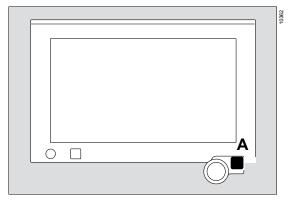

 Press the Audio paused key (A) below the screen, see page 19.

The A symbol and the remaining time for the suppressed alarm tone are displayed in the header bar.

After the alarm suppression has expired, the alarm tone immediately resumes if the cause of the alarm still exists.

While the alarm tone is suppressed, only those alarms with a higher alarm priority or internal priority rating are issued acoustically, see page 191.

Reactivating the alarm tone

• Press the Audio paused key again.


Downgrading and acknowledging alarm messages

Some alarms can be downgraded to low priority or can be cleared completely. The relevant alarms can be identified in the table "Alarm – Cause – Remedy" on page 191 by the following remedial messages:

Remedial message	Effect
Use "ALARM RESET" to downgrade alarm priority.	Alarm priority is changed to low.
Use "ALARM RESET" to acknowledge alarm.	Alarm is reset.

There are 2 options for downgrading or resetting the alarms:

Option 1:

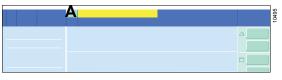

Or

Option 2:



Option 1	Option 2
Touch the ALARM RESET button (A) in the header bar and confirm.	In the <i>Alarms > Cur-</i> rent alarms dialog window (B), touch the Reset all button (C) and confirm.
All the alarms displayed in the alarm message field will be downgraded or reset.	All alarms will be down- graded or reset.

After an alarm is downgraded, the respective priority rating is maintained, see table "Alarm – Cause – Remedy" on page 191.

Opening the alarm history

The alarm history records all alarm messages for the current case in chronological sequence.

- 1 Open the *Alarms* dialog window.
- 2 Touch the Alarm history tab (A).




Use the rotary knob or the arrow buttons to scroll the cursor up (B) or down (C).

The alarm history is deleted when Perseus is shut down or a new case is started.


Adjusting the alarm limits

If an alarm is triggered because a lower limit or an upper limit is transgressed, it might be necessary to adjust the alarm limits. To do so, either set the alarm limits, see page 110, or change the alarm limits via the **Quick setup** window.

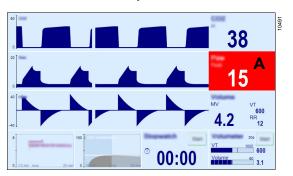

1 Touch the parameter field (A).

The **Quick setup** window opens, and the alarm limit (B) that was violated is already selected.

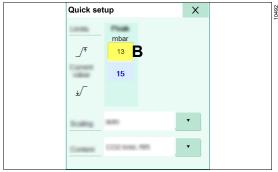

2 Adjust the value (B) and confirm.

Perseus can be configured so that the *Quick setup* window opens automatically in the event of an alarm, see page 155.

Adopting alarm settings when changing the ventilation mode

When the ventilation mode is changed, the alarm settings are also adapted.

Depending on the mode, alarm settings can either be adopted or set to Off.

Some modes can be configured whether or not the settings are adopted.

However, the settings can be adjusted at any time during operation.

Alarm or alarm	Mode				
limit	VC, VC - AF, PC, PC - APRV (optional) PSV (optional) with ΔPsupp ≥ 5 hPa	Pause, Man/Spon, PSV (optional) with ΔPsupp <5 hPa	Ext. FGO (optional)	CBM mode	
FiO ₂ low		Δlarm	settings are ad	onted	
inAgent high		Aldilli	settings are ad	opica	
Apnea (no CO2)		Alarm settings are adopted			
etCO2 high etCO2 low				Off	
inCO2 high					
FiO ₂ high	Is restored or remains	configurable, see page 156		configurable, see	
inAgent low	active			page 156	
xMAC low				Not activated	
Paw high		Alarm settings are		Alarm settings are	
Paw low		adopted		adopted	
MV high			Not mea- sured	6 11	
MV low		configurable, see page 156		configurable, see page 156	
Apnea (no flow) ¹⁾		1-3-		F-3	
Apnea (no pressure)	On	Off			

¹⁾ This alarm is only activated when the MV low alarm limit is also activated.

Activating the alarms related to volume

The upper alarm limit for **MV** is deactivated at the start of the ventilation, but can be set during the ventilation. The upper alarm limit for **VTi** is automatically set to 130 % of the set **VT** in volume-controlled ventilation modes. In pressure-controlled ventilation modes, it is deactivated at the start of the ventilation, but can be set during the ventilation.

The **MV low** alarm is delayed in certain cases and is issued as follows:

- No sooner than 90 seconds after a case starts,
- No sooner than 60 seconds after changing to a mode with greater respiratory support, see page 259.

 No sooner than 60 seconds after an Apnea (no flow) or Apnea (no pressure) alarm.

In volume-controlled ventilation modes, the alarm limit for the inspiratory tidal volume is automatically set to 130 % of the set tidal volume.

Resetting the Apnea (no CO₂) alarm

When changing to a ventilation mode with higher respiratory support, the *Apnea (no CO2)* alarm is reset. If the apnea situation persists, an alarm appears after the time specified in the table "Alarm delay and alarm escalation".

Alarm delay and alarm escalation

To prevent unnecessary alarms, some alarms are not displayed immediately after a limit violation, but after a delay. In addition, certain circumstances can cause the alarm priority to change.

Alarm	Priority		
	Low	Medium	High
inCO2 high etCO2 high etCO2 low FiO2 high Inspiratory N2O high		After two successive respiratory phases and 15 seconds	
inAgent low	After two successive respiratory phases and 15 seconds		
FiO ₂ low			After two successive respiratory phases and 15 seconds or after 30 seconds if no respiratory phases are detected

Alarm	Priority				
Alami	Low	Medium	High		
		After two successive respiratory phases and 15 seconds or after 30 seconds if no respiratory phases are detected	>165 seconds later		
Inspiratory xMAC high		insp. MAC ≥ 3 for more than 180 seconds	longer than 30 seconds: insp. MAC ≥ 3 and exp. MAC ≥ 2.5 or insp. MAC ≥ 5		
xMAC low	0 to 60 s	>60 s			
Apnea (no CO2) Apnea (no flow) Apnea (no pressure) Apnea		at the latest after 20 seconds (for RR ≥ 6) or At the latest after 35 seconds (for RR <6) or at the latest after 65 seconds for the Pause, Manual / Spontaneous and Ext. FGO modes	15 seconds later (for RR ≥ 6) or 30 seconds later (for RR <6)		
Apnea Ventilation	At the latest after 20 seconds (15 seconds for RRmin ≥ 4) (configurable, see page 155)				
No CO2 detected	>60 s				
Inspiratory tidal volume high Tidal volume not achieved		After 3 successive breaths			

Alone	Priority				
Alarm	Low	Medium	High		
Airway press. continuously high			> 15 seconds above the manually or auto- matically set limit		
Airway pressure negative			Pmean < -2 or Paw < -10		
PEEP/CPAP high		Airway pressure >(PEEP +5 hPa) during more than 10 consecutive breaths			
Airway pressure not achieved		>15 seconds (for RR			
Cardiac bypass mode still active?		If a minute volume of >50 % of the suggested value is measured after CBM mode has been activated for >60 seconds			
Fresh gas low or leakage	Breathing bag almost empty	Breathing bag empty	after 30 s or in the event of the additional alarm "Apnea (no flow)" or "Apnea (no pressure)" or in the event of the additional alarm "Emergency air inlet activated"		
O2 measurement not available N2O measurement not available Agent measurement not available Agent measurement temporarily inaccurate CO2 sensor accuracy low Agent measurement temporarily inaccurate N2O measurement temporarily inaccurate Measured gas concentrations temporarily inaccurate	>20 s				

Activation of alarms after breath detection

If no breaths have been detected even after leaving the *Standby* or *Pause* modes, the breathing gas is monitored with regard to an O2 concentration that is too low or an anesthetic gas concentration that is too high, regardless of the respiratory phase. At the same time, the message *Waiting for respiratory phases* is displayed in

the CO₂ waveform, and the symbol "Alarm temporarily inactive" is displayed beside the relevant parameters.

Once 2 breaths have been detected, the message and the symbol disappear. Only then, the respiratory-phase-dependent O2, CO2, N2O, and anesthetic gas alarms are completely activated.

Intelligent alarm behavior

Combined alarms

If multiple alarms occurring at the same time are caused by the same problem, they are combined into one alarm.

Problem	Alarms occurring at the same time	Combined alarm
Several causes of apnea are present.	Apnea (no flow) Apnea (no pressure) Apnea (no CO2)	Apnea
Faults in multiple components. This causes failure of a system function.	Example: Insp. press. sensor failure Exp. press. sensor failure	Ventilator failure

Suppressed alarms

Some low-priority alarms indicate a fault in a measurement function. If this measurement function monitors physiological parameters, alarms based on these parameters will not be generated.

Example:

Fault	Displayed alarm	Non-generated alarm
CO ₂ measurement is faulty	Sample line occluded	Apnea (no CO2)

Configuration

Device settings	149
Factory settings Start settings User-specific settings	149 149 149
Setting the date and time	149
Defining the start settings	150
Adjusting the settings System setup > Screen layout. System setup > Alarms System setup > Therapy System setup > Licenses/Options System setup > System. Resetting the start settings Resetting the consumptions General device information	150 150 152 157 163 164 170 170
Transferring device configurations	171
Importing the configurationExporting the configuration	171 171
Activating software options	172
Loading the activation code from a USB flash drive Entering the activation code Activating the licensed software option	172 172 172
Overview of configurable screen contents.	173
Waveforms and associated parameter fields	173 174

Device settings

Factory settings

Dräger delivers Perseus with factory settings that are used when starting the device for the first time. Service personnel can reset the device to the factory settings.

Start settings

Start settings take effect after every restart of the device or when starting a new case (touching one of the *New adult*, *New pediatric* or *New neonate* buttons). The start settings can be adjusted after the configuration password is entered.

If required, Perseus can be delivered with start settings that may differ from the factory settings.

User-specific settings

User-specific settings can be adjusted by the user without a configuration password. The settings take effect immediately but are discarded at the latest after a device restart.

Setting the date and time



Perseus can adopt the time from a network or from a device connected via MEDIBUS. The time synchronization takes place shortly after switch-on and at regular intervals thereafter.

If time synchronization is not set, the time can be changed manually:

Touch field (A).

Or

Set the date and time in the system setup.

The source for the time synchronization can also be set in the system setup. For more information, see page 164.

Defining the start settings

As soon as a vertical tab is selected, entry of the configuration password is required to access the settings in the **System setup** dialog window.

For information on the configuration password, see page 285.

Adjusting the settings

The following tables show all the setting possibilities in the **System setup** dialog window.

The respective factory settings are marked in **bold** format.

System setup > Screen layout

Vertical tab "General"

Headline/		Setting range			
Parameter	ŕ	Á	•	Description	
Color scheme	Day light; Day dark; Night			Set the color scheme, see page 43	
Screen brightness	10 to 100 80			Set the screen brightness.	

Vertical tab "Views"

Headline/	ine/ Setting range			
Parameter	Ť	Å	•	Description
Rename views	1 Standard; 2 Expert; 3 Expert			Specify names for screen layouts. Touch the button with the respective screen layout, enter the new name on the screen keyboard, and confirm with the ← button or with the rotary knob.
Number of wave- forms (view 3)		3; 4		
Default view	1 Stand	lard; 2 Expert;	3 Expert	Specify the standard view.
Save as system defaults	Current view (only available during operation)			Save current screen layout.
		All views		Save all screen layouts.

Vertical tab "Waveforms"

Headline/		Setting range		
Parameter	Ť	À	*	Description
Sweep speed [mm/s]	6.25 ; 12.5; 25	6.25 ; 12.5; 25	6.25; 12.5 ; 25	Specify sweep speed.
VT scale [mL]	Auto; 0 to 10; 0 to 50; 0 to 150; 0 to 500; 0 to 1000; 0 to 2000	Auto; 0 to 10; 0 to 50; 0 to 150; 0 to 500; 0 to 1000; 0 to 2000	Auto; 0 to 10; 0 to 50; 0 to 150; 0 to 500; 0 to 1000; 0 to 2000	Specify scale for volumeter.
Flow scale [L/min]	Auto -5 to 5 -10 to 10; -30 to 30; -60 to 60; -120 to 120	Auto -5 to 5 -10 to 10; -30 to 30; -60 to 60; -120 to 120	Auto -5 to 5 -10 to 10; -30 to 30; -60 to 60; -120 to 120	Specify scale for flow waveform.
O2 scale [%]	Auto; 0 to 100; 15 to 35; 25 to 45; 35 to 55; 45 to 65; 55 to 75; 65 to 85; 75 to 95; 85 to 105			Specify scale for O2 waveform.
CO2 scale		a]: Auto ; 0 to 6 : Auto ; 0 to 50;	Specify scale for CO2 waveform.	
Paw scale [mbar]; [hPa]; [cmH2O]	Auto ; -5 to 20;	-7.5 to 30; -10	Specify scale for the Paw wave- form.	
Agent prediction scale	Full range	0 to 2 xMAC; 0	Specify scale for anesthetic gas prediction.	
Flow-volume loop	ISC) standard ; Drä	ger	Specify coordinate axes for flow/volume loop display.

Vertical tab "Colors"

Headline/	Setting range			
Parameter	Ť	Á	*	Description
CO2; Paw; Flow, volume	Default color ; color palette with 7 additional colors			Specify parameter colors.
O2; Agent	Defa	ult color; ISO	color	

System setup > Alarms

Vertical tab "Alarm limits"

Headline/		Setting range		
Parameter	·[r	Á	*	Description
Alarm limits for ventil	ation and gases	for each patien	it category	
FiO2_/* [%]	19 to 99; Off	19 to 99; Off	19 to 99; Off 90	Inspiratory oxygen fraction
FiO2 <u>▼</u> / [%]	18 to 98 20	18 to 98 20	18 to 98 20	
etCO2 _/ ^a [%]; [kPa] [mmHg]	0.1 to 9.8; Off 7.0 1 to 75; Off 53	0.1 to 9.8; Off 7.0 1 to 75; Off 53	0.1 to 9.8; Off 7.0 1 to 75; Off 53	Expiratory CO2 concentration
etCO2 */ [%]; [kPa] [mmHg]	Off; 0.0 to 9.7; Off; 0 to 74	Off; 0.0 to 9.7; Off; 0 to 74	Off; 0.0 to 9.7; Off; 0 to 74	
inCO2 _/* [%]; [kPa] [mmHg]	0.1 to 1.4; Off 1.0 1 to 10; Off 8	0.1 to 1.4; Off 1.0 1 to 10; Off 8	0.1 to 1.4; Off 1.0 1 to 10; Off 8	Inspiratory CO2 concentration
Paw high _/* [mbar]; [hPa]; [cmH2O]	5 to 99 40	5 to 99 25	5 to 99 20	Airway pressure If the airway pressure is above the set value for Paw low for longer than 15 s, the alarm Airway press.
Paw low √ [mbar]; [hPa]; [cmH2O]	Auto ; 3 to 97	Auto ; 3 to 97	Auto ; 3 to 97	continuously high is issued. If the setting Auto is selected, the lower alarm limit is automatically adapted to the set value for PEEP.
inSev _/* [%]; [kPa]	0.10 to 9.95 4.40	0.10 to 9.95 5.10	0.10 to 9.95 6.70	Sevoflurane
inSev <u>▼</u> / [%]; [kPa]	Off ; 0.00 to 9.85	Off ; 0.00 to 9.85	Off ; 0.00 to 9.85	

Headline/	Setting range			
Parameter	Ť	Á	*	Description
Alarm limits for ventil	ation and gases	for each patien	t category	
inDes _/*	0.1 to 20.0	0.1 to 20.0	0.1 to 20.0	Desflurane
[%]; [kPa]	12.5	14.5	19.0	
inDes ▼/	Off ;	Off ;	Off ;	
[%]; [kPa]	0.0 to 19.9	0.0 to 19.9	0.0 to 19.9	
inEnf _/*	0.10 to 9.95	0.10 to 9.95	0.10 to 9.95	Enflurane
[%]; [kPa]	3.60	4.10	5.40	
inEnf ▼/	Off ;	Off ;	Off ;	
[%]; [kPa]	0.00 to 9.85	0.00 to 9.85	0.00 to 9.85	
inIso _/*	0.10 to 8.50	0.10 to 8.50	0.10 to 8.50	Isoflurane
[%]; [kPa]	2.40	2.80	3.70	
inIso ▼/	Off ;	Off ;	Off ;	
[%]; [kPa]	0.00 to 8.40	0.00 to 8.40	0.00 to 8.40	
inHal _/*	0.10 to 8.50	0.10 to 8.50	0.10 to 8.50	Halothane
[%]; [kPa]	1.60	1.90	2.40	
inHal ▼/	Off ;	Off ;	Off ;	
[%]; [kPa]	0.00 to 8.40	0.00 to 8.40	0.00 to 8.40	

Vertical tab "Alarm tone volume"

Headline/	Setting range			
Parameter	Ť	Ā	•	Description
Alarm tone volume		10 to 100 40		Set the alarm tone volume.
Minimum alarm tone volume		10 to 100 10		Set the minimum volume with which an alarm tone will be signaled.

Vertical tab "Autoset limits"

Headline/		Setting range					
Parameter	Ť	÷	*	Description			
Offset for "Autoset lin	Offset for "Autoset limits" function						
	Automatic adjustment of the parameters to current measured values, see page 111 By touching the Autoset button, the alarm limits are adjusted so that the upper alarm limit is above the current measured value by at least the percentage or value set here and the lower alarm value is correspondingly below it. Example: In the PC - CMV mode: measured MV: 5 L/min set offset: ±40 % new alarm limits: 7 and 3 L/min						
etCO2 ± [%]	20			In modes with low or no breathing support (Man / Spon, Ext. FGO (optional), CPAP / PSV (optional) and Pause), a further 20 percentage points are added to the configured value.			
Paw + [mbar]; [hPa]; [cmH2O]				The following parameters are taken into account when determining the Paw value: PIP, Pplat, Pinsp, PEEP, ΔPsupp and Phigh. In the Man / Spon and Pause modes, the new alarm limit is at least 25 hPa.			
MV ± [%]	Off; 20 to 80 40			In modes with low or no breathing support (Man / Spon, CPAP / PSV (optional) and Pause), a further 20 percentage points are added to the configured value.			
VTi + [%]		Off; 20 to 80 40		The setting for this parameter has no effect in volume-controlled ventilation modes.			

Vertical tab "Alarm config. 1"

Headline/	Setting range							
Parameter	Ť	Å	*	Description				
General alarm behav	General alarm behavior							
Open "Quick setup" if alarm occurs	On ; Off			Opens the Quick setup window automatically in the event of an alarm.				
"Second agent detected" alarm		On; Off		Issues an alarm when an anesthetic gas mixture is detected.				
"xMAC low" alarm		On; Off		Activate the xMAC low alarm.				
"FiO2 too high for neonates" thresh- old value [%]	50 ; 25 to 90			Set the FiO2 value that is considered critical for neonates. If this value is exceeded for a certain time, the FiO2 too high for neonates alarm will be triggered.				
"FiO2 too high for neonates" alarm after [h:mm]	Off; 0:10 to 9:50 0:15			Set the time after which the FiO2 too high for neonates alarm will be triggered.				
Priority of "Apnea ventilation" alarm	Medium; Low			Specify the alarm priority when the set minimum respiratory rate is not reached in PSV ventilation mode.				
Alarm behavior in "Pa	ause" mode							
Priority of "Pause time expired" alarm	High; Medium; Low			Specify the alarm priority for the alarm that is issued when the time set in the Pause mode has expired.				
Default value for "Timer" [mm:ss]	0:30 to 02:00 2:00	0:30 to 02:00 1:00	0:30 to 02:00 0:30	Specify default time for Pause.				

Vertical tab "Alarm config. 2"

Headline/		Setting range		
Parameter	Ť	Å	*	Description
Deactivate the alarm CPAP/PSV with ΔPs			ı, CPAP,	The CPAP and CPAP / PSV modes are optional.
FiO ₂ high	Yes; No	Yes; No	Yes; No	Specifies the alarm behavior when
MV low		Yes; No		changing to a different ventilation mode.
MV high		Yes; No		These settings apply only to a
xMAC low		Yes; No		change to a ventilation mode with lower or no breathing support (see
etCO2 low	Yes; No			page 259).
etCO2 high	Yes; No			The alarm behavior at the start of the therapy is defined by the con-
inCO2 high	Yes; No			figuration in the vertical tab Alarm
inAgent low	Yes; No			limits.
Deactivate the alarm	limit in cardiac	bypass mode (C	CBM)?	
FiO ₂ high	Yes; No			Specify alarm behavior in CBM
MV low	Yes; No			mode.
MV high	Yes ; No			1
inAgent low		Yes; No		

Vertical tab "Alarm config. 3"

Headline/		Setting range		
Parameter	Ť	Ã	Description	
Alarm limits for "Cylin	nder almost emp	ity"		
O2	[bar]; [kPax100]: Off; 15 to 50 20 [psi]: Off; 218 to 725 290			Specify alarm limits for supply pressure of connected gas cylinders.
Air	[bar]; [kPax100]: Off; 15 to 50 20 [psi]: Off; 218 to 725 290			
N2O	[bar]; [kPax100]: Off; 15 to 40 20 [psi]: Off; 218 to 580 290			

System setup > Therapy

Vertical tab "Vent. 1"

Headline/	Setting range			
Parameter	neter A B		Description	
Default ventilation mode	Buttons with	available ventil Man/Spon	ation modes	Specify the default ventilation mode at start of therapy.
VT and RR start setti	ngs			Specify the tidal volume and respi-
Based on	Patient car	tegory; Ideal bo	dy weight	ratory rate.
Selected: [Patient category]				Specify tidal volume and respiratory rate based on patient cate-
VT [mL]	20 to 2000 500	20 to 2000 150	20 to 2000 50	gory.
RR [1/min]	3 to 100 12	3 to 100 20	3 to 100 30	
Selected: [Ideal body weight] VT [mL] RR [1/min]	20 to 2000 100 kg (220 lbs): 700 75 kg (165 lbs): 520 15 kg (33 lbs): 110 5 kg (11 lbs): 35 3 to 100 100 kg (220 lbs): 10 75 kg (165 lbs): 12 15 kg (33 lbs): 26 5 kg (11 lbs): 32			Specify the tidal volume and respiratory rate based on ideal body weight. Set the tidal volume and the respiratory rate for the supporting points 5; 15; 75; 100 kg (11; 33; 165; 220 lbs). For calculated values for ideal body weight values that lie between these four supporting points, the start settings for tidal volume and respiratory rate are interpolated linearly. For ideal weight values lying outside these supporting points, calculation proceeds with the values of the highest or lowest supporting point.

The start settings for VT and RR influence the start values of the alarm limits for MV high, MV low, and VTi high:

MV high = $VT \times RR \times (1 + offset);$

minimal: 2.0 L/min

MV low = VT x RR x (1 – offset); minimal: 0.3 L/min

VTi high = VT x (1 + offset)

The "offset" value corresponds to the respective offset setting for automatic alarm adjustment. The "offset" value can be set in the **System setup** > **Alarms** > vertical tab **Autoset limits**.

The following applies in volume-controlled modes:

VTi high = $130 \% \times VT$

Vertical tab "Vent. 2"

Headline/	Setting range			D
Parameter	Ť	Å	*	Description
Start settings for ven	tilation			Specify start settings for the venti-
Pmax [mbar]; [hPa]; [cmH2O]	7 to 80 40	7 to 80 30	7 to 80 25	lation.
Pinsp [mbar]; [hPa]; [cmH2O]	3 to 80 15	3 to 80 15	3 to 80 15	
ΔPsupp [mbar]; [hPa]; [cmH2O]	Off; 1 to 80 10	Off; 1 to 80 10	Off; 1 to 80 10	
Insp. term. [%PIF]	5 to 80 25	5 to 80 25	5 to 80 25	
PEEP [mbar]; [hPa]; [cmH2O]	Off; 2 to 35 3	Off; 2 to 35 3	Off; 2 to 35 3	
Slope [s]	0 to 2 0.2	0 to 2 0.2	0 to 2 0.2	
RRmin [1/min]	Off; 3 to 25 6	Off; 3 to 25 10	Off; 3 to 25 15	
I:E	1:50 to 50:1 1:2.0	1:50 to 50:1 1:2.0	1:50 to 50:1 1:1.0	
%Tplat [%]	0 to 60 20	0 to 60 20	0 to 60 20	
Trigger [L/min]	0.3 to 15 4.0	0.3 to 15 2.0	0.3 to 15 1.0	
Sync.	On; Off	On; Off	On; Off	

Headline/		Setting range	
Parameter	Ť	Å	*
Start settings for APF	RV (optional)		
Phigh [mbar]; [hPa]; [cmH2O]	3 to 80 15	3 to 80 15	3 to 80 15
Plow	Off;	Off;	Off;
[mbar]; [hPa];	2 to 35	2 to 35	2 to 35
[cmH2O]	3	3	3
Tlow	0.2 to 10	0.2 to 10	0.2 to 10
[s]	4.0	2.0	1.0
Thigh [s]	0.2 to 10	0.2 to 10	0.2 to 10
	2.0	1.0	1.0
Slope	0 to 2	0 to 2	0 to 2
[s]	0.0	0.0	0.0

Vertical tab "Proced. 1"

Headline/	Setting range				
Parameter	Ť	Á	*	Description	
General settings					
Default maneuver		uitment; Multi-sto Insp./Exp. hold	Defines the default maneuver.		
Reminder [min]	Off , 10 to 180	Off , 10 to 180	Off , 10 to 180	If the value is not set to Variable, a reminder for a maneuver is given after the first switch to a controlled ventilation mode. Defines the time after which, following the end of a One-step recruitment or Multi-step recruitment maneuver, a reminder for a further maneuver is given.	
Layout					
Displayed parame- ter		Cdyn; VT		Defines which additional parameter is displayed in the One-step recruitment and Multi-step recruitment dialogs.	

Vertical tab "Proced. 2"

Headline/		Setting range		
Parameter	Ť	Á	•	Description
Default settings for o	ne-step recruitm			
Pressure [mbar]; [hPa]; [cmH2O]	PEEP + 1 to 80 30	PEEP + 1 to 80 25	PEEP + 1 to 80 20	Sets the pressure level for the maneuver.
Duration [s]	3 to 40 30	3 to 40 15	3 to 40 5	Sets the duration of the maneuver.
Default settings for m	nulti-step recruiti	ment		
PEEP max [mbar]; [hPa]; [cmH2O]	PEEP to 35 20	PEEP to 35 15	PEEP to 35 12	Sets the maximum PEEP pressure for the maneuver.
Pinsp max [mbar]; [hPa]; [cmH2O]	15 to 80 35	15 to 80 30	15 to 80 25	Sets the maximum inspiratory pressure for the maneuver.
ΔPressure [mbar]; [hPa]; [cmH2O]	5 to 30 10	5 to 30 10	5 to 30 10	Sets the pressure difference between Pinsp and PEEP with which pressure-controlled ventilation is carried out. If PEEP max or Pinsp max is reached in the course of the maneuver, the pressure difference will be reduced in steps until the other set value is also reached. The smallest possible value must be at least 3 above the set value for the pressure rise for each level.
Breaths/Step	1 to 20 3	1 to 20 4	1 to 20 5	Number of breaths at a pressure level during the increase or reduction
Breaths@Max	1 to 20 6	1 to 20 8	1 to 20 10	Number of breaths at the Pinsp max pressure level
Pressure rise per step [mbar]; [hPa]; [cmH2O]	2 to 10 5	2 to 10 4	2 to 10 3	Specifies the pressure by which PEEP and Pinsp will be increased in steps.
Pressure decrease per step when PEEP > 15 [mbar]; [hPa]; [cmH2O]	2 to 10 5	2 to 10 4	2 to 10 3	For PEEP >15: Specifies the pressure by which PEEP and Pinsp will be reduced in steps.

Headline/		Setting range		
Parameter	Ť	Á	*	Description
Pressure decrease per step when PEEP ≤ 15 [mbar]; [hPa]; [cmH2O]	1 to 10 2	1 to 10 2	1 to 10 2	For PEEP ≤15: Specifies the pressure by which PEEP and Pinsp will be reduced in steps.

Vertical tab "Fresh gas" (only with electronically controlled gas mixer)

Headline/	Setting range			
Parameter	Ť	Á	*	Description
Start settings for fres	h gas			Select start settings for the freshgas delivery.
FG O2 [%]	21 to 100 100	21 to 100 100	21 to 100 100	Set the O ₂ concentration.
FG flow [L/min]	0.20 to 15.00 2.00	0.20 to 15.00 2.00	0.20 to 15.00 2.00	Set the flow.
Minimal O2 flow (carrier gas: Air) [mL/min]	Off, 50 to 300 Off	Off, 50 to 300 Off	Off, 50 to 300 Off	Set the minimum O2 flow that is delivered when Air is used as carrier gas. Do not set this value too small; recommended is, e.g., 200 for adults, 100 for pediatric patients, and 50 for neonates.
Minimal O2 flow (carrier gas: N2O) [mL/min]	50 to 300 200	50 to 300 200	50 to 300 200	Set the minimal O ₂ flow that is delivered when N ₂ O is used as carrier gas. Do not set this value too small; recommended is, e.g., 200 for adults, 100 for pediatric patients, and 50 for neonates.
Carrier gas	Air; N2O	Air; N2O	Air; N2O	Set the carrier gas.

Vertical tab "Patient"

Headline/		Setting range		
Parameter	Ť	Å	*	Description
Default selection for "Start" dialog	New adult; New pediatric; New neonate			Load patient data: The button that is selected here is preselected the first time the Start dialog window is opened or after changes to the system setup.
Weight [kg]	30 to 300 80 67 to 661 176	5 to 50 25.0 12 to 110 55	0.4 to 10 3.0 0.9 to 22 6.6	Specify start settings for: - Patient weight - Patient height - Patient age
Height [cm] [in]	120 to 300 185 48 to 118 73	50 to 300 100 20 to 118 39	Off; 20 to 80 Off Off; 8 to 31 Off	If the Height parameter is configured to Off for the New neonate patient category, the corresponding setting control will not be displayed in the Start dialog window.
Age [years, Neo: months]	12 to 130 32	0 to 16 8	0 to 24 6	

Vertical tab "General"

Headline/	Setting range			
Parameter	Ť	*		Description
Pinsp changes with PEEP		On; Off		Changing PEEP automatically changes Pinsp, so that the difference between PEEP and Pinsp remains constant.
Ti changes with RR (I:E ratio is locked)		On ; Off		Changing the respiratory rate automatically changes Ti, so that the I:E ratio remains constant.
Breathing system warmer		On; Off		Switch breathing system warmer on or off.
Auto wake-up		On ; Off		The Start dialog window opens automatically when ventilation activity is detected (e.g., by squeezing the breathing bag several times).

Headline/		Setting range		
Parameter	Ť	Å	*	Description
Enable pause mode (deactivates moni- toring mode)		On; Off		Select whether Pause mode or Monitoring mode should be available. On: Pause mode is available. The mode has an adjustable timer, after which an alarm is generated. Pause can be activated from standby and from all ventilation modes. Off: Monitoring mode is available. Monitoring can be activated from Standby mode and from Man/Spon mode.

System setup > Licenses/Options

Vertical tab "Licenses/Options"

Headline/	Setting range			
Parameter	Ť	Ą	•	Description
Licenses for software	options			 Overview of available and active software options.
				 Activating software options. For more information, see page 172
Gas supply				
Disable N2O		On; Off		A device with connectors for nitrous oxide can be configured so that nitrous oxide is no longer displayed, nor can it be selected as a carrier gas. Prerequisite:
				 Current carrier gas is Air
				 Nitrous oxide is not connected or available To make the settings, restart Per- seus if necessary. On: Nitrous oxide cannot be deliv- ered. Off: Nitrous oxide can be deliv- ered.

System setup > System

Vertical tab "General"

Headline/		Setting range		D. a substitut
Parameter	Ť	Á	*	Description
Language		f available langi ish (United Sta		Select language. A flag symbol identifies the tabs that lead to the page with the language settings.
Time source	MEDIBUS 1; N	IEDIBUS 2; NT	P server; None	Select source for time synchronization. Prerequisite: Connected device supports this function.
IP address		4 numeric fields	3	Enter the IP address of the NTP server and confirm.
Date and time	day; month; year hour; minute			Set the date and time. The change is applied on leaving the "General" vertical tab.
Automatic switch to daylight savings time	On; Off			Activate or deactivate automatic switching to daylight saving time.
OR working hours	Hour : Minute to Hour : Minute 6:30 to 18:30			Set the working hours of the operating room. During this time, the gas measurement is kept in a preheated and calibrated state so measured values are available after only a short waiting period. However, this decreases the life span of the patient gas measurement module.
Device name	Rename Device name (up to 16 alphanumeric characters) Perseus			Change system name to, e.g., enter the installation site.
Configuration pass- word	Change password			Change the configuration password.
Reset all pages to		Factory defaults	3	Reset all the settings for the System setup dialog window to the factory settings.

Vertical tab "Units"

Headline/		Setting range		
Parameter	Ť	Å	•	Description
Weight	kg; lbs			Set units.
Height	cm ; in			
Airway pressure	mbar; hPa ; cmH2O			
Supply pressure	t	oar; kPa×100 ; ps		
CO ₂	%; kPa; mmHg			
Volatile agents		%; kPa		

Vertical tab "Auto On/Flush & Dry"

Headline/		Setting range		
Parameter	Å	Ā	•	Description
Auto On				Set the Auto On option. Perseus
Day and time	Monday to Sun Monday to Fri Saturday & Su	day: 6:30	nute	provides a calendar function which, when the device is shut down, suggests the next time the device should be ready for use,
When completed	Standby/ Shut down			depending on the weekday. To activate the automatic ready-foruse time for any desired day, select the corresponding day and set the time, see page 136. Pay attention to the correct setting of the gas supplies to be tested, see page 166.
Flush & Dry	•			
Flush duration between cases [min]		Off; 1 to 5 2		Set the flush cycle duration, e.g., at patient change. Off: Flush function is deactivated.
Flush & Dry duration before shutdown [min]		Off; 5 to 30 10		Set the duration of flushing and drying when the device is shut down. Off: Drying function is deactivated.

Vertical tab "System test"

Headline/	Setting range			
Parameter	Ť	ij	•	Description
General				
Always use walk- through mode		On; Off		When this function is switched on, system tests and leakage tests will always be executed in walk-through mode.
Sample line is con- nected during test	On ; Off			Specify whether the sample line is connected to the Y-piece or to the filter on the Y-piece during the automatic tests. When On is configured, the device can automatically check whether, e.g., breathing hoses or central supply hoses are incorrectly connected. When Off is configured, additional checks are required during the system test.
Test gas supply				
Central O2 supply		On; Off		Select which gas supplies are
Central Air supply		On; Off		tested during the automatic system test.
Central N2O supply		On; Off		If Auto On is configured, only
O2 cylinder		On; Off		check the gas supplies whose actual availability is guaranteed
Air cylinder		On; Off		after the automatic start-up. Set-
N2O cylinder		On; Off		ting both O2 sources to Off - in order to make the system test succeed - is not possible.
Verify O2 delivery		On ; Off		Select whether the system test will check that the O2 supply is actually delivering oxygen. The Sample line is connected during test parameter must also be set to On for this.

Headline/	Setting range			
Parameter	Ť	Á	•	Description
Test breathing circuit				
Test correct assembly of breathing hoses and Y-piece		On ; Off		Select whether the automatic test will check that the breathing circuit is correctly connected. During this test, a check is made as to whether the breathing gas can flow from the inspiratory port via the Y-piece to the expiratory port. If Sample line is connected during test is set to On, this test is carried out fully automatically. If Sample line is connected during test is set to Off, additional manual checks are required.

Vertical tab "Logbook"

Headline/		Setting range		
Parameter	Ť	† • •		Description
A logbook entry with	measured value	Create additional logbook entries with measured values.		
Every	1 min; 2 min; 5 min ; 10 min; 15 min			Create entries regularly.
For all high-priority alarms	On ; Off			Create entries for alarms.
For all medium-pri- ority alarms	On; Off			

Vertical tab "Sound volume"

Headline/	Setting range			
Parameter	† •		*	Description
Alarm tone volume	10 to 100 40			Set the alarm tone volume.
Minimum alarm tone volume	10 to 100 10			Set the minimum volume with which an alarm tone will be signaled.
Breathing sound volume (optional)	Off; 10 to 100 50			Set the breathing sound volume.

Vertical tab "Interfaces"

Headline/		Setting range)	
Parameter	Ť	Å ♣		Description
LAN		Configure the network.		
DHCP		On; Off		Make network settings.
IP address	XXX	(. XXX . XXX .	XXX	When using DHCP, consult with IT personnel to ensure that Perseus
Subnet mask	XXX	C. XXX . XXX .	XXX	is always assigned the same IP
Default gateway	XXX	(.XXX.XXX.	XXX	address by the DHCP server. Restart the device after each change to the network settings. The network settings are not affected by a reset to factory settings. Accept the changes to IP address Subnet mask, or Default gateway with the Apply button. The changes are only active after the device has been restarted.
MAC address				Displays the MAC address.
COM 1				Configure the COM interfaces.
Protocol	IV	IEDIBUS.X; No	ne	A baud rate of 19200 or 38400 is required for transmission of high-
Baud rate	1200; 2400; 48	00; 9600; 192 0	0 ; 38400	speed data, e.g., for waveforms.
COM 2				
Protocol	IV	IEDIBUS.X; No	ne	
Baud rate	1200; 2400; 48	9600; 192 0	0 ; 38400	
USB				Activate or deactivate the USB
USB interface		On; Off		interface.

Vertical tab "Infinity ID"

Headline/		Setting range		
Parameter	† •		*	Description
Monitoring of Infinity	ID accessories			
Breathing circuit		On; Off		Activate or deactivate the
Water trap		On; Off		Infinity ID functionality. On:
Flow sensors		On; Off		Generates a message when
CO2 absorber	On; Off			the maximum period of use is exceeded
				Generates a message when Infinity ID breathing hoses are incorrectly connected Off:
				 Messages are suppressed.
Exchange interval [da	ays]			
Breathing circuit		Off; 2 to 9 2		Set the exchange intervals for Infinity ID accessories.
Water trap		Off; 28		
Flow sensors		Off; 1 to 180 90		
CO2 absorber		Off; 1 to 28 7		

Vertical tab "Service"

Headline/	Setting range			
Parameter	Ť	Á	Description	
		The following functions are available after the appropriate credentials have been entered:		
				 Access to the service dialog

Resetting the start settings

Certain tabs in the **System setup** dialog window have a button for resetting the respective start settings to the factory settings.

Resetting the changes in a dialog window

- 1 Open the corresponding tab.
- 2 Touch the Factory defaults button and confirm.

Resetting the consumptions

The gas consumption can be reset in **Standby** > **System setup** > **System status** > **Consumption**.

Touch the Reset data button and confirm.

General device information

Further information is displayed in *Standby* > *System setup* > *System status*:



- A General information
 - Installed software version
 - Next maintenance date
- **B** QR code for further product information
- Scan the QR code with suitable equipment.


The QR code is decoded into an internet address which enables access to the stated information in a browser.

Transferring device configurations

Configurations can be exported and imported with the aid of a USB flash drive. Configurations from one device version (electronically or mechanically controlled gas mixer) can be transferred to other Perseus devices of the same version in this way.

Prerequisite: The USB flash drive is connected to the USB interface.

 Open System setup > Im/Export config. > Im/Export config. (A).

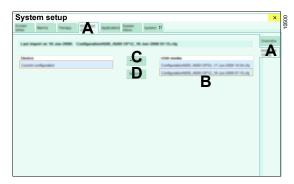

The configurations saved on the USB flash drive are displayed in a list (B). If not all of the configurations can be seen, delete all configurations from the USB flash drive that are not needed or move them to a subdirectory on the flash drive.

The following settings are neither imported nor exported:

- Device name
- Date and time
- IP address

Importing the configuration

- 1 Touch one of the configurations in the list (B).
- 2 Touch the *Import* button (C) and confirm.
- 3 Restart Perseus.

Exporting the configuration

 To export configurations, touch the Export button (D) and confirm.

Activating software options


The following software options require an activation code to be entered, followed by activation:

- Pressure Support
- APRV
- FiO₂ prediction
- Econometer trend
- Low-flow wizard
- Breathing sound
- Lung recruitment
- Lung recruitment reminder
- Inspiration hold / Expiration hold
- Logbook export

Trial licenses for these options are time-limited.

An activation code is linked with the serial number of the respective device and cannot be transferred. The activation codes can either be loaded from a USB flash drive or entered manually.

 Open page System setup > Licenses/Options > Licenses/Options (A).

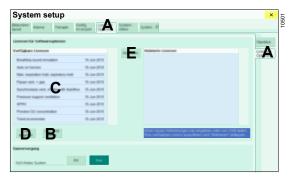

Loading the activation code from a USB flash drive

Prerequisite: A USB flash drive with valid licenses is connected to the USB interface.

• Touch the Load from USB button (B).

The activation codes are read and displayed in the list (C).

Entering the activation code

- 1 Touch the *Enter code* button (D).
- 2 Enter the activation code and confirm with OK.

The license is displayed in the list (C).

Activating the licensed software option

The licensed software options must be activated before they become available.

- Select the corresponding license from the list (C).
- 2 Touch the Activate button (E) and confirm.
- 3 After activating all desired licenses, restart Perseus.

Overview of configurable screen contents

During operation, waveforms and parameter fields are selected in the *Quick setup* window, see page 106.

Paw (4)

Parameters PIP, Pplat, Pmean, PEEP



Waveforms and associated parameter fields

Volume MV, VT, RR








Volume VT, MV, RR



CO₂ in/et, RR



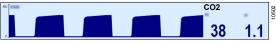

O2





Paw

O2, Δ





Paw (3)

Primary agent

Volume-controlled modes: Parameters PIP, Pplat, PEEP



All other modes: Parameters PIP, Pmean, PEEP


"Empty"






Parameter fields

Compliance trend

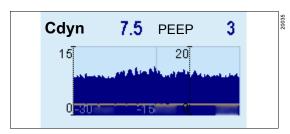


CO₂ in/et, RR

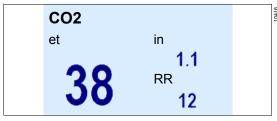


Paw

Compliance

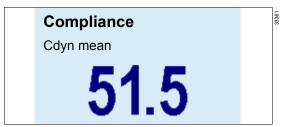



etCO₂

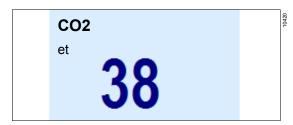

Paw (3)


Volume-controlled modes: Parameters PIP, Pplat, PEEP

All other modes:

Parameters PIP, Pmean, PEEP

Paw		10417
PIP	Pmean 7	
15	PEEP 3	

CO₂ in/et

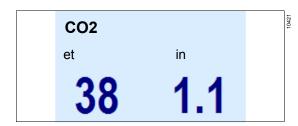

Paw (4)


Parameters PIP, Pplat, Pmean, PEEP:

Paw		104.59
PIP	Pplat	
15	15	
Pmean	PEEP	
7	3	


MVmand, spon

MV Mand Spon 9.0 0.0

ΔVΤ

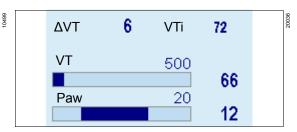

MV, VT, RR

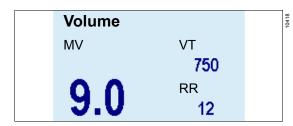

Volume MV, VT, RR



MV×CO₂ trend

WARNING

Risk of incorrect therapy settings


The numerical value displayed for the MV×CO2 parameter is not sufficiently accurate to enable therapy decisions to be made.


Do not make therapy decisions based solely on the displayed numerical value. Only the trend curve can be used for therapeutic decisions.





Ο2, Δ

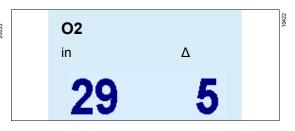

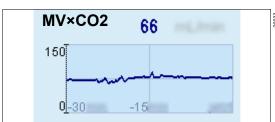


O₂ uptake trend

WARNING


Risk of incorrect therapy settings

The numerical value displayed for the O2 uptake parameter is not sufficiently accurate to enable therapy decisions to be made.


Do not make therapy decisions based solely on the displayed numerical value. Only the trend curve can be used for therapeutic decisions.

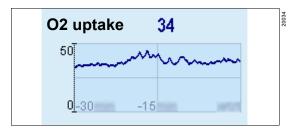

Econometer (optional)

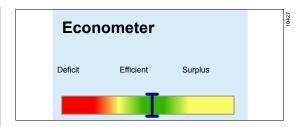

Low-flow wizard (optional)

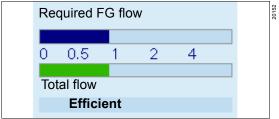

Primary agent



Gases in/et

Gase	s in/et		
	in	et	
02	30	27	
N2O	70	69	
Iso	0.85	0.65	


RR



Gas supply



Volumeter

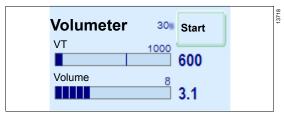

FiO₂ prediction (optional)



Stopwatch



Vaporizer setting (optional)




Timer




Empty




Econometer trend (optional)

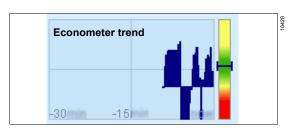

Agent prediction (optional)

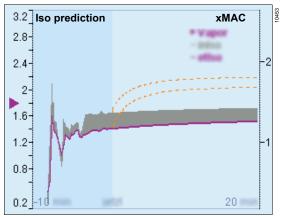

PV loop

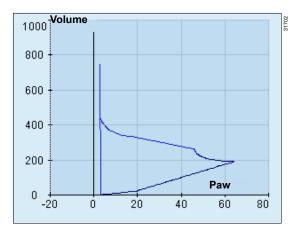

Flow tubes (electronically controlled gas mixer)

WARNING

Risk of mix-up

In some countries, the representation and order of the virtual flow tubes on the screen may differ from that illustrated here.


Pay attention to the labeling of the virtual flow tubes.

WARNING

Insufficient fresh-gas supply

The indication on the virtual flow tubes is intended only as additional information.

Do not use the virtual flow tubes alone when making therapeutic decisions.

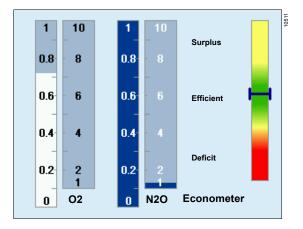


Flow tubes (mechanically controlled gas mixer)

WARNING

Risk of mix-up

In some countries, the representation and order of the virtual flow tubes on the screen may differ from that illustrated here.

Pay attention to the labeling of the virtual flow tubes.

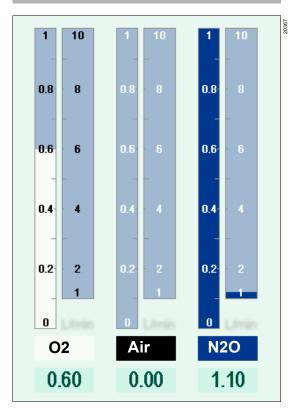

When the Agent prediction or PV loop parameter field is displayed, the flow tubes will be displayed at reduced size:



Troubleshooting

Leakage	181
Possible causes of leakageSystematic localization of leakages	181 181
Power supply failure	182
Mains power supply failureMains power supply failure and empty batteries	182
After power supply is restoredBlown fuses for the auxiliary power sockets	182 182 182
Emergency O2 delivery (electronically controlled gas mixer)	183
Failure of the gas supply	184
Failure of one gas supply Changing an empty gas cylinder Complete failure of the gas supply After the central gas supply is restored	184 184 184 185
Failure of fresh-gas delivery (electronically controlled gas mixture)	186
Ventilator failure	186
Gas measurement failure	187
Flow measurement failure	187
Screen fault/user interface failure	188
Complete failure	188
Problems with the anesthetic gas receiving system (AGS)	189
Replacing the anesthetic gas receiving system (AGS)	189
Problems with cylinder pressure reducers.	190
Problems with the Vapor View option	190
Support request	190
Alarm – Cause – Remedy	191

Leakage

Leakages may result in the system not being operational or being operational with limitations only.

WARNING

Risk due to leaks

As a result of leakage, ambient air may get into the breathing gas and cause a reduction in the depth of anesthesia, or ambient air may become enriched with anesthetic gas and consequently put the user at risk.

Perform the leakage test before using the device. Reduce leakage to a minimum.

WARNING

Risk of insufficient ventilation

Breathing gas may escape because of leakages, with the result that the applied volume is less than the set volume.

Perform the leakage test before using the device. Remedy all leakages.

Possible causes of leakage

- The CO₂ absorber or the CLIC adapter is not securely screwed to the breathing system.
- The APL valve is not correctly fitted to the breathing system or is not set to 30 hPa (cmH2O).
- The breathing bag, the breathing hoses, the Ypiece, or the microbial filter is incorrectly fitted or damaged.
- The flexible arm for the breathing bag (optional) is incorrectly fitted to the breathing system. The sealing ring is soiled or damaged.
- The water trap is not connected.
- The sample line is not connected, is kinked, or is leaking.
- The connections for the sample line are damaged.

- The O-rings on the inspiratory port or expiratory port are damaged, soiled, or missing.
- The flow sensors are incorrectly installed or damaged. The rear O-ring is missing.
- The upper part of the breathing system housing is incorrectly fitted or damaged.
- The valves or seals of the breathing system are damaged.
- The circuit plug is scratched or damaged.
- The filling or emptying connections on the vaporizer are leaking or are open. The vaporizer is incorrectly fitted. The O-ring is missing or damaged. The control dial is not set to the 0 position.

Systematic localization of leakages

To find causes of leakages, isolate individual components from the leakage test.

Component	Measure
Sample line	Remove the sample line and seal the Luer-Lock connector on the Y-piece.
Breathing hoses	Disconnect the breathing hoses. Connect the inspiratory port and expiratory port with a hose that is known to be without leakages. Connect the breathing bag directly to the breathing system.
Vaporizers	Remove the vaporizers.

- 1 Perform a leakage test, see page 123, and use the leakage assistant if necessary.
- 2 Contact service personnel if the leakages cannot be localized.

Power supply failure

Mains power supply failure

If mains power fails, Perseus automatically switches to the internal battery. A fully charged battery will maintain operation for at least 30 minutes.

Remaining battery charge is displayed on the status display.

The breathing system warmer is deactivated during battery operation. The peak inspiratory pressure *PIP* may be limited, but it will be at least 55 hPa (cmH2O).

CAUTION

Risk of device malfunction

If mains power fails, devices connected to the auxiliary power sockets are not supplied from the internal battery.

Ensure an alternative power supply for connected devices.

Mains power supply failure and empty batteries

If mains power fails and the batteries are empty, a signal tone is emitted. Manual ventilation and spontaneous breathing remain available. O2 and anesthetic agent can still be delivered using the emergency O2 delivery (with electronically controlled gas mixer) or the flow control valves (with mechanically controlled gas mixer) and connected vaporizers.

The following are not available:

- Ventilator
- Electronically controlled gas mixer
- Device monitoring and patient monitoring

WARNING

Risk of patient injury

If all power sources fail, the screen goes dark and automatic ventilation ends.

Manually ventilate the patient.

Further procedures:

- Check vaporizer setting.
- 2 Electronically controlled gas mixer: Use emergency O2 delivery.

Mechanically controlled gas mixer: Close the Air flow control valve and the N2O flow control valve and use only O2 as fresh gas.

- 3 Monitor the O2 flow on the O2 flowmeter (electronically controlled gas mixer) or on the total flow tube (mechanically controlled gas mixer).
- 4 Manually ventilate the patient.
- **5** Ensure corresponding substitute monitoring.

After power supply is restored

- 1 Restart the device, see page 80.
- 2 Charge empty batteries for at least 8 hours.
- 3 Check the displays for mains voltage and battery on the status display.

Blown fuses for the auxiliary power sockets

- 1 Disconnect mains power supply.
- 2 Remedy the malfunction.
- **3** For devices without isolation transformer: Replace the fuses.

For devices with isolation transformer:

- If necessary, wait 2 minutes for the circuit breaker to cool.
- Press the key on the circuit breaker back in or operate the switch on the isolation transformer.
- 4 Restore mains power supply.


Emergency O₂ delivery (electronically controlled gas mixer)

WARNING

Risk of increased anesthetic agent delivery

When the emergency O2 delivery (Add. O2) is in use, anesthetic agent continues to be delivered into the breathing system. When the emergency O2 delivery is used during low-flow anesthesia or minimal-flow anesthesia, an increased quantity of anesthetic agent may enter the breathing system. This may lead to an increased anesthetic gas concentration.

Carefully monitor the gas mixture.

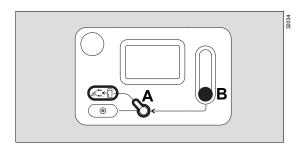

- 1 Check vaporizer setting.
- 2 Set the O₂ switch (A) upwards to the *Add. O*₂ position.
- 3 Open the flow control valve (B) on the O2 flowmeter and set the desired flow. This O2 flow flows through the vaporizer.

Failure of the gas supply

A failure of the central gas supply can result in simultaneous device malfunctions on all systems connected to it.

Perseus signals an alarm if the gas supply for the gases O₂, Air, or N₂O (optional) fails.

WARNING

Risk of contaminating the gas supply

When the central gas supply is connected, the smallest internal leakage can cause contamination of the supply gases.

If the central gas supply fails during operation, disconnect the hoses for the failed gas from the central supply.

- Open the corresponding gas cylinder.
- Restore central gas supply.

Electronically controlled gas mixer only: If the central gas supply for a gas fails and there is no sufficiently filled gas cylinder connected (see page 31 "Status display"), a substitute gas is used:

Failed gas	Substitute gas
O2	100 % Air
N2O	100 % O2
Air	100 % O2

The level of the fresh-gas flow remains constant.

Failure of one gas supply

Operation of fresh-gas delivery is still possible when supply of one gas fails. If, e.g., N2O fails, with the electronically controlled gas mixer, Air or 100 % O2 can be set as carrier gas. With the mechanically controlled gas mixer, open the corresponding flow control valve of the substitute gas.

Changing an empty gas cylinder

- 1 Close the valve of the empty gas cylinder.
- 2 Completely use up or completely vent any gas remaining in the pressure reducer and in the hose between Perseus and the gas cylinder. If there is no patient connected, venting can be performed as follows:
 - Disconnect the central O₂ supply.
 - Open the flow control valve of the O2 flowmeter. Wait until gas is no longer flowing.
 - Close the flow control valve of the O2 flowmeter again.
- 3 Unscrew the pressure reducer from the gas cylinder valve.
- 4 Replace the gas cylinder with a full gas cylinder.
- 5 Connect the pressure reducer to the new filled gas cylinder, see page 58.
- 6 Open the valve of the filled gas cylinder.

Complete failure of the gas supply

If the central gas supply for O2 and Air fails at the same time and there are no sufficiently filled gas cylinders connected, operation can be continued in automatic ventilation modes. This is possible because the ventilator does not require a drive gas

- 1 Remove the breathing bag.
- Continue automatic ventilation.

When the breathing bag is removed, the missing fresh-gas volume will automatically be filled by ambient air. This is likely to cause the alarm *Fresh gas low or leakage* to be triggered.

WARNING

Risk of patient recovering consciousness

If the gas supply fails completely, further operation takes place through gas supply of the anesthesia machine with ambient air. Anesthetic agents are no longer delivered and the inspiratory anesthetic gas concentration in the breathing gas decreases.

Monitor the gas mixture carefully and use intravenous anesthetic agents if need be.

CAUTION

Risk of increased anesthetic gas concentration in the ambient air

If the breathing bag is not connected, expiratory anesthetic gases may escape from the breathing system.

Ensure adequate circulation of the ambient air.

After the central gas supply is restored

- Connect the compressed gas hoses to the terminal units.
- 2 Close the gas cylinder valve on the corresponding gas cylinder again.

On devices that are equipped with Advanced Cylinder Support, the gas cylinder valves can remain open.

Failure of fresh-gas delivery (electronically controlled gas mixture)


If the fresh-gas delivery has failed, the emergency O2 delivery can be used to deliver oxygen and anesthetic agent. The current ventilation mode and the fresh-gas deficiency detection remain active.

WARNING

Risk of patient injury

If the gas mixer fails, no fresh gas is delivered.

Check vaporizer setting. Supply the patient with O2. Use emergency O2 delivery.

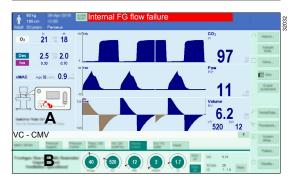

In the event of a fault, figures and instructions showing how to start the emergency O2 delivery are displayed in areas (A) and (B).

The emergency O2 delivery is started as follows:

- 1 Set the O2 switch (A) upwards to the Add. O2 position. (Pay attention to the figure on the screen.)
 - The *Internal FG flow failure* alarm will then be automatically downgraded.
- 2 Open the flow control valve on the O2 flowmeter and set the desired flow. This O2 flow flows through the vaporizer.
- 3 Check vaporizer setting.
- 4 Continuously monitor the O2 flow for the emergency delivery.
- If necessary, perform ventilation with ambient air, see "Complete failure of the gas supply" on page 184.
- If necessary, ventilate the patient with the manual resuscitator.

Ventilator failure

If the ventilator fails, only manual ventilation or spontaneous breathing remain possible. No other ventilation modes can be selected. The fresh-gas delivery remains ready for operation.

- 1 Switch to the *Man/Spon* ventilation mode.
- 2 Manually ventilate the patient.

Gas measurement failure

CAUTION

Risk due to gas measurement failure

If the gas measurement fails, the patient can no longer be adequately monitored.

- Ensure corresponding substitute monitoring.
- Check sample line and water trap for damage or blockage and resolve these as needed.
 Observe the exchange intervals.

 Arrange for appropriate substitute monitoring conforming to ISO 80601-2-55.

Flow measurement failure

If the flow measurement fails, normally the therapy can be continued. Thus, the flow sensors can be exchanged during the next **Standby**. There may be limitations in measured parameters or selection of therapy.

WARNING

Risk due to malfunction of the inspiratory flow measurement

If the inspiratory flow sensor malfunctions, the device automatically switches to pressure-controlled ventilation. If there is an additional malfunction of the inspiratory pressure sensor, the system automatically switches to the non-synchronized *PC - CMV* mode.

In either case, check the ventilation settings and adjust as needed.

- 1 Replace the flow sensor: To remove the flow sensors, follow the steps described in chapter "Inserting the flow sensors and the ports" in reverse order, see page 65. Then install the new flow sensors.
- 2 Perform the leakage test, see page 123.

Screen fault/user interface failure

The screen does not respond to operation. It has failed or the screen display is faulty.

- Switch Perseus to *Standby* mode: Press the ⁽⁾ button and confirm with the rotary knob.
- 2 Use emergency O2 delivery, see page 183.
- 3 Check vaporizer setting.
- 4 Manually ventilate the patient.
- **5** Ensure corresponding substitute monitoring.

Complete failure

The device no longer responds to operation.

- 1 Turn off the device with the main switch, see page 27.
- 2 Manually ventilate the patient.
- 3 Perform a start-up for operation in case of emergency, see page 82.

WARNING

Risk of device malfunction

If the breathing bag does not fill with fresh gas, the patient cannot be sufficiently ventilated.


- Check the oxygen supply and, if necessary, open the gas cylinder valves.
- If fresh gas still is not delivered or manual ventilation is not possible, close the flow control valve of the O2 flowmeter.
- Disconnect the patient from the device and use a replacement device!

Problems with the anesthetic gas receiving system (AGS)

Fault	Cause	Remedy
Flow indicator beneath the "restricted range"	The suction power of the ejector in the terminal unit of anesthetic gas scavenging system (AGSS) is insufficient.	Have the function of the AGSS terminal unit checked. Observe related instructions for use.
	Particle filter contaminated or blocked.	Replace the anesthetic gas receiving system (AGS) or have service personnel replace the particle filter.
Flow indicator above "normal range"	The suction power of the ejector in the AGSS terminal unit is too high.	Have the suction power of the ejector in the AGSS terminal unit adjusted to the working range of the AGS.
	Particle filter missing.	Replace the anesthetic gas receiving system (AGS) or have service personnel install the particle filter.

Replacing the anesthetic gas receiving system (AGS)

View from below:

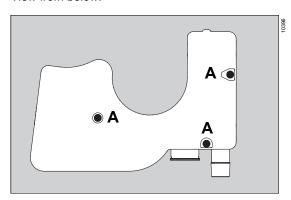

Disassembling the anesthetic gas receiving system (AGS)

- **1** Disconnect the scavenging hose.
- 2 If necessary, remove the transfer hose for the non-rebreathing system.
- 3 Unfasten the 3 screws (A)

4 Dispose of the AGS.

Assembling the anesthetic gas receiving system (AGS)

 Follow the above steps in reverse order to assemble the new AGS.

Problems with cylinder pressure reducers

Fault	Cause	Remedy
The connection between gas cylinder and pressure reducer leaks.		Replace the sealing ring.
The outlet pressure rises; the relief valve relieves the outlet of the pressure reducer.	The valve seat is soiled or damaged.	Close the gas cylinder valve. Repair by service personnel
Leakage in the housing area	The diaphragm is faulty.	Repair by service personnel

Problems with the Vapor View option

Fault	Cause	Remedy
Dräger Vapor 3000/ D- Vapor 3000 is not detected.	Vaporizer is not correctly mounted and locked.	Correctly mount and lock vaporizer.
	Sensor unit is faulty.	Have item repaired by service personnel.
Control dial setting position or filling level of a Dräger Vapor 3000/	Sensor unit or vaporizer is soiled.	Clean sensor unit and remove foreign matter.
D-Vapor 3000 is not detected.	Sensor unit or vaporizer is damaged.	Have item repaired by service personnel.
Dräger Vapor 3000/ D-	Control dial is in position T.	Turn control dial to position 0.
Vapor 3000 is not illuminated.	Illumination unit is soiled.	Clean illumination unit and remove foreign matter.
	Illumination unit is faulty.	Have item repaired by service personnel.
Illumination is turned on, even though no Dräger Vapor 3000/ D-Vapor 3000 is fitted.	Vapor View option is faulty.	Have item repaired by service personnel.

Support request

If the device is configured for remote maintenance, device information can be sent to Dräger in the event of a problem. Proceed as follows to send a support request to Dräger:

- 1 In the **Standby** mode, touch the **Tests...** button.
- 2 Touch the **Request support** button.

Alarm - Cause - Remedy

Alarm messages are displayed in hierarchal form in the alarm message field of the header bar, see page 140.

Different background colors indicate the priority levels of the alarms.

In the *Current alarms* table and *Alarm history* table, the priority of the alarm messages is also indicated by exclamation marks.

High	!!!	Red
Medium	!!	Yellow
Low	!	Cyan

In order to classify the alarms within an alarm priority, internal priority numbers are given in the table below. The most critical alarm is given the number 255. Lower numbers indicate a lower alarm priority.

The following table lists the alarm messages in alphabetical order. If an alarm occurs, the table helps to quickly identify causes and remedies. The possible causes and remedial measures should be looked through in the order they are listed until the alarm is resolved.

Some alarms appear in this table several times with different priorities because their priority may change under certain conditions, see page 144.

Pr	iority	Alarm	Cause	Remedy
!	100	Absorber disconnected?	Infinity ID CLIC absorber is	Check absorber.
			not correctly connected.	Use "ALARM RESET" to acknowledge alarm.
!!	100	"Add. O2" activated	O2 switch is set to "Add. O2".	Close the flow control valve of the O2 flowmeter. Set the O2 switch to "Aux. O2".
				Use "ALARM RESET" to downgrade alarm priority.
!	255	Agent measurement failed	There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			There is an internal failure.	Use an alternative gas measurement system.
				If the problem persists, call DrägerService.
!	255	Agent measurement not available	The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
			There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			Ambient temperature is too high.	Check ambient conditions.

Pr	iority	Alarm	Cause	Remedy
!	255	Agent measurement temporarily inaccurate	The sensor has not yet warmed up.	Wait for the automatic calibration.
			The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
!	220	Air cylinder almost empty	The cylinder is almost empty.	Replace the cylinder. Use the central supply.
!	255	Air cylinder empty	The cylinder is empty or closed.	Replace the cylinder. Use the central supply.
!	190	Air cylinder sensor?	Cylinder pressure sensor is not connected.	Check if cylinder pressure sensor is connected.
				Use "ALARM RESET" to acknowledge alarm.
!!	150	Air FG flow measurement	The measurement system for	Only use O2 as fresh gas.
		failed	the Air fresh-gas flow has failed.	Use the total flow tube to check the fresh-gas flow. Adjust the fresh-gas flow so that it equals or exceeds the minute volume.
!!!	110	Air supply low	Central supply pressure and cylinder pressure are low.	Check central Air supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	255	Airway press. continuously	Airway pressure has been	Check ventilation settings.
		high	continuously high.	Check breathing hoses, breathing system, and anes- thetic gas scavenging sys- tem.
				In Man/Spon mode, check the APL valve setting.
				Check lower alarm limit for airway pressure.
!!!	255	Airway pressure high	irway pressure high The upper alarm limit for the airway pressure has been exceeded. The applied inspiratory pressure is higher than the set value.	Check patient condition.
				Check ventilation settings.
				Check alarm limit.
			Breathing hoses are blocked or the tube is kinked.	Check breathing circuit and tube.

Pri	ority	Alarm	Cause	Remedy
!!!	255	Airway pressure negative	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			Suction maneuver during ventilation.	Check the bronchial suction system.
			Failure of the anesthetic gas scavenging system.	Check anesthetic gas scavenging system.
!!	10	Airway pressure not achieved	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
!!!	220	Apnea	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
!!	255	Apnea	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
!!!	220	Apnea (no CO ₂)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
			Sample line is not connected.	Connect sample line to breathing circuit.
!!	255	Apnea (no CO ₂)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
			Sample line is not connected.	Connect sample line to breathing circuit.
!!!	220	Apnea (no flow)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
			Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			The breathing hoses are blocked or leaking.	Check breathing circuit and tube.

Pr	iority	Alarm	Cause	Remedy
!!	255	Apnea (no flow)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
				Check spontaneous breathing ability of the patient.
			Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			The breathing hoses are blocked or leaking.	Check breathing circuit and tube.
!!!	220	Apnea (no pressure)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
			Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			The breathing hoses are blocked or leaking.	Check breathing circuit and tube.
!!	255	Apnea (no pressure)	No breathing or ventilation.	Start manual ventilation!
				Check ventilation settings.
			Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			The breathing hoses are blocked or leaking.	Check breathing circuit and tube.
!!	0	Apnea Ventilation	No inspiratory effort of the patient detected.	Check spontaneous breathing ability of the patient.
				Adjust the setting for "Trig- ger".
				Change to pressure-controlled or volume-controlled ventilation mode.
				Use "ALARM RESET" to downgrade alarm priority.
!	0	Apnea Ventilation	No inspiratory effort of the patient detected.	Check spontaneous breathing ability of the patient.
				Adjust the setting for "Trig- ger".
				Change to pressure-controlled or volume-controlled ventilation mode.
!	100	"Audio paused" key stuck	Key is stuck or was pressed for more than 10 seconds.	Ventilation performance is not affected.
				If the problem persists, call DrägerService.

Pr	iority	Alarm	Cause	Remedy
!!	100	Backup speaker failure	The backup speaker for	Call DrägerService.
			alarm tones is faulty.	Use "ALARM RESET" to acknowledge alarm.
!	170	Bag pressure sensor failure	Sensor calibration failed.	Manually check the filling level of the breathing bag.
				Perform the system test.
!	180	Battery charge low	The battery charge is low and	Restore mains power supply.
			the mains power supply is not available.	The breathing system warmer has been switched off. Check the breathing circuit for condensate. Increase the freshgas flow if necessary.
!!	30	Battery charge very low	The battery charge is critical and the mains power supply is not available. The device	Make sure that the mains power supply is correctly connected.
		will shut down in the next 5 minutes.		The breathing system warmer has been switched off. Check the breathing circuit for condensate. Increase the freshgas flow if necessary.
				Once the battery has been depleted, ventilate the patient manually.
!!	170	Battery failure	The battery is faulty. If the mains power supply fails, the device will switch off immediately.	Call DrägerService.
!	100	Battery temperature high	The battery temperature is high. Charging of the battery has been suspended to pro-	Ensure that the system is connected to the mains power supply.
			tect it from damage.	Check the ambient temperature.
!	135	Breathing bag too small?	The breathing bag is too small. The tidal volume cannot be delivered.	Use correct breathing bag.
			The breathing bag hose is kinked, too long, or too thin.	Check the hose and replace if necessary.
				Increase the fresh-gas flow.
				Use "ALARM RESET" to acknowledge alarm.
!	100	Breathing circuit expired	Accessory has been used too long.	Replace the accessory if necessary.
				Use "ALARM RESET" to acknowledge alarm.

Pri	ority	Alarm	Cause	Remedy
!!!	255	Breathing system temp. high?	The breathing system warmer is faulty.	Check the inspiratory breathing gas temperature as close to the Y-piece as possible.
				Use longer inspiratory hose.
				Remove breathing system cover.
				Turn off the breathing system warmer.
				Turn off the device with the main switch. Use "Add. O2". Ventilate the patient manually.
				If the breathing gas temperature is too high, ventilate the patient with a manual resuscitator.
!	100	Breathing system warmer failure	The breathing system warmer is faulty.	Check breathing circuit for condensation. Increase freshgas flow if necessary.
				Call DrägerService.
			An internal temperature sensor is faulty.	Call DrägerService.
!	100	Cardiac bypass mode still active?	A significant minute volume was measured during cardiac bypass mode.	Deactivate the cardiac bypass mode.
				Use "ALARM RESET" to acknowledge alarm.
!	100	Central Air supply high	The central supply pressure is high. The gas delivery might fail.	Check central supply.
!	255	Central Air supply low	Central supply pressure is low.	Check central supply.
!	100	Central N2O supply high	The central supply pressure is high. The gas delivery might fail.	Check central supply.
!	255	Central N2O supply low	Central supply pressure is low.	Check central supply.
!	100	Central O2 supply high	The central supply pressure is high. The gas delivery might fail.	Check central supply.
!!!	210	Central O2 supply low	Central supply pressure is low.	Check central O2 supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!	255	Central O2 supply low	Central supply pressure is low.	Check central supply.

Pr	iority	Alarm	Cause	Remedy
!	100	CO2 absorber expired	Accessory has been used too long.	essary.
				Use "ALARM RESET" to acknowledge alarm.
!	255	CO2 measurement failed	There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			There is an internal failure.	Use an alternative gas measurement system.
				If the problem persists, call DrägerService.
!	255	CO2 sensor accuracy low	The sensor has not yet warmed up.	Wait for the automatic calibration.
			The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
!	0	Communication failure	The network connection	Re-establish the connection.
			could not be established.	Use "ALARM RESET" to acknowledge alarm.
				Check the connections from/to Connectivity Converter CC300.
			A network certificate has expired or Connectivity Converter CC300 is faulty.	Call DrägerService.
			In the network, there is another anesthesia machine with the same integrated system ID.	Check the integrated system ID.
!!	135	Cooling fan failure	An internal fan for evacuating gases is faulty.	To prevent potential damage, switch off system at your earliest convenience. Increased risk of fire.
				Call DrägerService.
!!	140	CPAP changed to "Off"	The set CPAP pressure could not be achieved due to leakage. The system has	Check the breathing circuit for tight connections and leakages.
			changed the CPAP setting to	Increase fresh-gas flow.
			"Off".	Reapply the CPAP setting.
				Use "ALARM RESET" to acknowledge alarm.

150 Emergency air inlet activated There was not enough gas to ventilate the patient. To maintain a minimum ventilation, the device has started to use ambient air. Increase fresh-gas flow Check the breathing big with 02 flush. Increase fresh-gas flow Check the breathing big tight connections and le ages. Use correct breathing big oright connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big oright connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big tight connections and le ages. Use correct breathing big with calles to the pages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use correct breathing big with connections and le ages. Use lock the breathing big with connections and le ages. Use lock the breathing big with connections and le ages. Use lock the breathing big with connections and le ages. Use lock the breathing big with connections and le ages. Use lock the breathing big with connections and le ages. Use lock the breathing big with connections and leads. Use lock the breathing big with connections and leads. Use lock the breathing big with conn	
the device has started to use ambient air. Check the breathing cirr tight connections and leages. Use correct breathing be check ventilation. Fasure that a suitable stute monitoring is availate to accessory essary. Use "ALARM RESET" to acknowledge alarm. Perform the leakage test. The breathing system has been used too accessory essary. Use "Outendary to accessory essary. It acknowledge alarm. Perform the settendary to accessory essary. Use "Outendary to accessory essary. Use "ALARM RESET" to acknowledge alarm. Perform the leakage test. It acknowledge alarm. Perform the settendary to accessory essary. Use "ALARM reset" t	, e.g.,
ambient air. Check the breathing circ tight connections and le ages.	<i>i</i> .
!! 135 etCO2 high etCO2 has exceeded the upper alarm limit. Check ventilation. !! 135 etCO2 low etCO2 is below the lower alarm limit. Check ventilation. !! 100 Exp. press. sensor failure Sensor calibration failed. Ensure that a suitable stute monitoring is availad Perform the system test ! 100 Expiratory flow sensor expired Accessory has been used too long. Replace the accessory essary. ! 190 Expiratory flow sensor not calibrated. The breathing system has been replaced or disconnected since last calibration. Perform the leakage test. !! 190 Expiratory flow sensor not calibrated. The breathing system has been replaced or disconnected since last calibration. Perform the leakage test. !!! 200 External fresh-gas outlet failure? Failure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode. Use "O2+" button to det mine flow direction of freathing tem or breathing bag fill tem or breathing bag fill external fresh-gas outlet to another ventilation mode.	
upper alarm limit. !! 135 etCO2 low etCO2 is below the lower alarm limit. !! 100 Exp. press. sensor failure Sensor calibration failed. !! 100 Expiratory flow sensor expired	ag.
100 Exp. press. sensor failure Sensor calibration failed. Ensure that a suitable stute monitoring is availad Perform the system test. 100 Expiratory flow sensor expired Accessory has been used too long. Replace the accessory essary. 190 Expiratory flow sensor not calibrated The sensor is not calibrated. The breathing system has been replaced or disconnected since last calibration. Perform the leakage test 190 Expiratory flow sensor not calibrated. The breathing system has been replaced or disconnected since last calibration. Replace the flow sensor form the leakage test. 190 External fresh-gas outlet failure? Failure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode. Sensor calibration 100 Expiratory flow sensor not calibrated. The breathing system has been replaced or disconnected since last calibration. 100 External fresh-gas outlet failure? Failure when switching to external fresh-gas outlet. Great a suitable state to the monitoring is availad to the failure monitoring is availad to the monitoring is availated.	
tute monitoring is availated Perform the system test Replace the accessory essary. I 190 Expiratory flow sensor not calibrated calibrated Expiratory flow sensor not calibrated. The sensor is not calibrated. The breathing system has been replaced or disconnected since last calibration. Failure of the flow sensor. External fresh-gas outlet failure? External fresh-gas outlet failure when switching to external fresh-gas outlet to another ventilation mode. Expiratory flow sensor not calibrated. The sensor is not calibrated. The breathing system has been replaced or disconnected since last calibration. Failure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode.	
! 100 Expiratory flow sensor expired Accessory has been used too long. Replace the accessory essary. ! 190 Expiratory flow sensor not calibrated. The breathing system has been replaced or disconnected since last calibration. Perform the leakage test. !!! 200 External fresh-gas outlet failure? Failure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode. Use "O2+" button to det mine flow direction of free gas: -If internal breathing to external fresh-gas outlet to another ventilation mode.	
essary. 190 Expiratory flow sensor not calibrated The sensor is not calibrated The breathing system has been replaced or disconnected since last calibration. 200 External fresh-gas outlet failure? Failure when switching to external fresh-gas outlet to another ventilation mode. Sarary. Use "ALARM RESET" to acknowledge alarm. Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perform the leakage test Perfor	t.
190 Expiratory flow sensor not calibrated The sensor is not calibrated. The breathing system has been replaced or disconnected since last calibration. Replace the flow sensor form the leakage test.	if nec-
calibrated The breathing system has been replaced or disconnected since last calibration. Failure of the flow sensor. Replace the flow sensor form the leakage test. Pailure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode.	:О
form the leakage test. !!! 200 External fresh-gas outlet failure? External fresh-gas outlet external fresh-gas outlet. Failure when switching to external fresh-gas outlet. Failure when switching from external fresh-gas outlet to another ventilation mode. form the leakage test. Use "O2+" button to det mine flow direction of from gas: -If internal breathing external fresh-gas outlet to external fresh-gas outlet.	st.
failure? external fresh-gas outlet. mine flow direction of fr Failure when switching from gas: -If internal breathin external fresh-gas outlet to tem or breathing bag fill another ventilation mode. external fresh-gas outle	r. Per-
not availableIf gas flow of the external fresh-gas let, external fresh-gas of can be used. Internal browning system can only be when the breathing bag connected (ventilation wambient air only).	resh ng sys- l, et is ws out- soutlet reath- used j is not with
Check fresh-gas setting Call DrägerService.	<u>jo.</u>
!! 50 Filling level of vaporizer low The filling level of the active vaporizer is low. Check filling level. Refill necessary.	l if
Use "ALARM RESET" to downgrade alarm priorities.	
! 100 Filling level of vaporizer low The filling level of an inactive Check filling level. Refill vaporizer is low.	

Pri	ority	Alarm	Cause	Remedy
!!	10	FiO2 high	FiO ₂ has exceeded the upper alarm limit.	Check FG O2.
!!! 255	255	55 FiO2 low	FiO2 has fallen below the lower alarm limit.	Check O2 concentration and fresh-gas settings.
				Check the breathing system for high leakages.
				Check O2 supply.
!!	135	FiO ₂ too high for neonates	FiO2 has exceeded the threshold value for longer	Check O2 concentration and fresh-gas settings.
			than the time configured in the system setup.	Use "ALARM RESET" to acknowledge alarm.
!!	30	Flow control valve still open	At least one flow control valve is still open.	Close all flow control valves.
!	80	Flow sensor calibration required	The flow sensors have not been calibrated for more than 24 hours or since the device was turned on.	Perform the leakage test.
!!!	100	Fresh gas low or leakage	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Refill the breathing system immediately, e.g., with O2 flush.
				Check fresh-gas settings and position of breathing bag.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
				Check tube or mask.
!!	150	Fresh gas low or leakage	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
				Check tube or mask.
				Use "ALARM RESET" to downgrade alarm priority.
!	170	Fresh gas low or leakage	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
				Check tube or mask.
!!	100	Fresh-gas flow high	The total fresh-gas flow is greater than 15 L/min.	Reduce fresh-gas flow.

Pr	iority	Alarm	Cause	Remedy
!!	60	Fresh-gas flow inaccurate	The delivered fresh-gas flow differs from the set fresh-gas flow.	Make sure that sufficient fresh-gas and anesthetic agent are delivered.
				Check the measured gas concentrations.
				Use "ALARM RESET" to downgrade alarm priority.
!!	50	Fresh-gas flow inaccurate	The accuracy of the fresh-gas flow measurement is	Use the total flow tube to verify the current fresh-gas flow.
			reduced.	Check the measured gas concentrations.
				Use "ALARM RESET" to downgrade alarm priority.
			Perform the system test.	
				If the problem persists, call DrägerService.
!!	100	00 Gas sensor failure	Gas sensor failure The patient-gas measurement module has failed.	Use an alternative gas measurement system.
				Use "ALARM RESET" to downgrade alarm priority.
				Call DrägerService.
!!	255	Hose connected to wrong port	A breathing hose is not correctly connected.	Connect breathing hoses correctly.
!!	100	Hose does not fit to pat.	The detected breathing hose	Use compatible accessory.
		category	is not suitable for the selected patient category.	Use "ALARM RESET" to acknowledge alarm.
!	100	Hose does not fit to pat.	The detected breathing hose	Use compatible accessory.
		category	is not suitable for the selected patient category.	Use "ALARM RESET" to acknowledge alarm.
!!	150	inCO2 high	Soda lime is depleted.	Check soda lime.
				Increase fresh-gas flow.
				Check fresh-gas settings.
			There is internal leakage in the breathing system or in the coaxial breathing hose.	Replace the breathing system or the coaxial breathing hose.
			Gas measurement is inaccurate due to high respiratory rate.	Adjust alarm limits if necessary.
			Large dead space.	Check the ventilation settings and the breathing circuit.
!	100	Infinity ID breathing circuit	An incompatible accessory is	Check accessory.
		not compatible	connected.	Use "ALARM RESET" to acknowledge alarm.

Pr	iority	Alarm	Cause	Remedy
!	100	Infinity ID CO2 absorber	An incompatible accessory is	Check accessory.
		not compatible	connected.	Use "ALARM RESET" to acknowledge alarm.
!	100	Infinity ID water trap not	An incompatible accessory is	Check accessory.
		compatible	connected.	Use "ALARM RESET" to acknowledge alarm.
!!	100	Insp. press. sensor failure	Sensor calibration failed.	Change to Man/Spon mode and ventilate manually.
				Perform the system test.
!!!	255	Inspiratory desflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!!	255	Inspiratory desflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!	75	Inspiratory desflurane low	gas concentration is below the lower alarm limit.	Check vaporizer and freshgas settings.
				Refill the vaporizer.
				Check the breathing system for high leakages.
			The soda lime has dried out.	Replace the soda lime.
!!!	255	Inspiratory enflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!!	255	Inspiratory enflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!	75	Inspiratory enflurane low	gas concentration is below	Check vaporizer and freshgas settings.
			the lower alarm limit.	Refill the vaporizer.
				Check the breathing system for high leakages.
			The soda lime has dried out.	Replace the soda lime.
!	100	Inspiratory flow sensor expired	Accessory has been used too long.	Replace the accessory if necessary.
				Use "ALARM RESET" to acknowledge alarm.
!	190	Inspiratory flow sensor not calibrated	The sensor is not calibrated. The breathing system has been replaced or discon- nected since last calibration.	Perform the leakage test.
			Failure of the flow sensor.	Replace the flow sensor. Perform the leakage test.

Pri	ority	Alarm	Cause	Remedy
!!!	255	Inspiratory halothane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
#!	255	Inspiratory halothane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!	75	Inspiratory halothane low	The inspiratory anesthetic gas concentration is below	Check vaporizer and freshgas settings.
			the lower alarm limit.	Refill the vaporizer.
				Check the breathing system for high leakages.
			The soda lime has dried out.	Replace the soda lime.
!!!	255	Inspiratory isoflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!!	255	Inspiratory isoflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!	75	Inspiratory isoflurane low	The inspiratory anesthetic gas concentration is below the lower alarm limit.	Check vaporizer and freshgas settings.
				Refill the vaporizer.
				Check the breathing system for high leakages.
			The soda lime has dried out.	Replace the soda lime.
!!	10	Inspiratory N2O high	Inspiratory N2O exceeds 82	Check fresh-gas composition.
			%.	Press the O2+ button to flush the breathing system.
!!!	255	Inspiratory sevoflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!!	255	Inspiratory sevoflurane high	Inspiratory anesthetic gas concentration has exceeded the upper alarm limit.	Check vaporizer and freshgas settings.
!	75	Inspiratory sevoflurane low	The inspiratory anesthetic gas concentration is below the lower alarm limit.	Check vaporizer and freshgas settings.
				Refill the vaporizer.
				Check the breathing system for high leakages.
			The soda lime has dried out.	Replace the soda lime.

Pri	ority	Alarm	Cause	Remedy
!!	255	Inspiratory tidal volume	The delivered inspiratory tidal	Check ventilation settings.
		high	volume is higher than the set value.	Check patient compliance. Check if the patient is breathing spontaneously.
				Check the breathing circuit for tight connections and leakages.
!!	255	Inspiratory tidal volume	The delivered inspiratory tidal	Check ventilation settings.
		high	volume exceeds the upper alarm limit.	Check patient compliance. Check if the patient is breathing spontaneously.
				Check the breathing circuit for tight connections and leakages.
				Check alarm limit.
!!	50	Inspiratory tidal volume high	The ventilation settings are not adequate for AutoFlow.	Change the ventilation mode or increase the tidal volume.
!!!	255	Inspiratory xMAC high	The inspiratory anesthetic gas concentration has exceeded 3 xMAC for more than 30 seconds. The expiratory anesthetic gas concentration has exceeded 2.5 xMAC for more than 30 seconds.	Check vaporizer and freshgas settings.
			The inspiratory anesthetic gas concentration has exceeded 5 xMAC.	Check vaporizer and freshgas settings.
!!	255	Inspiratory xMAC high	The inspiratory anesthetic gas concentration has exceeded 3 xMAC for more than 180 seconds.	Check vaporizer and freshgas settings.
!!	100	Internal device temperature high	A ventilation slot at the rear of the device is blocked.	Check the ventilation slots. Ensure air flow at the rear of the device.
			Ambient temperature is too high.	Check ambient conditions.
			A fan is faulty.	Call DrägerService.
			Excessive ventilation settings are applied (e.g., high respiratory rate, high inspiratory pressure, short slopes).	Check ventilation settings.

Pri	iority	Alarm	Cause	Remedy
!	255	Internal device temperature high	A ventilation slot at the rear of the device is blocked.	Check the ventilation slots. Ensure air flow at the rear of the device.
			Ambient temperature is too high.	Check ambient conditions.
			Excessive ventilation settings are applied (e.g., high respiratory rate, high inspiratory pressure, short slopes).	Check ventilation settings.
			A fan is faulty.	Call DrägerService.
!!!	150	Internal FG flow failure	The internal gas delivery system is not operational. A system test may be able to resolve the issue.	Deliver O2: 1. Set the O2 switch to "Add. O2". 2. Set the O2 flowmeter to the desired flow. Check the vaporizer setting. Make sure that fresh gas is reaching the patient.
				Perform the system test after finishing the case.
				If the problem persists, call DrägerService.
!	150	Internal FG flow failure	The internal gas delivery system is not operational. A system test may be able to resolve the issue.	Set the O2 flowmeter to the desired flow. Check the vaporizer setting. Make sure that fresh gas is reaching the patient.
				Perform the system test after finishing the case.
				If the problem persists, call DrägerService.
!	60	License expired	next startup, some functions will no longer be available.	Use "ALARM RESET" to acknowledge alarm.
				To order a permanent license, call Dräger.
!	50	License will expire soon	A trial license will expire within the next 14 days.	Use "ALARM RESET" to acknowledge alarm.
				To order a permanent license, call Dräger.
!!!	0	Loss of data	An internal memory failure has occurred. System data	Check current settings and default settings.
			and system settings are lost.	Call DrägerService.
!	50	Maintenance will be due	Maintenance will be due	Call DrägerService.
		soon	within the next 30 days.	Use "ALARM RESET" to acknowledge alarm.

Pr	iority	Alarm	Cause	Remedy
!		Measured gas concentra-	The measured values are out	Check patient condition.
		tions out of range	of the measurement range.	Check the vaporizer setting and the fresh-gas settings.
				Use an alternative gas measurement system.
				Perform the system test.
!	255	Measured gas concentra- tions temporarily inaccu-	The sensor has not yet warmed up.	Wait for the automatic calibration.
		rate	The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
!	0	MEDIBUS COM 1 failure	Communication via the corre-	Re-establish the connection.
			sponding COM port is inter- rupted.	Use "ALARM RESET" to acknowledge alarm.
			The configured baud rate is not sufficient for the amount of data to be transferred.	Increase the baud rate. Check the configuration of the external device.
!	0	MEDIBUS COM 2 failure	Communication via the corresponding COM port is interrupted.	Re-establish the connection.
				Use "ALARM RESET" to acknowledge alarm.
			The configured baud rate is not sufficient for the amount of data to be transferred.	Increase the baud rate. Check the configuration of the external device.
!!	30	Minute volume high Upper alarm limit for the minute volume has been exceeded.		Check spontaneous breathing.
			Check ventilation settings (e.g., VT, Pinsp, RR).	
				In Pressure Support, correct the trigger threshold if necessary.
				Check alarm limit.
			Flow measurement is inaccurate.	Replace the expiratory flow sensor. Perform the leakage test.

Pr	iority	Alarm	Cause	Remedy
!!	10	Minute volume low	The minute volume is below	Check patient condition.
			the lower alarm limit.	Check tube or mask.
				Check ventilation settings.
				Check alarm limit.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
			Flow measurement is inaccurate.	Replace the expiratory flow sensor. Perform the leakage test.
!	220	N2O cylinder almost empty	The cylinder is almost empty.	Replace the cylinder. Use the central supply.
!	255	N2O cylinder empty	The cylinder is empty or closed.	Replace the cylinder. Use the central supply.
!	190	N2O cylinder sensor?		Check if cylinder pressure sensor is connected.
				Use "ALARM RESET" to acknowledge alarm.
!!	150	N2O FG flow measurem.	The measurement system for	Only use O ₂ as fresh gas.
		failed	the N2O fresh-gas flow has failed.	Use the total flow tube to check the fresh-gas flow. Adjust the fresh-gas flow so that it equals or exceeds the minute volume.
!	255	N2O measurement failed	There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			_	Use an alternative gas measurement system.
				If the problem persists, call DrägerService.
!	255	N2O measurement not available	The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
			There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			Ambient temperature is too high.	Check ambient conditions.

Pri	ority	Alarm	Cause	Remedy
!	255	N2O measurement temporarily inaccurate	The sensor has not yet warmed up.	Wait for the automatic calibration.
			The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
!!!	110	N2O supply low	Central supply pressure and cylinder pressure are low.	Check central N2O supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	110	No Air delivery	Air is not available. Gas mixer is using 100 % O2 instead.	Check central Air supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!	100	No CO2 detected	Ventilation was started, but	Check patient condition.
			no exhaled CO2 was detected for more than 60 seconds.	Check sample line, water trap, and patient-side filter.
				Use "ALARM RESET" to acknowledge alarm.
!!!	200	No fresh-gas flow	No fresh-gas flow is set.	Open the flow control valves.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	110	No N2O delivery	N2O is not available. Gas mixer is using O2 instead.	Check central N2O supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	210	No O2 delivery	O2 is not available. Gas mixer is using Air instead.	Check central O2 supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	120	O No ventilation An internal error occurred. Ventilation has stopped.	Change the ventilation mode to Man/Spon and ventilate the patient manually, or change back to automatic ventilation.	
			Check the upper alarm limit for the airway pressure.	
				If the problem persists, call DrägerService.
!	220	O2 cylinder almost empty	The cylinder is almost empty.	Replace the cylinder. Use the central supply.
!!!	210	O2 cylinder empty	The cylinder is empty or closed.	Replace the cylinder. Use the central supply.
				Use "ALARM RESET" to downgrade alarm priority.

Pr	ority	Alarm	Cause	Remedy
!	255	O2 cylinder empty	The cylinder is empty or closed.	Replace the cylinder. Use the central supply.
!	190	O2 cylinder sensor?	Cylinder pressure sensor is not connected.	Check if cylinder pressure sensor is connected.
				Use "ALARM RESET" to acknowledge alarm.
!!	150	O ₂ FG flow measurement	The measurement system for	Only use O2 as fresh gas.
		failed	the O ₂ fresh-gas flow has failed.	Use the total flow tube to check the fresh-gas flow. Adjust the fresh-gas flow so that it equals or exceeds the minute volume.
!!	255	O2 measurement not available	The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
			There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
			Ambient temperature is too high.	Check ambient conditions.
!	255	O2 measurement temporarily inaccurate	The sensor has not yet warmed up.	Wait for the automatic calibration.
			The ambient air used for calibrating the sensor was impure.	Position the device in an environment with clean ambient air.
				Wait for the automatic calibration.
!!	255	O2 sensor failure The O2 sensor in the patient- gas measurement module is	Use alternative O2 measuring system.	
			faulty.	Use "ALARM RESET" to downgrade alarm priority.
				Perform the system test.
				Call DrägerService.
!!!	210	O2 supply low	Central supply pressure and cylinder pressure are low.	Check central O2 supply or use cylinder.
				Use "ALARM RESET" to downgrade alarm priority.
!	100	On/Standby key stuck	Key is stuck or was pressed for more than 10 seconds.	Ventilation performance is not affected.
				If the problem persists, call DrägerService.
!	75	Patient-gas measurement module is calibrating	The automatic calibration of the patient-gas measurement module is in progress.	Wait for the calibration to complete.

Priority		Alarm	Cause	Remedy
!!!	100	Pause time expired	Ventilation and gas delivery have been paused longer than the set pause time.	Resume ventilation or adjust timer.
!!	100	Pause time expired	Ventilation and gas delivery have been paused longer than the set pause time.	Resume ventilation or adjust timer.
!	100	Pause time expired	Ventilation and gas delivery have been paused longer than the set pause time.	Resume ventilation or adjust timer.
!!	60	PEEP not applied	Failure of the PEEP valve.	Check the breathing circuit for tight connections and leakages.
				Perform the system test.
!!	50	PEEP/CPAP high	The expiratory limb is blocked.	Check expiratory breathing hose and breathing system.
				Perform the leakage test.
			Failure of the anesthetic gas scavenging system.	Check anesthetic gas scavenging system.
!!	100	PEEP/CPAP low	Fresh-gas flow is insufficient, breathing bag is blocked or positioned incorrectly.	Check fresh-gas settings and position of breathing bag.
			Leakage or disconnection.	Check the breathing circuit for tight connections and leakages.
			Failure of the anesthetic gas scavenging system.	Check anesthetic gas scavenging system.
!	170	Power failure	Mains power supply is not available. The device has switched to battery operation.	Restore mains power supply.
!!	100	Power supply failure	Internal fault in the power supply.	Operation of the device can be continued.
				Call DrägerService.
!!	100	Pressure sensor failure	There is condensate in the breathing hoses.	Check the breathing hoses.
			Sensor calibration failed.	Perform the system test.
				If the problem persists, call DrägerService.
!	255	Pressure-relief valve opened	Pressure in the breathing system is too high.	Check APL valve and freshgas settings.

!!! 2	200	Determitionals atrials		
		Rotary knob stuck	Key is stuck or was pressed for more than 10 seconds.	Ventilation performance is not affected.
				Press and turn rotary knob repeatedly. If alarm condition persists, settings cannot be adjusted anymore.
				Disconnect the patient and ventilate with a manual resuscitator.
				Turn off the device with the main switch. Use "Add. O2". Ventilate the patient manually.
! 1	170	Sample line disconnected?	The sample line or the water trap has been disconnected.	Check sample line and water trap.
! 1	170	Sample line occluded	Sample line is occluded.	Check sample line, water trap, and patient-side filter.
! 7	75	Second agent detected	A second anesthetic agent has been detected.	Flush the system if necessary.
				Check fresh-gas settings.
				Wait for transition phase to end.
! 6	60	Service date reached	Maintenance is due.	Use "ALARM RESET" to
				acknowledge alarm.
!!! 2	200	Set O2 switch to "Add. O2"	The internal gas delivery sys-	Call DrägerService. Set the O2 switch to "Add.
::: 2	200	Set 02 switch to Add. 02	tem has failed and the O2 switch is still set to "Aux. O2".	O2". Set the O2 flowmeter to the desired flow.
!! C)	Speaker failure	The loudspeaker is faulty.	Call DrägerService.
				Use "ALARM RESET" to downgrade alarm priority.
!! 1	100	Synchronized ventilation	Inspiratory efforts of the	Check ventilation settings.
		failure	patient cannot be detected due to faulty sensors.	Change to a non-synchronized ventilation mode.
				Change to Man/Spon mode and ventilate manually.
!!! 1	100	System failure	Failure of internal components. System is no longer	Ventilate the patient with a manual resuscitator.
			operational.	Call DrägerService.
! 2	255	Temperature of TurboVent 2 blower module high	Temperature of TurboVent 2 blower module is too high.	Reduce airway pressure or respiratory rate.

Pri	ority	Alarm	Cause	Remedy
!!!	0	Therapy settings not applied	The last changes to the therapy settings were not applied.	Change to Man/Spon and restart the ventilation mode. Use "ALARM RESET" to
				acknowledge alarm.
				If the problem persists, call DrägerService.
!!	100	Third agent detected	A mixture of more than 2 anesthetic agents has been detected.	Flush the system if necessary.
				Check fresh-gas settings.
				Wait for transition phase to end.
			There is electromagnetic interference.	Check for electromagnetic radiation in the vicinity.
!!	30	Tidal volume not achieved	The delivered inspiratory tidal	Check ventilation settings.
			volume is lower than the set value.	Check Pmax setting.
			value.	Check patient compliance. Check if the patient is breathing spontaneously.
!!!	100	TurboVent 2 failure	The TurboVent 2 blower module is not operational or overheated.	Start manual ventilation!
				Replace TurboVent 2 blower module.
!!!	100	TurboVent 2 not connected	The TurboVent 2 blower module is not connected.	module.
				Use "ALARM RESET" to downgrade alarm priority.
!!!	100	Unknown alarm	Error in the internal alarm system.	Check patient condition.
				Call DrägerService.
!	0	USB write error	The USB flash drive is full, faulty, write-protected, or not compatible. The USB flash drive is not correctly connected. The USB flash drive is not correctly formatted.	Check the USB flash drive.
				Use "ALARM RESET" to acknowledge alarm.
!	100	Vapor View failure (left slot)	Failure of Vapor View hardware.	If the problem persists, call DrägerService.
				Use "ALARM RESET" to acknowledge alarm.
!	100	Vapor View failure (middle slot)	Failure of Vapor View hardware.	If the problem persists, call DrägerService.
				Use "ALARM RESET" to acknowledge alarm.
!	100	Vapor View failure (right slot)	Failure of Vapor View hardware.	If the problem persists, call DrägerService.
				Use "ALARM RESET" to acknowledge alarm.

Priority		Alarm	Cause	Remedy
!	100	Vapor View not operational	Failure of Vapor View hard-	Call DrägerService.
			ware.	Use "ALARM RESET" to acknowledge alarm.
!!	100	Vaporizer open	Control dial is not in position "0".	Close vaporizer to prevent accumulation of anesthetic agent in the breathing system.
!!!	120	Ventilator failure	Failure of the pressure sen-	Start manual ventilation!
			sors.	Perform the system test.
				If the problem persists, call DrägerService.
!!	100	Vol. Ctrl. ventilation failure	The inspiratory flow sensor is faulty.	Change to a pressure-controlled ventilation mode.
				Replace the flow sensor. Perform the leakage test.
!	100	Water trap disconnected?	Infinity ID water trap is not correctly connected.	Check water trap.
!	100	Water trap expired	Accessory has been used too long.	Replace the accessory if necessary.
				Use "ALARM RESET" to acknowledge alarm.
!	170	Water trap full	The water trap of the gas measurement is full.	Check water trap.
			Sample line is occluded.	Check sample line, water trap, and patient-side filter.
!!	80	xMAC low	Inspiratory and expiratory gas concentrations are lower than the automatically calculated	Check patient condition.
				Check filling level. Refill if necessary.
limit.	iimit.	Check vaporizer setting.		
				Check the breathing system and the breathing bag for leakages.
				If the current xMAC is acceptable, use "ALARM RESET" to acknowledge the alarm.

Priority		Alarm	Cause	Remedy
!	170 xMAC low Inspiratory and expiratory gas	Check patient condition.		
		concentrations are lower that the automatically calculated limit.	Check filling level. Refill if necessary.	
			iiiiii.	Check vaporizer setting.
				Check the breathing system and the breathing bag for leakages.
			If the current xMAC is acceptable, use "ALARM RESET" to acknowledge the alarm.	
!!	100	Y-piece is open or not con- nected	The Y-piece is open or not connected during flush or Flush & Dry procedure.	Make sure that the Y-piece is sealed.

Service

Safety information	215
Transport for maintenance purposes	216
Additional information concerning the ceiling-mounted version	216
Definition of service terminology	216
Inspection	217
Remote service	
Safety checks	217
Maintenance	219
Repair	220

Safety information

WARNING

Risk due to inappropriately reprocessed products

The product may be contaminated with infectious agents.

Before maintenance and before returning the product for repair, reprocess the product. Carry out the reprocessing as described in the reprocessing instructions supplied with the product.

WARNING

Risk if service is not performed regularly

Wear and material fatigue of the components may lead to device failure and malfunctions.

Perform service at the specified intervals.

WARNING

Risk if service is not performed properly

Personal injury and property damage may occur if service is not performed properly.

Service must be performed by those user groups that are assigned to the particular measure.

WARNING

Risk if maintenance is not performed properly

If the device is connected to the power supply or the gas supply during maintenance, there is a risk of personal injury and property damage.

Before performing maintenance, disconnect all electrical connections from the power supply and all gas connections from the gas supply.

WARNING

Risk when the housing is being opened

Under the housing, there are live electrical components, which may cause an electric shock.

The housing may only be opened by those user groups that are assigned to that particular measure.

WARNING

Risk of patient injury

If maintenance activities are carried out during ventilation, the patient will be put at

Only carry out maintenance activities when there is no patient connected to the device.

Transport for maintenance purposes

Additional information concerning the ceiling-mounted version

WARNING

Risk of injury or damage to the device

If the medical device is not correctly docked or has been undocked from the ceiling supply unit, it may fall from the device mounting interface or tip over.

The medical device may only be docked to the ceiling supply unit or be undocked from it by service personnel.

WARNING

Risk of injury or damage to the device

If the medical device was placed on the service cart, it does not fulfill the tipping stability requirements as per IEC 80601-2-13 and IEC 60601-1, and there is an increased risk of tipping over.

- Only service personnel may move the medical device on the service cart.
- Preparation for transport and actual transport of the medical device on the service cart must only be carried out if there is no patient in the immediate vicinity.
- For docking to the ceiling supply unit, undocking from the ceiling supply unit, and for transporting on the service cart, contact service personnel.

Definition of service terminology

Concept	Definition
Service	All measures (inspection, mainte- nance, repair) intended to maintain or restore the functional integrity of a product
Inspection	Measures intended to determine and assess the current state of a product
Mainte- nance	Regular specified measures intended to maintain the functional integrity of a product
Repair	Measures intended to restore the functional integrity of a product after a failure

A service contract with Dräger is recommended.

Inspection

Measure	Interval	User group
Inspection and safety check	Every 12 months	Service personnel

Remote service

Perseus supports the following Remote Service functions:

- Help ticket
- Remote Device Check
- Proactive call

Contact the responsible DrägerService organization for further information on the Remote Service function.

Safety checks

Safety checks are no substitute for preventive maintenance measures (including preventive replacement of wear parts) as identified by the manufacturer.

Performing the safety checks

- Check whether the corresponding instructions for use are available in the correct version.
- 2 Perform a functional test of the following functions in accordance with the instructions for use:
 - Emergency O2 delivery
- 3 Check that the product is in good condition:
 - All labels are complete and legible
 - No visible damage to the following components:
 - Trolley and castors
 - Housing parts
 - Brakes
 - Auxiliary power sockets
 - Vaporizer mount
 - Water trap mount
 - O-rings of the water trap mount

- Screen
- Gas inlets
- Status display
- Breathing system
- AGS and AGS valves
- Hoses and cables
- Strain reliefs for AGS, compressed gas hoses, cables
- Optional: Holding arms
- Optional: Flexibility trolley, castors, and locking mechanisms
- Optional: Pressure reducers and their sensor lines
- Fuses which are accessible from the outside are in compliance with the specified values
- Check that the country-specific labeling of the gas type matches the screen display.
- Check the end stop of the holding arms for patient monitors.
- 4 Check the electrical safety in accordance with the IEC 62353 standard.
- 5 Check the safety equipment:
 - Functional integrity of optical and acoustic alarm generators
 - O2 switch on gas mixing unit (electronically controlled gas mixer)
 - Internal battery
 - Optional: Check the function of the minimum O2 delivery (mechanically controlled gas mixer).

- **6** Check the accuracy of the gas measurement based on a certified test gas concentration:
 - Anesthetic gas measurement: Isoflurane, 1 Vol% Sevoflurane, 1 Vol% Accuracy ±0.35 Vol%
 - N2O measurement, 70 Vol% Accuracy ±7.6 Vol%
 - CO2 measurement, 5 Vol% Accuracy ±0.83 Vol%
- 7 Check the accuracy of the O₂ measurement:
 - Ambient air 21 Vol%
 Accuracy ±3 Vol%
 - 100 Vol%Accuracy –5 Vol%
- 8 Check the sample gas flow of the patient-gas measurement module:
 - Accuracy 200 ±20 mL/min
- 9 Check the patient gas measurement module for leakage:
 - Leakage at –200 hPa (cmH2O)
 - <20 hPa/min (cmH2O/min)</p>
- **10** Check the non-return valve of the central gas supply for leakage:
 - Leakage ≤20 mL/min
- 11 Check the operational readiness by means of a system test.
- **12** Optional: Check the accuracy of the pressure gauge for the internal breathing system:
 - Accuracy 30 hPa (cmH2O) ±10 hPa (cmH2O)
- 13 Optional: Check the external cylinder pressure reducers
 - Annually: Relief valve
 7.5 to 8 kPa x 100
 - Every 6 years: output pressure 5 to 6 kPa x 100
- **14** Optional: Check the functional integrity of the safety sensor.
- 15 Optional: Check the functional integrity of the service cart.

Maintenance

Component	Interval	Measure	User group
CO2 absorber/soda lime with disposable dust filter	If colored violet or according to the config- ured Infinity ID exchange interval	Replace, see page 76	Users
Water trap	If required, if soiled, or according to the configured Infinity ID exchange interval	Replace, see page 77	Users
Flow sensors	If required, if calibration is no longer possible, or according to the config- ured Infinity ID exchange interval	Replace, see page 187	Users
AGS	If required, if filter is soiled, or flow is no longer achieved	Replace, see page 189	Users
AGS filter	If required, if filter is soiled, or flow is no longer achieved	Replace	Service personnel
O-rings on the water trap mount	Every 2 years	Replace	Service personnel
Sealing, sieve insert, and sintered filter on the cylinder connection of the external cylinder pressure reducer (optional)	Every 3 years	Replace	Specialized service personnel
Air filter mat	Every 3 years	Replace	Service personnel
 Patient gas measurement module 			
 Power supply unit 			
CLIC adapter (optional)	Every 4 years	Replace	Users
Lead-gel battery (2 pieces)	Every 2 years	Replace	Service personnel

Component	Interval	Measure	User group
Pressure reducer (optional)			
Cylinder pressure reducer, type MP:			
complete	Every 6 years	Replace	Service personnel
Cylinder pressure reducer, type MK:			
Sealing	Every 3 years		
 Sieve insert 	Every 3 years		
 Sintered metal filter on the cylinder con- nection 	Every 3 years		
 Diaphragm unit 	Every 6 years	Replace	Specialized service
 Valve cone 	Every 6 years		personnel
 Blow-off valve 	Every 6 years		
 Safety valve spring 	Every 6 years		

It is recommended that only original parts from Dräger are used and that the parts are replaced by Dräger.

Repair

Repairs may be performed only by specialized service personnel.

It is recommended that only original parts from Dräger are used and that the parts are replaced by Dräger.

Disposal

Safety information	222
Disposing of the device	222
Disposing of accessories	222

Safety information

WARNING

Risk due to inappropriately reprocessed products

The product may be contaminated with infectious agents.

Before disposal, reprocess the product. Carry out the reprocessing as described in the reprocessing instructions supplied with the product.

Disposing of the device

The disposal of electrical and electronic devices is subject to special guidelines. This device must be disposed of in accordance with national regulations. In countries of the European Union, Dräger will organize the return of the device. Additional information is available at www.draeger.com (search term: WEEE).

Disposing of accessories

When disposing of the following accessory parts, observe the hospital hygiene regulations and the respective instructions for use:

- Flow sensors
- Breathing hoses
- Filter, HME, HMEF
- Breathing bag
- Masks
- Water trap
- CLIC absorber, Infinity ID CLIC absorber
- Soda lime

Dispose on the following articles according to hospital hygiene regulations:

- Sample line
- Disposable dust filter
- AGS

Technical data

General information	224
Ambient conditions	224
Fresh-gas delivery	225
Fresh-gas delivery (electronically controlled gas mixer)	225
Fresh-gas delivery (mechanically controlled gas mixer)	226
Ventilator	226
Breathing system	229
External fresh-gas outlet	231
Anesthetic gas receiving system	231
Measuring systems and displays	231
Display of calculated values	237
Operating characteristics	238
Device outputs	244
Relevant standards	247
Diagrams Pressure characteristics and flow characteristics of the breathing system without breathing hoses and filters (conforming	249
to ISO 8835-2 and ISO 80601-2-13) Response times in event of concentration changes	
Device combinations	250
EMC declaration	251 251
Emission of high-frequency energy	252

Connections to IT networks	253
LAN interface	253
RS232 interface	254
USB interface	255
Open-source software	256

General information

Units of measurement for pressure

User's operating location

All specified tolerances apply for 20 °C (68 °F), 60 % relative humidity, and 1013 hPa (760 mmHg).

The accuracies indicated below change according to ambient pressure, temperature, and relative humidity. If one of the ambient conditions is changed up to the permissible limit, the accuracy of the corresponding value can change by up to 50 %. If more than one of the ambient conditions are changed, the accuracy may change by up to 100 %.

Example: Accuracy of a measured pressure value: ±4 % under standard conditions. At 10 °C, the accuracy changes to ±6 %; at 10 °C and 20 % relative humidity, to ±8 %.

All patient-related volumes and flow values are standardized to the conditions in the lung. (BTPS)

1 hPa = 1 mbar = 1 cmH2O

100 kPa = 0.1 MPa = 1 bar = 1 kPa x 100

At the front at a distance of 1 m (39 in) and a

height of 1.5 m (59 in)

Ambient conditions

During operation

Temperature 10 to 40 °C

(50 to 104 °F)

Ambient pressure 620 to 1060 hPa (9.0 to 15.3 psi)

20 to 95 %, non-condensing

CO2 concentration 300 to 1000 ppm

Height Up to 4000 m (13123 ft)

During storage and transport

Relative humidity

Temperature

Device without battery -20 °C to 60 °C

(-4 °F to 140 °F)

Battery -15 °C to 40 °C

(5 °F to 104 °F)

Ambient conditions (continued)

For storage longer than 12 months —15 °C to 25 °C

(5 °F to 77 °F)

Maximum storage duration without

recharging

180 days

Ambient pressure 500 to 1060 hPa (7.3 to 15.3 psi)

Relative humidity 10 to 95 %, without condensate formation

CO2 concentration Not relevant

The permissible ambient conditions depend on the accessories used. Follow the corresponding instructions for use

sponding instructions for use.

Fresh-gas delivery

All data are standardized to STPD conditions.

O2 flush 25 to 75 L/min at 2.7 to 6.9 kPa x 100 (39 to

100 psi; 0.27 to 0.69 MPa) supply pressure

O2 flow for Aux. O2 and Add. O2

Range Off; 2 to at least 10 L/min at 2.7 kPa x 100

(39 psi or 0.27 MPa) supply pressure

Accuracy $\pm 10 \%$ of the set value for flows >2.0 L/min

Resolution of displayed value 1.0 L/min (up to 10 L/min)

5.0 L/min (above 10 L/min)

Fresh-gas delivery (electronically controlled gas mixer)

All data are standardized to STAPD conditions.

O₂ concentration FG O₂

Setting range 21 to 100 Vol% (carrier gas: Air)

25 to 100 Vol% (carrier gas: N2O)

Accuracy ±5 % or ±2 Vol% (the larger value applies)

Fresh-gas flow FG flow

Setting range Off; 0.2 to 15 L/min

Accuracy ±10 % or ±50 mL/min (the larger value applies)

Fresh-gas delivery (mechanically controlled gas mixer)

All data are standardized to STAPD conditions.

O2 concentration 21 to 100 Vol%

Setting range for fresh-gas flow 0 to at least 12 L/min (O2, Air, and N2O)

Electronic measurement of fresh-gas flow

Range 0 to 15 L/min (O₂, Air, and N₂O)

Accuracy ±10 % or ±0.12 L/min (the larger value applies)

Resolution of the value displayed on the

screen

0.02 L/min (from 0.2 to 0.5 L/min) 0.05 L/min (from 0.5 to 1 L/min) 0.1 L/min (from 1 to 15 L/min)

0.01 L/min (from 0 to 0.2 L/min)

Resolution of the value displayed in the sta-

tus display

0.1 L/min

Total flow tube All data are standardized to STPD conditions.

Range 0 to 10 L/min

Accuracy ±10 % of the set value at 100 % O2 and for

flows >1 L/min

Resolution of displayed value 0.5 L/min

Ventilator

Electronically driven ventilator, fresh-gas decoupled

Time-based settings:

Respiratory rate RR 3 to 100 /min

±10 % of the set value or ±1/min (the larger

value applies)

Minimum respiratory rate in PSV RRmin

mode

Off, 3 to 25 /min

±10 % of the set value or ±1/min (the larger

value applies)

Inspiratory time Ti 0.2 to 10 s

±10 % of the set value or ±100 ms (the larger

value applies)

Maximum inspiratory time for supported

breaths (fixed setting)

Patient category "adult" 4 s or 1 / (2 x RRmin) (the smaller value applies)

Patient categories "pediatric patient" and

"neonate"

1.5 s or 1 / (2 x RRmin) (the smaller value

applies)

Ventilator (continued)

Pressure rise time (Slope) 0 to 2 seconds;

±20 % of the set value or ±200 ms (the larger

value applies)

Ratio of inspiratory time to expiratory time

1:50 to 50:1

Time at upper pressure level (Thigh) in PC - APRV mode

0.2 to 10 s \pm 10 % of the set value or \pm 100 ms (the larger

value applies)

Time at lower pressure level (Tlow) in PC - APRV mode

0.2 to 10 s

 ± 10 % of the set value or ± 100 ms (the larger value applies)

Inspiration termination criterion (Insp. term.)

5 to 80 %PIF

±20 % of the set value or ±2.5 L/min (the larger value applies)

Percentage ratio of plateau time to inspiratory time of the mandatory breaths in the VC - CMV , VC - SIMV, and

0 to 60 %

VC - SIMV / PS

±20 % of the set value or ±200 ms (the larger value applies)

VC - SIMV / PS (%Tplat)

(%Tplat)
Inflation time (Duration) including rising and

3 to 40 s

falling slopes in the One-step recruitment maneuver

±10 % of the set value

Number of breaths at each level (Breaths/Step), with the exception of the highest level in the Multi-step recruitment maneuver

1 to 20

Number of breaths at the highest level (Breaths@Max) in the Multi-step recruitment maneuver

1 to 20

Reminder for repeated One-step recruitment or Multi-step recruitment maneuver

Off, 10 to 180 min ±1 min

Volume-based and flow-based settings:

Tidal volume VT

20 to 2000 mL, standardized to BTPS condi-

tions

±10 % of the set value or ±15 mL (the larger

value applies)

The applied VT is corrected by the determined breathing hose compliance. The applied VT is also corrected by the sampling flow of the built-in patient-gas measurement module as soon as

CO₂ respiratory phases are detected.

Trigger sensitivity Trigger

0.3 to 15 L/min

±20 % of the set value or ±1 L/min (the larger

value applies)

Ventilator (continued)

Inspiratory flow Flow for adults Minimum 1 L/min, maximum >180 L/min Results from the VT / Pinsp and Ti settings

Inspiratory flow Flow for children and neo-Minimum 1 L/min, maximum 60 L/min nates Results from the VT / Pinsp and Ti settings

Pressure-related settings:

Inspiratory pressure Pinsp PEEP +1 to 80 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Pressure limitation Pmax PEEP +5 to 80 hPa (cmH2O);

±10 % of the set value or ±3 hPa (cmH2O) (the

larger value applies)

Lower pressure level Plow in APRV mode Off, 2 to 35 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Upper pressure level Phigh in APRV mode Plow +1 to 80 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Pressure amplitude above PEEP: ΔPsupp

In the CPAP/PSV and PC - BIPAP / PS

modes

Off, 1 to (80 - PEEP) hPa (cmH2O),

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

In the VC - SIMV / PS and VC - SIMV /

PS / AutoFlow modes

Off, 1 to (Pmax - PEEP) hPa (cmH2O),

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Positive end-expiratory pressure PEEP

Off, 2 to 35 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Continuous positive pressure CPAP

Off. 2 to 35 hPa (cmH2O):

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Inspiratory pressure Pressure in the One-

step recruitment maneuver

3 to 80 hPa (cmH2O) and greater than PEEP:: ±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Highest inspiratory pressure Pinsp max in the Multi-step recruitment maneuver

15 to 80 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Highest PEEP pressure PEEP max in the Multi-step recruitment maneuver

2 to 35 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O) (the

larger value applies)

Pressure amplitude above PEEP: ΔPressure in the Multi-step recruitment maneuver 5 to 30 hPa (cmH2O);

±10 % of the set value or ±2 hPa (cmH2O)

(the larger value applies)

Ventilator (continued)

Minimum pressure limit as per ISO 8835-5

and ISO 80601-2-13

0 hPa (cmH2O)

Breathing system

Total volume including the disposable CO2

absorber CLIC Absorber Free

In the Man / Spon mode Typically 2.18 L In automatic ventilation modes Typically 2.18 L

Compliance (without breathing circuit)

In the Man / Spon mode Typically 2.18 mL/hPa (mL/cmH2O)

corresponds to 65.4 mL at 30 hPa (cmH2O)

In automatic ventilation modes Typically 0.28 mL/hPa (mL/cmH2O)

corresponds to 8.4 mL at 30 hPa (cmH2O)

Soda lime filling volume for absorber container

Reusable CO₂ absorber, filled 1500 ml 1300 ml

Disposable CO₂ absorber CLIC absorber +008

Disposable CO₂ absorber CLIC absorber 1200 ml

Free

Flexible arm for breathing bag

Volume 0.13 L

Compliance 0.13 mL/hPa (mL/cmH2O)

corresponds to 3.9 mL at 30 hPa (cmH2O)

<150 mL/min at 30 hPa (cmH2O) standardized Total system leakage

to BTPS conditions

API valve

Setting range Open, 5 to 70 hPa (cmH2O)

Accuracy (in range from 15 to 30 L/min): ±20 % of the set value or ±3 hPa (the larger

value applies), but not more than +10 hPa

(cmH2O)

Pressure drop at 30 L/min (ATPD), fully

opened

Dry: 2.1 hPa (cmH2O) Wet: 2.2 hPa (cmH2O)

Breathing system (continued)

Breathing system resistance at 60 L/min (reusable CO2 absorber or disposable CO2 absorber, with or without flexible arm for breathing bag, normal operation, filled with Drägersorb 800 Plus)

According to ISO 8835-2, dry, with adult

breathing hose set M30146

According to ISO 8835-2, dry, without hoses

Inspiratory: -3.5 hPa (cmH2O) Expiratory: 4.1 hPa (cmH2O)

Inspiratory: -1.1 hPa (cmH2O) Expiratory: 3.4 hPa (cmH2O)

Recommendation for breathing hoses

All compliances and volumes indicated include inspiratory and expiratory filters

VT <50 mL Maximum compliance: 2.0 mL/hPa (mL/cmH₂O)

VT =50 to 300 mL Maximum compliance: 4.0 mL/hPa (mL/cmH2O)

VT >300 mL Maximum compliance: 6.0 mL/hPa (mL/cmH2O)

Maximum length 200 cm (78.7 in)

350 cm (137.8 in) (with restrictions for compliance correction and pressure metering accu-

racy)

Maximum hose volume 1.6 L

2.7 L (with restrictions for compliance correction

and pressure metering accuracy)

Recommendation for a breathing bag hose if no

flexible arm is in use.

Maximum resistance 1.5 hPa (cmH2O) at 60 L/min

3 hPa (cmH2O) at 60 L/min with increased

fresh-gas consumption

Maximum length 180 cm (70.9 in)

350 cm (137.8 in) with increased fresh-gas con-

sumption

Maximum hose volume 0.7

0.9 L with increased fresh-gas consumption

Recommendation for the bag size

Volume 0.5 L to 5 L (should have at least double the

applied tidal volume)

External fresh-gas outlet

Connection 22 mm outer taper / 15 mm inner taper (ISO)

Delivery See "Fresh-gas delivery"

Pressure limitation Not pressure-limited

Anesthetic gas receiving system

Suction flow

Normal range 32 to 50 L/min
At lower end of restricted range 14 L/min

Maximum fresh-gas flow to prevent contaminat-

ing ambient air

For external breathing systems (normal

range)

For external breathing systems (restricted

range)

For internal breathing systems (restricted

range)

16 I /min

5 L/min

7 I /min

Measuring systems and displays

As Perseus uses more accurate values for internal calculation and alarming, small deviations due to rounding may occur.

Airway pressure

Airway pressure Paw
Plateau pressure Pplat
Positive end-expiratory pressure PEEP
Peak inspiratory pressure PIP
Mean airway pressure Pmean

Range -20 to +120 hPa (cmH2O)

Accuracy ±4 % of the measured value or ±2 hPa (cmH2O)

(the larger value applies)

Resolution of displayed value 1 hPa (cmH2O)

Pressure gauge for indicating the pressure in the internal breathing system

Range -20 to +80 hPa (cmH2O)

Accuracy ±5 % of the measured value or ±2 hPa (cmH2O)

(the larger value applies)

Resolution of displayed value 5 hPa (cmH2O)

Volume

Measured volumes are corrected by the determined breathing hose compliance. As soon as CO2 respiratory phases are detected, measured volumes are also corrected by the suction flow of the integral patient-gas measurement module and according to the determined gas composition.

Tidal volume

Inspiratory VTi Expiratory VT

Range 0 to 2500 mL

Accuracy ±8 % of measured value or ±15 mL (the larger

value applies)

Resolution of displayed value 1 mL

ΔVΤ

Range 0 to 2500 mL

Accuracy ±16 % of measured value or ±30 mL (the larger

value applies)

Resolution of displayed value 1 mL

Minute volume

Total MV

Mandatory MVmand
Spontaneous MVspon
Range 0 to 40 L/min

Accuracy ±8 % of measured value or ±0.2 L/min (the

larger value applies)

Resolution of displayed value 0.1 or 0.01 L/min
To...90 <45 s (RR ≥6 /min)

<105 s (RR <6 /min)

Low-flow wizard

Range 0 to 8 L/min

Accuracy ±25 % of measured value or ±100 mL/min (the

larger value applies)

Respiratory rate

Total RR

Spontaneous RRspon
Mandatory RRmand
Range 0 to 150 /min

Accuracy ±1 /min or 10 % (the larger value applies)

Resolution of displayed value 1 /min

To...90 <45 s (RR ≥6 /min)

<105 s (RR <6 /min)

Dynamic compliance

Compliance Cdyn

Mean compliance Cdyn mean

Range 0 to 200 mL/hPa (mL/cmH₂O)

Accuracy ±30 % or ±3 mL/hPa (mL/cmH2O) (the larger

value applies)

If during automatic ventilation there is still an inspiratory flow at the end of inspiration, the device will determine how much volume could not be applied. This volume will then be

included with the normal volume in the compli-

ance calculation.

With increasing spontaneous breathing activity, the compliance values can be distorted. In this

case the measurement accuracy may be

reduced.

Resolution of displayed value 0.1 mL/hPa (mL/cmH2O)

Resistance R

Range 0 to 100 hPa / (L/s) (cmH2O / (L/s))

Accuracy $\pm 3 \text{ hPa} / (\text{L/s}) \text{ (cmH2O} / (\text{L/s})) \text{ or } \pm 30 \text{ % of the}$

measured value (the larger value applies) (With increasing spontaneous breathing activity, the R values can be severely distorted. Hence maintaining measurement accuracy cannot be guaranteed with spontaneous breathing).

.

Resolution of displayed value 1 hPa / (L/s) (cmH2O / (L/s))

Ε Elastance

0.005 to 10 hPa/mL (cmH2O/mL) Range

Resolution of displayed value 0.001 hPa/mL (cmH2O/mL)

Gas measurement Sidestream gas measurement (The sample gas is fed back into the breathing system and included in calculations for measurements and

delivery);

All values are measured under ATPS condi-

Sample gas flow standardized to STPD condi-

tions.

The measurement is corrected for ambient pres-

Due to the T_{10...90} time and the sampling rate, the accuracies of the measured values for O2, CO₂, N₂O₂, and anesthetic agent may deviate at respiratory rates of 75 /min or higher and an I:E ratio of 1:2. The influence of respiratory rate and the I:E ratio on the accuracy has been verified in a simulated breathing system using a rectangular waveform for the gas concentration.

End-tidal measured values are calculated for each breath from the local maxima and minima of the real-time measurements during expira-

200 mL/min ±10 % Sample gas flow

Maximum time until emptying of the water

trap is necessary

System response time

Sensor sampling rate <50 ms

Time after switch-on until the specified accuracy is attained

Time until CO2 measured values are dis-

played with reduced accuracy

Typical delays

Cross sensitivity

41 h (sample gas under BTPS conditions at 23 °C ambient temperature)

The system response time results from the typical delay and the gas type specific T_{10...90} time.

<500 s

<90 s

<5 s

None with respect to alcohol (<3000 ppm in

Acetone (<1000 ppm), methane, water vapor,

NO and CO

Drift Compensated by automated cyclic zeroing

Ambient air is fed into the breathing system

during automatic zeroing

O2

Range 0 to 100 Vol%

Accuracy $\pm (2.5 \text{ Vol\%} + 2.5 \% \text{ relative})$

Resolution of displayed value 1 Vol%
T10...90 <500 ms

CO₂

Range 0 to 13.6 Vol%

0 to 13.6 kPa 0 to 102 mmHg

Accuracy $\pm (0.43 \text{ Vol}\% + 8 \% \text{ relative})$

±(0.43 kPa + 8 % relative) ±(3.3 mmHg + 8 % relative)

Resolution of displayed value 0.1 Vol% 0.1 kPa

1 mmHg

T10...90 <350 ms

N₂O

Range 0 to 100 Vol%

Accuracy $\pm (2 \text{ Vol}\% + 8 \% \text{ relative})$

Resolution of displayed value 1 Vol%
T10...90 <500 ms

Anesthetic gases

Range

 Halothane
 0 to 8.5 Vol% (kPa)

 Isoflurane
 0 to 8.5 Vol% (kPa)

 Enflurane
 0 to 10 Vol% (kPa)

 Sevoflurane
 0 to 10 Vol% (kPa)

 Desflurane
 0 to 20 Vol% (kPa)

Accuracy $\pm (0.2 \text{ Vol\% } +15 \text{ % relative})$ $\pm (0.2 \text{ kPa} +15 \text{ % relative})$

Resolution of displayed value 0.1 Vol% (kPa) for desflurane

0.01 Vol% (kPa) for all other anesthetic gases

T10...90 <500 ms

Detection Automatic

Primary gas At the latest at 0.3 Vol%

Secondary gas At the latest at 0.4 Vol% or 0.1 xMAC (the larger

value applies)

With a desflurane concentration greater than 4 Vol%, mixture detection occurs at the latest when the concentration of the second anesthetic gas rises above 10 % of the desflurane

concentration.

The secondary gas becomes the primary gas when the expiratory xMAC value is more than 0.2 xMAC above that of the primary gas.

Minimum displayed concentration The specified detection thresholds refer to rising

anesthetic gas concentrations (e.g., at the start of an operation). If the anesthetic gas concentration falls, a concentration of down to 0.05 Vol% will be measured, based on the last anesthetic agent detected. Below this concen-

anesthetic agent detected. Below this concertration, a value of 0 Vol% will be displayed.

xMAC Based on patient age, anesthetic agent, and

nitrous oxide concentration

(the xMAC is corrected for ambient pressure)

Range 0 to 9.9

Accuracy Corresponds to that of the gas measurement

Resolution of displayed value 0.1

Measurement of supply pressures

Central supply

Range 0 to 9.8 kPa x 100

0 to 140 psi 0 to 0.98 MPa

Accuracy ±4 % or ±0.2 kPa x 100 (the larger value

applies)

±4 % or ±3 psi (the larger value applies) ±4 % or ±0.02 MPa (the larger value applies)

Resolution of displayed value 0.1 kPa x 100

1 psi 0.01 MPa

Gas cylinders

Range 0 to 250 kPa x 100

0 to 3600 psi 0 to 25 MPa

Accuracy ±4 % or ±6 kPa x 100 (the larger value applies)

±4 % or ±87 psi (the larger value applies) ±4 % or ±0.6 MPa (the larger value applies)

Resolution of displayed value 1 kPa x 100

1 psi 0.1 MPa

Test results

Total leakage

Range 10 to 5000 mL/min

Resolution of displayed value 1 mL/min

Leakage in automatic ventilation

Range 10 to 5000 mL/min

Resolution of displayed value 1 mL/min
Hose compliance Compliance

Range 0 to 9.9 mL/hPa (mL/cmH2O)
Resolution of displayed value 0.1 mL/hPa (mL/cmH2O)

Display of calculated values

Measurement of consumption and production

CO₂ production of the patient

Range 0 to 9999 mL/min

Accuracy ±25 % or 100 mL/min (the larger value applies)

Resolution of displayed value 1 mL/min

O2 uptake of the patient

Range 0 to 9999 mL/min

Accuracy ±25 % or 100 mL/min (the larger value applies)

Resolution of displayed value 1 mL/min

Fresh-gas consumption

Range 0 to 99999 L, only gas delivered by the gas

mixer is measured

Accuracy $\pm 15 \%$ Resolution of displayed value 1 L

Display of calculated values (continued)

Anesthetic agent

Range 0 to 999.9 mL liquid

Accuracy ±25 % Resolution of displayed value 0.1 mL

Waveforms O2 concentration

Primary anesthetic agent concentration

CO2 concentration Airway pressure Volume (only for loops)

Flow

Sweep speed 6.25; 12.5; 25 mm/s

Scale

Airway pressure Paw (t) —20 to 80 hPa (cmH2O)

Flow (t) —120 to 120 L/min

Volume (t) 0 to 2000 mL O2 (t) 0 to 100 Vol%

CO₂ (t) 0 to 100 mmHg (0 to 12 Vol%, 0 to 12 kPa)

Anesthetic agent (t)

 Halothane
 0 to 5 Vol% (kPa)

 Enflurane
 0 to 6 Vol% (kPa)

 Isoflurane
 0 to 5 Vol% (kPa)

 Sevoflurane
 0 to 10 Vol% (kPa)

 Desflurane
 0 to 20 Vol% (kPa)

Loops Pressure-Volume

Flow-Volume

Operating characteristics

Mains power supply

Electrical mains connection for the trolley

version

Without isolation transformer 100 to 240 $V\sim 50/60$ Hz With isolation transformer 100 to 127 $V\sim 50/60$ Hz

Or

220 to 240 V~ 50/60 Hz

Maximum power consumption 12 A

Electrical mains connection for the ceiling-

mounted version

Ceiling-mounted version (without isola- 100 to 240 V~ 50/60 Hz

tion transformer)

Maximum power consumption 4 A

Power cable

Maximum length 5 m (16.4 ft)

Protective ground resistance Maximum 0.1 ohm

12 A

0.65 A

12 A

Operating voltage ≥250 V
Operating current ≥16 A

Power consumption

At 230 V~

Standby 0.30 A
Typical (without charging the internal 0.35 A

battery)

Maximum (with auxiliary power sock-

ets)

Maximum (ceiling-mounted version) 4 A

For 110 V~

Standby 0.55 A

Typical (without charging the internal

battery)

Maximum (with auxiliary power sock-

ets)

Maximum (ceiling-mounted version) 4 A

Power consumption

Standby 55 W
Typical 70 W
Maximum (trolley version) 2.2 kW
Maximum (ceiling-mounted version) 400 W

Inrush current (trolley version)

With isolation transformer Approx. 26 to 48 A peak Approx. 18 to 34 A quasi-RMS

Without isolation transformer Approx. 8 to 14 A peak

Approx. 6 to 10 A quasi-RMS

Inrush current (ceiling-mounted version)

Without isolation transformer Approx. 8 to 14 A peak

Approx. 6 to 10 A quasi-RMS

Internal battery

Type Lead-gel battery

Sealed, maintenance-free

Capacity 7.2 Ah Voltage 24 V

Fuse F15A 80V UL248-14, breaking capacity 1000 A,

size 19.7 mm * 19 mm * 5 mm

Current Maximum 15 A

Backup time with new and fully charged bat-

tery

Minimum 30 minutes
Typical 150 minutes
Charging time (to reach full power) At least 8 hours
Charging power Max. 50 W

Gas supply

Gas quality

Oil content <0.1 mg/m3

Dew point 5 °C (41 °F) at ambient temperature

Particle size Dust-free air (filtered with pore size <1 μm)

Supply pressure for O₂, Air, N₂O 2.7 to 6.9 kPa x 100

39 to 100 psi 0.27 to 0.69 MPa

Maximum short-term peak inlet flows at 6.9 kPa x 100 (100 psi or 0.69 MPa) supply

pressure

O2 135 L/min (applies only when there is no distri-

bution piece for the central O₂ supply)

Air

Without suction unit 50 L/min Including a directly connected Dräger 130 L/min

ejector suction unit

N2O 40 L/min
Drive gas Not needed

Gas supply connection Depending on configuration: DIN, NIST, DISS,

Air Liquide, SIS

Gas cylinders (dimensions)

Diameter 100 to 140 mm (3.94 to 5.51 in)

102 to 106 mm (4.00 to 4.18 in) for versions with hanger yokes for gas cylinders with pin-index

connection

Maximum height 830 mm (32.68 in)

750 mm (29.52 in) for versions with hanger yokes for gas cylinders with pin-index connec-

tion

Pressure reducers

Version Single-stage pressure reducers

Permissible inlet pressure range (Pv)

Air, O2 11 to 200 kPa x 100 N2O 11 to 60 kPa x 100

Nominal outlet pressure (PA) Fixed setting

See "Flow characteristics"

Due to the design, the outlet pressure rises as

the cylinder pressure falls.

Air, O₂

With inlet pressure Pv = 11 kPa x 100and nominal flow QN = 80 L/min

al flow Qn = 80 L/min

With inlet pressure Pv = 200 kPa x 100 and nominal flow Qn = 80 L/min

4.55 kPa x 100 ±50 kPa

4.75 kPa x 100 ±50 kPa

N₂O

With inlet pressure PV = 11 kPa x 100and nominal flow QN = 65 L/min

iominal flow Qn = 65 L/min

With inlet pressure Pv = 60 kPa x 100 and nominal flow Qn = 65 L/min

4.95 kPa x 100 ±50 kPa

5.2 kPa x 100 ±50 kPa

Nominal flow

Air, O2 QN = 80 L/minN2O QN = 65 L/min

Performance data See "Flow characteristics"

Relief valve

Opening pressure (7.5+0.5) kPa x 100
Minimum flow 220 L/min Air

Input connections Complying with either

EN 850, ISO 407 DIN 477 T1 NF E 29-650 BS 341; Part 1 NBN 226 ISO 5145

Noise emissions from device Free field measurements complying with

ISO 3744

Average sound pressure level Leq(A) during

ventilation with typical settings

≤42 dB(A)

Sound pressure level L(A) of the alarm tones at the workstation, measured according to

IEC 60601-1-8 Edition 2.1: Alarm tone sequence

Alarm tone volumes (all priorities)

Settable from >50 dB(A) to <75 dB(A)

Secondary acoustic alarm signal ≥55 dB(A) and ≤75 dB(A)

Mains power supply failure alarm ≥55 dB(A) and ≤75 dB(A)

Dimensions of trolley version (can deviate with accessory equipment)

 Width
 1150 mm (45.23 in)

 Height
 1480 mm (58.27 in)

 Depth
 790 mm (31.1 in)

Dimensions of ceiling-mounted version without flexibility trolley (can deviate with accessory equipment)

 Width
 1150 mm (45.23 in)

 Height
 1045 mm (41.14 in)

 Depth
 790 mm (31.1 in)

Dimensions of ceiling-mounted version with flexibility trolley (can deviate with accessory equipment)

 Width
 1150 mm (45.23 in)

 Height
 1480 mm (58.27 in)

 Depth
 790 mm (31.1 in)

Work surface dimensions

Width approx. 850 mm (33.46 in)

Depth approx. 350 to 500 mm (13.78 to 19.69 in)

Suitable for paper size DIN A3

Weight of trolley version

cylinder holders)

Nominal configuration consisting of trolley version with isolation transformer, electronic 3-gas gas mixer, screen mounted centrally on the column, plug-in connections for 2 vaporizers, Vapor View option, breathing system, CLIC adapter and CLIC absorber, breathing hoses, central supply hoses (5 m (16.4 ft)), scavenging hose (5 m (16.4 ft))

Various attached parts (e.g., baskets, flexible breathing bag holder, handles, vaporizer parking holders, cylinder pressure reducers,

Hanger yokes with pin-index connection

Plug-in connection for 3 vaporizers with Vapor View option

Mechanically controlled gas mixer

Drawer module (standard) Drawer module (large)

Bronchial suction with swivel arm and

accessories

Permissible total weight

Weight of ceiling-mounted version

Nominal configuration (no flexibility trolley)

hoses.

hose (1.5 m (59 in))

Nominal configuration of the ceilingmounted version on a flexibility trolley

Maximum configuration

Permissible total weight

Approx. 160 kg (353 lbs)

Approx. 10 kg (22 lbs)

+ approx. 3 kg (7 lbs)

+ approx. 3 kg (7 lbs)

+ approx. 3 kg (7 lbs)

+ approx. 12 kg (26 lbs)

+ approx. 16 kg (35 lbs)

+ approx. 6 kg (13 lbs)

330 kg (728 lbs)

approx. 100 kg (220.0 lbs)

Nominal configuration consisting of ceilingmounted version with electronic 3-gas gas mixer, screen mounted centrally on the column, plug-in connections for 2 vaporizers, lockable drawer.

CLIC adapter and CLIC absorber, breathing

central supply hoses (1.5 m (59 in)), scavenging

165 kg (364 lbs)

approx. 200 kg (440.9 lbs)

Refer to the instructions for use of the ceiling supply unit and observe?the maximum load

capacity.

RFID system

Operating frequency 13.56 MHz \pm 50 ppm (broadband)

Transmitter power \leq 42 dB μ A/m (200 mW \pm 1 dB)

Modulation ASK (Amplitude-Shift Keying)

Electromagnetic compatibility Tested in compliance with IEC 60601-1-2

Protection classes

Device I, in compliance with IEC 60601-1

Applied parts (connections for breathing TYPE BF

hoses)

Penetration of liquids IP20 according to IEC 60529, meets IEC 60601-2-13

meets inc occur

Classification in compliance with Directive

93/42/EEC, Annex IX

UMDNS Code Universal Medical Device 10-134 Nomenclature System – nomenclature system

for medical devices

Use of latex Perseus is not made with natural rubber latex.

II b

Device outputs

Serial interfaces COM 1 and COM 2

Only connect devices that meet the requirements of IEC 60950-1/IEC 62368-1 for ungrounded SELV circuits and the requirements of IEC 60601-1 (as of the 2nd edition) for touch-

able secondary circuits with a maximum nomi-

nal voltage of 24 V DC.

Typically <2 s

Protocol MEDIBUS.X

Alarm delay time (measured from time of

query)

Connector 9-way Sub-D

Baud rate 1200, 2400, 4800, 9600, 19200, 38400 baud

Data bits8ParityEvenStop bits1

Device outputs (continued)

Pin assignment

 Pin 1
 n/c

 Pin 2
 RXD

 Pin 3
 TXD

Pins 4, 6 Pins 4 and 6 are connected internally

Pin 5 SHLD-GND

Pins 7, 8 Pin 7 and Pin 8 are connected internally

Pin 9 n/c

Housing SHLD-GND

USB interface Only connect USB storage media that do not

have their own power supply. Do not connect

any charging cables.

Type USB type A connector; USB 1.1

Supported devices USB flash drive formatted with FAT16 or FAT32

Dräger recommends storage media with FIPS 140-2 compatible encryption and their own

encryption facility.

LAN interface Only for Dräger Remote Service

Only connect devices and/or networks that meet the requirements of IEC 60950-1/IEC 62368-1 for ungrounded SELV circuits and the requirements of IEC 60601-1 (as of the 2nd edition) for touchable secondary circuits with a maximum

Only for workstation lights approved by Dräger,

nominal voltage of 24 V DC.

Type RJ45 plug

Transfer speed 10Base-T, IEEE 802.3 Clause 14. Requires at least a CAT3 cable.

Connection for external workstation light

(optional)

Distribution piece for central O2 supply

(optional)

Supply pressure 2.7 to 6.9 kPa x 100

39 to 100 psi 0.27 to 0.69 MPa

see list of accessories.

Maximum permissible extraction flow rate 40 L/min

Device outputs (continued)

Observe the maximum current for each power Mains power sockets

socket, the total current for all power sockets, and the total permissible leakage current. Mains power sockets have no battery back-up and are independent of the switch-on state of

the device.

Power sockets (version without isolation transformer)

IEC Power socket type

Power sockets with protective ground 2 (maximum 3.15 A per power socket) Power sockets without protective ground 2 (maximum 1 A per power socket)

Fuse type Fuses, T3,15AH250V IEC 60127-2/V and

T1AH250V IEC 60127-2/V, size 5 mm * 20 mm

Power sockets

(version with isolation transformer)

Power socket type Country-specific

Fuse type Electronic automatic circuit breakers: 4 x 3 A. slow-blow, for voltages up to 250 V and break-

ing capacity 2000 A.

And 1 x 8 A, slow-blow, for voltages up to 250 V and breaking capacity 2000 A, in each case

conforming to IEC 60934

Power sockets with protective ground

Total current

4 (maximum 3 A per power socket)

Max. 8 A

Main switch for isolation transformer and all

power sockets

Relevant standards

In addition to the standards listed here, this medical device also meets various other standards, e.g., standards concerning special national requirements.

IEC 60601-1 (2nd edition) Part 1:

Medical electrical equipment Requirements for safety

IEC 60601-1-2 Part 1-2:

Medical electrical equipment General requirements for safety; Collateral stan-

dard: Electromagnetic compatibility; Require-

ments and tests

IEC 60601-1-8 Parts 1-8:

Medical electrical equipment General requirements for safety: Collateral stan-

dard: General requirements, tests and guidance for alarm systems in medical electrical systems

IEC 60601-2-13 Parts 2-13:

Medical electrical equipment Particular requirements for the safety of

anaesthetic systems

ISO 8835-2 Part 2:

Systems for inhalational anesthesia Anesthesia breathing systems

ISO 8835-3 Part 3:

Systems for inhalational anesthesia Transfer and receiving systems of active

anaesthetic gas scavenging systems

Part 4: ISO 8835-4

Anesthetic vaporizers

ISO 8835-5 Part 5:

Anesthesia ventilators Systems for inhalational anesthesia

ISO 21647 Particular requirements for the basic safety and Medical electrical equipment

essential performance of respiratory gas moni-

tors

Relevant standards (continued)

For devices manufactured since July 2014, the following also apply:

IEC 60601-1 (3rd edition) Medical electrical equipment

ISO 80601-2-13 Medical electrical equipment

ISO 80601-2-55:2011 (devices with production date July 2014 - April 2017)
ISO 80601-2-55:2018 (devices with produc-

ISO 80601-2-55:2018 (devices with production date from May 2017 onwards)

IEC 60601-1-2 Medical electrical equipment

IEC 60601-1-8

Part 1:

General requirements for basic safety and essential performance

Parts 2-13:

Particular requirements for basic safety and essential performance of an anaesthetic work-station

Parts 2-55:

Particular requirements for basic safety and essential performance of respiratory gas monitors

Part 1-2:


General requirements for safety - Collateral standard: Electromagnetic compatibility - Requirements and tests

Parts 1-8:

General requirements for basic safety and essential performance - Collateral standard Alarm systems - General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems

Diagrams

Pressure characteristics and flow characteristics of the breathing system without breathing hoses and filters (conforming to ISO 8835-2 and ISO 80601-2-13)

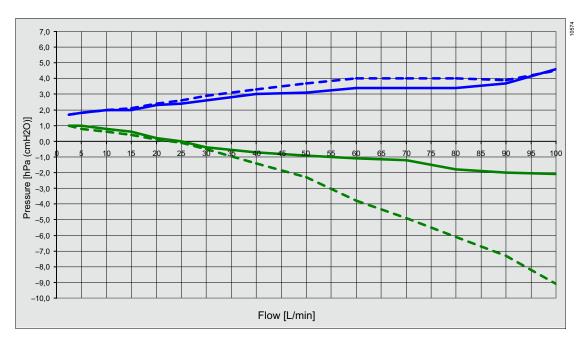

Pexsp [hPa (cmH2O)], ventilation drive inactive

Pexsp [hPa (cmH2O)] Man / Spon

Pinsp [hPa (cmH2O)], ventilation drive inactive

Pinsp [hPa (cmH2O)] Man / Spon

Breathing system, dry, with filled reusable		Resistance [hPa (cmH2O)]		
CO2 absorber and soda lime dust filter	Peak flow in use [L/min]	Man / Spon	Man / Spon	
MX50115		Inspiratory	Expiratory	
No breathing circuit or inspiratory filter	60	-1.1	3.4	
	30	-0.4	2.6	
	15	0.6	2.0	
	2.5	1.0	1.7	
With breathing circuit for adults M30146, inspiratory filter MP01730	60	-3.5	4.1	
	30	-1.1	2.9	

I(:()) ahearhar and eada lime dust tiltar	Peak flow in use [L/min]	Resistance [hPa (cmH2O)]		
		Man / Spon		
		Inspiratory	Expiratory	
With breathing circuit for pediatric patients M27542, inspiratory filter MP01735	15	-0.7	3.0	
With breathing circuit for neonates MP00333, inspiratory filter MP01735	2.5	1.1	1.7	

Response times in event of concentration changes

Typical response times (To..90) for an oxygen concentration change from 21 Vol% to 100 Vol% at the following fresh-gas flows:

	2 L/min	4 L/min	8 L/min	O2 flush
Test lung for adults (MP02400), breathing circuit (MP00300), breathing bag 2 L (MP00222) VT=500 mL, RR=10 /min, I:E=1:2	712 s	174 s	32 s	9 s
Test lung for neonates (8410079), breathing circuit (MP00333), breathing bag 1 L (MP00383) VT=30 mL, RR=30 /min, I:E=1:1	91 s	64 s	46 s	7 s

Device combinations

This device can be operated in combination with other Dräger devices or with devices from other manufacturers. Follow the accompanying documents of the individual devices.

If a device combination is not approved by Dräger, the safety and the correct functioning of the individual devices can be compromised. The operating organization must ensure that the device combination complies with the applicable editions of the relevant standards for medical devices.

Device combinations approved by Dräger meet the requirements of the following standards:

- IEC 60601-1 (electrical safety, mechanical safety, software)
- IEC 60601-1-2 (EMC)
- IEC 60601-1-8 (alarm systems)

EMC declaration

General information

This device was tested for electromagnetic compatibility using accessories from the list of accessories. Other accessories may only be used if they do not compromise the electromagnetic compatibility. The use of non-compliant accessories may result in increased electromagnetic emissions or decreased electromagnetic immunity of the device.

This device may be used in the direct vicinity of other devices only if Dräger has approved this device arrangement. If no approval has been given by Dräger, it must be ensured that this device functions correctly in the desired arrangement before use. The instructions for use for the other devices must be followed.

Electromagnetic environment

This device may be used only in environments specified in section "Environments of use" on page 16.

Emissions	Compliance
Radiated emissions	Class A, group 1 (30 MHz to 1 GHz)
Conducted emissions	Class A, Group 1 (150 kHz to 30 MHz)

NOTE

The emissions characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 class A). If it is used in a residential environment (for which CISPR 11 class B is normally required), this equipment might not offer adequate protection to radio-frequency communication services. The user might need to take mitigation measures, such as relocating or re-orienting the equipment.

Immunity against	Test level and required electromagnetic environment	
Electrostatic discharge (ESD) (IEC 61000-4-2)	Contact discharge: ±8 kV	
	Air discharge: ±15 kV	
Fast transient electrical disturbances (bursts)	Power cable: ±2 kV	
(IEC 61000-4-4)	Longer signal input lines/output lines: ±1 kV	
Impulse voltages (surges) (IEC 61000-4-5)	Voltage, external conductor – external conductor: ±1 kV	
	Voltage, external conductor – protective ground conductor: ±2 kV	
Magnetic fields at mains frequency (IEC 61000-4-8)	50 Hz: 30 A/m	
Voltage dips and short interruptions in the supply voltage (IEC 61000-4-11)	Voltage dips of 30 % to 100 %, 8.3 ms to 5 s, different phase angles	
Radiated high-frequency disturbances (IEC 61000-4-3)	80 MHz to 2.7 GHz: 3 V/m	

Immunity against	Test level and required electromagnetic environment
Conducted high-frequency disturbances (IEC 61000-4-6)	150 kHz to 80 MHz: 3 V, ISM bands: 6 V
,	Various frequencies from 385 MHz to 5785 MHz: 9 V/m to 28 V/m

Recommended separation distances from wireless communication devices

To ensure that the full functional integrity of this device is not compromised, there must be a separation distance of at least 1.0 m (3.3 ft) between this device and wireless high-frequency communication equipment.

Emission of high-frequency energy

This medical device is equipped with an RFID module (radio frequency identification) to enable wireless communication with Infinity ID accessories.

This medical device has been designed and manufactured to comply with emission limit values for high-frequency energy. These limit values are incorporated in international safety standards such as IEC 60601-1-2 and standards for radio equipment such as EN 300330 and have been defined by regulatory authorities.

The RFID system of this medical device complies with Part 15 of the FCC regulations and the license-free RSS regulations of Industry Canada. Operation is subject to the following 2 conditions:

- This medical device does not cause any harmful interference.
- 2 The medical device is not liable to damage caused by the reception of interference, including interference causing undesired operating conditions.

Changes and modifications that have not been expressly approved by Dräger may result in the user no longer being permitted to operate the device.

Dräger hereby declares that this medical device, including its radio equipment, is in compliance with Directive 2014/53/EU.

The complete EU Declaration of Conformity can be viewed at the following internet address: http://www.draeger.com/doc-radio

Connections to IT networks

In an IT network, data can be exchanged by means of wired or wireless technologies. An IT network can be any data interface (e.g., RS232, LAN, USB, printer interface) that is described in standards and conventions.

During operation, this device can exchange information with other devices by means of IT networks. The following functions are supported:

- Display of waveforms and parameter data
- Transmission of existing alarm conditions
- Transfer of device settings and patient data
- Service mode, access to logbooks

Connecting this device to a network that incorporates other devices or making subsequent changes to that network can lead to new risks for patients, users, and third parties. Before the device is connected to the network or the network is changed, these risks must be identified, analyzed, and evaluated by the hospital IT representative in accordance with the IEC 80001-1 standard (risk management for medical IT networks). Appropriate measures must be taken on the basis of the results.

Examples of subsequent changes to the network:

- Changing the network configuration
- Removing devices from the network
- Adding new devices to the network
- Performing upgrades or updates on devices that are connected to the network

LAN interface

Service

In conjunction with the Dräger SCG (ServiceConnect Gateway) or a DrägerService computer, the LAN interface enables the following functions:

- Using the SNMP protocol:
 Monitoring the service status of the device, querying the service status, support during the installation of device software and during software download, configuration support
- Using the FTP protocol (as a client):
 Querying the device status, support during the installation of device software and during software download, configuration support

The following personal data are transmitted unencrypted over the interface:

Logbook with details of age, weight, and height of the patient

Time synchronization

The LAN interface allows synchronization with an NTP (Network Time Protocol) server using the NTP protocol.

Required characteristics

The LAN must ensure the connection between the device and the following destinations:

- Service Connect Gateway or DrägerService Computer
- NTP server

Connections bet	Connections between Perseus and the DrägerService computer or the NTP server					
Function	Protocol	Perseus port	Direction	Remote port	Remote partner	Connection
SNMP V3	UDP	161	\longleftrightarrow	>1023	SCG	
TCP SNMP V3 (trap)	UDP	>1023	→	162	SCG	
FTP (command)	TCP	>1023	→	21	SCG	New, estab- lished
FTP (command)	TCP	>1023	+	21	SCG	Established
FTP (data)	TCP	>1023	→	>1023	SCG	New, estab- lished
FTP (data)	TCP	>1023	—	>1023	SCG	Established
SNTP	UDP	>1023	\longleftrightarrow	123	NTP server	
DHCP	UDP	67	\leftrightarrow	67	DHCP server	

Typical data volume:

- Update of the device firmware: typically 50 MB
- Help ticket (system log for service purposes): typically 3 MB

While using the service functions, the device typically causes an average network load of up to 150 KB/s. During normal use, the bandwidth used is negligibly low.

Hazardous situations

NOTE

Risk when using the device in an overloaded network

An overload of the device due to high network loading (e.g., caused by denial-of-service attacks) may lead to a shut-down of the device-side network interface. The interface will not be available again until the device is restarted. This function is used to protect the primary functions of the device against unauthorized intrusion.

RS232 interface

The RS232 interface supports the MEDIBUS.X protocol. MEDIBUS.X is a communication protocol for data exchange between Perseus and, e.g., the following external medical or non-medical equipment:

- Hemodynamic monitor
- Data management system
- Computer
- RS232 to Ethernet converter

The transferred data include the following information:

- Settings
- Measured values
- Waveforms
- Text messages
- Alarm status

Take note of the documentation for the following communication protocols before transferring the data:

- MEDIBUS.X, Rules and Standards for Implementation (9052607)
- MEDIBUS.X, Profile Definition for Data Communication V1.n (9052608)

The documents are only available in English.

The following personal data are transmitted unencrypted over the interface:

 Therapy data with details of age, weight, and height of the patient

WARNING

Risk due to incompletely transferred data

Data (e.g., measured values, alarms) which the anesthesia machine transfers to other systems such as patient monitors or EMR systems may be displayed there incompletely or incorrectly.

- Only use these data for information purposes.
- Do not use these data for patient monitoring or device monitoring.
- Do not make diagnostic or therapeutic decisions on the basis of these data.

With the aim of improving the clinical process, the data can be used to construct a distributed alarm system with unacknowledged alarm transmission conforming to IEC 60601-1-8. However, the data must not be used as a substitute for Perseus as the primary alarm source.

Alarm conditions are only transmitted in the network if they are present at the time of an active query by a connected device. An alarm history is not transmitted.

Note: Each connected device within the distributed alarm system must be marked with the following warning: No guaranteed reception of alarm signals.

Required characteristics

The RS232 interface is a point-to-point connection. A connected device must prevent access by unauthorized users to the data that are sent over the RS232 interface and must itself be protected from infections by malware and computer viruses.

USB interface

The USB interface supports the transfer of data to external storage media. Existing data on external storage media can be deleted during this operation.

The following personal data are transmitted unencrypted over the interface:

- User logbook with details of age, weight, and height of the patient
- Screenshots potentially showing information about the age, weight, and height of the patient

Required characteristics

A connected device must conform to the mass storage medium USB device class (e.g., connecting devices to charge the battery is not intended.). Dräger recommends the use of FIPS 140-2 compatible storage media with hardware encryption.

Hazardous situations

Connecting active devices to the Perseus USB interface can cause Perseus to restart.

Open-source software

Dräger devices that use software may use opensource software, depending on their setup. Opensource software may be subject to different terms of license. Additional information regarding the open-source software used in this device is available at the following web page:

www.draeger.com/opensource

Principles of operation

Description of the ventilation modes	258
Meaning and function of the therapy controls Degree of respiratory support	
Ventilation modes and effective parameters Pressure-supported ventilation	
Pressure-controlled ventilation	263
Volume-controlled ventilation	
Minimum O2 delivery	268
Influence of patient category, weight, and	
age on device behavior	269
Influence of patient categoryInfluence of the ideal patient weight and patient	269
height	269
Influence of the patient age	269
Infinity ID accessory support	269
Infinity ID functionality	269
Schematic diagram of alarm tones	271
Tone sequence for various alarm priorities Tone signals during operation	
Combination lock on the drawer (optional).	272
Opening the combination lock	272
Changing the combination	272
Loss of the combination	273

Description of the ventilation modes

Meaning and function of the therapy controls

Therapy controls	Meaning / Function		
%Tplat	Plateau time as a percentage of the inspiratory time Ti in VC - CMV mode		
RR	Respiratory rate		
RRmin	Minimum respiratory rate at which suppose Support mode.	orted breaths are applied in Pressure	
Insp. term.	When the flow falls below this flow value supported breath is interrupted.	e (in % of the measured peak flow), a	
PEEP/CPAP	Positive end-expiratory pressure / Conti Pressure that is always maintained.	nuous Positive Airway Pressure	
ΔΡѕирр	Pressure difference of a pressure-supported breath between PEEP level and inspiratory pressure. This pressure support is only available if synchronization of spontaneous breathing (Sync. on) is switched on. When pressure support is switched on, the naming of the following ventilation modes changes:		
	Without pressure support	With pressure support	
	PC - BIPAP	PC - BIPAP / PS	
	VC - SIMV VC - SIMV / PS		
	VC - SIMV / AutoFlow	VC - SIMV / PS / AutoFlow	
Phigh	Upper pressure level in APRV	•	
Pinsp	Inspiratory pressure		
Plow	Lower pressure level in APRV		
Pmax	Upper pressure limit in volume-controlled ventilation. When this pressure is reached, the breath is held at this level until the set inspiratory time Ti is reached.		
Trigger	Flow that, when exceeded, triggers a su	pported breath	
Thigh	Period of time in APRV during which the upper pressure level is maintained.		
Ti	Inspiratory time		
Tlow	Period of time in APRV during which the lower pressure level is maintained.		
Slope	Period of time during which a pressure rise from the PEEP or CPAP pressure to the inspiratory pressure or PSV pressure takes place. This time determines the steepness of the rise in pressure from the lower to the upper level.		
VT	Tidal volume		

Therapy controls	Meaning / Function		
Sync. on/Sync. off	Switching spontaneous respiratory support on / off Switching the synchronization on or off causes the following change to the ventilation mode:		
		Pressure-controlled	Volume-controlled
	Sync. on	PC - BIPAP	VC - SIMV
	Sync. off PC - CMV VC - CMV		
	the patient's by adapting the expiratory phatory breath candult) or 1.5 sinspiratory floger window, and if no spontant	preathing effort. In doing so, the mandatory breaths and the case, an inspiratory trigger windown be initiated prematurely by useconds (patient categories Pew reaches the set value of the a premature mandatory breath is	in the inspiratory trigger window, a

Degree of respiratory support

Respiratory support	Ventilation mode
None Standby, Pause, Ext. FGO	
Low	Man/Spon, CPAP, CPAP / PSV with ΔPsupp <5 hPa (cmH2O)
Medium	CPAP / PSV with ΔPsupp ≥5 hPa (cmH2O)
High	Volume-controlled modes Pressure-controlled modes

Ventilation modes and effective parameters


Group	Tab	Ventilation mode	Base parameters (normal therapy bar)	Additional parameters (expanded therapy bar)
Manual ventila- tion / Sponta- neous breathing	Man/Spon	Manual / Spontaneous	CPAP ¹⁾	
Pressure-sup- ported ventila- tion	PSV ¹⁾	CPAP / PSV	Trigger ΔPsupp RRmin PEEP Slope	Insp. term.
Pressure-controlled ventilation	PC	PC - CMV	Pinsp RR PEEP Ti Sync. off	Slope
		PC - BIPAP	Pinsp ΔPsupp ¹⁾ = Off RR PEEP Ti Sync. on	Trigger Slope
		PC - BIPAP / PS ¹⁾	Pinsp ΔPsupp ¹⁾ >0 RR PEEP Ti Sync. on	Trigger Insp. term. ¹⁾ Slope
	APRV ²⁾	PC - APRV	Phigh Thigh Slope Plow Tlow	

Group	Tab	Ventilation mode	Base parameters (normal therapy bar)	Additional parameters (expanded therapy bar)
Volume-con- trolled ventila- tion	VC - AF	VC - CMV / AutoFlow	Pmax VT RR PEEP Ti Sync. off	Slope
		VC - SIMV / AutoFlow	Pmax VT RR PEEP Ti Sync. on	Trigger ΔPsupp ¹⁾ = Off Slope
VC		VC - SIMV / PS / Auto- Flow ¹⁾	Pmax VT RR PEEP Ti Sync. on	Trigger ΔPsupp ¹⁾ >0 Insp. term. ¹⁾ Slope
	VC	VC - CMV	Pmax VT RR PEEP Ti Sync. off	%Tplat
		VC - SIMV	Pmax VT RR PEEP Ti Sync. on	Trigger ΔPsupp ¹⁾ = Off %Tplat
		VC - SIMV / PS ¹⁾	Pmax VT RR PEEP Ti Sync. on	Trigger ΔPsupp ¹⁾ >0 Insp. term. ¹⁾ Slope %Tplat

Requires software option PSV
 Requires software option APRV

Pressure-supported ventilation

PSV (optional)

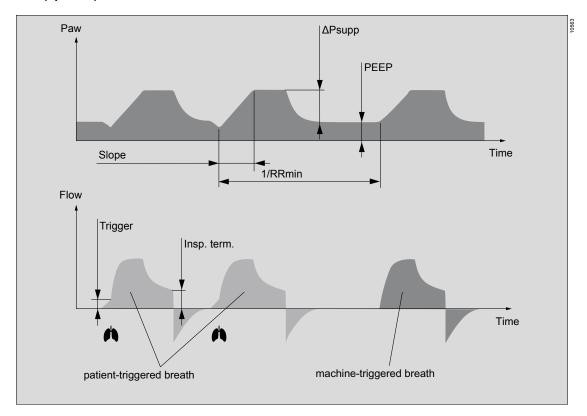

CPAP / PSV (optional)

- Spontaneous breathing
- Spontaneous breathing with continuous positive pressure level with or without pressure support


Each detected inspiratory effort at CPAP level induces a patient-triggered, flow-controlled, and pressure-supported breath. Point in time, number, and duration of pressure-supported breaths are determined by the patient. When no inspiratory effort is detected, mechanical breaths are delivered at the set minimum respiratory rate RRmin and pressure support ΔPsupp.

Patient-triggered breaths are ended as soon as the inspiratory flow falls below the flow defined by the Insp. term. setting. The duration of a machine-triggered breath is additionally determined by the patient category and the set minimum respiratory rate RRmin.

Purely spontaneous breathing at CPAP level is achieved by setting ΔPsupp=Off. In this case the patient does not receive any more machine-triggered breaths.

Pressure-controlled ventilation

PC

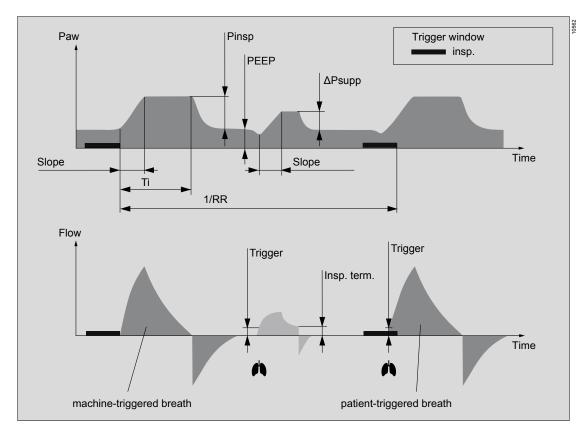

PC - CMV

- Pressure-controlled
- Time-controlled
- Machine-triggered

The mandatory breaths are machine-triggered and are not triggered by the patient.

PC - BIPAP


- Pressure-controlled
- Time-controlled
- Machine-triggered or patient-triggered
- Synchronized inspiration

In PC - BIPAP, the patient can breathe spontaneously at any time, while the number of mandatory breaths is predefined. When synchronization is switched on, the breaths are adapted to the spontaneous breathing efforts of the patient. If spontaneous breathing effort by the patient is detected during the inspiratory trigger window, a patient-triggered breath will be initiated.

PC - BIPAP / PS (optional)

This mode is similar to PC - BIPAP, except that the patient's spontaneous breathing at the PEEP level during the expiratory phase is pressure-supported with Δ Psupp when outside the trigger window.

APRV (optional)

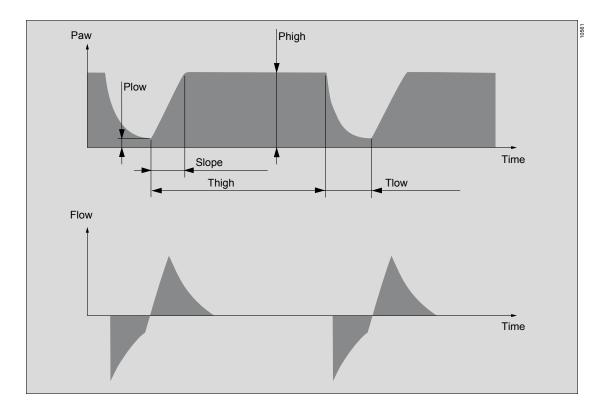

PC - APRV

- Pressure-controlled
- Time-controlled
- Machine-triggered
- Spontaneous breathing at continuous positive airway pressure with short pressure releases

In PC - APRV, the patient's spontaneous breathing occurs at the upper pressure level Phigh. This pressure level is maintained for the duration of Thigh.

The number of pressure releases is determined by the Thigh and Tlow settings. The releases are time-controlled and are not triggered by the patient. The duration is determined by Tlow.


Volume-controlled ventilation

Compliance correction

The applied VT is corrected by the determined breathing hose compliance, i.e., an additional volume is delivered in order to ensure the

application of the volume to the patient. The applied VT is also corrected by the sampling flow of the built-in patient-gas measurement module as soon as CO₂ respiratory phases are detected.

VC

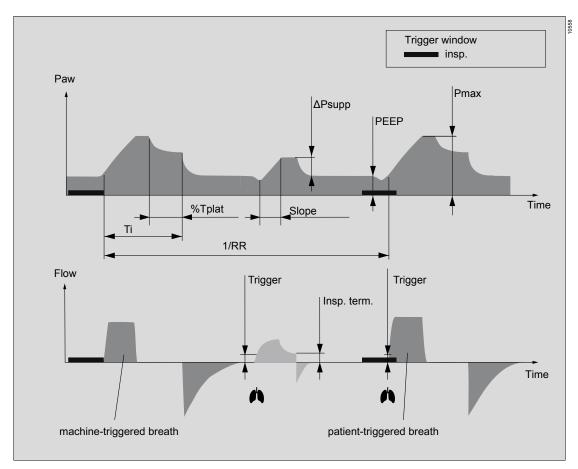

VC - CMV

- Volume-controlled
- Pressure-limited
- Time-controlled

- Machine-triggered
- Constant inspiratory flow

In this volume-controlled ventilation mode, the patient receives the set tidal volume VT with each mandatory breath.

VC - SIMV


- Volume-controlled
- Pressure-limited
- Time-controlled
- Machine-triggered or patient-triggered
- Constant inspiratory flow
- Synchronized inspiration

In VC - SIMV, the patient can breathe spontaneously at any time, while the number of mandatory breaths is predefined. When synchronization is switched on, the breaths are adapted to the spontaneous breathing of the patient. If inspiratory effort by the patient is detected during the inspiratory trigger window, a patient-triggered breath will be initiated.

VC - SIMV / PS

This mode is similar to VC - SIMV, except that the patient's spontaneous breathing at the PEEP level during the expiratory phase is pressure-supported with ΔPsupp when outside the trigger window.

VC - AF

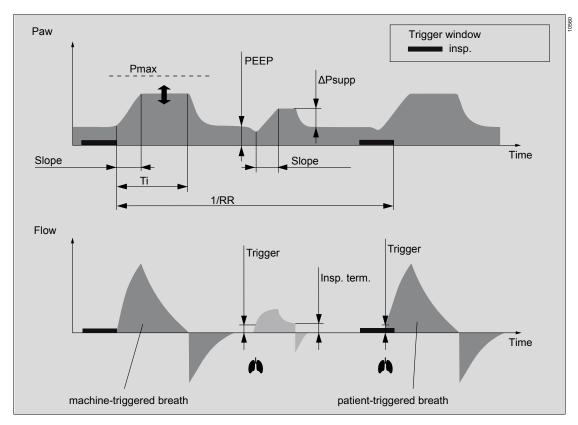

With AutoFlow, the set tidal volume VT is applied for all mandatory volume-controlled breaths with the lowest required pressure. The patient can breathe spontaneously throughout the entire respiratory cycle, during both inspiration and expiration. The pressure patterns and flow patterns of the mechanical inspiratory breaths correspond to those of pressure-controlled ventilation.

Due to the patient's spontaneous breathing efforts or compliance changes in the lungs, the tidal volume in an individual breath may deviate from the set tidal volume VT. However, averaged over time a tidal volume corresponding to the set volume VT is applied.

If no automatic ventilation has previously taken place, a volume-controlled test breath with constant inspiratory flow is performed when starting a ventilation mode with AutoFlow in order to estimate the lung parameters. The inspiratory pressure required at the start is determined from this test breath. Each additional breath-related readjustment of the inspiratory pressure is limited to ±3 hPa (cmH2O). The pressure difference (inspiratory pressure - PEEP) is at least 5 hPa (cmH2O) and the upper inspiratory pressure limit is set by Pmax. If the set value for VT is reduced, the inspiratory pressure will be reduced by a greater amount if necessary.

VC - CMV / AutoFlow

- Volume-controlled
- Pressure-limited
- Time-controlled
- Machine-triggered
- Decelerating inspiratory flow

The mandatory breaths are machine-triggered and are not triggered by the patient.

VC - SIMV / AutoFlow

- Volume-controlled
- Pressure-limited
- Time-controlled
- Machine-triggered or patient-triggered

- Decelerating inspiratory flow
- Synchronized inspiration

In VC - SIMV / AutoFlow, the patient can breathe spontaneously at any time, while the number of mandatory breaths is predefined. When synchronization is switched on, the breaths are adapted to the spontaneous breathing efforts of the patient. If spontaneous breathing effort by the patient is detected during the inspiratory trigger window, a patient-triggered breath will be initiated.

VC - SIMV / PS / AutoFlow

This mode is similar to VC - SIMV / AutoFlow, except that the patient's spontaneous breathing at the PEEP level during the expiratory phase is pressure-supported with Δ Psupp when outside the trigger window.

Minimum O₂ delivery

Gas mixer	Minimum FG-O2 concentration		Minimum O2 flow	
Gas Illixei	Carrier gas Air	Carrier gas N2O	William O2 now	
Electronically controlled	21 %	25 %	Configurable for each patient category, see chapter "Vertical tab "Fresh gas" (only with electronically controlled gas mixer)" When the minimum O2 delivery switches on, the FG O2 therapy control is selected in addition to the selected therapy control. When the active set value is changed, FG O2 changes automatically with it.	
			Continuously adjustable with flow control valves.	
			The minimum O2 delivery interrupts the N2O flow in the following cases:	
Mechanically controlled	21 %	21 %	 N2O flow control valve open and O2 flow control valve closed 	
			 O2 flow less than 200 mL/min O2 will continue to be delivered in the event of an N2O failure. 	

Influence of patient category, weight, and age on device behavior

Influence of patient category

- Alarm limits and start settings for therapy
- Volumeter scale
- Flow measurement and software algorithms to suppress artifacts
- Maximum duration of a pressure-supported breath

Influence of the ideal patient weight and patient height

Ideal patient weight describes the portion of the body that is relevant to setting the ventilation parameters (patient body weight minus assumed excess fat).

In the **Adult** and **Ped** patient categories, the ideal patient weight is calculated from the entered patient height.

In the **Neo** patient category, the ideal patient weight is equal to the entered patient weight.

The calculated ideal patient weight affects:

- Start settings for tidal volume VT
- Start settings for respiratory rate RR
- Start settings for VT and MV alarm limits
- Flow trigger

VT and RR are only dependent on the ideal patient weight when the Ideal body weight function is selected in System setup > Therapy, see page 157.

Changing the patient weight during an automatic ventilation has no effect on the current ventilation settings.

Influence of the patient age

During operation, the set age influences:

Calculation of MAC value

Infinity ID accessory support

Perseus can be operated with accessories with Infinity ID functionality.

- Infinity ID breathing circuit
- Infinity ID WaterLock 2 water trap
- Infinity ID flow sensors
- Infinity ID CLIC absorber

The Infinity ID functionality can be configured, see page 169.

Infinity ID functionality

WARNING

If no Infinity ID accessories are used, the additional functions such as exchange interval monitoring and anti-interchange security are not available.

CAUTION

Risk of incorrect values for resistance and compliance

If an unused Infinity ID accessory is in the immediate vicinity of the device, values may be transferred inadvertently from this accessory.

Do not keep unused Infinity ID accessories in the vicinity of the device.

CAUTION

Risk of inappropriate compliance values

When accessories are added to an Infinity ID breathing circuit, the values for compliance and leakage may deviate from those saved on the breathing circuit.

Always perform the leakage test before starting therapy in order to determine the actual values for compliance and resistance. If the test cannot be performed because the patient is already connected, particular attention is required during ventilation.

CAUTION

Malfunctions of the Infinity ID functionality

Particular EMC situations or defects in Infinity ID components can cause permanent alarms.

To prevent distracting the user under these circumstances, contact service personnel to deactivate the Infinity ID alarms.

Exchange interval monitoring

An automatic monitoring of the period of life is available for Infinity ID WaterLock 2 water traps, Infinity ID CLIC absorbers, Infinity ID breathing circuits, and Infinity ID flow sensors.

An exceeded period of use is signaled during the system test.

The exchange interval for the connected Infinity ID accessories can be adjusted.

This interval must be specified in accordance with the applicable hygiene regulations or the requirements stated in the instructions for use of the corresponding accessory.

WARNING

Risk of inappropriate operating life

Exchange monitoring only considers the absolute period of use and not the current condition of the Infinity ID accessory and therefore does not excuse the user from periodic checks of the accessory.

The exchange interval which can be set for the exchange monitoring does not represent a guarantee for the maximum period of use of the accessory.

Anti-interchange security

WARNING

Risk of interchanged or incorrect breathing hoses

The Infinity ID function to prevent interchanging the breathing hoses does not excuse the user from checking the accessories.

The Infinity ID function to prevent interchanging does not represent a guarantee that the hoses are correctly connected.

When Infinity ID breathing hoses and Infinity ID breathing bags are used, the incorrect connection of breathing hoses and the breathing bag is detected and reported. Hoses that are incorrectly connected with the breathing system trigger an automatic alarm.

Schematic diagram of alarm tones

Tone sequence for various alarm priorities

Alarm priority	Standard (according to IEC 60601-1-8)	Repetitive
High	Depending on the overall alarm situation, this tone sequence may be played as a 5-tone sequence due to the timing of the individual alarms.	Yes
Medium		Yes
Low		No

The described acoustic alarm signals are handled by a backup loudspeaker if the primary loudspeaker fails. This plays the tone sequences of the high and medium alarm priority with constant tone frequency and constant volume.


Tone signals during operation

When	Signal
Therapy start or change of ventilation mode	-
Timeout	


Combination lock on the drawer (optional)

The lower drawer on the large drawer module is provided with a combination lock. When delivered, the combination on the lock is set to **3333**.


Opening the combination lock

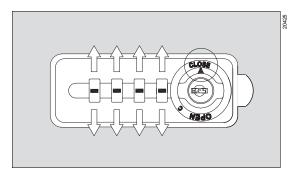

1 Set the combination with the digit wheels.

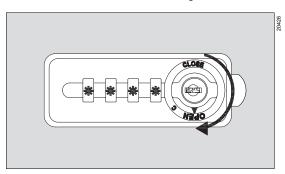


2 To open the combination lock, turn the knob to **OPEN**.

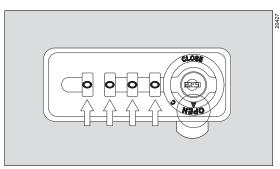


The combination lock is opened. The digit wheels return automatically to **0000**.

Alternatively, the combination lock can be opened with the supplied key. For information on the procedure, see "Loss of the combination" on page 273.


Changing the combination

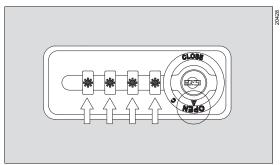

While the combination lock is open, set the current combination.

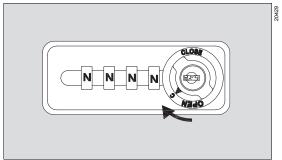


- 2 Turn the knob to position C.
- 3 Set the new combination.

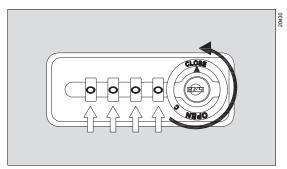


4 Turn the knob to **CLOSE**.

The new combination is stored.

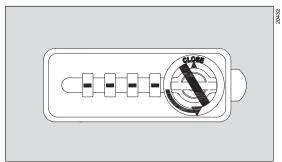


2 Turn the key counterclockwise to the final position stop.

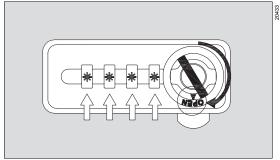

Turn the knob to **OPEN**. The digit wheels show the current combination.

Loss of the combination

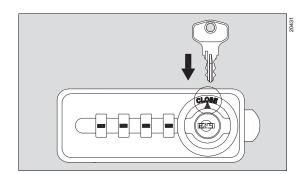

1 Insert the key.

Annex

Abbreviations	275
Symbols	277
Product labels	280
Overview of the manu structure	282

Abbreviations

Abbreviation	Explanation
%, Vol %	Percentage gas ratio, related to total volume
A	Ampere
Add. O2	Emergency O ₂ delivery
AGS	Anesthetic gas receiving system
AGSS	Anesthetic gas scavenging system
Air	Medical compressed air
APL	Adjustable Pressure Limitation, adjustable pressure limitation
APRV	Airway Pressure Release Ventilation
ASA	American Society of Anesthesi- ologists, american society of anesthesiologists
ATPS	Ambient Temperature and Pressure, Saturated Ambient temperature and ambient pressure, 100 % relative humidity
Aux. O2	O2 insufflation
BIPAP	Biphasic Positive Airway Pres- sure, spontaneous breathing under continuous positive airway pressure with 2 different pres- sure levels
BMI	Body mass index
BTPS	Body Temperature and Pressure, Saturated 37 °C (98.6 °F), ambient pressure, 100 % relative humidity
CAL	Display when a measurement value is calibrated.
CBM mode	Cardiac bypass mode
Cdyn	Dynamic compliance (patient)
CISPR	Comité International Spécial des Perturbations Radioélectriques International special committee on radio interference

A la la una di addi a un	F		
Abbreviation	Explanation		
cmH2O	Centimeters of water		
CMV	Controlled Mandatory Ventilation		
CO	Carbon monoxide		
CO ₂	Carbon dioxide		
COM	Serial port		
CPAP	Continuous Positive Airway Pressure, continuous positive airway pressure		
CSA	Canadian Standards Agency		
dB(A)	Sound pressure level, A-weighted		
Des	Desflurane		
ΔΟ2	Difference between inspiratory and expiratory O2 concentration		
ΔPsupp	Pressure support above PEEP		
EMC	Electromagnetic compatibility		
Enf	Enflurane		
ERR	Display when a measured value cannot be determined.		
ESD	Electrostatic Discharge Electrostatic Discharge		
FG	Fresh gas		
FiO ₂	Inspiratory oxygen fraction		
FTP	File Transfer Protocol		
GPL	General Public Licence		
Hal	Halothane		
HF	High-frequency		
HME	Heat and moisture exchanger		
HMEF	HME filter		
hPa	Hectopascal		
Hz	Hertz		
I:E	Ratio of inspiratory time to expiratory time		
ID	Identification		

A b b way i ati a ia	Evalenation		
Abbreviation	Explanation		
Insp. term.	Inspiration termination criterion in % based on peak inspiratory flow		
Iso	Isoflurane		
kg	Kilogram		
kPa	Kilopascal		
L	Liter		
LAN	Local area network		
lbs	Pound; unit of mass		
LED	Light-emitting diode		
LGPL	Lesser General Public Licence		
MAC	Minimum Alveolar Concentration		
Man/Spon Manual / Spon- taneous	Manual ventilation / Sponta- neous breathing		
mbar	Millibar		
MEDIBUS.X	Communication protocol for medical devices with uniform data definition for all devices		
min	Minute		
mL	Milliliter		
mmHg	Millimeter of mercury		
MPa	Megapascal		
MV	Minute volume		
N2O	Nitrous oxide		
NTP	Network Time Protocol, standard for synchronizing clocks		
O2	Oxygen		
O2+	O2 flush		
Pa	Pascal; unit of pressure		
Paw	Airway pressure		
PEEP	Positive end-expiratory pressure		
Pinsp	Inspiratory pressure		
PIP	Peak inspiratory pressure		
Pmax	Maximum pressure		
Pmean	Mean pressure		

Abbreviation	Explanation		
png	Graphics format		
Pplat	Plateau pressure		
ppm	Parts per million		
R	Resistance		
RFID	Radio Frequency Identification, radio frequency identification		
RR	Respiratory rate		
RRmin	Minimum respiratory rate		
Sev	Sevoflurane		
Slope	Pressure rise time		
SNMP	Simple Network Management Protocol		
STAPD	Standard Temperature, Ambient Pressure, Dry 20 °C (68 °F), ambient pres- sure, dry gas		
STPD	Standard Temperature and Pressure, Dry 20 °C (68 °F), 1013 hPa, dry gas		
TC	Time constant		
Ti	Inspiratory time		
UMDNS	Universal Medical Device Nomenclature System Nomenclature for medical devices		
USB	Universal Serial Bus		
V	Volt		
VT	Tidal volume		
xMAC	Accumulated multiple of the MAC values of anesthetic agents and N2O		

Symbols

Symbol	Explanation			
•••	Manufacturer			
×xxx ×xx	Date of manufacture			
	WEEE marking			
[]	Observe the instructions for use			
(3)	Warning! Strictly follow these instructions for use			
\triangle	Caution! Follow the accompanying documentation (Symbol)			
	Attention! (safety sign)			
	Group <i>Views</i>			
u	Group <i>Trends/Data</i>			
\triangle	Group <i>Alarms</i>			
Eml	Group Procedures			
F	Group System setup			
\bigcirc	Group Start/Standby Device on/ Standby			
I	Main switch on			
Ö	Main switch off			
.	Patient category Neo			
- Á	Patient category Ped			

Symbol	Explanation		
Ť	Patient category Adult		
Audio paused	Temporarily suppress acoustic alarm signal		
	Audio paused Acoustic alarm signal is temporarily suppressed.		
\bigotimes	Alarm inactive		
Ŵ	Alarm temporarily inactive		
- D-	Mains power		
×	Mains power unavailable		
	Battery completely charged		
	Battery empty		
	Battery is charging.		
-	Central gas supply connected and pressure within specified range		
	Central gas supply not con- nected or pressure not within specified range		
	Gas cylinder full		
Û	Gas cylinder empty or gas cylinder valve turned off		
Ä	Gas cylinder pressure sensor not connected		
T O	Symbol for programmed Auto On		

Symbol	Explanation	Symbol	Explanation
<u>*</u>	Key for switching on and off and dimming the workplace illumination	\rightarrow	Expiration Identification mark on the breathing system and the breathing system cover
Add. O ₂	Emergency O2 delivery (Add. O2)		Breathing bag
*	Applied part of type BF (body floating)		Vaporizer plug-in system, "fixed" position
\Diamond	Potential equalization connector		Auto Exclusion Plug-in connec-
×	Closes the dialog window	Dräger-Vapor ® Auto Exclusion Operating Instructions	tion
/ *	Upper alarm limit	REF	Part number
	Lower alarm limit	SN	Serial number
	No alarm limit	LOT	Lot number
^	Spontaneous breathing activity by the patient		Use by: YYYY-MM-DD Expiration date
†	In lists: One line up	*	Keep away from sunlight
+	In lists: One line down		Storage temperature
†	In lists: One page up	<u>1</u>	
†	In lists: One page down	<u> </u>	Relative humidity
^		€	Ambient pressure
	Risk of crushing		Do not use if package damaged
	ESD warning label, follow the information on electromagnetic compatibility.	8	Do not reuse
<u> </u>	Locked	∜ SPARE PART	Spare part
3	Unlocked	뫔	LAN connection
	Inspiration Identification mark on the breath-	•	USB port
	ing system and the breathing system cover	☼	Identifies the interface for the workstation light
		₹	External fresh-gas outlet

Symbol	Explanation		
	CO2 absorber bypass		
← J	Enter key		
ECD	Connection for Embedded Control Display		
*	Indicates a changed view which has not yet been saved		
p	Identifies the tabs that lead to the page with the language settings.		
$\overline{}$	Read the flow at the center of the float.		
ACS	Advanced Cylinder Support		
N ₂ O	Identifies N2O cylinders. The color code conforms to the locally applicable standard.		
O ₂	Identifies O ₂ cylinders. The color code conforms to the locally applicable standard.		
AIR	Identifies Air cylinders. The color code conforms to the locally applicable standard.		
□ >	Gas outlet		
-	Gas inlet		
Total	Display on the total flow tube indicating the cumulative value of the individual flows		
MR	MR unsafe Do not use this device near MRI scanners.		
	Exemplary representation of the weight distribution of the nominal weight and the maximum total weight, see "Technical data".		
	Key for releasing the locking brake		
	Key for releasing the locking brake		

Symbol	Explanation
\uparrow	Up arrow key for lifting the ceiling supply unit
↓	Down arrow key for lowering the ceiling supply unit
	Other Dräger devices
	Ceiling-mounted version of the Perseus
	Docking or undocking of Perseus

Product labels

Product label	Explanation
□	When connecting auxiliary devices, be aware of the leakage current. Observe chapters "Connecting additional devices to the auxiliary power sockets" and "Technical data".
INCENTRATION IN THE PARTY OF TH	Transport instructions, "Intrahospital transport" on page 56
WINDSPIRE SECTION	Ensure that the control dial of the vaporizer is correctly positioned. Do not leave the control dial in the "T" position while the vaporizer is connected to the medical device.
≥ 200 mm (≥ 8 in)	Maintain the correct minimum distance of 200 mm (8 in) between the electrical connections and the gas cylinder.
max Max	Observe the correct flow of the anesthetic gas receiving system, see 132.

Product label	Explanation
	Observe the weight of the nominal configuration and the permissible total weight, see "Technical data".
△	Observe the weight of the nominal configuration and the permissible total weight, see "Technical data".
凸 nom. 65 kg	Observe the weight of the nominal configuration for the ceiling-mounted version on the flexibility trolley and the total permissible weight, see "Technical data".
mounting / dismounting only by Dräger service or authorized service personnel	Assembly and disassembly should only be performed by DrägerService or authorized service personnel.
	If a non-rebreathing system is used, make sure that the <i>Ext. FGO</i> operating mode is used. For additional information, see "Special forms of therapy" in chapter "Operation".

Overview of the menu structure

The following tables list the grouped buttons of the main menu bar with the resulting dialog windows of the same name and the tabs. For information on operation, see "Operating concept" on page 39.

Group \triangle

Button in main menu bar	Horizontal tab	Vertical tab	Description
Alarms	Limits		Displaying or changing alarm limits
	Current alarms		Displaying information on active alarms
	Alarm history		Viewing the alarm history
	Settings		Setting the alarm tone volume Activating or deactivating CO ₂ alarms ¹⁾ Switching CBM mode on or off ¹⁾
CO ₂ alarms off ^{1), 2)}			Deactivating CO2 alarms
Autoset limits 1), 3)			Automatically adapts alarm limits to current measured or set values
Exit CBM 1), 4)			Exit CBM mode

- 1) Only during operation, not in Standby mode
- 2) Only in the modes: Manual / Spontaneous, Ext. FGO, Pause
- 3) Only in the modes: PSV, PC, PC APRV (optional), VC CMV / AutoFlow, VC CMV
- 4) Only in CBM mode

Group ___

Button in main menu bar	Horizontal tab	Vertical tab	Description
Views ¹⁾			Switching to other configured views Resetting current view to start settings Displaying alarm limits, units, mini- trends, and loops
ili View ¹⁾ ili View ¹⁾ ili View ¹⁾			Switching between the 3 configured views.
Export screenshot			Exporting screenshots to USB flash drive

¹⁾ Only during operation, not in Standby mode

Group 🔎

Button in main menu bar	Horizontal tab	Vertical tab	Description
Trends/Data	Graphical trends	Overview	Displaying trends of measured values
		Vent. 1	in graphic form
		Vent. 2	
		Anesthesia	
	Tabular trends	Overview	Displaying trends of measured values
		Vent. 1	in table form
		Vent. 2	
		Anesthesia	
	Values	Ventilation 1)	Displaying overview of current measured values
		Gases 1)	sureu values
		System	
	Logbook		Displaying the logbook
	Export 2)		Exporting data to USB flash drive

- Only during operation, not in *Standby* mode
 Only in *Standby* mode

Group 🐚

Button in main menu bar	Horizontal tab	Vertical tab	Description
Procedures ¹⁾	Insp./Exp. hold		Starts a maneuver during ventilation.
	One-step recruit- ment		
	Multi-step recruit- ment		
Flush ²⁾			Flushing the breathing system

- Only during operation, not in *Standby* mode
 Only in *Standby* mode

Group [₹

Button in main menu bar	Horizontal tab	Vertical tab	Description
System setup			Configuration of device functions and start settings, see page 150
Patient 1)			Setting patient data
Tests ²⁾			Displaying test results Testing the system Flushing the breathing system

- 1) Only during operation, not in Standby mode
- 2) Only in **Standby** mode

Group (1)

Button in main menu bar	Horizontal tab	Vertical tab	Description
Start 1)			Beginning or continuing a case
Standby ²⁾			Ending the case

- 1) Only in **Standby** mode
- 2) Only during operation, not in Standby mode

Password

Configuration password for Perseus A500 Software 2.0n

Cut-out from the instructions for use Perseus A500 Software 2.0n

To prevent unauthorized alteration, the start settings for Perseus A500 are protected by the following configuration password:

0000



Information for the configuration password

To prevent unauthorized alteration, the start settings for Perseus A500 are protected by a password with 0 to 8 digits. Information on the start settings, see page 150.

The configuration password appears on this page of the instructions for use. Cut out the area with the password and keep in a place which is safe from access by unauthorized persons.

The configuration password can only be reset by specialized service personnel.

This page has been left blank intentionally.

Index

A	APL_valve
Abbreviations	Fitting
Accessories	Set
Connecting	APRV
Accessories and consumptions	Audio paused
Activation code	Auto On
Additional devices	AutoFlow
Connecting	Auxiliary power sockets
Advanced Cylinder Support	
Alarm tone volume	В
Set153, 167	- "
Setting110	Battery
Alarms	Charging
Acknowledging alarm messages 141	Blower module
Adjusting alarm limits	Installing 63
Alarm settings when changing ventilation	Breathing bag
	Connecting
modes	Breathing hoses
Breath detection	Fitting
Causes	Breathing sound
CBM mode	Setting the volume
	Breathing system
Colors	Assembling 63
Details	Checking
Displaying alors limits	Flushing
Displaying alarm limits	Inserting66
Downgrading alarm messages	Pressure/flow characteristics 249
Dräger Vapor 3000	Warmer
Gas cylinders	Breaths
MV alarms	Synchronizing
Priorities	
Quick setup	C
Remedy	
Setting alarm limits automatically 111, 154	Cable channel
Setting alarm response	Cardiac bypass mode 93
Setting alarm tone volume	CBM mode
Setting limits	Ceiling-mounted version 47
Setting the volume	Central brake
Specifying alarm limits	Central gas supply
Suppressing the alarm tone	Connecting
Viewing the alarm history	Connections
Anesthetic gas concentration	Change of patient
Prediction	Checklist
Anesthetic gas receiving system (AGS)18, 21, 132	System test
Connecting	Circuit plug
Replacing	CLIC adapter
Anti-interchange security	Mounting
	Clinical users

CO2 absorber	Emergency start-up
Filling and mounting	Environments of use
CO2 alarms	Expiratory port
Adjusting	Exporting configurations
Color scheme	External fresh-gas outlet
Setting108	
Colors for parameters	F
Specifying	
Colors of the control elements	Factory settings
COM 1, COM 227, 168	Resetting to factory settings 164
Combination lock	Failure
Changing the combination	Flow measurement
Loss of the combination	Fresh-gas delivery
Opening	Gas measurement
Complete failure	Gas supply
Compliance correction	Failure of the gas supply
Configurable screen contents	Filter
Confirm	Fitting
Consumables	FiO2 concentration
Connecting	Prediction
Consumptions	Flexible arm
Displaying	Mounting
Resetting	Flow measurement faulty
Contraindications	Flow sensors
Cylinder pressure reducer problems 190	Inserting
	Flow waveform scale
D	Specifying
Data	Flush & Dry
Data Experting 113	Flush function
Exporting	Front Coiling mounted version
Date Set	Ceiling-mounted version
	Trolley version
Daylight savings time	Fuses
	ruses
Specifying	
Device combinations	G
Transferring171	Gas cylinders
Device name	Advanced Cylinder Support
Changing	Changing
Diagrams	Checking
Display for current time and Auto On 31, 32	Connecting
Drying function	Connections
brying function	Gas flow plan
=	Gas inlets
E	Gas measurement faulty
Econometer (optional)	Gas mixer (electronically controlled)
Electromagnetic compatibility11	Gas mixer (mechanically controlled)
EMC	Total flow tube 32, 38, 82, 89
EMC declaration	Gas mixtures for the carrier gas
Emergency O2 delivery	3. 2 3 3 3 3
Checking	
OHECKING	

Gas scavenging	M
Ensuring	MAC
Gas supply	Main screen
Checking128	Main switch
Establishing57	Mains power supply
	Establishing
H	Mains power supply failure
Hardware 40	Maintenance
Hardware	Maneuvers (optional)
Hose configurations	Manual ventilation
Fitting	Measured values
Titung	Viewing
	Menu structure
1	Minimum O ₂ delivery
Illumination	Minimum O ₂ flow
Importing configurations	Mini-trends
Indications	Minute volume measurement
Infinity ID	Monitoring
Accessory exchange interval monitoring 270	
Accessory replacement intervals 169	N
Accessory support	
Activating the function	N ₂ O
InfinityID	Operation without N2O
Inspection	Name
Inspiratory port	Changing
Interfaces	Names for screen layouts
IT networks	Configuring
	Network interface
L	New case
T	Defining
Language	Non-rebreathing system
Select	Connecting
Latex	Using
Leakage Causes	3
Current values	0
Locating	O
Leakage assistant	O2 flowmeter
Leakage test	O2 flush
Linked therapy controls	Checking
Logbook	O2 insufflation
Creating entries	Using
Viewing	O2 switch
Loops	O2 waveform scale
Configuring	Specifying
Using106	O2+
Low-flow wizard	Checking
Lung recruitment	Using 91 Open-source software 256
	Operation 256
	Safety information 85
	Galoty information

P	Screen
Parameter fields	Contents
Parameters	Failure
Setting43	Saving the layout
Password	Setting brightness 108, 150
Patient category	Screen brightness
Patient data	Setting
	Screenshot
Changing113	Exporting
Checking	Service
Loading87	Service personnel
Pause	Setting procedure
Potential equalization	Canceling
Establishing54	Settings
Potential equalization pin	Resetting to starting settings
Power sockets	Soda lime
Power supply failure	Checking
Prediction	Software options
Anesthetic gas concentration	Activating
FiO2 concentration	Displaying
Pressure-controlled ventilation	Sound volume
Pressure-supported ventilation260, 262	Setting alarm tones
r ressure-supported vertiliationzoo, zoz	Setting breathing sound
_	Spontaneous breathing 90, 260
Q	Start dialog
Quick setup	Opening
Quick Scrap	Start settings
_	Gas delivery
R	Patient category
Rating plate	Resetting
Rear	Ventilation
Ceiling-mounted version	Start values 41
Trolley version	Stopwatch
Recruitment95	Using
Remote control	Storage
Remote service	Symbols
Repair220	Synchronized ventilation
Reprocessing personnel	System test
Resistance	Checklist
Respiratory support	Configuring
Rotary knob	Test types
riotary knob	
0	Т
S	Technical data
Safety checks	Technical data
Safety sensor	
Sample line	Testing the system
Connecting	Therapy 116
Scale	Ending
Adjusting	Setting and starting
-,	Therapy controls
	Set

Tipping stability Increasing	Volumeter Specifying the scale 151
Trademarks	Specifying the scale
Transport	Volume-controlled ventilation 260, 265
Trends	volume controlled ventuation 200, 200
Displaying	•••
Exporting	W
MVxCO2	Water trap
Oxygen uptake	Checking
	Emptying or replacing
Device	Waveforms 151
Turn on	Configuring
Device	Overview
Turning on the device	Saving view
Automatic	Specifying the scale
Device80	Working hours
	Set
U	
Units	X
	xMAC
Displaying	XIVIAC 100, 112
Set	
USB interface	
User group requirements	
User groups	
Users	
User-specific settings	
V	
Valves	
Inserting	
Vapor View	
Vaporizer	
Checking127	
Fitting61	
Illumination	
Using	
Ventilation mode	
Changing	
Ventilation modes	
Ventilation parameters	
Linked setting	
Ventilator failure	
View	
Available views	
Changing	
Customizing	
Volume	
Setting the breathing sound volume153, 167	

These instructions for use are provided for customer information only and are only updated or exchanged upon customer request.

CE 2123

((

Directive 2014/53/EU concerning radio equipment

Manufacturer

Drägerwerk AG & Co. KGaA

Moislinger Allee 53-55 23542 Lübeck Germany

+49 451 8 82-0

FAX +49 451 8 82-2080

http://www.draeger.com

9510595 – 6500.510 en © Drägerwerk AG & Co. KGaA Edition: 3 – 2020-01

(Edition: 1 – 2018-10)

Dräger reserves the right to make modifications to the medical device without prior notice.

