

Bovie® ELECTROSURGICAL GENERATOR



USER'S GUIDE



USER'S GUIDE

This manual and the equipment it describes are for use only by qualified medical professionals trained in the particular technique and surgical procedure to be performed. It is intended as a guide for using the Bovie® IDS-300 only.

Additional technical information is available in the *Bovie® IDS-300 Service Guide*.

Equipment Covered in this Manual

Bovie® IDS-300

Reference No.: IDS-300

For Information Contact

Bovie Medical Corporation • 5115 Ulmerton Road, Clearwater, Florida 33760 USA

U.S. Phone 1-800-537-2790 Fax 1-800-323-1640 • International Phone +1-727-384-2323 Fax +1-727-347-9144

www.boviemed.com • sales@boviemed.com



EC REP

Peter J. Smith Medical Products Marketing
18 Yeates Close
Thame OX9 3AR, UK



Made in USA
Printed in USA

©2011 Bovie Medical Corporation. All rights reserved. Contents of this publication may not be reproduced without the written permission of Bovie Medical Corporation.

Bovie® Part Number MC-55-059-001 Rev. 5

CONVENTIONS USED IN THIS GUIDE

WARNING:

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION:

Indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE:

Indicates an operating tip, a maintenance suggestion, or a hazard that may result in product damage.

TABLE OF CONTENTS

- Equipment Covered in this Manualiii
- For Information Contactiii
- Conventions Used in this Guideiii
- Introducing the Bovie® IDS-3001-1**
 - Key Features1-2
 - Components and Accessories1-3
 - Safety1-4
- Controls, Indicators, and Receptacles2-1**
 - Front Panel2-2
 - Symbols on the Front Panel2-3
 - Preset Controls2-4
 - Cut and Blend Controls2-5
 - Coag Controls2-6
 - Bipolar Controls2-7
 - Indicators2-8
 - Power Switch and Receptacles2-9
 - Rear Panel2-10
 - Symbols on the Rear Panel2-10
- Getting Started3-1**
 - Initial Inspection3-2
 - Installation3-2
 - Function Checks3-2
 - Setting Up the Unit3-2
 - Checking the Return Electrode Alarm3-2
 - Confirming Modes3-3
 - Checking Bipolar Mode (with bipolar footswitch)3-3
 - Checking Monopolar Mode (with monopolar footswitch)3-3
 - Checking Monopolar Mode (with handswitch)3-3
 - Performance Checks3-3
- Using the Bovie® IDS-3004-1**
 - Inspecting the Generator and Accessories4-2
 - Setup Safety4-2
 - Setting Up4-3
 - Preparing for Monopolar Surgery4-4
 - Applying the Return Electrode4-4
 - Connecting Accessories4-4
 - Preparing for Bipolar Surgery4-5
 - Setting and Recalling Memory Presets4-5
 - Memory4-5
 - Memory Function Overview4-5
 - Setting Your Presets4-6
 - Memory Feature (Last Selected Preset)4-6
 - Activating the Unit4-7
 - Activation Safety4-8

Maintaining the Bovie® IDS-300	5-1
Cleaning	5-2
Periodic Inspection	5-2
Fuse Replacement	5-2
Troubleshooting	6-1
Repair Policy and Procedures	7-1
Responsibility of the Manufacturer	7-2
Returning the Generator for Service	7-2
Step 1 – Obtain a Returned Goods Authorization Number	7-2
Step 2 – Clean the Generator	7-2
Step 3 – Ship the Generator	7-2
Technical Specifications	A-1
Performance Characteristics	A-2
Input Power	A-2
Duty Cycle	A-2
Dimensions and Weight	A-2
Operating Parameters	A-2
Transport and Storage	A-2
Audio Volume	A-3
Return Electrode Sensing	A-3
Low Frequency (50-60 Hz) Leakage Current	A-3
High Frequency (RF) Leakage Current	A-4
Standards and IEC Classifications	A-4
Class I Equipment (IEC 60601-1)	A-4
Type CF Equipment (IEC 60601-1) / Defibrillator Proof	A-4
Drip Proof (IEC 60601-2-2)	A-4
Electromagnetic Interference	A-4
Electromagnetic Compatibility (IEC 60601-1-2 and IEC 60601-2-2)	A-4
Voltage Transients (Emergency Generator Mains Transfer)	A-4
Output Characteristics	A-5
Maximum Output for Monopolar and Bipolar Modes	A-5
EMC Compliance	A-5
Output Power Curves	A-9
Warranty	B-1

LIST OF FIGURES

Figure 2 – 1	Layout of controls, indicators, and receptacles on the front panel.....	2-2
Figure 2 – 2	Controls for setting and recalling presets.....	2-4
Figure 2 – 3	Controls for the Cut and Blend modes.....	2-5
Figure 2 – 4	Controls for the Coag mode.....	2-6
Figure 2 – 5	Controls for the Bipolar mode.....	2-7
Figure 2 – 6	Indicators for power, return electrodes, and footswitch control.....	2-8
Figure 2 – 7	Location of the unit power switch and front panel receptacles.....	2-9
Figure 2 – 8	Layout of connectors and controls on the rear panel.....	2-10
Figure 5 – 1	Fuse holder.....	5-2
Figure A – 1	Display setting vs. output power in watts at rated load.....	A-9
Figure A – 2	Output voltage (V_{peak}) versus power setting at rated load.....	A-9
Figure A – 3	Output power vs impedance for Cut I mode.....	A-10
Figure A – 4	Output power vs impedance for Cut II mode.....	A-10
Figure A – 5	Output power versus impedance for Blend I (minimum) mode.....	A-11
Figure A – 6	Output power versus impedance for Blend II (maximum) mode.....	A-11
Figure A – 7	Output power vs impedance for Pinpoint mode.....	A-12
Figure A – 8	Output power vs impedance for Spray mode.....	A-12
Figure A – 9	Output power vs impedance for Bipolar mode.....	A-13

INTRODUCING THE BOVIE® IDS-300

This section includes the following information:

- *Key Features*
- *Components and Accessories*
- *Safety*

.....
CAUTIONS:

Read all warnings, cautions, and instructions provided with this generator before using.

.....
Read the instructions, warnings, and cautions provided with electro-surgical accessories before using. Specific instructions are not included in this manual.
.....

KEY FEATURES

The Bovie® IDS-300 includes the latest technology. This unit offers unsurpassed performance, flexibility, reliability, and convenience.

It includes the following features:

- **Two Cut Modes, Cut I & Cut II**

Two cut modes give the surgeon flexibility to cut all types of tissue without losing performance.

Cut I generates constant output power over a wide range of impedances. Refer to Figure A-1 in the *Technical Specifications* section of this guide.

Cut II is a softer cut that generates constant output power over a small range of impedances. Refer to Figure A-2 in the *Technical Specifications* section of this guide.

- **Blend with 10 settings**

The Blend mode is a combination of cutting and hemostasis. The IDS-300™ gives the surgeon freedom to adjust the desired level of hemostasis. A setting of 1 is minimal blend with maximum cutting effect. A setting of 10 is maximum hemostasis (blend) with minimal cutting effect. This adjustment is easily achieved by an incremental adjustment. Refer to Section 2, *Controls, Indicators, and Receptacles, Cut and Blend Controls*. The Blend mode improves the rate of targeted tissue desiccation without increasing the power delivered by the generator.

- **Presets**

The surgeon can store 10 user-defined presets for easy recall of frequently used settings.

- **Two levels of coagulation: Pinpoint and Spray**

Pinpoint provides precise control of bleeding in localized areas.

Spray provides greater control of bleeding in highly vascular tissue over broad surface areas.

- **Return electrode sensing and contact quality monitoring**

The IDS-300™ incorporates a return electrode contact quality monitoring system (Bovie NEM™). This system detects the type of return electrode: solid or split. The system also continually monitors the contact quality between the patient and the split return electrode. This feature is designed to minimize patient burns at the return electrode site.

- **FDFS™ (Fast Digital Feedback System)**

The FDFS™ (Fast Digital Feedback System) measures voltage and current at 5,000 times a second and immediately adjusts the power to varying impedance during the electrosurgical procedure. The unit's digital technology senses and responds to changes in tissue and density. Unlike analog, this feature reduces the need to adjust power settings manually.

NOTICE:

The Bovie NEM™ system recommends that you use a split return electrode.

- **Memory**

The unit automatically powers up to the last selected preset settings.

- **Isolated RF output**

This minimizes the potential of alternate site burns.

- **Standard connectors**

These connectors accept the latest monopolar and bipolar instruments. Refer to Section 2, *Controls, Indicators, and Receptacles* to learn more.

- **Self diagnostics**

These diagnostics continually monitor the unit to ensure proper performance.

COMPONENTS AND ACCESSORIES

You should receive the following components with your generator:

- Bovie® IDS-300
- Hospital-grade power cord (110 VAC and 220 VAC)
- User's Guide
- Service Guide

SAFETY

The safe and effective use of electrosurgery depends to a large degree on factors solely under the control of the operator. There is no substitute for a properly trained and vigilant medical staff. It is important that they read, understand, and follow the operating instructions supplied with this electrosurgical equipment.

Physicians have used electrosurgical equipment safely in numerous procedures. Before starting any surgical procedure, the surgeon should be familiar with the medical literature, complications, and hazards of using electrosurgery in that procedure.

To promote the safe use of the Bovie® IDS-300, this section presents the warnings and cautions that appear throughout this user's guide. It is important that you read, understand, and follow the instructions in these warnings and cautions so that you can operate this equipment with maximum safety. It is also important that you read, understand, and follow the instructions for use in this user's guide.

WARNINGS:

Hazardous Electrical Output - This equipment is for use only by trained, licensed physicians.

Danger: Fire / Explosion Hazard - Do not use the Bovie® IDS-300 in the presence of flammable materials.

Fire / Explosion Hazard - The following substances will contribute to increased fire and explosion hazards in the operating room:

- Flammable substances (such as alcohol based skin prepping agents and tinctures)
- Naturally occurring flammable gases which may accumulate in body cavities such as the bowel
- Oxygen enriched atmospheres
- Oxidizing agents (such as nitrous oxide [N₂O] atmospheres).

The sparking and heating associated with electrosurgery can provide an ignition source. Observe fire precautions at all times. When using electrosurgery in the same room with any of these substances or gases, prevent their accumulation or pooling under surgical drapes, or within the area where electrosurgery is performed.

Connect the power cord to a properly polarized and grounded power source with the frequency and voltage characteristics that match those listed on the back of the unit.

Electric Shock Hazard - Connect the generator power cord to a properly grounded receptacle. Do not use power plug adapters.

Electric Shock Hazard - Always turn off and unplug the generator before cleaning.

Fire Hazard - Do not use extension cords.

Patient Safety - Use the generator only if the self-test has been completed as described. Otherwise, inaccurate power outputs may result.

Failure of the high frequency electrosurgical equipment could result in an unintended increase of output power.

The instrument receptacles on this generator are designed to accept only one instrument at a time. Do not attempt to connect more than one instrument at a time into a given receptacle. Doing so will cause simultaneous activation of the instruments.

Use the lowest output setting necessary to achieve the desired surgical effect. Use the active electrode only for the minimum time necessary in order to lessen the possibility of unintended burn injury. Pediatric applications and/or procedures performed on small anatomic structures may require reduced power settings. The higher the current flow, and the longer the current is applied, the greater the possibility of unintended thermal damage to tissue, especially during use on small structures.

WARNINGS:

Use electrosurgery with caution in the presence of internal or external devices such as pacemakers or pulse generators. Interference produced by the use of electrosurgical devices can cause devices such as pacemakers to enter an asynchronous mode or can block the pacemaker effect entirely. Consult the device manufacturer or hospital Cardiology Department for further information when use of electrosurgical appliances is planned for patients with cardiac pacemakers or other implantable devices.

If the patient has an Implantable Cardioverter Defibrillator (ICD), contact the ICD manufacturer for instructions before performing an electrosurgical procedure. Electrosurgery may cause multiple activation of ICDs.

Do not use electrosurgical equipment unless properly trained to use it in the specific procedure being undertaken. Use by physicians without such training has resulted in serious, unintended patient injury, including bowel perforation and unintended, irreversible tissue necrosis.

For surgical procedures where the high frequency current could flow through parts of the body having a relatively small cross-sectional area, the use of bipolar techniques may be desirable to avoid unwanted coagulation.

In some circumstances, potential exists for alternate site burns at points of skin contact (e.g., between the arm and the side of the body). This occurs when electrosurgical current seeks a path to the return electrode that includes the skin-to-skin contact point. Current passing through small skin-to-skin contact points is concentrated and may cause a burn. This is true for grounded, ground referenced, and isolated output generators.

To reduce the potential for alternate site burns, do one or more of the following:

- Avoid skin-to-skin contact points, such as fingers touching leg, when positioning the patient.
- Place 5 to 8 cm (2 to 3 in.) of dry gauze between contact points to ensure that contact does not occur.
- Position the return electrode to provide a direct current route between the surgical site and the return electrode which avoids skin-to-skin contact areas.
- In addition, place patient return electrodes according to the manufacturer's instructions.

Potential for alternate site burns increases if the return electrode is compromised. Bovie Medical Corporation recommends the use of split return electrodes and Bovie® generators with a contact quality monitoring system.

Do not wrap the accessory cords or return electrode cords around metal objects. This may induce currents that could lead to shocks, fires, or injury to the patient or surgical team.

.....
CAUTIONS:

At no time should you touch the active electrode or bipolar forceps. A burn could result.

Do not stack equipment on top of the generator or place the generator on top of electrical equipment. These configurations are unstable and/or do not allow adequate cooling.

Provide as much distance as possible between the electrosurgical generator and other electronic equipment (such as monitors). An activated electrosurgical generator may cause interference with them.

Non-function of the generator may cause interruption of surgery. A backup generator should be available for use.

Do not turn the activation tone down to an inaudible level. The activation tone alerts the surgical team when an accessory is active.

When using a smoke evacuator in conjunction with the electrosurgical generator, place the smoke evacuator a distance from the generator and set the generator volume control at a level that ensures that the activation tones can be heard.

The use of high frequency current can interfere with the function of other electromagnetic equipment.

When high frequency surgical equipment and physiological monitoring equipment are used simultaneously on the same patient, place any monitoring electrodes as far as possible from the surgical electrodes. Monitoring systems incorporating high frequency current-limiting devices are recommended.

Do not use needles as monitoring electrodes during electrosurgical procedures. Inadvertent electrosurgical burns may result.

To avoid the possibility of an electrosurgical burn to either the patient or the physicians, do not allow the patient to come in contact with a grounded metal object during activation. When activating the unit, do not allow direct skin contact between the patient and the physician.

Remove any loose fitting jewelry from the patient before activation.

Examine all accessories and connections to the electrosurgical generator before use. Ensure that the accessories function as intended. Improper connection may result in arcs, sparks, accessory malfunction, or unintended surgical effects.

When not using active accessories, place them in a holster or in a clean, dry, non-conductive, and highly visible area not in contact with the patient. Inadvertent contact with the patient may result in burns.

Studies have shown that smoke generated during electrosurgical procedures can be potentially harmful to patients and the surgical team. These studies recommend adequately ventilating the smoke by using a surgical smoke evacuator or other means.¹

.....
1. U.S. Department of Health and Human Services. National Institute for Occupational Safety and Health (NIOSH). Control of Smoke from Laser / Electric Surgical Procedures. HAZARD CONTROLS, Publication No. 96-128, September, 1996.

NOTICES:

If required by local codes, connect the generator to the hospital equalization connector with an equipotential cable.

Do not clean the generator with abrasive cleaning or disinfectant compounds, solvents, or other materials that could scratch the panels or damage the generator.

CONTROLS, INDICATORS, AND RECEPTACLES

This section describes:

- The Front and Rear Panels*
- Controls, Indicators, Receptacles, and Ports*











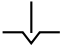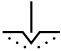


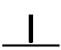



FRONT PANEL

Figure 2 – 1 Layout of controls, indicators, and receptacles on the front panel




Symbols on the Front Panel

Refer to the following table for descriptions of symbols found on the front panel of the Bovie® IDS-300.

SYMBOLS	DESCRIPTION
<i>Cut Controls</i>	
	Cut Mode
	Blend Mode
<i>Coag Controls</i>	
	Pinpoint Mode
	Spray Mode
<i>Bipolar Controls</i>	
	Bipolar Mode
<i>Indicators</i>	
	Split Return Electrode
	Solid Return Electrode
<i>Regulatory Symbology</i>	
	Read instructions before use.
	Defibrillator Proof Type CF Equipment
	RF Isolated – patient connections are isolated from earth at high frequency.
<i>Power Switch and Handpiece Connectors</i>	
	Return Electrode Receptacle
	Caution High Voltage
	Cut Mode
	Coag Mode
	Monopolar Handpiece Receptacle
	Bipolar Mode
	Bipolar Handpiece Receptacle

PRESET CONTROLS

Figure 2 – 2 Controls for setting and recalling presets

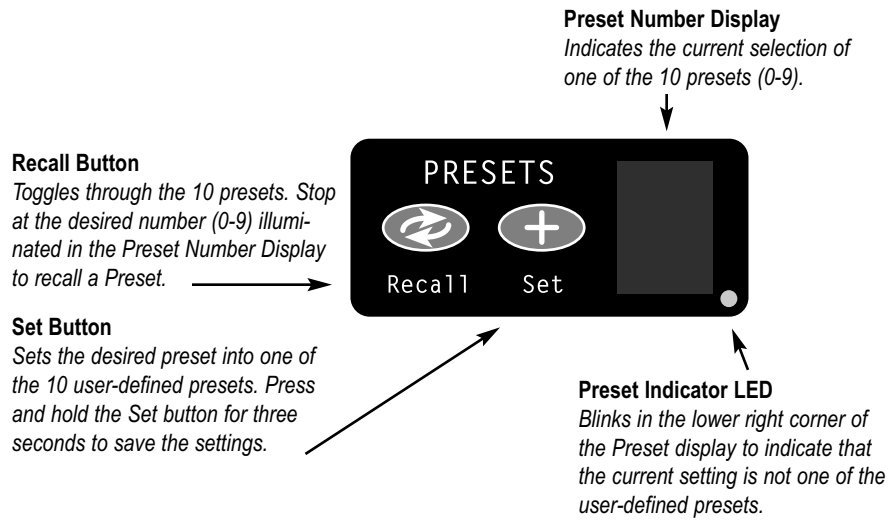


NOTICES:

The Bovie® IDS-300 incorporates 10 factory-set presets that are all set to zero and can be reset to your preferred settings.

Set and Recall are disabled while the unit is activated.

CUT AND BLEND CONTROLS

Figure 2 – 3 Controls for the Cut and Blend modes

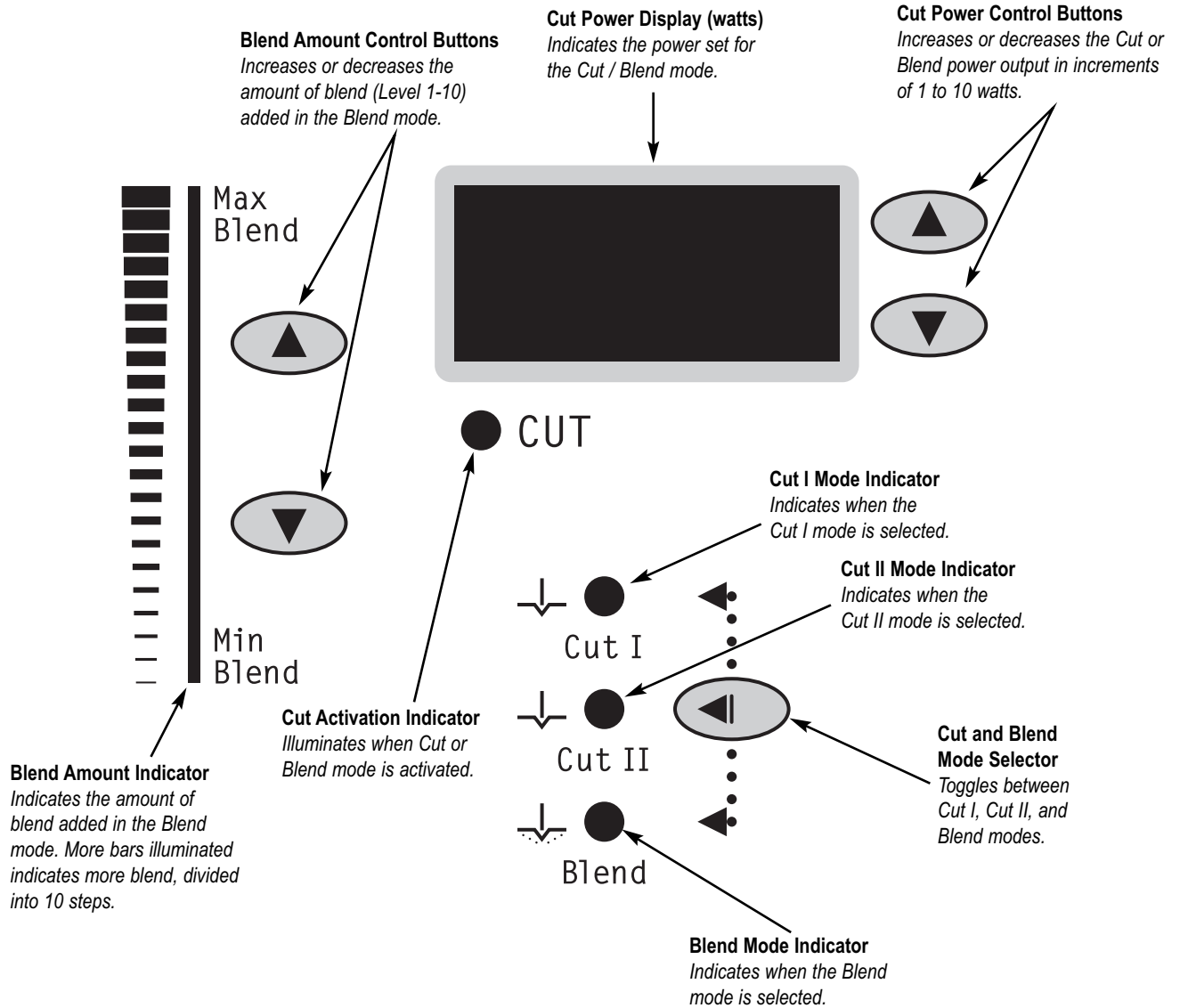


NOTICE:

When selecting the Blend mode, the unit defaults to a setting of minimum blend (only the first bar is illuminated).


COAG CONTROLS

Figure 2 – 4 Controls for the Coag mode

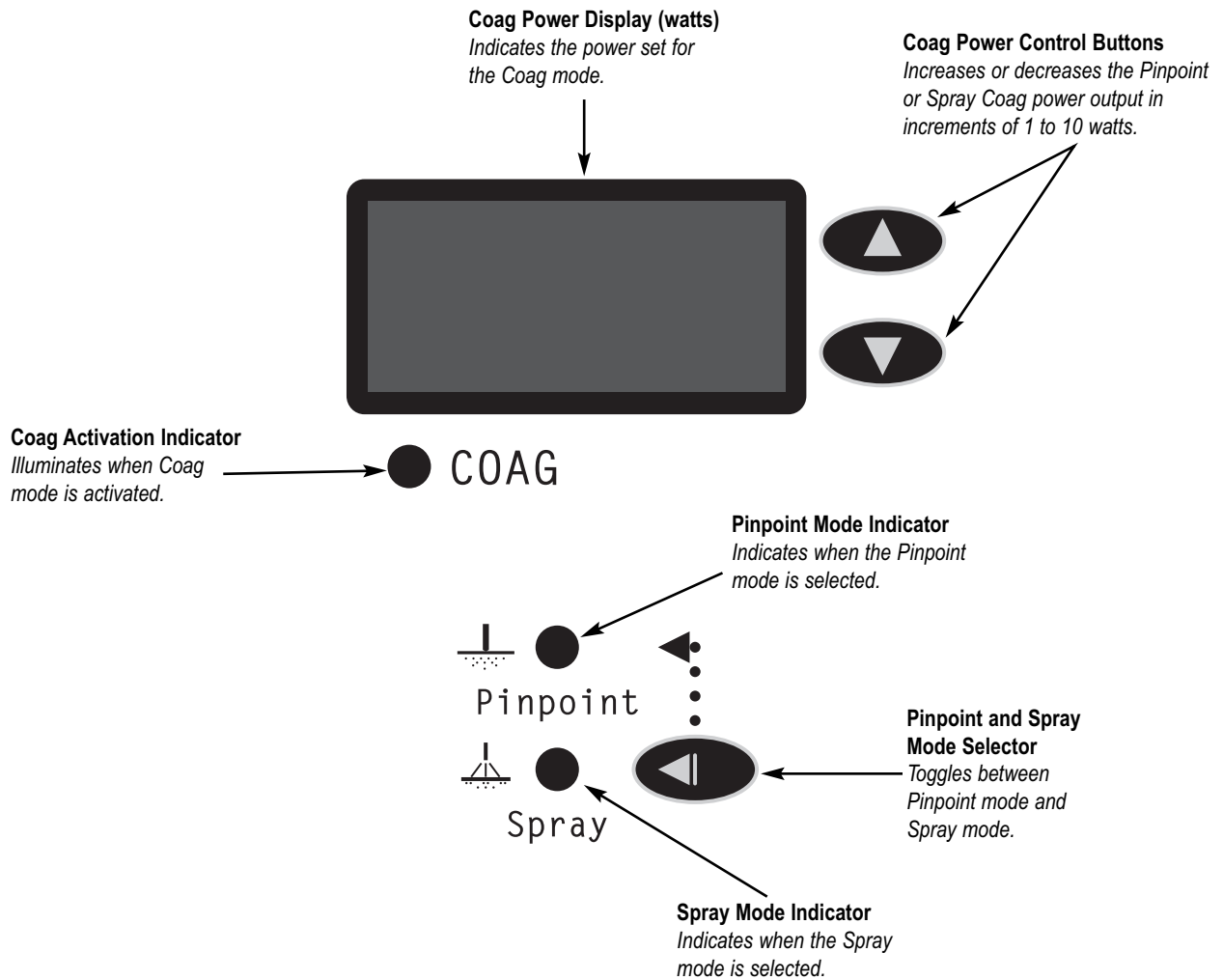


BIPOLAR CONTROLS

Figure 2 – 5 Controls for the Bipolar mode

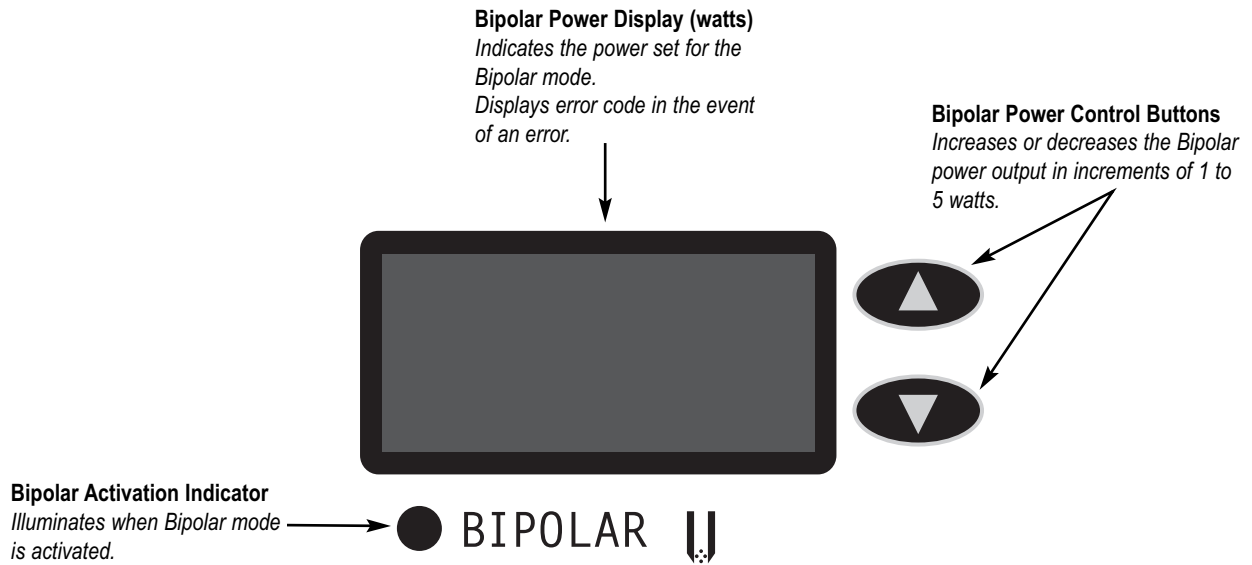


INDICATORS

Figure 2 – 6 Indicators for power, return electrodes, and footswitch control

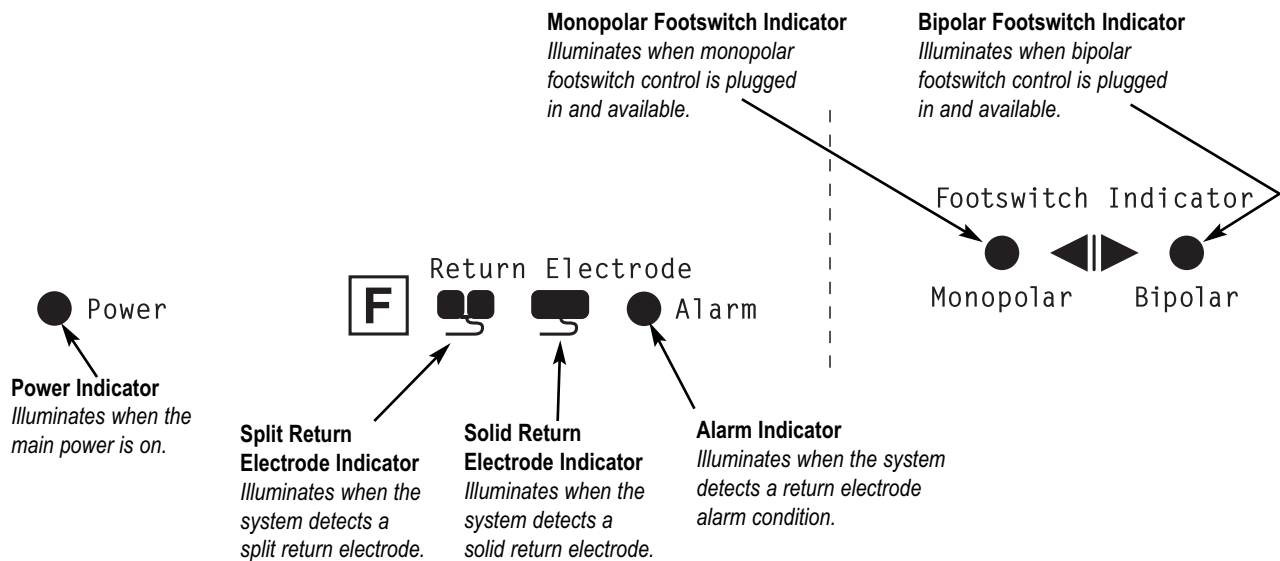


POWER SWITCH AND RECEPTACLES

Figure 2 – 7 Location of the unit power switch and front panel receptacles

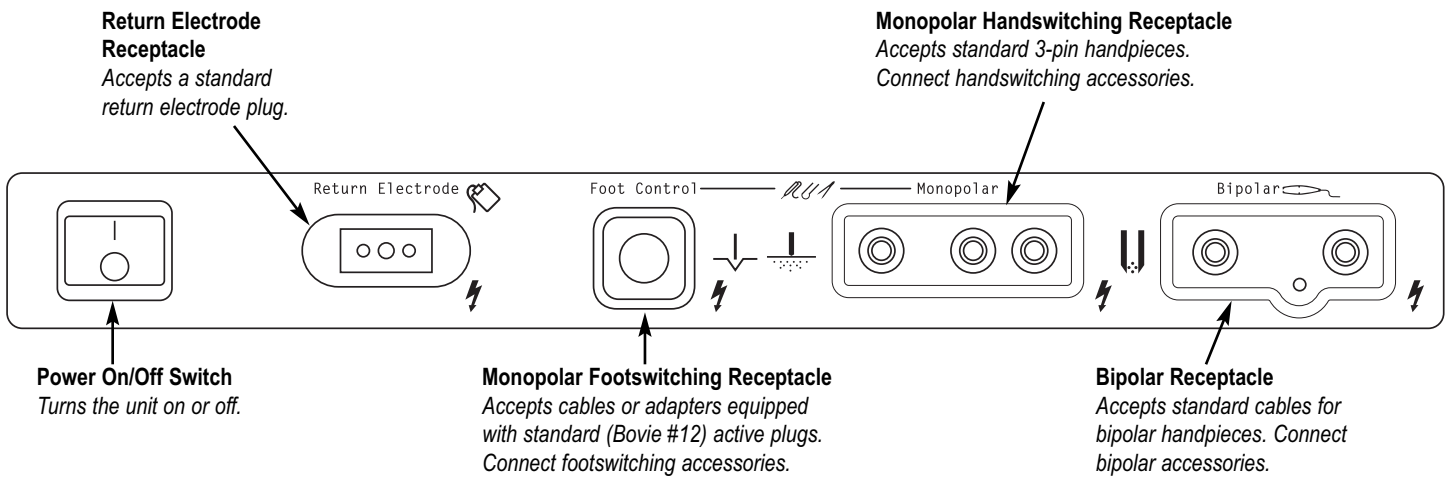

REAR PANEL

Figure 2 – 8 Layout of connectors and controls on the rear panel

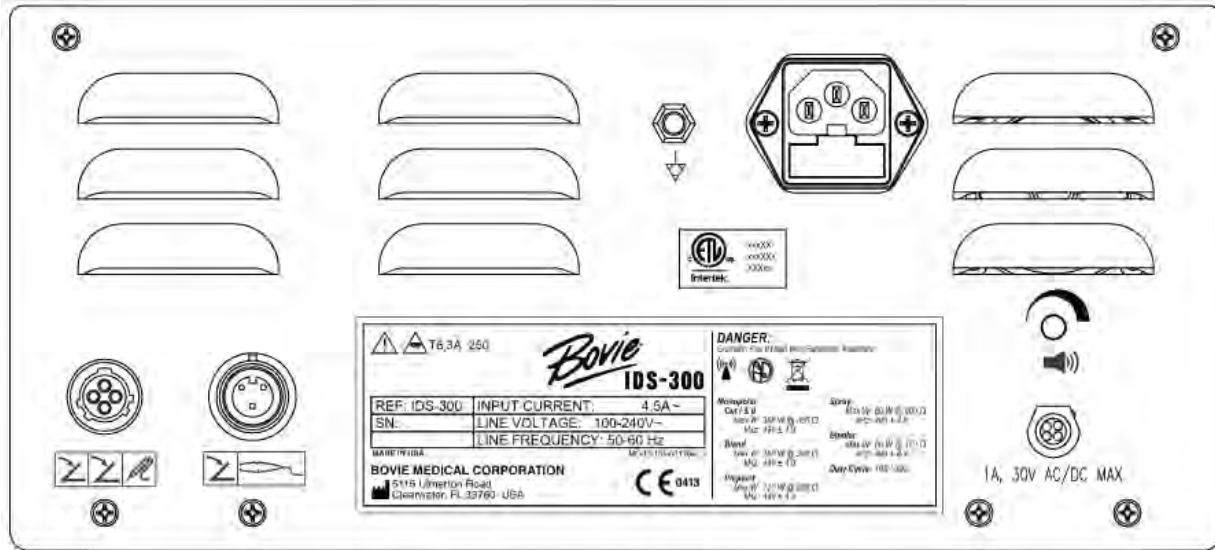

Symbols on the Rear Panel

Refer to the following table for descriptions of symbols found on the rear panel of the Bovie® IDS-300.

SYMBOLS	DESCRIPTION
	Equipotential Ground Stud
	Non-ionizing Radiation
	Volume Control
	Danger - Explosion Risk If Used With Flammable Anesthetics.
	Fuse Enclosed
	Relay Connector
	Monopolar Footswitch Input Jack
	Bipolar Footswitch Input Jack
	Manufacturer
	Caution, Consult Accompanying Documents
	Do Not Dispose of Unit in Municipal Waste Stream.

NOTICE:

* Please note that infected medical devices must be disposed of as medical/biohazard waste and cannot be included in used electronic equipment disposal/recycling programs. In addition, certain electronic products must be returned directly to Bovie Medical Corporation. Contact your Bovie® sales representative for return instructions.

GETTING STARTED

This section includes the following information:

- *Initial Inspection*
- *Installation*
- *Function Checks*
- *Performance Checks*

INITIAL INSPECTION

When you first unpack your Bovie® IDS-300, inspect it visually:

- Look for any signs of damage.
- Verify that the shipping package contains all items listed on the packing list.

If the unit or any accessories are damaged, notify Bovie Medical Corporation's Customer Service immediately. Do not use any damaged equipment.

INSTALLATION

Place the Bovie® IDS-300 on any flat surface with a tilt angle not more than 10°. The unit relies on natural convection cooling. Do not block its bottom or rear vents. Ensure that air flows freely on all sides of the unit.

WARNING:

Connect the power cord to a properly polarized and grounded power source with the frequency and voltage characteristics that match those listed on the back of the unit.

FUNCTION CHECKS

Upon initial installation of the unit, perform the tests listed below. Refer to the figures in the previous chapter for the location of connectors and controls.

WARNING:

At no time should you touch the active electrode or bipolar forceps. A burn could result.

Setting Up the Unit

1. Verify that the Power Switch is in the Off (O) position and that no accessories are connected to the unit.
2. Connect a hospital grade power cable to the AC power cable receptacle on the back of the unit, then to a properly grounded wall outlet.
3. Connect a two-button monopolar pencil to the appropriate receptacle. The use of Bovie® pencils is recommended.
4. Do not connect a patient return electrode at this time.
5. Turn the unit on by switching the power switch to the On (I) position.

Checking the Return Electrode Alarm

1. Adjust the power settings for each mode (Cut, Coag, Bipolar) to one watt.
2. Press the Coag button of the pencil. Verify that an alarm sounds for three seconds and the patient return electrode sensing alarm indicator light illuminates, indicating that no return electrode is connected to the unit.
3. Verify that adjusting the volume control on the back of the unit while the alarm is sounding does not change the alarm volume.

Confirming Modes

Confirm that you can select each mode and adjust the power up and down.

Checking Bipolar Mode (with bipolar footswitch)

1. Plug in the Bipolar footswitch. Verify that the Bipolar footswitch indicator illuminates.
2. Press the pedal on the Bipolar footswitch. Verify that the Bipolar mode activation indicator illuminates and that the system generates the Bipolar activation tone.
3. While activating the Bipolar mode, rotate the volume control over the full range to verify that the sound is audible throughout the range.
4. Confirm that releasing the pedal returns the unit to an idle state.

Checking Monopolar Mode (with monopolar footswitch)

1. Plug in the Monopolar footswitch. Verify that the monopolar footswitch indicator illuminates.
2. Connect a solid return electrode to the return electrode receptacle. Verify that the green solid return electrode indicator illuminates.
3. Press the Cut pedal (yellow) on the footswitch. Verify that the Cut mode activation indicator illuminates and that the system generates the Cut activation tone.
4. While activating the Cut mode, rotate the volume control over the full range to verify that the sound is audible throughout the range.
5. Press the Coag pedal (blue) on the footswitch. Verify that the Coag mode activation indicator illuminates and that the system generates the Coag activation tone.
6. While activating the Coag mode, rotate the volume control over the full range to verify that the sound is audible throughout the range.

Checking Monopolar Mode (with handswitch)

1. Connect a handswitching handpiece to the Monopolar handpiece receptacle.
2. Connect a solid return electrode to the return electrode receptacle. Verify that the green solid return electrode indicator illuminates.
3. Activate, one at a time, the Cut and Coag handswitching controls. Verify that each control causes the correct indicator and tone to sound.

PERFORMANCE CHECKS

After the unit has passed the preliminary functional test, it is ready for performance testing. A qualified biomedical engineer who is thoroughly familiar with electrosurgical devices should conduct this testing. The testing should include checking all modes of operation for proper function and power output.

USING THE BOVIE® IDS-300

This section contains the following procedures:

- *Inspecting the Generator and Accessories*
- *Setup Safety*
- *Setting Up*
- *Preparing for Monopolar Surgery*
- *Preparing for Bipolar Surgery*
- *Setting and Recalling Memory Presets*
- *Activating the Unit*
- *Activation Safety*

.....
CAUTIONS:

Read all warnings, cautions, and instructions provided with this generator before use.

.....
Read the instructions, warnings, and cautions provided with electro-surgical accessories before use. Specific instructions are not included in this manual.
.....

INSPECTING THE GENERATOR AND ACCESSORIES

Before each use of the Bovie® IDS-300, verify that the unit and all accessories are in good working order:

- Inspect for damage to the Electrosurgical Generator and all its connections.
- Verify that the appropriate accessories and adapters are present.
- Inspect all cords and connectors for signs of wear, damage, and abrasion.
- Verify that no errors occur when you turn on the unit.

SETUP SAFETY

WARNINGS:

Hazardous Electrical Output - This equipment is for use only by trained, licensed physicians.

Electric Shock Hazard - Connect the generator power cord to a properly grounded receptacle. Do not use power plug adapters.

Connect the power cord to a properly polarized and grounded power source with the frequency and voltage characteristics that match those listed on the back of the unit.

Fire Hazard - Do not use extension cords.

Patient Safety - Use the generator only if the self-test has been completed as described. Otherwise, inaccurate power outputs may result.

The instrument receptacles on this generator are designed to accept only one instrument at a time. Do not attempt to connect more than one instrument at a time into a given receptacle. Doing so will cause simultaneous activation of the instruments.

Failure of the high frequency electrosurgical equipment could result in an unintended increase of output power.

Do not use electrosurgical equipment unless properly trained to use it in the specific procedure being undertaken. Use by physicians without such training has resulted in serious, unintended patient injury, including bowel perforation and unintended, irreversible tissue necrosis.

For surgical procedures where the high frequency current could flow through parts of the body having a relatively small cross-sectional area, the use of bipolar techniques may be desirable to avoid unwanted coagulation.

If the patient has an Implantable Cardioverter Defibrillator (ICD), contact the ICD manufacturer for instructions before performing an electrosurgical procedure. Electrosurgery may cause multiple activation of ICDs.

In some circumstances, potential exists for alternate site burns at points of skin contact (e.g., between the arm and the side of the body). This occurs when electrosurgical current seeks a path to the patient return electrode that includes the skin-to-skin contact point. Current passing through small skin-to-skin contact points is concentrated and may cause a burn. This is true for grounded, ground referenced, and isolated output generators.

To reduce the potential for alternate site burns, do one or more of the following:

- Avoid skin-to-skin contact points, such as fingers touching leg, when positioning the patient.
- Place 5 to 8 cm (2 to 3 in.) of dry gauze between contact points to ensure that contact does not occur.
- Position the return electrode to provide a direct current route between the surgical site and the return electrode which avoids skin-to-skin contact areas.
- In addition, place return electrodes according to the manufacturer's instructions.

Potential for alternate site burns increases if the return electrode is compromised. Bovie Medical Corporation recommends the use of split return electrodes and Bovie® generators with a contact quality monitoring system.

.....
CAUTIONS:

Do not stack equipment on top of the generator or place the generator on top of electrical equipment. These configurations are unstable and/or do not allow adequate cooling.

.....
Provide as much distance as possible between the electrosurgical generator and other electronic equipment (such as monitors). An activated electrosurgical generator may cause interference with them.

.....
Non-function of the generator may cause interruption of surgery. A backup generator should be available for use.

.....
Do not turn the activation tone down to an inaudible level. The activation tone alerts the surgical team when an accessory is active.

.....
When using a smoke evacuator in conjunction with the electrosurgical generator, place the smoke evacuator a distance from the generator and set the generator volume control at a level that ensures that the activation tones can be heard.
.....

NOTICE:

If required by local codes, connect the generator to the hospital equalization connector with an equipotential cable.

SETTING UP

1. Verify that the generator is Off by pressing the power switch Off (O).
2. Place the generator on a stable flat surface, such as a table, platform, or medical cart. Carts with conductive wheels are recommended. For details, refer to the procedures for your institution or to local codes. Provide at least 10 to 15 cm (4 to 6 in.) of space from the sides and top of the generator for cooling. Normally, the top, sides, and rear panel are warm when you use the generator continuously for extended periods of time.
3. Plug the generator power cord into the AC Power Cable Receptacle on the rear panel.
4. Plug the generator power cord into a grounded receptacle.
5. Turn on the generator by pressing the power switch On (I). Verify the following:
 - All visual indicators and displays on the front panel illuminate.
 - Activation tones sound to verify that the speaker is working properly.
6. If the self-test is successful, a tone sounds. Verify the following:
 - A Cut mode is selected; a Coag mode is selected.
 - Each display shows a power setting. The unit automatically powers up to the last selected preset settings.
 - The Patient Return Electrode Alarm Indicator illuminates red.

If the self-test is not successful, an alarm tone sounds. An error code will appear in the Bipolar display, in most cases, the generator is disabled. Note the error code and refer to Section 6, *Troubleshooting*.

Once the self-test is successful, connect the accessories and set the generator controls. Refer to *Preparing for Monopolar Surgery* or *Preparing for Bipolar Surgery* later in this section.

PREPARING FOR MONOPOLAR SURGERY

Monopolar surgery requires a return electrode.

Applying the Return Electrode

To maximize patient safety, Bovie Medical Corporation recommends using a split return electrode and a Bovie® generator with a contact quality monitoring system (Bovie NEM™).

NOTICE:

The Bovie NEM™ system recommends that you use a split return electrode.

Refer to the manufacturer's instructions for application site and placement procedures. When using metal plate return electrodes, use a conductive gel specifically designed for electrosurgery. Select a return electrode site with good blood flow. While a properly applied electrode results in minimal tissue heating beneath the electrode, a good blood flow helps carry heat away from the site.

1. Connect the cable to the Return Electrode receptacle on the front of the unit.
The unit will automatically sense the presence of a split or solid return electrode and, if a split return electrode is used, will constantly monitor the resistance at the contact between the electrode and the patient.
2. Adjust the Blend setting to the desired amount of hemostasis (Level 1 - 10). Adjustment is preformed by pressing the up or down buttons next to the Blend setting indicator.

Select the desired power settings for Cutting. Adjustment is preformed by pressing the up or down buttons next to the Cut display.

Select the mode of operation for Coagulation, either Pinpoint or Spray.

Select the desired power setting for Coagulation. Adjustment is preformed by pressing the up or down buttons next to the Coag display.

Connecting Accessories

1. Connect a 3-pin monopolar device into the monopolar receptacle on the front of the unit.


If footswitching control capabilities are preferred, connect the Bovie® monopolar footswitch to the appropriate footswitch connecting socket on the rear of the unit.

If you are using...	Connect it to...
Standard 3-pin handswitching pencil	Monopolar handswitching receptacle
Footswitching pencil	Monopolar footswitching receptacle

To activate the Monopolar mode, depress the cut or coag button on the monopolar handpiece or the cut or coag pedal on the monopolar footswitch.

Blend Controls

Blend settings can be adjusted to a desired amount of hemostasis (Level 1-10). Ascending illuminated bars indicate increased hemostasis. Increase and decrease the amount of blend added to the Blend mode by pressing the Blend amount control

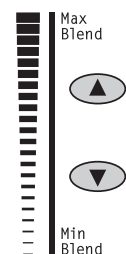

arrowed buttons.

NOTICE:

There are 10 levels of blend available in the Blend Mode.

When selecting the Blend mode, the unit defaults to a setting of minimum blend (only the first bar is illuminated).

PREPARING FOR BIPOLAR SURGERY

1. Connect a Bipolar cable to the Bipolar receptacle on the front of the unit.
2. Connect a forceps instrument to the bipolar cable.
3. Connect the bipolar footswitch to the bipolar footswitch connecting socket located on the rear of the unit.

To activate the Bipolar mode, depress the pedal on the bipolar footswitch.

SETTING AND RECALLING MEMORY PRESETS

The Bovie® IDS-300 incorporates 10 user-defined memory preset settings for easy recall of frequently used settings in all three modes.

Memory

The Memory feature allows the Bovie® IDS-300 (unit) to display the last selected Preset when the generator is turned on. When activated by the handpiece or footswitch, the unit will operate in that particular mode and power setting.

The small red blinking dot in the lower right hand corner of the Preset display lets the user know that the Preset values have been adjusted.

All **new settings must be saved** as a Preset to be available at startup or as a Preset selection (0 through 9) when using the unit.

Memory Function Overview

- The unit powers up with the last selected preset (0-9).
- Mode (Cut and Coag) membrane switches are disabled during activation.
- Blend amount control buttons are disabled during activation.
- Recall and Set membrane switches are disabled during activation.
- During activation, the activated mode can be adjusted up and or down a maximum of four steps. Refer to the following table for power increments.

POWER SETTINGS	INCREMENTS	Example 1
1-50 Watts	1 Watts	<i>The unit is activated using the same preset values as described in Example 2 of this section.</i> While activated, the Cut 1 power output of 30 watts can be adjusted 4 steps down to 26 watts or 4 steps up to 34 watts. The Pinpoint and Bipolar can be adjusted to display a different setting but can not be saved during activation.
50-100 Watts	2 Watts	
100-200 Watts	5 Watts	
200-300 Watts	10 Watts	

- While operating the unit outside of a user-defined preset (small red dot will be blinking in lower right corner of the Preset display as an indicator), the unit temporarily stores the power setting for the activated mode (Cut, Coag, or Bipolar). This temporary power setting is available until either the unit is reset, a preset is selected, or the power setting for the mode in use is adjusted and the unit is again activated.
- Presets only store one Cut mode (Cut I or Cut II, or Blend) and power setting, one Blend level (if applicable), one Coag mode (Pinpoint or Spray) and power setting, and Bipolar power setting. When storing, only the information displayed in the display windows will be saved to the unit's memory.

Setting Your Presets



Select the desired preset (0-9) by pressing the recall button.

Select the desired modes to be stored by pressing the mode membrane switches (Cut and Coag).

If presetting the Blend mode, select the desired level of hemostasis (Blend Bar 1-10) by pressing the Blend amount control button.

Select the desired power (Cut, Coag, and Bipolar) to be stored by using the power output up and down membrane switches.

Once all of the settings are selected, depress and hold the Set button for three seconds. To indicate the settings have been stored, the Preset Memory Number (0-9) will blink.

To recall a Preset, repeatedly press the Recall button to toggle through all of the presets.

NOTICES:

The Bovie® IDS-300 incorporates 10 factory-set presets that are all set to zero and can be reset to your preferred settings.

A small red dot blinking in the lower right corner of the Preset indicator display indicates that the unit is not presently set to a user-defined preset.

Set and Recall buttons are disabled while the unit is activated.

Presets only store one Cut mode (Cut I or Cut II, or Blend) and power setting, one Blend level (if applicable), one Coag mode (Pinpoint or Spray) and power setting, and Bipolar power setting. When storing, only the information displayed in the display windows will be saved to the unit's memory.

Memory Feature (Last Selected Preset)

The Memory feature allows the unit to display the last selected power preset when the generator is turned on.

NOTICE:

To have a setting selection available at startup or to be one of the 10 user-defined presets, the adjustment to the mode and/or power settings must be saved by pressing the Set button on the Preset display panel.

Examples

Examples 2 through 6 explain how the Memory and temporary memory features work and what happens when the power and/or mode is adjusted but not saved as one of the 10 Preset selections. Example 5 explains what happens when the power and/or mode is adjusted and saved as a new Preset setting:

#2. The physician performs a surgical procedure using Preset 2. The Preset has been stored with the following mode and power:

- The mode is set to Cut I
- The power setting for Cut I is 30 watts
- The power setting for Pinpoint is 15 watts
- The power setting for Bipolar is 20 watts.

The procedure is completed and the unit is switched off.

The next time the unit is switched on, the number 2 Preset will be displayed and available when the unit is activated. The number 2 Preset will be the same as the modes and settings indicated above.

#3. The physician performs a surgical procedure using Preset 2 (same as Example #1 Preset values).

He adjusts the power settings for each mode but does not store the new settings into the Preset.

The next time the unit is switched on, the number 2 Preset will be displayed and available when the unit is activated. The number 2 Preset will be the same as the modes and settings indicated in Example #1.

#4. The physician performs a surgical procedure using Preset 2 (same as Example #1 Preset values).

He changes the settings by selecting the Cut II mode. The displayed power will remain at 30 watts. The physician then adjusts the power to 100 watts. He resumes the procedure now using Cut II at 100 watts. He then switches the mode back to Cut I. The power output returns to 30 watts as stored in the # 2 Preset. The physician switches again to the Cut II mode and the output power returns to the temporary memory of 100 watts as previously selected. The procedure is completed without saving any modes or power settings. The next time the unit is switched on, the number 2 Preset will be displayed and available when the unit is activated. The number 2 Preset will be the same as the modes and settings indicated in Example #1.

#5. The physician performs a surgical procedure using Preset 2 (same as Example #1 Preset values).

As required in the procedure, he selects the Blend mode (the Blend Amount Indicator illuminates to one bar indicating the Blend mode can be increased to the preferred amount of blend). He adjusts the hemostasis level up to a 30% blend but does not store the new settings into the Preset. The next time the unit is switched on, the number 2 Preset will be displayed and available when the unit is activated. The number 2 Preset will be the same as the modes and settings indicated in Example #1.

#6. The physician performs a surgical procedure using Preset 2 (same as Example #1 Preset values). He adjusts the power settings for a Cut mode, a Coag mode, and a Bipolar mode and presses the Store button for three seconds to save the new settings as Preset number 2. The next time the unit is switched on, the number 2 Preset will be displayed and available when the unit is activated. *The number 2 Preset will now be the last saved Preset settings for Preset 2.*

ACTIVATING THE UNIT

NOTICE:

Review Activation Safety on page 6 of this section before activating the unit. When you turn on your unit remember the following feature:

The Bovie® IDS-300 will power up to the modes and settings displayed when the unit was last activated. For example, if you set Cut I mode at 50 watts and activate the unit, then turn the unit off, it will automatically return to Cut I mode at 50 watts when you turn it on again. Similarly, if you set Pinpoint mode at 40 watts and activate the unit before you turn it off, it will return to Pinpoint mode at 40 watts when you turn it on again.

1. Monopolar Cut - select the mode of operation for Cut: Cut I, Cut II, or Blend then select the desired Cut power settings by pressing the up and down buttons next to the Cut power output display.
2. If using Blend, vary the Blend setting by pressing the up and down buttons next to the blend amount indicator graph.
3. Monopolar Coag - select the mode of operation for coagulation: Pinpoint or Spray, then select the coagulation power settings by pressing the up and down buttons next to the Coag power output display.
4. Bipolar - adjust the Bipolar power settings by pressing the up and down buttons next to the Bipolar power output display.
5. Activate the generator by pressing the appropriate button on the handpiece or pedal on the footswitch.

NOTICE:

Monopolar and bipolar footswitching operations are controlled by independent foot controls.

ACTIVATION SAFETY

WARNINGS:

Do not wrap the accessory cords or patient return electrode cords around metal objects. This may induce currents that could lead to shocks, fires, or injury to the patient or surgical team.

Danger: Fire / Explosion Hazard - Do not use the Bovie® IDS-300 in the presence of flammable anesthetics.

Fire / Explosion Hazard - The following substances will contribute to increased fire and explosion hazards in the operating room:

- Flammable substances (such as alcohol based skin prepping agents and tinctures)
- Naturally occurring flammable gases that may accumulate in body cavities such as the bowel
- Oxygen enriched atmospheres
- Oxidizing agents (such as nitrous oxide [N₂O] atmospheres).

The sparking and heating associated with electrosurgery can provide an ignition source. Observe fire precautions at all times. When using electrosurgery in the same room with any of these substances or gases, prevent their accumulation or pooling under surgical drapes, or within the area where electrosurgery is performed.

Use the lowest output setting necessary to achieve the desired surgical effect. Use the active electrode only for the minimum time necessary in order to lessen the possibility of unintended burn injury. Pediatric applications and/or procedures performed on small anatomic structures may require reduced power settings. The higher the current flow, and the longer the current is applied, the greater the possibility of unintended thermal damage to tissue, especially during use on small structures.

Use electrosurgery with caution in the presence of internal or external devices such as pacemakers or pulse generators. Interference produced by the use of electrosurgical devices can cause devices such as pacemakers to enter an asynchronous mode or can block the pacemaker effect entirely. Consult the device manufacturer or hospital Cardiology Department for further information when use of electrosurgical appliances is planned for patients with cardiac pacemakers or other implantable devices.

CAUTIONS:

The use of high frequency current can interfere with the function of other electromagnetic equipment.

When high frequency surgical equipment and physiological monitoring equipment are used simultaneously on the same patient, place any monitoring electrodes as far as possible from the surgical electrodes.

Do not use needles as monitoring electrodes during electrosurgical procedures. Inadvertent electrosurgical burns may result.

To avoid the possibility of an electrosurgical burn to either the patient or the physicians, do not allow the patient to come in contact with a grounded metal object during activation. When activating the unit, do not allow direct skin contact between the patient and the physician.

Remove any jewelry from the patient before activation.

Studies have shown that smoke generated during electrosurgical procedures can be potentially harmful to patients and the surgical team. These studies recommend adequately ventilating the smoke by using a surgical smoke evacuator or other means.¹

Examine all accessories and connections to the electrosurgical generator before use. Ensure that the accessories function as intended. Improper connection may result in arcs, sparks, accessory malfunction, or unintended surgical effects.

When not using active accessories, place them in a holster or in a clean, dry, non-conductive, and highly visible area not in contact with the patient. Inadvertent contact with the patient may result in burns.

1. U.S. Department of Health and Human Services. National Institute for Occupational Safety and Health (NIOSH). Control of Smoke from Laser / Electric Surgical Procedures. HAZARD CONTROLS, Publication No. 96-128, September, 1996.

MAINTAINING THE BOVIE® IDS-300

This section covers the following topics:

- *Cleaning*
- *Periodic Inspection*
- *Fuse Replacement*

Bovie Medical Corporation recommends that you complete periodic inspection and performance testing. Perform inspections and performance testing every six months. A qualified biomedical technician should conduct this testing to ensure that the unit is operating effectively and safely.

CLEANING

After each use, clean the unit.

WARNING:

Electric Shock Hazard - Always turn off and unplug the generator before cleaning.

NOTICE:

Do not clean the generator with abrasive cleaning or disinfectant compounds, solvents, or other materials that could scratch the panels or damage the generator.

1. Turn off the generator, and unplug the power cord from the wall outlet.
2. Thoroughly wipe all surfaces of the generator and power cord with a mild cleaning solution or disinfectant and a damp cloth. Follow the procedures approved by your institution or use a validated infection control procedure. Do not allow fluids to enter the chassis. Do not sterilize the generator.

PERIODIC INSPECTION

Every six months, visually inspect the Bovie® IDS-300 for signs of wear or damage.

In particular, look for any of the following problems:

- Damage to the power cord
- Damage to the power cable receptacle
- Obvious damage to the unit
- Damage to any receptacle
- Accumulation of lint or debris in or around the unit

FUSE REPLACEMENT

Fuses for the unit reside directly below the Power Cable Receptacle on the rear of the unit.

To replace the fuses, follow this procedure:

1. Unplug the power cord from the wall outlet.
2. Remove the power cord from the Power Cable Receptacle on the rear panel.
3. To release the fuse drawer, insert a small flathead screwdriver into the slot on the drawer below the power cord receptacle. Then, slide the drawer out.
4. Remove the two fuses and replace them with new fuses with the same values.
5. Insert the fuse holder into the Power Cable Receptacle.

NOTICE:

If the unit does not display an error and does not power on, check fuses.

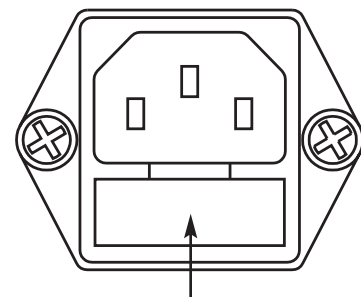


Figure 5 – 1 Fuse holder

TROUBLESHOOTING

This section includes Error Code Descriptions and actions to take to resolve them.

The Bovie® IDS-300 includes automatic self-diagnostics. If the diagnostics detect an error, the system displays an error code, sounds an audible tone, and deactivates the unit output power.

Most error codes result from faults in accessories attached to the unit. The following table lists the error codes, describes the errors, and recommends actions to take to resolve the errors.

All error codes are displayed in the Bipolar display. If the unit displays any other error code, it requires service.

Error Code	Description	Recommended Action
F1	Cut handpiece button may be stuck	<ol style="list-style-type: none"> 1. Turn off, then turn on the generator. Do not press buttons or activate footpedals during the self-test. 2. If the error code reappears, disconnect all accessories. Turn off, then turn on the generator again. 3. If the problem persists, disconnect the handpiece or footswitch and repeat the restart. 4. If the error code reappears, record the number and call Bovie Medical Corporation customer service.
F2	Coag handpiece button may be stuck	
F3	Cut footswitch pedal may be stuck	
F4	Coag footswitch pedal may be stuck	
F5	Bipolar footswitch pedal may be stuck	
F6	Simultaneous activation error	<p>The unit does not allow simultaneous activation of the cut and coagulation modes. The activation mode is “first come, first serve.” This means that whichever mode is selected first will be the function the unit is activated to dispense. An example of this functionality includes, when the handpiece Cut button is pressed, the unit is activated for Cut. If a footswitch is simultaneously pressed for Coag, the unit will continue in the Cut mode as long as the handpiece Cut button is pressed. If the Cut button is released, the unit will sense an error and both functions will be disabled.</p> <ol style="list-style-type: none"> 1. Release either the cut or coag button on the handpiece, or the cut or coag pedal on the footswitch. 2. If the error code reappears, record the number and contact Bovie Medical Corporation customer service.
E1	Output current out of specification	<ol style="list-style-type: none"> 1. Turn the unit off. 2. Turn the unit on. 3. If the error code reappears, record the number and contact Bovie Medical Corporation customer service.
E2	Output current sensors delta error	
E3	Output voltage sensors delta error	
E4	System power supply voltages error	
E5	Internal temperature of a section of the unit exceeded the limit.	<ol style="list-style-type: none"> 1. Turn the unit off. 2. Allow the unit to cool for 20 minutes. 3. Turn the unit on. 4. If the error code reappears, record the number and contact Bovie Medical Corporation customer service.
E6		
E7		
E8	NEM circuit error	<ol style="list-style-type: none"> 1. Turn the unit off. 2. Turn the unit on. 3. If the error code reappears, record the number and contact Bovie Medical Corporation customer service.

NOTICE:

If the unit does not power on to display an error, check fuses as described in Section 5 of this guide.

REPAIR POLICY AND PROCEDURES

Refer to this section for information on:

- Responsibility of the Manufacturer*
- Returning the Generator for Service*

RESPONSIBILITY OF THE MANUFACTURER

Bovie® is responsible for the safety, reliability, and performance of the generator only under the following circumstances:

- The user has followed the Installation and Setup Procedures in this User's Guide.
- Persons authorized by Bovie Medical Corporation performed assembly operation, readjustments, modifications, or repairs.
- The electrical installation of the relevant room complies with local codes and regulatory requirements, such as IEC and BSI.
- Equipment use is in accordance with the Bovie Medical Corporation instructions for use.

Please note that infected medical devices must be disposed of as medical/biohazard waste and cannot be included in used electronic equipment disposal/recycling programs. In addition, certain electronic products must be returned directly to Bovie Medical Corporation. Contact your Bovie Medical Corporation representative for return instructions.

For warranty information, refer to Appendix B - *Warranty*.

RETURNING THE GENERATOR FOR SERVICE

Before you return the generator, call your Bovie Medical Corporation representative for assistance. If instructed to send the generator to Bovie Medical Corporation, first obtain a Returned Goods Authorization Number. Then, clean the Generator and package securely to ensure proper protection of the unit. So as to aid in the processing of the unit, please be sure to include a reference to the Bovie® Return Goods Authorization Number on the outside of the box and ship directly to Bovie Medical Corporation.

Step 1 – Obtain a Returned Goods Authorization Number

Call the Bovie Medical Corporation Customer Service Center to obtain a Returned Goods Authorization Number. Have the following information ready when you call:

- Hospital / clinic name / customer number
- Telephone number/fax number
- Department / address, city, state, and zip code
- Model number
- Description of the problem
- Type of repair to be done
- P.O. number

Step 2 – Clean the Generator

WARNING:

Electric Shock Hazard - Always turn off and unplug the generator before cleaning.

NOTICE:

Do not clean the generator with abrasive cleaning or disinfectant compounds, solvents, or other materials that could scratch the panels or damage the generator.

- A. Turn off the generator, and unplug the power cord from the wall outlet.
- B. Thoroughly wipe all surfaces of the generator and power cord with a mild cleaning solution or disinfectant and a damp cloth. Follow the procedures approved by your institution or use a validated infection control procedure. Do not allow fluids to enter the chassis. You cannot sterilize the generator.

Step 3 – Ship the Generator

- A. Attach a tag to the generator that includes the Returned Goods Authorization Number and the information (hospital, phone number, etc.) listed in *Step 1 – Obtain a Returned Goods Authorization Number*.
- B. Be sure the generator is completely dry before you pack it for shipment. Although the preference is to have the Generator repackaged using its original packaging, Bovie understands that this may not always be possible. If necessary, contact Customer Service for the proper packaging to ship the unit. Please be sure to include a reference of the Bovie Return Goods Authorization Number on the outside of the box/container.
- C. Ship the generator, prepaid, to the address given to you by the Bovie Medical Corporation Service Center.

TECHNICAL SPECIFICATIONS

All specifications are nominal and subject to change without notice. A specification referred to as “typical” is within $\pm 20\%$ of a stated value at room temperature (25° C / 77° F) and a nominal input power voltage.

PERFORMANCE CHARACTERISTICS

Input Power

Input Voltage	100-240 ~ VAC \pm 10%
Mains line frequency range (nominal):	50 – 60 Hz
Power consumption:	560 VA
Fuses (two):	6.3 A (slow blow)

Duty Cycle

Under maximum power settings and rated load conditions (Pure Cut, 300 watt @ 300 ohm load), the generator is suitable for activation times of 10 seconds ON followed by 30 seconds OFF for one hour.

The internal temperature of the unit is continuously monitored. If the temperature rises above 85⁰ C, the alarm will sound and output power will be deactivated.

Dimensions and Weight

Width	31.1 cm (12.25 in.)	Depth	41.3 cm (16.25 in.)
Height	15.3 cm (6.00 in.)	Weight	< 8.75 kg (< 19 lbs)

Operating Parameters

Ambient temperature range	10° to 40° C
Relative humidity	30% to 75%, non-condensing
Atmospheric pressure	70kPa to 106kPa
Warm-up time	If transported or stored at temperatures outside the operating temperature range, allow one hour for the generator to reach room temperature before use.

Transport and Storage

Generator should fit on all standard Carts for monopolar generators. The device should be stored and used in a room temperature of approximately 77⁰ F/25⁰ C.

Ambient temperature range	-40° to +70° C
Relative humidity	10% to 100%, including condensation
Atmospheric pressure	50kPa to 106kPa

Audio Volume

The audio levels stated below are for activation tones (cut, coag, and bipolar) and alarm tones (return electrode and system alarms) at a distance of one meter. Alarm tones meet the requirements for IEC 60601-2-2.

Activation Tone

Volume (adjustable)	40 to 65 dB
Frequency	Cut I: 610 Hz Cut II: 610 Hz Blend: 610 Hz Pinpoint: 910 Hz Spray: 910 Hz Bipolar: 910 Hz
Duration	Continuous while the generator is activated

Alarm Tone

Volume (not adjustable)	70 dB \pm 5 dB
Frequency	2 kHz $\frac{1}{2}$ seconds / 1 kHz $\frac{1}{2}$ seconds
Duration	2 seconds

Return Electrode Sensing

The system presents audible and visible alarms when it senses no return electrode.

Solid	Trip resistance: 0 Ω to 5 Ω \pm 3 Ω Continuous measurement: Once the system establishes the solid return electrode resistance, an increase to 20 Ω \pm 5 Ω in resistance will cause an alarm. When the alarm condition exists, the system deactivates output power.
Split	Trip resistance: 10 Ω \pm 5 Ω to 135 Ω \pm 10 Ω Continuous measurement: Once the system establishes the split return electrode resistance, an increase of 40% in resistance will cause an alarm. When the alarm condition exists, the system deactivates output power.

Low Frequency (50-60 Hz) Leakage Current

Enclosure source current, ground open	< 500 μ A
Source current, patient leads, all outputs	Normal polarity, intact ground: < 10 μ A Normal polarity, ground open: < 10 μ A Reverse polarity, ground open: < 10 μ A
Sink current at high line, all inputs	< 10 μ A

High Frequency (RF) Leakage Current

Bipolar RF leakage current	< 63 mA _{rms} at 80 watts
Monopolar RF leakage current (additional tolerance)	< 150 mA _{rms}

STANDARDS AND IEC CLASSIFICATIONS

Class I Equipment (IEC 60601-1)

Accessible conductive parts cannot become live in the event of a basic insulation failure because of the way in which they are connected to the protective earth conductor.

Type CF Equipment (IEC 60601-1) / Defibrillator Proof



The Bovie® IDS-300 provides a high degree of protection against electric shock, particularly regarding allowable leakage currents. It is type CF equipment. Patient connections are isolated from earth and resist the effects of defibrillator discharge.

Drip Proof (IEC 60601-2-2)

The generator enclosure is constructed so that liquid spillage in normal use does not wet electrical insulation or other components which, when wet, are likely to affect adversely the safety of the generator.

Electromagnetic Interference

When other equipment is placed on or beneath a Bovie® IDS-300, the unit can be activated without interference. The generator minimizes electromagnetic interference to video equipment used in the operating room.

Electromagnetic Compatibility (IEC 60601-1-2 and IEC 60601-2-2)

The Bovie® IDS-300 complies with the appropriate IEC 60601-1-2 and IEC 60601-2-2 specifications regarding electromagnetic compatibility.

Voltage Transients (Emergency Generator Mains Transfer)

The Bovie® IDS-300 operates in a safe manner when the transfer is made between line AC and an emergency generator voltage source.

OUTPUT CHARACTERISTICS

Maximum Output for Monopolar and Bipolar Modes

Power readouts agree with actual power into rated load to within 20% or 5 watts, whichever is greater.

Mode	Output Power	Output Frequency	Repetition Rate	Vpeak max	Crest Factor* (Rated Load)
Cut I	300 W @ 300 Ω	490 kHz \pm 4.9 kHz	N / A	1000V	1.6 \pm 20%
Cut II	300 W @ 300 Ω	490 kHz \pm 4.9 kHz	N / A	1000V	1.6 \pm 20%
Blend (Max)	200 W @ 300 Ω	490 kHz \pm 4.9 kHz	30 kHz \pm 5 kHz	2000V	3.5 \pm 20%
Pinpoint	120 W @ 500 Ω	490 kHz \pm 4.9 kHz	30 kHz \pm 5 kHz	2400V	4.5 \pm 20%
Spray	80 W @ 500 Ω	490 kHz \pm 4.9 kHz	30 kHz \pm 5 kHz	4000V	6.5 \pm 20%
Bipolar	80 W @ 150 Ω	490 kHz \pm 4.9 kHz	N/A	450V	1.6 \pm 20%

- an indication of a waveform's ability to coagulate bleeders without a cutting effect.

EMC COMPLIANCE

Special precautions should be taken regarding the Bovie® IDS-300. Medical Electrical Equipment needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in this manual.


Understand that only the Accessories supplied with or ordered from Bovie® should be used with your device. The use of accessories, transducers, and cables other than those specified, may result in increased Emissions or decreased Immunity of the IDS-300. The Bovie® IDS-300 and its accessories are not suitable for interconnection with other equipment.

Portable and mobile RF communications equipment can affect Medical Electrical Equipment. The Bovie® IDS-300 should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the IDS-300 should be observed to verify normal operation in the configuration in which it will be used.

Recommended separation distances between portable and mobile RF communications equipment and the IDS-300			
The IDS-300 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the IDS-300 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the IDS-300 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	separation distance according to frequency of transmitter		
	150 kHz to 80 MHz $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$	80 MHz to 800 MHz $d = \left[\frac{3.5}{E_1} \right] \sqrt{P}$	800 MHz to 2.5 GHz $d = \left[\frac{7}{E_1} \right] \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Guidance and manufacturer' s declaration – electromagnetic emissions		
The IDS-300 is intended for use in the electromagnetic environment listed below. The customer or the user of the IDS-300 should assure that is is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment - guidance
RF Emissions CISPR 11	Group 2	The IDS-300 must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
RF Emissions CISPR 11	Class A	The IDS-300 is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used in domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

Guidance and manufacturer' s declaration – electromagnetic immunity			
The IDS–300 is intended for use in the electromagnetic environment listed below. The customer or the user of the IDS–300 should assure that is is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000–4–2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000–4–4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000–4–5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000–4–11	<5 % U_t (<95 % dip in U_t) for 0.5 cycle <40 % U_t (<60 % dip in U_t) for 5 cycles 70 % U_t (<30 % dip in U_t) for 25 cycles <5 % U_t (>95 % dip in U_t) for 5 sec	<5 % U_t (<95 % dip in U_t) for 0.5 cycle <40 % U_t (<60 % dip in U_t) for 5 cycles 70 % U_t (<30 % dip in U_t) for 25 cycles <5 % U_t (>95 % dip in U_t) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the IDS–300 requires continued operation during power mains interruptions, it is recommended that the IDS–300 be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000–4–8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE U_t is the a.c. mains voltage prior to application of the test level.			

Guidance and manufacturer' s declaration – electromagnetic immunity continued...			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000–4–6	3 Vrms 150 kHz to 80 MHz	3 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the IDS–300, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$
Radiated RF IEC 61000–4–3	3 V/m 80 MHz to 2.5 GHz	3 V/m	$d = \left[\frac{3.5}{E_1} \right] \sqrt{P}$ <p>80 MHz to 800 MHz</p> $d = \left[\frac{7}{E_1} \right] \sqrt{P}$ <p>800 MHz to 2.5 GHz</p> <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m)</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol. </p>
<p>NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.</p> <p>NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p>^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicated theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location which the IDS–300 is used exceeds the applicable RF compliance level above, the IDS–300 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the IDS–300.</p> <p>^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than $[V_1]$ V/m.</p>			

OUTPUT POWER CURVES

Figure A-1 illustrates output power delivered to rated load for all available modes. Figure A-2 illustrates the maximum peak voltage available at a given power setting and output mode. Figures A-3 through A-9 illustrate specific output power delivered to a range of load resistances for each mode.

Figure A – 1 Display setting vs. output power in watts at rated load

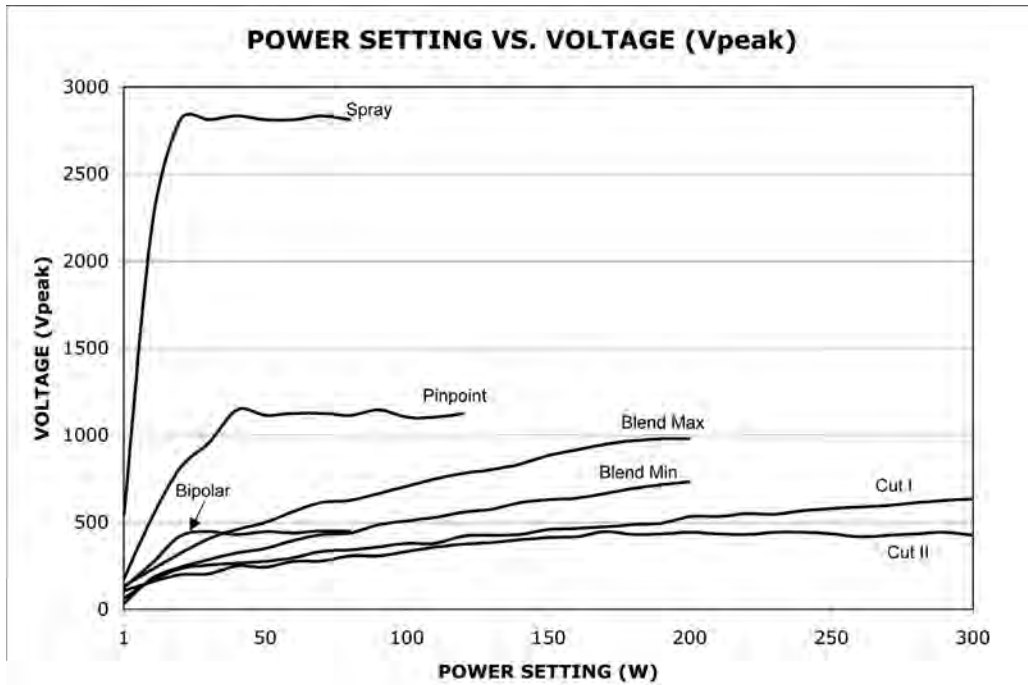


Figure A – 2 Output voltage (Vpeak) versus power setting at rated load

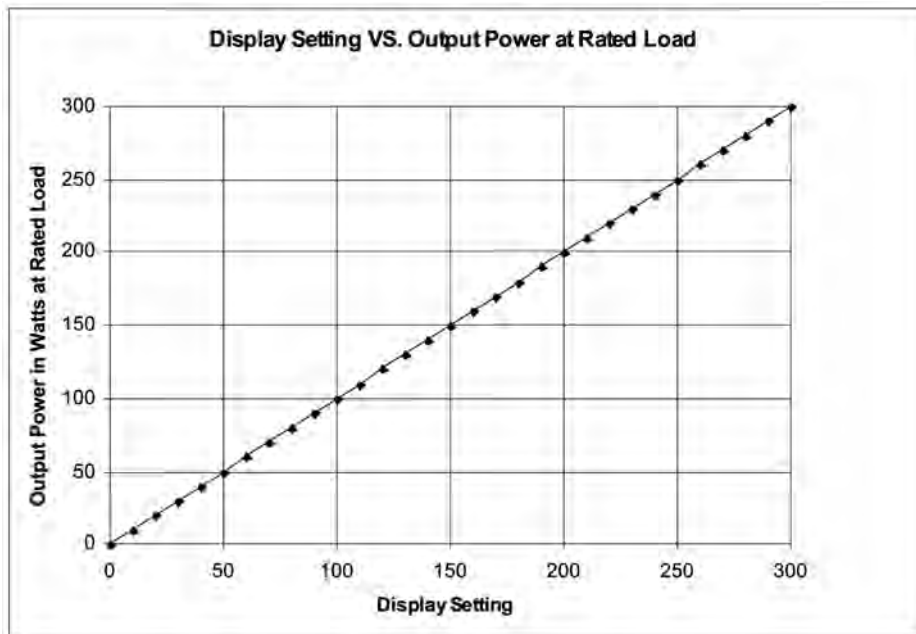


Figure A – 3 Output power versus impedance for Cut I mode

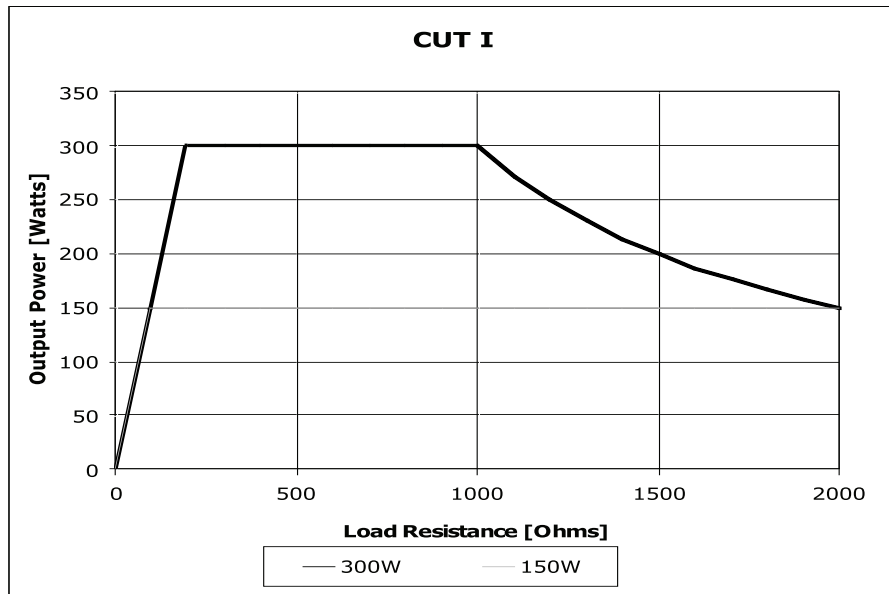


Figure A – 4 Output power versus impedance for Cut II mode

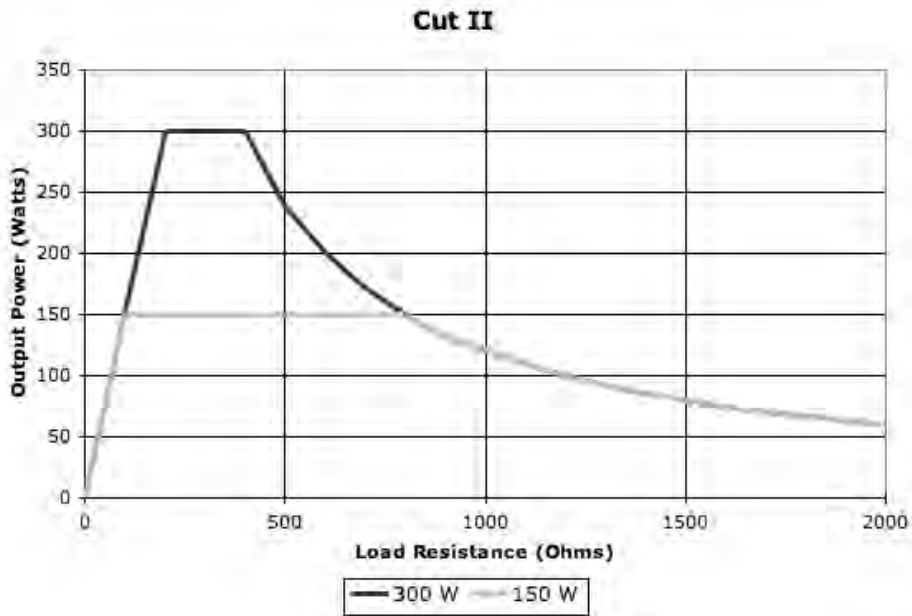


Figure A – 5 Output power versus impedance for Blend Min mode



Figure A – 6 Output power vs impedance for Blend Max mode

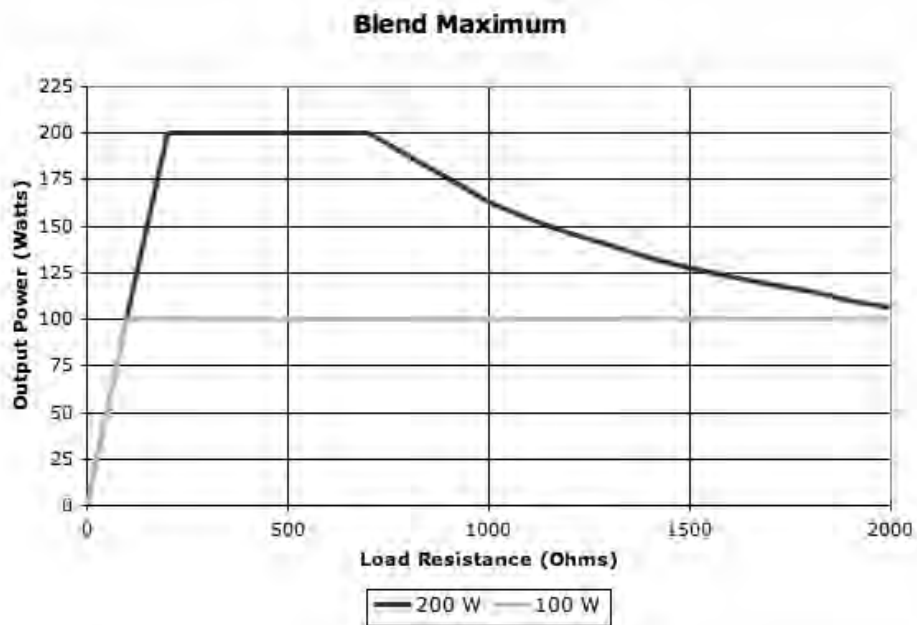


Figure A- 7 Output power vs impedance for Pinpoint mode

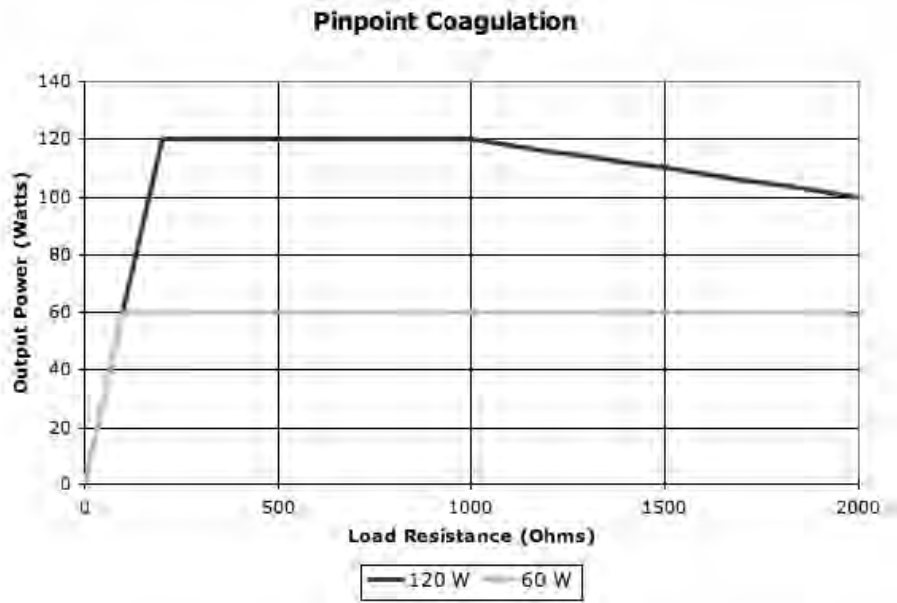


Figure A - 8 Output power vs impedance for Spray mode


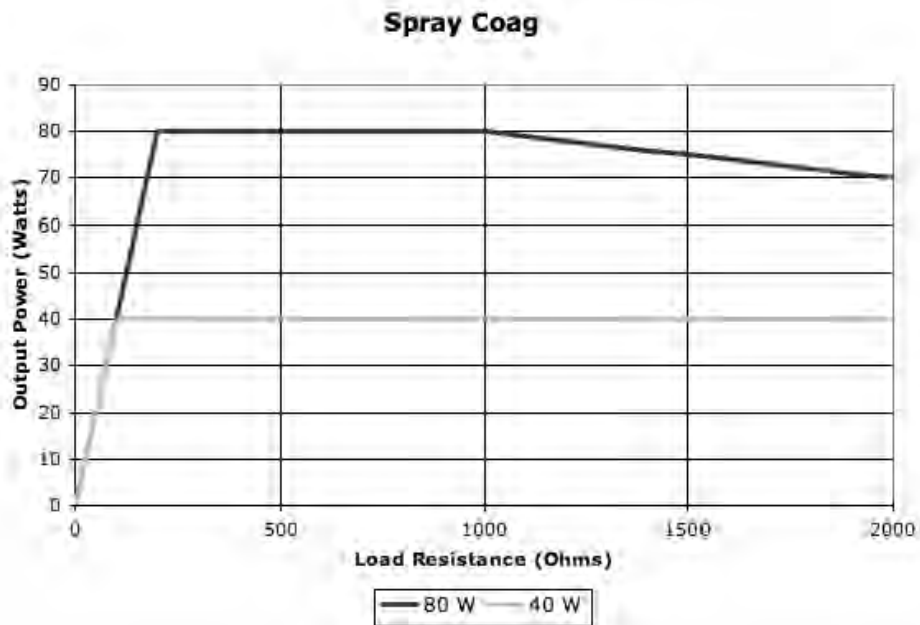


Figure A-9 Output power vs impedance for Bipolar mode

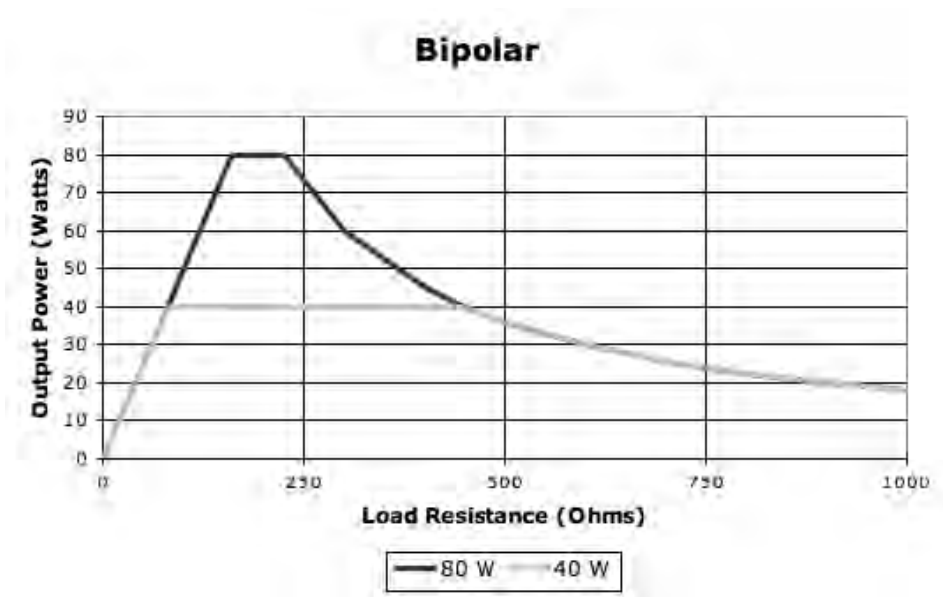

WARRANTY

Bovie Medical Corporation, warrants each product manufactured by it to be free from defects in material and workmanship under normal use and service for the period(s) set forth below.

Bovie Medical Corporation's obligation under this warranty is limited to the repair or replacement, at its sole option, of any product, or part thereof, which has been returned to it or its Distributor within the applicable time period shown below after delivery of the product to the original purchaser, and which examination discloses, to Bovie Medical Corporation's satisfaction, that the product is indeed, defective.

This warranty does not apply to any product, or part thereof, which has been repaired or altered outside Bovie Medical Corporation's factory in a way so as, in Bovie Medical Corporation's judgment, to affect its stability or reliability, or which has been subjected to misuse, neglect, or accident.

The warranty periods for Bovie Medical Corporation products are as follows:

- Electrosurgical Generators: Two years from date of shipment
- Mounting Fixtures (all models): Two years from date of shipment
- Footswitches (all models): Ninety days from date of shipment
- Patient Return Electrodes: Shelf life only as stated on packaging
- Sterile Single Use Accessories: Only as stated on packaging
- Handpiece: Only as stated on packaging

This warranty is in lieu of all other warranties, express or implied, including without limitation, the warranties of merchantability and fitness for a particular purpose, and of all other obligations or liabilities on the part of Bovie Medical Corporation.

Bovie Medical Corporation neither assumes nor authorizes any other person to assume for it any other liability in connection with the sale or use of any of Bovie Medical Corporation's products.

Notwithstanding any other provision herein or in any other document or communication, Bovie Medical Corporation's liability with respect to this agreement and products sold hereunder shall be limited to the aggregate purchase price for the goods sold by Bovie Medical Corporation to the customer.

Bovie Medical Corporation disclaims any liability hereunder or elsewhere in connection with the sale of this product, for indirect or consequential damages.

This warranty and the rights and obligations hereunder shall be construed under and governed by the laws of the State of Florida, USA.

The sole forum for resolving disputes arising under or relating in any way to this warranty is the District Court of the County of Pinellas, State of Florida, USA.

Bovie Medical Corporation, its dealers, and representatives reserve the right to make changes in equipment built and/or sold by them at any time without incurring any obligation to make the same or similar changes on equipment previously built and/or sold by them.

Bovie Medical Corporation



**5115 Ulmerton Road
Clearwater, FL 33760**

U.S. Phone 1-800-537-2790 • Fax 1-800-323-1640

Int'l. Phone +1-727-384-2323 • Fax +1-727-347-9144

www.boviemed.com • sales@boviemed.com