


ELx405™ Select Deep Well Washer

BioTek's deep well microplate washer aspirates and dispenses into 96- and 384-well plates ≤50 mm in height. ANSI/SBS standard microplates are also accommodated, resulting in a total solution for all of a laboratory's solid bottom microplate washing needs.

The ELx405™ Select Deep Well Microplate Washer is a robot compatible, full plate washer incorporating BioTek's patented Dual-Action™ manifold with independent filling and evacuation control for precise overfill washing and overflow protection.

Available low-flow rates and angled dispensing make the ELx405 particularly useful in cell-based assays. With optional patented Ultrasonic Advantage™, BioTek washers are the only products available today with built in ultrasonic cleaners able to automatically prime with cleaning solution, sonicate tubes to remove any protein or salt crystal build up and re-prime with rinse or wash buffer leaving the instrument clean and ready for its next wash.

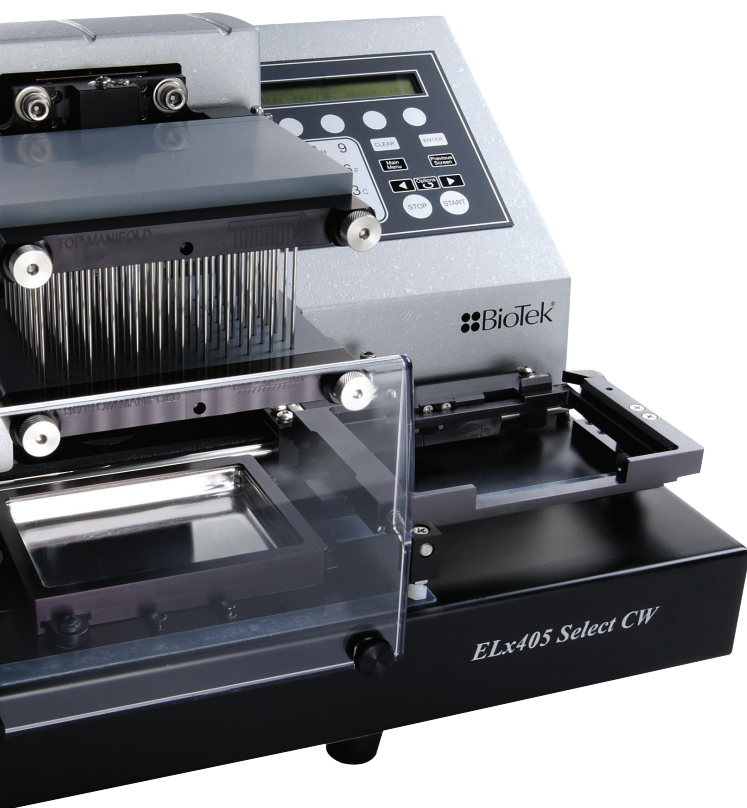

Features:

- Washes microplates ≤50 mm in height
- 96- and 384-well formats
- Optional cell washing with gentle, low-flow rates and 20° angled dispense tubes
- Dual-Action™ manifold for independent dispense and aspiration control
- Available automatic 1-to-4 wash buffer switching
- Optional patented Ultrasonic Advantage™ – built in ultrasonic cleaner for manifold self-maintenance
- Biomagnetic separation ready
- Software control from onboard keypad or Liquid Handling Control™ PC Software
- 21 CFR Part 11 compliance with LHC™ Secure Software
- BioStack™ Microplate Stacker compatibility for standard height plates


Typical Applications:

Deep Well Plates

- Plasmid DNA purification
- Serum/plasma sample preparation
- Cell signaling – phospho flow setup for flow cytometry
- Hematopoietic progenitor cell enumeration
- Cell culture – bacterial, yeast and mammalian cells grown in suspension

Standard Height Plates

- ELISAs
- MSD assays
- HCS immune cytochemistry
- Cell-based assays
- Magnetic bead assays – such as gene expression and cytokine panels
- ELISPOT assays

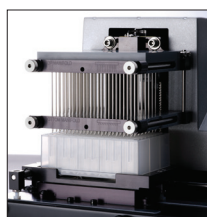

The ELx405™ Select Deep Well Washer aspirates and dispenses into full size 96- and 384-well blocks with unattended operation.

Configurations:

| Configurations | Part # | 96- /384- well | Deep well | Cell wash | Wash buffer switch | Ultrasonic Advan- tage | Biomag sep |
|-------------------|-------------|----------------------|--------------|--------------|--------------------------|------------------------------|---------------|
| ELx405™ Select | ELX405UD | • | • | | | | • |
| | ELX405USD | • | • | | | • | • |
| | ELX405UVD | • | • | | • | | • |
| | ELX405UVSD | • | • | | • | • | • |
| | ELX405UCWD | • | • | • | | | • |
| | ELX405UCWSD | • | • | • | | • | • |
| | ELX405UCWVD | • | • | • | • | | • |
| | ELX405UCWSD | • | • | • | • | • | • |

Optional Accessories:

- Dispense/Waste Systems – choice of 4L or 10L bottles and standard or high flow* vacuum pumps
- Direct Drain Waste System
- Product Qualification Package
- BioStack™ Microplate Stacker
- Liquid Handling Control™ PC Software
- Magnets – choice of 96- or 384-well formats

Technical Details:

| | |
|---------------------|---|
| Assays: | ELISA Cell-based assays (<i>configuration dependent</i>) Magnetic bead <ul style="list-style-type: none"> • Multiplex assays • Bead-based ELISA |
| Microplate types: | 96- and 384-well Low profile and standard height Deep well up to 50 mm |
| Magnet: | Choice of high strength, flat or ring designs 96- and 384-well formats |
| Onboard software: | 2x24 character backlit LCD display |
| Software: | Optional LHC™ PC Software for wash protocol programming and execution |
| Manifold type: | 96-/384-well washing - Dual-Action™ 96-tube (8x12) |
| Washing speed: | 96-well (96-tube manifold) - <30 sec: 3 asp./disp. cycles, 300 µL/well |
| Fluid delivery: | Internal positive displacement pump |
| Volume range: | 50 – 3,000 µL/well |
| Buffer selection: | Optional switching of up to 4 wash buffers |
| Wash cycles: | 1 – 10 |
| Dispense precision: | ≤3% CV |
| Residual volume: | ≤2 µL/well |
| Shaking: | User-programmable speeds and timing |
| Soak time: | 1 – 600 seconds |
| Power: | 100 – 240 Volts AC. 50/60 Hz |
| Dimensions: | Depth: 17" (43.2 cm) Width: 14" (35.6 cm) Height: 11" (27.9 cm) |
| Weight: | 30 lbs (13.5 kg) |

Regulatory

In Vitro Diagnostic use configurations available. CE and TUV marked, RoHS compliant.

*High flow vacuum pump recommended for 384-well washing with buffers not containing surfactant.



The ELx405 is Luminex® xMAP® approved.
xMAP® is a registered trademark of Luminex Corporation.



BioTek Instruments, Inc.
Highland Park, P.O. Box 998
Winooski, Vermont 05404-0998, USA

Phone: 802-655-4040 • Toll-Free: 888-451-5171
Outside the USA: 802-655-4740
www.biotek.com

Technical details are subject to change.